



LINHA ESTÉRIL
NÃO REUTILIZÁVEL

INSTRUÇÃO
DE USO

Instrumental não articulado cortante

CADASTRO ANVISA nº 80482490003

*NOTA: Visando facilidade de acesso as informações, conforme IN nº 4 / 2012 estabelecida pela ANVISA, a **BIOMETAL** disponibiliza todas as INSTRUÇÕES DE USO de seus produtos através do endereço eletrônico <http://www.biometal.med.br>; e também através de consulta pública no site da ANVISA = <https://consultas.anvisa.gov.br/#/saude/25351404414201687/?numeroRegistro=80482490025>.*

É importante verificar a correlação adequada da revisão da INSTRUÇÃO DE USO, conforme orientação contidas no impresso “ALERTA DE INSTRUÇÃO DE USO” – que segue junto ao produto para saúde adquirido, e está replicado na última deste documento.

Produtos Fabricados e Distribuídos por:



BIOMETAL INDÚSTRIA E COMÉRCIO EIRELI

CNPJ: 07.165.928/0001-75

Endereço: V LUIZ VARGA 2010 - GLEBA BEATRIZ - LIMEIRA/SP

CEP: 13.486-606

Fones: (19) 3495-2485 / (19) 3495 3296

RESPONSÁVEL TÉCNICO: RENATA VICELLI – CREA/SP 5063105420

CARACTERÍSTICAS E ESPECIFICAÇÕES TÉCNICAS


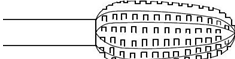
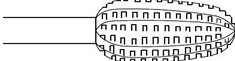
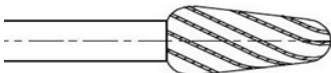
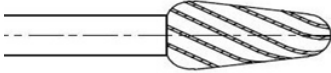
NOME TÉCNICO: Brocas Cirúrgicas

CLASSE DE RISCO ANVISA: II – médio baixo

NOME COMERCIAL do CADASTRO ANVISA: Brocas Cortantes, Fresas e Lâminas Estéreis

MODELOS COMERCIAIS NESTE CADASTRO:



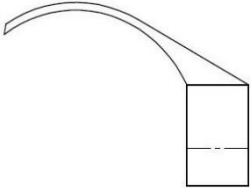

TABELA 1 – MODELOS COMERCIAIS

IMAGEM ILUSTRATIVA	MODELO COMERCIAL	MEDIDAS		ENGATES
		Ø (mm)	Comprimentos (mm)	
	BROCA CORTANTE	Ø 0,60 / 1,00 / 1,40 / 1,60 / 1,8 / 2,30 / 2,70 / 3,10 / 3,50 / 4,00 / 4,50 / 5,00 / 5,50 / 6,0 / 7,0 mm	50 / 55 / 60 / 65 / 70 / 75 / 80 / 85 / 90 / 95 / 100 / 105 / 110 / 115 / 120 / 125 / 130 / 135 /140 / 145 / 150 mm	C6 / C8 / C10 / C12 / C14 / C16
	BROCA CORTANTE CÍLINDRICA CROSS-CUT	Ø 4,0 Ø 5,0 Ø 6,0 Ø 7,0 Ø 8,0	50 / 55 / 60 / 65 / 70 / 75 / 80 / 85 / 90 / 95 / 100 / 105 / 110 / 115 / 120 / 125 / 130 / 135 /140 / 145 / 150 mm	C10 / C12 / C14 / C16
	BROCA CORTANTE CÍLINDRICA	Ø 4,0 Ø 5,0 Ø 6,0 Ø 7,0 Ø 8,0	50 / 55 / 60 / 65 / 70 / 75 / 80 / 85 / 90 / 95 / 100 / 105 / 110 / 115 / 120 / 125 / 130 / 135 /140 / 145 / 150 mm	C10 / C12 / C14 / C16
	BROCA CORTANTE BUD CROSS-CUT	Ø 4,0 Ø 5,0 Ø 6,0 Ø 7,0 Ø 8,0	50 / 55 / 60 / 65 / 70 / 75 / 80 / 85 / 90 / 95 / 100 / 105 / 110 / 115 / 120 / 125 / 130 / 135 /140 / 145 / 150 mm	C10 / C12 / C14 / C16
	BROCA CORTANTE BUD	Ø 4,0 Ø 5,0 Ø 6,0 Ø 7,0 Ø 8,0	50 / 55 / 60 / 65 / 70 / 75 / 80 / 85 / 90 / 95 / 100 / 105 / 110 / 115 / 120 / 125 / 130 / 135 /140 / 145 / 150 mm	C10 / C12 / C14 / C16

INSTRUÇÃO DE USO

Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003



	BROCA AUTOMÁTICA PARA TREPANAÇÃO	\varnothing 14.00/11.00 \varnothing 11.00/8.00 \varnothing 9.00/6.00	único	HUDSON SMITH
	FRESA CRANIÓTOMO	\varnothing 2,00 mm \varnothing 2,30 mm \varnothing 2,50 mm \varnothing 2,70 mm \varnothing 3,00 mm \varnothing 3,20 mm \varnothing 3,70 mm	único	-----
	LÂMINA REMOVEDORA DE ACETABULO	\varnothing 40,0 – 80,0	Inicial e Final	-----
	BROCA BUCO MAXILO	\varnothing 0,9 / 1,0 / 1,1 / 1,2 / 1,3 / 1,4 / 1,5 / 1,6 / 1,7 / 1,8 / 1,9 / 2,0 / 2,1 /	50 / 55 / 60 / 65 / 70 / 75 / 80 / 85 / 90 / 95 / 100 / 105 / 110 / 115 / 120 / 125 / 130 / 135 / 140 / 145 / 150 mm	Com stop de 4 / 5 / 6 / 7 / 8 / 9 / 10 / 11 / 12 / 13 / 14 / 15 / 16 / 17 / 18 / 19 / 20 / 21 / 22 / 23 / 24 / 25 / 30 / 35 / mm

PRODUTO ESTÉRIL

MÉTODO DE ESTERILIZAÇÃO:

- ✓ Óxido de Etileno, ou
- ✓ Radiação GAMA

**ver indicação do método aplicado na rotulagem do produto adquirido.*

INSTRUÇÃO DE USO

Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003















VALIDADE: 05 anos - após a data da esterilização.

PRODUTO DE USO ÚNICO

REPROCESSAMENTO PROIBIDO

TABELA 2 – LEGENDA DOS SÍMBOLOS UTILIZADOS NAS EMBALAGENS

 REF	Código do produto		Não utilizar se a embalagem estiver danificada
 LOT	Código de lote		Manter ao abrigo do sol
	Data de Validade		Manter seco
	Data de Fabricação		Frágil, manusear com cuidado
	Consultar instruções de uso		Produto Estéril – Esterilizado utilizando óxido de etileno
	Produto de uso único		Esterilizado utilizando irradiação

APRESENTAÇÃO

Todos os materiais que integram o CADASTRO ANVISA nº **80482490003** são instrumentos de corte sem articulação, cuja finalidade é auxiliar na perfuração, corte ou fresagem dos tecidos ósseos durante diversas especialidades cirúrgicas.

Esses instrumentos devem ser utilizados em motores Cirúrgicos, conforme encaixe normatizados nos modelos indicados na descrição de cada produto.

É importante consultar-se esta INSTRUÇÃO DE USO, o rótulo do produto e a técnica operatória reserva ao médico para obter informações quanto as indicações e contraindicações de utilização, bem como aos riscos associados, modo de uso e alertas importantes.

Todos os materiais que integram o CADASTRO ANVISA nº **80482490003** são produtos invasivos cirurgicamente e não reutilizáveis.

COMPOSIÇÃO

Os materiais utilizados para fabricação dos instrumentos cirúrgicos foram selecionados em conformidade com as normas regulamentares vigentes.

São utilizados na fabricação dos instrumentais os aços inoxidáveis de acordo com a norma ASTM F899 - “Standard Specification for Wrought Stainless Steel for Surgical Instruments”: 420B.

O processo de fabricação preserva as propriedades do metal, a estabilidade, e a resistência ao trabalho e a corrosão.

INDICAÇÃO E FINALIDADE

Todos os materiais que integram o CADASTRO ANVISA nº **80482490003** são classificadas tecnicamente como instrumentos cirúrgicos, e podem ser empregadas para PERFURAÇÃO, CORTE ÓSSEO, DESBASTE ÓSSEO, ACABAMENTO DE ORIFÍCIOS ou POLIR tecidos ósseos e/ou outras estruturas de materiais (por exemplo: compósitos, cimentos cirúrgicos).

CONTRAINDICAÇÃO

Os instrumentais cirúrgicos não devem ser utilizados para finalidades diferentes daquelas para as quais foram projetados, de modo que quaisquer outras utilizações são consideradas contraindicadas ou sem substrato científico.

FORMA DE APRESENTAÇÃO

Os produtos são comercializados na condição estéril, embaladas unitariamente em sistema de DUPLO BLISTER SELADO COM PAPEL TYVEK®, e acondicionadas em uma CARTONAGEM – este sistema de embalagem promove proteção da condição estéril, bem como do produto até seu uso clínico.

Toda linha de brocas, fresas e lâminas fabricados na BIOMETAL contêm informações gravadas a laser para fins de rastreabilidade.

Sobre o BLISTER EXTERNO e sobre a CARTONAGEM de cada unidade é fixado um rótulo, que fornece os seguintes dados do produto: LOTE, DESCRIÇÃO completa, Nº CADASTRO ANVISA, NOME TÉCNICO E COMERCIAL, DATA DE FABRICAÇÃO, DATA DE VALIDADE e MÉTODO DE ESTERILIZAÇÃO empregado.

Dentro de cada embalagem é disponibilizado de forma adicional:


- ➔ 3 VIAS da ETIQUETA DE RASTREABILIDADE = a serem utilizadas para fins de rastreabilidade e;
- ➔ 1 VIA do impresso “ALERTA DE INSTRUÇÃO DE USO” = contendo as informações necessárias que permite a posterior consulta de sua IU – INSTRUÇÃO DE USO.

MODELOS COMERCIAIS

A seguir serão apresentados todos os modelos disponíveis contemplados no atual CADASTRO ANVISA nº 80482490003:

ANEXO – IMAGENS GRÁFICAS

INSTRUMENTAL ÉSTERIL NÃO ARTICULADO CORTANTE – BIOMETAL

BROCA CORTANTE ESTÉRIL	
	
CODIGO	DESCRIÇÃO
06.5.1006.006E	Broca Cortante Ø 0,6x 50,0x C6
06.5.1106.006E	Broca Cortante Ø 0,6x 55,0x C6
06.5.1206.006E	Broca Cortante Ø 0,6x 60,0x C6
06.5.1306.006E	Broca Cortante Ø 0,6x 65,0x C6
06.5.1406.006E	Broca Cortante Ø 0,6x 70,0x C6
06.5.1506.006E	Broca Cortante Ø 0,6x 75,0x C6
06.5.1606.006E	Broca Cortante Ø 0,6x 80,0x C6
06.5.1706.006E	Broca Cortante Ø 0,6x 85,0x C6
06.5.1806.006E	Broca Cortante Ø 0,6x 90,0x C6
06.5.1906.006E	Broca Cortante Ø 0,6x 95,0x C6
06.5.2006.006E	Broca Cortante Ø 0,6x 100,0x C6
06.5.2106.006E	Broca Cortante Ø 0,6x 105,0x C6

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2206.006E	Broca Cortante Ø 0,6x 110,0x C6
06.5.2306.006E	Broca Cortante Ø 0,6x 115,0x C6
06.5.2406.006E	Broca Cortante Ø 0,6x 120,0x C6
06.5.2506.006E	Broca Cortante Ø 0,6x 125,0x C6
06.5.2606.006E	Broca Cortante Ø 0,6x 130,0x C6
06.5.2706.006E	Broca Cortante Ø 0,6x 135,0x C6
06.5.2806.006E	Broca Cortante Ø 0,6x 140,0x C6
06.5.2906.006E	Broca Cortante Ø 0,6x 145,0x C6
06.5.3006.006E	Broca Cortante Ø 0,6x 150,0x C6
06.5.1006.010E	Broca Cortante Ø 1,0x 50,0x C6
06.5.1106.010E	Broca Cortante Ø 1,0x 55,0x C6
06.5.1206.010E	Broca Cortante Ø 1,0x 60,0x C6
06.5.1306.010E	Broca Cortante Ø 1,0x 65,0x C6
06.5.1406.010E	Broca Cortante Ø 1,0x 70,0x C6
06.5.1506.010E	Broca Cortante Ø 1,0x 75,0x C6
06.5.1606.010E	Broca Cortante Ø 1,0x 80,0x C6
06.5.1706.010E	Broca Cortante Ø 1,0x 85,0x C6
06.5.1806.010E	Broca Cortante Ø 1,0x 90,0x C6
06.5.1906.010E	Broca Cortante Ø 1,0x 95,0x C6
06.5.2006.010E	Broca Cortante Ø 1,0x 100,0x C6
06.5.2106.010E	Broca Cortante Ø 1,0x 105,0x C6
06.5.2206.010E	Broca Cortante Ø 1,0x 110,0x C6
06.5.2306.010E	Broca Cortante Ø 1,0x 115,0x C6
06.5.2406.010E	Broca Cortante Ø 1,0x 120,0x C6
06.5.2506.010E	Broca Cortante Ø 1,0x 125,0x C6
06.5.2606.010E	Broca Cortante Ø 1,0x 130,0x C6
06.5.2706.010E	Broca Cortante Ø 1,0x 135,0x C6
06.5.2806.010E	Broca Cortante Ø 1,0x 140,0x C6
06.5.2906.010E	Broca Cortante Ø 1,0x 145,0x C6
06.5.3006.010E	Broca Cortante Ø 1,0x 150,0x C6
06.5.1008.010E	Broca Cortante Ø 1,0x 50,0x C8
06.5.1108.010E	Broca Cortante Ø 1,0x 55,0x C8
06.5.1208.010E	Broca Cortante Ø 1,0x 60,0x C8
06.5.1308.010E	Broca Cortante Ø 1,0x 65,0x C8
06.5.1408.010E	Broca Cortante Ø 1,0x 70,0x C8
06.5.1508.010E	Broca Cortante Ø 1,0x 75,0x C8
06.5.1608.010E	Broca Cortante Ø 1,0x 80,0x C8
06.5.1708.010E	Broca Cortante Ø 1,0x 85,0x C8
06.5.1808.010E	Broca Cortante Ø 1,0x 90,0x C8
06.5.1908.010E	Broca Cortante Ø 1,0x 95,0x C8
06.5.2008.010E	Broca Cortante Ø 1,0x 100,0x C8
06.5.2108.010E	Broca Cortante Ø 1,0x 105,0x C8
06.5.2208.010E	Broca Cortante Ø 1,0x 110,0x C8
06.5.2308.010E	Broca Cortante Ø 1,0x 115,0x C8
06.5.2408.010E	Broca Cortante Ø 1,0x 120,0x C8
06.5.2508.010E	Broca Cortante Ø 1,0x 125,0x C8
06.5.2608.010E	Broca Cortante Ø 1,0x 130,0x C8
06.5.2708.010E	Broca Cortante Ø 1,0x 135,0x C8
06.5.2808.010E	Broca Cortante Ø 1,0x 140,0x C8

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2908.010E	Broca Cortante Ø 1,0x 145,0x C8
06.5.3008.010E	Broca Cortante Ø 1,0x 150,0x C8
06.5.1008.014E	Broca Cortante Ø 1,4x 50,0x C8
06.5.1108.014E	Broca Cortante Ø 1,4x 55,0x C8
06.5.1208.014E	Broca Cortante Ø 1,4x 60,0x C8
06.5.1308.014E	Broca Cortante Ø 1,4x 65,0x C8
06.5.1408.014E	Broca Cortante Ø 1,4x 70,0x C8
06.5.1508.014E	Broca Cortante Ø 1,4x 75,0x C8
06.5.1608.014E	Broca Cortante Ø 1,4x 80,0x C8
06.5.1708.014E	Broca Cortante Ø 1,4x 85,0x C8
06.5.1808.014E	Broca Cortante Ø 1,4x 90,0x C8
06.5.1908.014E	Broca Cortante Ø 1,4x 95,0x C8
06.5.2008.014E	Broca Cortante Ø 1,4x 100,0x C8
06.5.2108.014E	Broca Cortante Ø 1,4x 105,0x C8
06.5.2208.014E	Broca Cortante Ø 1,4x 110,0x C8
06.5.2308.014E	Broca Cortante Ø 1,4x 115,0x C8
06.5.2408.014E	Broca Cortante Ø 1,4x 120,0x C8
06.5.2508.014E	Broca Cortante Ø 1,4x 125,0x C8
06.5.2608.014E	Broca Cortante Ø 1,4x 130,0x C8
06.5.2708.014E	Broca Cortante Ø 1,4x 135,0x C8
06.5.2808.014E	Broca Cortante Ø 1,4x 140,0x C8
06.5.2908.014E	Broca Cortante Ø 1,4x 145,0x C8
06.5.3008.014E	Broca Cortante Ø 1,4x 150,0x C8
06.5.1006.014E	Broca Cortante Ø 1,4x 50,0x C6
06.5.1106.014E	Broca Cortante Ø 1,4x 55,0x C6
06.5.1206.014E	Broca Cortante Ø 1,4x 60,0x C6
06.5.1306.014E	Broca Cortante Ø 1,4x 65,0x C6
06.5.1406.014E	Broca Cortante Ø 1,4x 70,0x C6
06.5.1506.014E	Broca Cortante Ø 1,4x 75,0x C6
06.5.1606.014E	Broca Cortante Ø 1,4x 80,0x C6
06.5.1706.014E	Broca Cortante Ø 1,4x 85,0x C6
06.5.1806.014E	Broca Cortante Ø 1,4x 90,0x C6
06.5.1906.014E	Broca Cortante Ø 1,4x 95,0x C6
06.5.2006.014E	Broca Cortante Ø 1,4x 100,0x C6
06.5.2106.014E	Broca Cortante Ø 1,4x 105,0x C6
06.5.2206.014E	Broca Cortante Ø 1,4x 110,0x C6
06.5.2306.014E	Broca Cortante Ø 1,4x 115,0x C6
06.5.2406.014E	Broca Cortante Ø 1,4x 120,0x C6
06.5.2506.014E	Broca Cortante Ø 1,4x 125,0x C6
06.5.2606.014E	Broca Cortante Ø 1,4x 130,0x C6
06.5.2706.014E	Broca Cortante Ø 1,4x 135,0x C6
06.5.2806.014E	Broca Cortante Ø 1,4x 140,0x C6
06.5.2906.014E	Broca Cortante Ø 1,4x 145,0x C6
06.5.3006.014E	Broca Cortante Ø 1,4x 150,0x C6
06.5.1006.016E	Broca Cortante Ø 1,6x 50,0x C6
06.5.1106.016E	Broca Cortante Ø 1,6x 55,0x C6
06.5.1206.016E	Broca Cortante Ø 1,6x 60,0x C6
06.5.1306.016E	Broca Cortante Ø 1,6x 65,0x C6
06.5.1406.016E	Broca Cortante Ø 1,6x 70,0x C6

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.1506.016E	Broca Cortante Ø 1,6x 75,0x C6
06.5.1606.016E	Broca Cortante Ø 1,6x 80,0x C6
06.5.1706.016E	Broca Cortante Ø 1,6x 85,0x C6
06.5.1806.016E	Broca Cortante Ø 1,6x 90,0x C6
06.5.1906.016E	Broca Cortante Ø 1,6x 95,0x C6
06.5.2006.016E	Broca Cortante Ø 1,6x 100,0x C6
06.5.2106.016E	Broca Cortante Ø 1,6x 105,0x C6
06.5.2206.016E	Broca Cortante Ø 1,6x 110,0x C6
06.5.2306.016E	Broca Cortante Ø 1,6x 115,0x C6
06.5.2406.016E	Broca Cortante Ø 1,6x 120,0x C6
06.5.2506.016E	Broca Cortante Ø 1,6x 125,0x C6
06.5.2606.016E	Broca Cortante Ø 1,6x 130,0x C6
06.5.2706.016E	Broca Cortante Ø 1,6x 135,0x C6
06.5.2806.016E	Broca Cortante Ø 1,6x 140,0x C6
06.5.2906.016E	Broca Cortante Ø 1,6x 145,0x C6
06.5.3006.016E	Broca Cortante Ø 1,6x 150,0x C6
06.5.1008.018E	Broca Cortante Ø 1,8x 50,0x C8
06.5.1108.018E	Broca Cortante Ø 1,8x 55,0x C8
06.5.1208.018E	Broca Cortante Ø 1,8x 60,0x C8
06.5.1308.018E	Broca Cortante Ø 1,8x 65,0x C8
06.5.1408.018E	Broca Cortante Ø 1,8x 70,0x C8
06.5.1508.018E	Broca Cortante Ø 1,8x 75,0x C8
06.5.1608.018E	Broca Cortante Ø 1,8x 80,0x C8
06.5.1708.018E	Broca Cortante Ø 1,8x 85,0x C8
06.5.1808.018E	Broca Cortante Ø 1,8x 90,0x C8
06.5.1908.018E	Broca Cortante Ø 1,8x 95,0x C8
06.5.2008.018E	Broca Cortante Ø 1,8x 100,0x C8
06.5.2108.018E	Broca Cortante Ø 1,8x 105,0x C8
06.5.2208.018E	Broca Cortante Ø 1,8x 110,0x C8
06.5.2308.018E	Broca Cortante Ø 1,8x 115,0x C8
06.5.2408.018E	Broca Cortante Ø 1,8x 120,0x C8
06.5.2508.018E	Broca Cortante Ø 1,8x 125,0x C8
06.5.2608.018E	Broca Cortante Ø 1,8x 130,0x C8
06.5.2708.018E	Broca Cortante Ø 1,8x 135,0x C8
06.5.2808.018E	Broca Cortante Ø 1,8x 140,0x C8
06.5.2908.018E	Broca Cortante Ø 1,8x 145,0x C8
06.5.3008.018E	Broca Cortante Ø 1,8x 150,0x C8
06.5.1008.023E	Broca Cortante Ø 2,3x 50,0x C8
06.5.1108.023E	Broca Cortante Ø 2,3x 55,0x C8
06.5.1208.023E	Broca Cortante Ø 2,3x 60,0x C8
06.5.1308.023E	Broca Cortante Ø 2,3x 65,0x C8
06.5.1408.023E	Broca Cortante Ø 2,3x 70,0x C8
06.5.1508.023E	Broca Cortante Ø 2,3x 75,0x C8
06.5.1608.023E	Broca Cortante Ø 2,3x 80,0x C8
06.5.1708.023E	Broca Cortante Ø 2,3x 85,0x C8
06.5.1808.023E	Broca Cortante Ø 2,3x 90,0x C8
06.5.1908.023E	Broca Cortante Ø 2,3x 95,0x C8
06.5.2008.023E	Broca Cortante Ø 2,3x 100,0x C8
06.5.2108.023E	Broca Cortante Ø 2,3x 105,0x C8

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2208.023E	Broca Cortante Ø 2,3x 110,0x C8
06.5.2308.023E	Broca Cortante Ø 2,3x 115,0x C8
06.5.2408.023E	Broca Cortante Ø 2,3x 120,0x C8
06.5.2508.023E	Broca Cortante Ø 2,3x 125,0x C8
06.5.2608.023E	Broca Cortante Ø 2,3x 130,0x C8
06.5.2708.023E	Broca Cortante Ø 2,3x 135,0x C8
06.5.2808.023E	Broca Cortante Ø 2,3x 140,0x C8
06.5.2908.023E	Broca Cortante Ø 2,3x 145,0x C8
06.5.3008.023E	Broca Cortante Ø 2,3x 150,0x C8
06.5.1006.023E	Broca Cortante Ø 2,3x 50,0x C6
06.5.1106.023E	Broca Cortante Ø 2,3x 55,0x C6
06.5.1206.023E	Broca Cortante Ø 2,3x 60,0x C6
06.5.1306.023E	Broca Cortante Ø 2,3x 65,0x C6
06.5.1406.023E	Broca Cortante Ø 2,3x 70,0x C6
06.5.1506.023E	Broca Cortante Ø 2,3x 75,0x C6
06.5.1606.023E	Broca Cortante Ø 2,3x 80,0x C6
06.5.1706.023E	Broca Cortante Ø 2,3x 85,0x C6
06.5.1806.023E	Broca Cortante Ø 2,3x 90,0x C6
06.5.1906.023E	Broca Cortante Ø 2,3x 95,0x C6
06.5.2006.023E	Broca Cortante Ø 2,3x 100,0x C6
06.5.2106.023E	Broca Cortante Ø 2,3x 105,0x C6
06.5.2206.023E	Broca Cortante Ø 2,3x 110,0x C6
06.5.2306.023E	Broca Cortante Ø 2,3x 115,0x C6
06.5.2406.023E	Broca Cortante Ø 2,3x 120,0x C6
06.5.2506.023E	Broca Cortante Ø 2,3x 125,0x C6
06.5.2606.023E	Broca Cortante Ø 2,3x 130,0x C6
06.5.2706.023E	Broca Cortante Ø 2,3x 135,0x C6
06.5.2806.023E	Broca Cortante Ø 2,3x 140,0x C6
06.5.2906.023E	Broca Cortante Ø 2,3x 145,0x C6
06.5.3006.023E	Broca Cortante Ø 2,3x 150,0x C6
06.5.1010.027E	Broca Cortante Ø 2,7x 50,0x C10
06.5.1110.027E	Broca Cortante Ø 2,7x 55,0x C10
06.5.1210.027E	Broca Cortante Ø 2,7x 60,0x C10
06.5.1310.027E	Broca Cortante Ø 2,7x 65,0x C10
06.5.1410.027E	Broca Cortante Ø 2,7x 70,0x C10
06.5.1510.027E	Broca Cortante Ø 2,7x 75,0x C10
06.5.1610.027E	Broca Cortante Ø 2,7x 80,0x C10
06.5.1710.027E	Broca Cortante Ø 2,7x 85,0x C10
06.5.1810.027E	Broca Cortante Ø 2,7x 90,0x C10
06.5.1910.027E	Broca Cortante Ø 2,7x 95,0x C10
06.5.2010.027E	Broca Cortante Ø 2,7x 100,0x C10
06.5.2110.027E	Broca Cortante Ø 2,7x 105,0x C10
06.5.2210.027E	Broca Cortante Ø 2,7x 110,0x C10
06.5.2310.027E	Broca Cortante Ø 2,7x 115,0x C10
06.5.2410.027E	Broca Cortante Ø 2,7x 120,0x C10
06.5.2510.027E	Broca Cortante Ø 2,7x 125,0x C10
06.5.2610.027E	Broca Cortante Ø 2,7x 130,0x C10
06.5.2710.027E	Broca Cortante Ø 2,7x 135,0x C10
06.5.2810.027E	Broca Cortante Ø 2,7x 140,0x C10

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2910.027E	Broca Cortante Ø 2,7x 145,0x C10
06.5.3010.027E	Broca Cortante Ø 2,7x 150,0x C10
06.5.1006.027E	Broca Cortante Ø 2,7x 50,0x C6
06.5.1106.027E	Broca Cortante Ø 2,7x 55,0x C6
06.5.1206.027E	Broca Cortante Ø 2,7x 60,0x C6
06.5.1306.027E	Broca Cortante Ø 2,7x 65,0x C6
06.5.1406.027E	Broca Cortante Ø 2,7x 70,0x C6
06.5.1506.027E	Broca Cortante Ø 2,7x 75,0x C6
06.5.1606.027E	Broca Cortante Ø 2,7x 80,0x C6
06.5.1706.027E	Broca Cortante Ø 2,7x 85,0x C6
06.5.1806.027E	Broca Cortante Ø 2,7x 90,0x C6
06.5.1906.027E	Broca Cortante Ø 2,7x 95,0x C6
06.5.2006.027E	Broca Cortante Ø 2,7x 100,0x C6
06.5.2106.027E	Broca Cortante Ø 2,7x 105,0x C6
06.5.2206.027E	Broca Cortante Ø 2,7x 110,0x C6
06.5.2306.027E	Broca Cortante Ø 2,7x 115,0x C6
06.5.2406.027E	Broca Cortante Ø 2,7x 120,0x C6
06.5.2506.027E	Broca Cortante Ø 2,7x 125,0x C6
06.5.2606.027E	Broca Cortante Ø 2,7x 130,0x C6
06.5.2706.027E	Broca Cortante Ø 2,7x 135,0x C6
06.5.2806.027E	Broca Cortante Ø 2,7x 140,0x C6
06.5.2906.027E	Broca Cortante Ø 2,7x 145,0x C6
06.5.3006.027E	Broca Cortante Ø 2,7x 150,0x C6
06.5.1010.031E	Broca Cortante Ø 3,1x 50,0x C10
06.5.1110.031E	Broca Cortante Ø 3,1x 55,0x C10
06.5.1210.031E	Broca Cortante Ø 3,1x 60,0x C10
06.5.1310.031E	Broca Cortante Ø 3,1x 65,0x C10
06.5.1410.031E	Broca Cortante Ø 3,1x 70,0x C10
06.5.1510.031E	Broca Cortante Ø 3,1x 75,0x C10
06.5.1610.031E	Broca Cortante Ø 3,1x 80,0x C10
06.5.1710.031E	Broca Cortante Ø 3,1x 85,0x C10
06.5.1810.031E	Broca Cortante Ø 3,1x 90,0x C10
06.5.1910.031E	Broca Cortante Ø 3,1x 95,0x C10
06.5.2010.031E	Broca Cortante Ø 3,1x 100,0x C10
06.5.2110.031E	Broca Cortante Ø 3,1x 105,0x C10
06.5.2210.031E	Broca Cortante Ø 3,1x 110,0x C10
06.5.2310.031E	Broca Cortante Ø 3,1x 115,0x C10
06.5.2410.031E	Broca Cortante Ø 3,1x 120,0x C10
06.5.2510.031E	Broca Cortante Ø 3,1x 125,0x C10
06.5.2610.031E	Broca Cortante Ø 3,1x 130,0x C10
06.5.2710.031E	Broca Cortante Ø 3,1x 135,0x C10
06.5.2810.031E	Broca Cortante Ø 3,1x 140,0x C10
06.5.2910.031E	Broca Cortante Ø 3,1x 145,0x C10
06.5.3010.031E	Broca Cortante Ø 3,1x 150,0x C10
06.5.1006.031E	Broca Cortante Ø 3,1x 50,0x C6
06.5.1106.031E	Broca Cortante Ø 3,1x 55,0x C6
06.5.1206.031E	Broca Cortante Ø 3,1x 60,0x C6
06.5.1306.031E	Broca Cortante Ø 3,1x 65,0x C6
06.5.1406.031E	Broca Cortante Ø 3,1x 70,0x C6

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.1506.031E	Broca Cortante Ø 3,1x 75,0x C6
06.5.1606.031E	Broca Cortante Ø 3,1x 80,0x C6
06.5.1706.031E	Broca Cortante Ø 3,1x 85,0x C6
06.5.1806.031E	Broca Cortante Ø 3,1x 90,0x C6
06.5.1906.031E	Broca Cortante Ø 3,1x 95,0x C6
06.5.2006.031E	Broca Cortante Ø 3,1x 100,0x C6
06.5.2106.031E	Broca Cortante Ø 3,1x 105,0x C6
06.5.2206.031E	Broca Cortante Ø 3,1x 110,0x C6
06.5.2306.031E	Broca Cortante Ø 3,1x 115,0x C6
06.5.2406.031E	Broca Cortante Ø 3,1x 120,0x C6
06.5.2506.031E	Broca Cortante Ø 3,1x 125,0x C6
06.5.2606.031E	Broca Cortante Ø 3,1x 130,0x C6
06.5.2706.031E	Broca Cortante Ø 3,1x 135,0x C6
06.5.2806.031E	Broca Cortante Ø 3,1x 140,0x C6
06.5.2906.031E	Broca Cortante Ø 3,1x 145,0x C6
06.5.3006.031E	Broca Cortante Ø 3,1x 150,0x C6
06.5.1012.035E	Broca Cortante Ø 3,5x 50,0x C12
06.5.1112.035E	Broca Cortante Ø 3,5x 55,0x C12
06.5.1212.035E	Broca Cortante Ø 3,5x 60,0x C12
06.5.1312.035E	Broca Cortante Ø 3,5x 65,0x C12
06.5.1412.035E	Broca Cortante Ø 3,5x 70,0x C12
06.5.1512.035E	Broca Cortante Ø 3,5x 75,0x C12
06.5.1612.035E	Broca Cortante Ø 3,5x 80,0x C12
06.5.1712.035E	Broca Cortante Ø 3,5x 85,0x C12
06.5.1812.035E	Broca Cortante Ø 3,5x 90,0x C12
06.5.1912.035E	Broca Cortante Ø 3,5x 95,0x C12
06.5.2012.035E	Broca Cortante Ø 3,5x 100,0x C12
06.5.2112.035E	Broca Cortante Ø 3,5x 105,0x C12
06.5.2212.035E	Broca Cortante Ø 3,5x 110,0x C12
06.5.2312.035E	Broca Cortante Ø 3,5x 115,0x C12
06.5.2412.035E	Broca Cortante Ø 3,5x 120,0x C12
06.5.2512.035E	Broca Cortante Ø 3,5x 125,0x C12
06.5.2612.035E	Broca Cortante Ø 3,5x 130,0x C12
06.5.2712.035E	Broca Cortante Ø 3,5x 135,0x C12
06.5.2812.035E	Broca Cortante Ø 3,5x 140,0x C12
06.5.2912.035E	Broca Cortante Ø 3,5x 145,0x C12
06.5.3012.035E	Broca Cortante Ø 3,5x 150,0x C12
06.5.1006.035E	Broca Cortante Ø 3,5x 50,0x C6
06.5.1106.035E	Broca Cortante Ø 3,5x 55,0x C6
06.5.1206.035E	Broca Cortante Ø 3,5x 60,0x C6
06.5.1306.035E	Broca Cortante Ø 3,5x 65,0x C6
06.5.1406.035E	Broca Cortante Ø 3,5x 70,0x C6
06.5.1506.035E	Broca Cortante Ø 3,5x 75,0x C6
06.5.1606.035E	Broca Cortante Ø 3,5x 80,0x C6
06.5.1706.035E	Broca Cortante Ø 3,5x 85,0x C6
06.5.1806.035E	Broca Cortante Ø 3,5x 90,0x C6
06.5.1906.035E	Broca Cortante Ø 3,5x 95,0x C6
06.5.2006.035E	Broca Cortante Ø 3,5x 100,0x C6
06.5.2106.035E	Broca Cortante Ø 3,5x 105,0x C6

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2206.035E	Broca Cortante Ø 3,5x 110,0x C6
06.5.2306.035E	Broca Cortante Ø 3,5x 115,0x C6
06.5.2406.035E	Broca Cortante Ø 3,5x 120,0x C6
06.5.2506.035E	Broca Cortante Ø 3,5x 125,0x C6
06.5.2606.035E	Broca Cortante Ø 3,5x 130,0x C6
06.5.2706.035E	Broca Cortante Ø 3,5x 135,0x C6
06.5.2806.035E	Broca Cortante Ø 3,5x 140,0x C6
06.5.2906.035E	Broca Cortante Ø 3,5x 145,0x C6
06.5.3006.035E	Broca Cortante Ø 3,5x 150,0x C6
06.5.1012.040E	Broca Cortante Ø 4,0x 50,0x C12
06.5.1112.040E	Broca Cortante Ø 4,0x 55,0x C12
06.5.1212.040E	Broca Cortante Ø 4,0x 60,0x C12
06.5.1312.040E	Broca Cortante Ø 4,0x 65,0x C12
06.5.1412.040E	Broca Cortante Ø 4,0x 70,0x C12
06.5.1512.040E	Broca Cortante Ø 4,0x 75,0x C12
06.5.1612.040E	Broca Cortante Ø 4,0x 80,0x C12
06.5.1712.040E	Broca Cortante Ø 4,0x 85,0x C12
06.5.1812.040E	Broca Cortante Ø 4,0x 90,0x C12
06.5.1912.040E	Broca Cortante Ø 4,0x 95,0x C12
06.5.2012.040E	Broca Cortante Ø 4,0x 100,0x C12
06.5.2112.040E	Broca Cortante Ø 4,0x 105,0x C12
06.5.2212.040E	Broca Cortante Ø 4,0x 110,0x C12
06.5.2312.040E	Broca Cortante Ø 4,0x 115,0x C12
06.5.2412.040E	Broca Cortante Ø 4,0x 120,0x C12
06.5.2512.040E	Broca Cortante Ø 4,0x 125,0x C12
06.5.2612.040E	Broca Cortante Ø 4,0x 130,0x C12
06.5.2712.040E	Broca Cortante Ø 4,0x 135,0x C12
06.5.2812.040E	Broca Cortante Ø 4,0x 140,0x C12
06.5.2912.040E	Broca Cortante Ø 4,0x 145,0x C12
06.5.3012.040E	Broca Cortante Ø 4,0x 150,0x C12
06.5.1006.040E	Broca Cortante Ø 4,0x 50,0x C6
06.5.1106.040E	Broca Cortante Ø 4,0x 55,0x C6
06.5.1206.040E	Broca Cortante Ø 4,0x 60,0x C6
06.5.1306.040E	Broca Cortante Ø 4,0x 65,0x C6
06.5.1406.040E	Broca Cortante Ø 4,0x 70,0x C6
06.5.1506.040E	Broca Cortante Ø 4,0x 75,0x C6
06.5.1606.040E	Broca Cortante Ø 4,0x 80,0x C6
06.5.1706.040E	Broca Cortante Ø 4,0x 85,0x C6
06.5.1806.040E	Broca Cortante Ø 4,0x 90,0x C6
06.5.1906.040E	Broca Cortante Ø 4,0x 95,0x C6
06.5.2006.040E	Broca Cortante Ø 4,0x 100,0x C6
06.5.2106.040E	Broca Cortante Ø 4,0x 105,0x C6
06.5.2206.040E	Broca Cortante Ø 4,0x 110,0x C6
06.5.2306.040E	Broca Cortante Ø 4,0x 115,0x C6
06.5.2406.040E	Broca Cortante Ø 4,0x 120,0x C6
06.5.2506.040E	Broca Cortante Ø 4,0x 125,0x C6
06.5.2606.040E	Broca Cortante Ø 4,0x 130,0x C6
06.5.2706.040E	Broca Cortante Ø 4,0x 135,0x C6
06.5.2806.040E	Broca Cortante Ø 4,0x 140,0x C6

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2906.040E	Broca Cortante Ø 4,0x 145,0x C6
06.5.3006.040E	Broca Cortante Ø 4,0x 150,0x C6
06.5.1012.045E	Broca Cortante Ø 4,5x 50,0x C12
06.5.1112.045E	Broca Cortante Ø 4,5x 55,0x C12
06.5.1212.045E	Broca Cortante Ø 4,5x 60,0x C12
06.5.1312.045E	Broca Cortante Ø 4,5x 65,0x C12
06.5.1412.045E	Broca Cortante Ø 4,5x 70,0x C12
06.5.1512.045E	Broca Cortante Ø 4,5x 75,0x C12
06.5.1612.045E	Broca Cortante Ø 4,5x 80,0x C12
06.5.1712.045E	Broca Cortante Ø 4,5x 85,0x C12
06.5.1812.045E	Broca Cortante Ø 4,5x 90,0x C12
06.5.1912.045E	Broca Cortante Ø 4,5x 95,0x C12
06.5.2012.045E	Broca Cortante Ø 4,5x 100,0x C12
06.5.2112.045E	Broca Cortante Ø 4,5x 105,0x C12
06.5.2212.045E	Broca Cortante Ø 4,5x 110,0x C12
06.5.2312.045E	Broca Cortante Ø 4,5x 115,0x C12
06.5.2412.045E	Broca Cortante Ø 4,5x 120,0x C12
06.5.2512.045E	Broca Cortante Ø 4,5x 125,0x C12
06.5.2612.045E	Broca Cortante Ø 4,5x 130,0x C12
06.5.2712.045E	Broca Cortante Ø 4,5x 135,0x C12
06.5.2812.045E	Broca Cortante Ø 4,5x 140,0x C12
06.5.2912.045E	Broca Cortante Ø 4,5x 145,0x C12
06.5.3012.045E	Broca Cortante Ø 4,5x 150,0x C12
06.5.1006.045E	Broca Cortante Ø 4,5x 50,0x C6
06.5.1106.045E	Broca Cortante Ø 4,5x 55,0x C6
06.5.1206.045E	Broca Cortante Ø 4,5x 60,0x C6
06.5.1306.045E	Broca Cortante Ø 4,5x 65,0x C6
06.5.1406.045E	Broca Cortante Ø 4,5x 70,0x C6
06.5.1506.045E	Broca Cortante Ø 4,5x 75,0x C6
06.5.1606.045E	Broca Cortante Ø 4,5x 80,0x C6
06.5.1706.045E	Broca Cortante Ø 4,5x 85,0x C6
06.5.1806.045E	Broca Cortante Ø 4,5x 90,0x C6
06.5.1906.045E	Broca Cortante Ø 4,5x 95,0x C6
06.5.2006.045E	Broca Cortante Ø 4,5x 100,0x C6
06.5.2106.045E	Broca Cortante Ø 4,5x 105,0x C6
06.5.2206.045E	Broca Cortante Ø 4,5x 110,0x C6
06.5.2306.045E	Broca Cortante Ø 4,5x 115,0x C6
06.5.2406.045E	Broca Cortante Ø 4,5x 120,0x C6
06.5.2506.045E	Broca Cortante Ø 4,5x 125,0x C6
06.5.2606.045E	Broca Cortante Ø 4,5x 130,0x C6
06.5.2706.045E	Broca Cortante Ø 4,5x 135,0x C6
06.5.2806.045E	Broca Cortante Ø 4,5x 140,0x C6
06.5.2906.045E	Broca Cortante Ø 4,5x 145,0x C6
06.5.3006.045E	Broca Cortante Ø 4,5x 150,0x C6
06.5.1012.050E	Broca Cortante Ø 5,0x 50,0x C12
06.5.1112.050E	Broca Cortante Ø 5,0x 55,0x C12
06.5.1212.050E	Broca Cortante Ø 5,0x 60,0x C12
06.5.1312.050E	Broca Cortante Ø 5,0x 65,0x C12
06.5.1412.050E	Broca Cortante Ø 5,0x 70,0x C12

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.1512.050E	Broca Cortante Ø 5,0x 75,0x C12
06.5.1612.050E	Broca Cortante Ø 5,0x 80,0x C12
06.5.1712.050E	Broca Cortante Ø 5,0x 85,0x C12
06.5.1812.050E	Broca Cortante Ø 5,0x 90,0x C12
06.5.1912.050E	Broca Cortante Ø 5,0x 95,0x C12
06.5.2012.050E	Broca Cortante Ø 5,0x 100,0x C12
06.5.2112.050E	Broca Cortante Ø 5,0x 105,0x C12
06.5.2212.050E	Broca Cortante Ø 5,0x 110,0x C12
06.5.2312.050E	Broca Cortante Ø 5,0x 115,0x C12
06.5.2412.050E	Broca Cortante Ø 5,0x 120,0x C12
06.5.2512.050E	Broca Cortante Ø 5,0x 125,0x C12
06.5.2612.050E	Broca Cortante Ø 5,0x 130,0x C12
06.5.2712.050E	Broca Cortante Ø 5,0x 135,0x C12
06.5.2812.050E	Broca Cortante Ø 5,0x 140,0x C12
06.5.2912.050E	Broca Cortante Ø 5,0x 145,0x C12
06.5.3012.050E	Broca Cortante Ø 5,0x 150,0x C12
06.5.1008.050E	Broca Cortante Ø 5,0x 50,0x C8
06.5.1108.050E	Broca Cortante Ø 5,0x 55,0x C8
06.5.1208.050E	Broca Cortante Ø 5,0x 60,0x C8
06.5.1308.050E	Broca Cortante Ø 5,0x 65,0x C8
06.5.1408.050E	Broca Cortante Ø 5,0x 70,0x C8
06.5.1508.050E	Broca Cortante Ø 5,0x 75,0x C8
06.5.1608.050E	Broca Cortante Ø 5,0x 80,0x C8
06.5.1708.050E	Broca Cortante Ø 5,0x 85,0x C8
06.5.1808.050E	Broca Cortante Ø 5,0x 90,0x C8
06.5.1908.050E	Broca Cortante Ø 5,0x 95,0x C8
06.5.2008.050E	Broca Cortante Ø 5,0x 100,0x C8
06.5.2108.050E	Broca Cortante Ø 5,0x 105,0x C8
06.5.2208.050E	Broca Cortante Ø 5,0x 110,0x C8
06.5.2308.050E	Broca Cortante Ø 5,0x 115,0x C8
06.5.2408.050E	Broca Cortante Ø 5,0x 120,0x C8
06.5.2508.050E	Broca Cortante Ø 5,0x 125,0x C8
06.5.2608.050E	Broca Cortante Ø 5,0x 130,0x C8
06.5.2708.050E	Broca Cortante Ø 5,0x 135,0x C8
06.5.2808.050E	Broca Cortante Ø 5,0x 140,0x C8
06.5.2908.050E	Broca Cortante Ø 5,0x 145,0x C8
06.5.3008.050E	Broca Cortante Ø 5,0x 150,0x C8
06.5.1012.055E	Broca Cortante Ø 5,5x 50,0x C12
06.5.1112.055E	Broca Cortante Ø 5,5x 55,0x C12
06.5.1212.055E	Broca Cortante Ø 5,5x 60,0x C12
06.5.1312.055E	Broca Cortante Ø 5,5x 65,0x C12
06.5.1412.055E	Broca Cortante Ø 5,5x 70,0x C12
06.5.1512.055E	Broca Cortante Ø 5,5x 75,0x C12
06.5.1612.055E	Broca Cortante Ø 5,5x 80,0x C12
06.5.1712.055E	Broca Cortante Ø 5,5x 85,0x C12
06.5.1812.055E	Broca Cortante Ø 5,5x 90,0x C12
06.5.1912.055E	Broca Cortante Ø 5,5x 95,0x C12
06.5.2012.055E	Broca Cortante Ø 5,5x 100,0x C12
06.5.2112.055E	Broca Cortante Ø 5,5x 105,0x C12

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2212.055E	Broca Cortante Ø 5,5x 110,0x C12
06.5.2312.055E	Broca Cortante Ø 5,5x 115,0x C12
06.5.2412.055E	Broca Cortante Ø 5,5x 120,0x C12
06.5.2512.055E	Broca Cortante Ø 5,5x 125,0x C12
06.5.2612.055E	Broca Cortante Ø 5,5x 130,0x C12
06.5.2712.055E	Broca Cortante Ø 5,5x 135,0x C12
06.5.2812.055E	Broca Cortante Ø 5,5x 140,0x C12
06.5.2912.055E	Broca Cortante Ø 5,5x 145,0x C12
06.5.3012.055E	Broca Cortante Ø 5,5x 150,0x C12
06.5.1008.055E	Broca Cortante Ø 5,5x 50,0x C8
06.5.1108.055E	Broca Cortante Ø 5,5x 55,0x C8
06.5.1208.055E	Broca Cortante Ø 5,5x 60,0x C8
06.5.1308.055E	Broca Cortante Ø 5,5x 65,0x C8
06.5.1408.055E	Broca Cortante Ø 5,5x 70,0x C8
06.5.1508.055E	Broca Cortante Ø 5,5x 75,0x C8
06.5.1608.055E	Broca Cortante Ø 5,5x 80,0x C8
06.5.1708.055E	Broca Cortante Ø 5,5x 85,0x C8
06.5.1808.055E	Broca Cortante Ø 5,5x 90,0x C8
06.5.1908.055E	Broca Cortante Ø 5,5x 95,0x C8
06.5.2008.055E	Broca Cortante Ø 5,5x 100,0x C8
06.5.2108.055E	Broca Cortante Ø 5,5x 105,0x C8
06.5.2208.055E	Broca Cortante Ø 5,5x 110,0x C8
06.5.2308.055E	Broca Cortante Ø 5,5x 115,0x C8
06.5.2408.055E	Broca Cortante Ø 5,5x 120,0x C8
06.5.2508.055E	Broca Cortante Ø 5,5x 125,0x C8
06.5.2608.055E	Broca Cortante Ø 5,5x 130,0x C8
06.5.2708.055E	Broca Cortante Ø 5,5x 135,0x C8
06.5.2808.055E	Broca Cortante Ø 5,5x 140,0x C8
06.5.2908.055E	Broca Cortante Ø 5,5x 145,0x C8
06.5.3008.055E	Broca Cortante Ø 5,5x 150,0x C8
06.5.1014.060E	Broca Cortante Ø 6,0x 50,0x C14
06.5.1114.060E	Broca Cortante Ø 6,0x 55,0x C14
06.5.1214.060E	Broca Cortante Ø 6,0x 60,0x C14
06.5.1314.060E	Broca Cortante Ø 6,0x 65,0x C14
06.5.1414.060E	Broca Cortante Ø 6,0x 70,0x C14
06.5.1514.060E	Broca Cortante Ø 6,0x 75,0x C14
06.5.1614.060E	Broca Cortante Ø 6,0x 80,0x C14
06.5.1714.060E	Broca Cortante Ø 6,0x 85,0x C14
06.5.1814.060E	Broca Cortante Ø 6,0x 90,0x C14
06.5.1914.060E	Broca Cortante Ø 6,0x 95,0x C14
06.5.2014.060E	Broca Cortante Ø 6,0x 100,0x C14
06.5.2114.060E	Broca Cortante Ø 6,0x 105,0x C14
06.5.2214.060E	Broca Cortante Ø 6,0x 110,0x C14
06.5.2314.060E	Broca Cortante Ø 6,0x 115,0x C14
06.5.2414.060E	Broca Cortante Ø 6,0x 120,0x C14
06.5.2514.060E	Broca Cortante Ø 6,0x 125,0x C14
06.5.2614.060E	Broca Cortante Ø 6,0x 130,0x C14
06.5.2714.060E	Broca Cortante Ø 6,0x 135,0x C14
06.5.2814.060E	Broca Cortante Ø 6,0x 140,0x C14

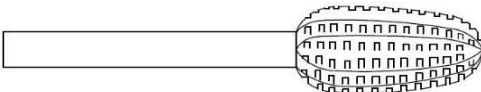
INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

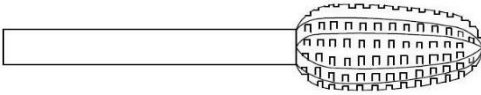
06.5.2914.060E	Broca Cortante Ø 6,0x 145,0x C14
06.5.3014.060E	Broca Cortante Ø 6,0x 150,0x C14
06.5.1010.060E	Broca Cortante Ø 6,0x 50,0x C10
06.5.1110.060E	Broca Cortante Ø 6,0x 55,0x C10
06.5.1210.060E	Broca Cortante Ø 6,0x 60,0x C10
06.5.1310.060E	Broca Cortante Ø 6,0x 65,0x C10
06.5.1410.060E	Broca Cortante Ø 6,0x 70,0x C10
06.5.1510.060E	Broca Cortante Ø 6,0x 75,0x C10
06.5.1610.060E	Broca Cortante Ø 6,0x 80,0x C10
06.5.1710.060E	Broca Cortante Ø 6,0x 85,0x C10
06.5.1810.060E	Broca Cortante Ø 6,0x 90,0x C10
06.5.1910.060E	Broca Cortante Ø 6,0x 95,0x C10
06.5.2010.060E	Broca Cortante Ø 6,0x 100,0x C10
06.5.2110.060E	Broca Cortante Ø 6,0x 105,0x C10
06.5.2210.060E	Broca Cortante Ø 6,0x 110,0x C10
06.5.2310.060E	Broca Cortante Ø 6,0x 115,0x C10
06.5.2410.060E	Broca Cortante Ø 6,0x 120,0x C10
06.5.2510.060E	Broca Cortante Ø 6,0x 125,0x C10
06.5.2610.060E	Broca Cortante Ø 6,0x 130,0x C10
06.5.2710.060E	Broca Cortante Ø 6,0x 135,0x C10
06.5.2810.060E	Broca Cortante Ø 6,0x 140,0x C10
06.5.2910.060E	Broca Cortante Ø 6,0x 145,0x C10
06.5.3010.060E	Broca Cortante Ø 6,0x 150,0x C10
06.5.1016.070E	Broca Cortante Ø 7,0x 50,0x C16
06.5.1116.070E	Broca Cortante Ø 7,0x 55,0x C16
06.5.1216.070E	Broca Cortante Ø 7,0x 60,0x C16
06.5.1316.070E	Broca Cortante Ø 7,0x 65,0x C16
06.5.1416.070E	Broca Cortante Ø 7,0x 70,0x C16
06.5.1516.070E	Broca Cortante Ø 7,0x 75,0x C16
06.5.1616.070E	Broca Cortante Ø 7,0x 80,0x C16
06.5.1716.070E	Broca Cortante Ø 7,0x 85,0x C16
06.5.1816.070E	Broca Cortante Ø 7,0x 90,0x C16
06.5.1916.070E	Broca Cortante Ø 7,0x 95,0x C16
06.5.2016.070E	Broca Cortante Ø 7,0x 100,0x C16
06.5.2116.070E	Broca Cortante Ø 7,0x 105,0x C16
06.5.2216.070E	Broca Cortante Ø 7,0x 110,0x C16
06.5.2316.070E	Broca Cortante Ø 7,0x 115,0x C16
06.5.2416.070E	Broca Cortante Ø 7,0x 120,0x C16
06.5.2516.070E	Broca Cortante Ø 7,0x 125,0x C16
06.5.2616.070E	Broca Cortante Ø 7,0x 130,0x C16
06.5.2716.070E	Broca Cortante Ø 7,0x 135,0x C16
06.5.2816.070E	Broca Cortante Ø 7,0x 140,0x C16
06.5.2916.070E	Broca Cortante Ø 7,0x 145,0x C16
06.5.3016.070E	Broca Cortante Ø 7,0x 150,0x C16
06.5.1010.070E	Broca Cortante Ø 7,0x 50,0x C10
06.5.1110.070E	Broca Cortante Ø 7,0x 55,0x C10
06.5.1210.070E	Broca Cortante Ø 7,0x 60,0x C10
06.5.1310.070E	Broca Cortante Ø 7,0x 65,0x C10
06.5.1410.070E	Broca Cortante Ø 7,0x 70,0x C10

06.5.1510.070E	Broca Cortante Ø 7,0x 75,0x C10
06.5.1610.070E	Broca Cortante Ø 7,0x 80,0x C10
06.5.1710.070E	Broca Cortante Ø 7,0x 85,0x C10
06.5.1810.070E	Broca Cortante Ø 7,0x 90,0x C10
06.5.1910.070E	Broca Cortante Ø 7,0x 95,0x C10
06.5.2010.070E	Broca Cortante Ø 7,0x 100,0x C10
06.5.2110.070E	Broca Cortante Ø 7,0x 105,0x C10
06.5.2210.070E	Broca Cortante Ø 7,0x 110,0x C10
06.5.2310.070E	Broca Cortante Ø 7,0x 115,0x C10
06.5.2410.070E	Broca Cortante Ø 7,0x 120,0x C10
06.5.2510.070E	Broca Cortante Ø 7,0x 125,0x C10
06.5.2610.070E	Broca Cortante Ø 7,0x 130,0x C10
06.5.2710.070E	Broca Cortante Ø 7,0x 135,0x C10
06.5.2810.070E	Broca Cortante Ø 7,0x 140,0x C10
06.5.2910.070E	Broca Cortante Ø 7,0x 145,0x C10
06.5.3010.070E	Broca Cortante Ø 7,0x 150,0x C10

BROCA CORTANTE CILINDRICA CROSS-CUT ESTÉRIL	
	
CODIGO	DESCRIÇÃO
06.5.1010.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 50,0x C10
06.5.1110.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 55,0x C10
06.5.1210.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 60,0x C10
06.5.1310.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 65,0x C10
06.5.1410.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 70,0x C10
06.5.1510.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 75,0x C10
06.5.1610.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 80,0x C10
06.5.1710.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 85,0x C10
06.5.1810.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 90,0x C10
06.5.1910.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 95,0x C10
06.5.2010.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 100,0x C10
06.5.2110.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 105,0x C10
06.5.2210.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 110,0x C10
06.5.2310.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 115,0x C10
06.5.2410.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 120,0x C10
06.5.2510.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 125,0x C10
06.5.2610.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 130,0x C10
06.5.2710.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 135,0x C10
06.5.2810.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 140,0x C10
06.5.2910.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 145,0x C10
06.5.3010.140E	Broca Cilíndrica Cross-Cut Ø 4,0x 150,0x C10
06.5.1012.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 50,0x C12
06.5.1112.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 55,0x C12

06.5.1212.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 60,0x C12
06.5.1312.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 65,0x C12
06.5.1412.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 70,0x C12
06.5.1512.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 75,0x C12
06.5.1612.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 80,0x C12
06.5.1712.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 85,0x C12
06.5.1812.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 90,0x C12
06.5.1912.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 95,0x C12
06.5.2012.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 100,0x C12
06.5.2112.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 105,0x C12
06.5.2212.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 110,0x C12
06.5.2312.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 115,0x C12
06.5.2412.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 120,0x C12
06.5.2512.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 125,0x C12
06.5.2612.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 130,0x C12
06.5.2712.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 135,0x C12
06.5.2812.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 140,0x C12
06.5.2912.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 145,0x C12
06.5.3012.150E	Broca Cilíndrica Cross-Cut Ø 5,0x 150,0x C12
06.5.1014.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 50,0x C14
06.5.1114.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 55,0x C14
06.5.1214.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 60,0x C14
06.5.1314.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 65,0x C14
06.5.1414.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 70,0x C14
06.5.1514.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 75,0x C14
06.5.1614.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 80,0x C14
06.5.1714.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 85,0x C14
06.5.1814.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 90,0x C14
06.5.1914.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 95,0x C14
06.5.2014.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 100,0x C14
06.5.2114.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 105,0x C14
06.5.2214.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 110,0x C14
06.5.2314.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 115,0x C14
06.5.2414.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 120,0x C14
06.5.2514.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 125,0x C14
06.5.2614.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 130,0x C14
06.5.2714.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 135,0x C14
06.5.2814.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 140,0x C14
06.5.2914.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 145,0x C14
06.5.3014.160E	Broca Cilíndrica Cross-Cut Ø 6,0x 150,0x C14
06.5.1016.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 50,0x C16
06.5.1116.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 55,0x C16
06.5.1216.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 60,0x C16
06.5.1316.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 65,0x C16
06.5.1416.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 70,0x C16
06.5.1516.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 75,0x C16
06.5.1616.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 80,0x C16
06.5.1716.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 85,0x C16
06.5.1816.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 90,0x C16

06.5.1916.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 95,0x C16
06.5.2016.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 100,0x C16
06.5.2116.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 105,0x C16
06.5.2216.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 110,0x C16
06.5.2316.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 115,0x C16
06.5.2416.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 120,0x C16
06.5.2516.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 125,0x C16
06.5.2616.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 130,0x C16
06.5.2716.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 135,0x C16
06.5.2816.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 140,0x C16
06.5.2916.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 145,0x C16
06.5.3016.170E	Broca Cilíndrica Cross-Cut Ø 7,0x 150,0x C16
06.5.1016.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 50,0x C16
06.5.1116.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 55,0x C16
06.5.1216.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 60,0x C16
06.5.1316.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 65,0x C16
06.5.1416.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 70,0x C16
06.5.1516.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 75,0x C16
06.5.1616.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 80,0x C16
06.5.1716.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 85,0x C16
06.5.1816.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 90,0x C16
06.5.1916.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 95,0x C16
06.5.2016.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 100,0x C16
06.5.2116.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 105,0x C16
06.5.2216.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 110,0x C16
06.5.2316.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 115,0x C16
06.5.2416.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 120,0x C16
06.5.2516.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 125,0x C16
06.5.2616.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 130,0x C16
06.5.2716.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 135,0x C16
06.5.2816.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 140,0x C16
06.5.2916.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 145,0x C16
06.5.3016.180E	Broca Cilíndrica Cross-Cut Ø 8,0x 150,0x C16

BROCA CORTANTE CILÍNDRICA ESTÉRIL	
	
CODIGO	DESCRIÇÃO
06.5.1010.240E	Broca Cilíndrica Ø 4,0x 50,0x C10
06.5.1110.240E	Broca Cilíndrica Ø 4,0x 55,0x C10
06.5.1210.240E	Broca Cilíndrica Ø 4,0x 60,0x C10
06.5.1310.240E	Broca Cilíndrica Ø 4,0x 65,0x C10
06.5.1410.240E	Broca Cilíndrica Ø 4,0x 70,0x C10
06.5.1510.240E	Broca Cilíndrica Ø 4,0x 75,0x C10
06.5.1610.240E	Broca Cilíndrica Ø 4,0x 80,0x C10
06.5.1710.240E	Broca Cilíndrica Ø 4,0x 85,0x C10
06.5.1810.240E	Broca Cilíndrica Ø 4,0x 90,0x C10

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.1910.240E	Broca Cilíndrica Ø 4,0x 95,0x C10
06.5.2010.240E	Broca Cilíndrica Ø 4,0x 100,0x C10
06.5.2110.240E	Broca Cilíndrica Ø 4,0x 105,0x C10
06.5.2210.240E	Broca Cilíndrica Ø 4,0x 110,0x C10
06.5.2310.240E	Broca Cilíndrica Ø 4,0x 115,0x C10
06.5.2410.240E	Broca Cilíndrica Ø 4,0x 120,0x C10
06.5.2510.240E	Broca Cilíndrica Ø 4,0x 125,0x C10
06.5.2610.240E	Broca Cilíndrica Ø 4,0x 130,0x C10
06.5.2710.240E	Broca Cilíndrica Ø 4,0x 135,0x C10
06.5.2810.240E	Broca Cilíndrica Ø 4,0x 140,0x C10
06.5.2910.240E	Broca Cilíndrica Ø 4,0x 145,0x C10
06.5.3010.240E	Broca Cilíndrica Ø 4,0x 150,0x C10
06.5.1012.250E	Broca Cilíndrica Ø 5,0x 50,0x C12
06.5.1112.250E	Broca Cilíndrica Ø 5,0x 55,0x C12
06.5.1212.250E	Broca Cilíndrica Ø 5,0x 60,0x C12
06.5.1312.250E	Broca Cilíndrica Ø 5,0x 65,0x C12
06.5.1412.250E	Broca Cilíndrica Ø 5,0x 70,0x C12
06.5.1512.250E	Broca Cilíndrica Ø 5,0x 75,0x C12
06.5.1612.250E	Broca Cilíndrica Ø 5,0x 80,0x C12
06.5.1712.250E	Broca Cilíndrica Ø 5,0x 85,0x C12
06.5.1812.250E	Broca Cilíndrica Ø 5,0x 90,0x C12
06.5.1912.250E	Broca Cilíndrica Ø 5,0x 95,0x C12
06.5.2012.250E	Broca Cilíndrica Ø 5,0x 100,0x C12
06.5.2112.250E	Broca Cilíndrica Ø 5,0x 105,0x C12
06.5.2212.250E	Broca Cilíndrica Ø 5,0x 110,0x C12
06.5.2312.250E	Broca Cilíndrica Ø 5,0x 115,0x C12
06.5.2412.250E	Broca Cilíndrica Ø 5,0x 120,0x C12
06.5.2512.250E	Broca Cilíndrica Ø 5,0x 125,0x C12
06.5.2612.250E	Broca Cilíndrica Ø 5,0x 130,0x C12
06.5.2712.250E	Broca Cilíndrica Ø 5,0x 135,0x C12
06.5.2812.250E	Broca Cilíndrica Ø 5,0x 140,0x C12
06.5.2912.250E	Broca Cilíndrica Ø 5,0x 145,0x C12
06.5.3012.250E	Broca Cilíndrica Ø 5,0x 150,0x C12
06.5.1014.260E	Broca Cilíndrica Ø 6,0x 50,0x C14
06.5.1114.260E	Broca Cilíndrica Ø 6,0x 55,0x C14
06.5.1214.260E	Broca Cilíndrica Ø 6,0x 60,0x C14
06.5.1314.260E	Broca Cilíndrica Ø 6,0x 65,0x C14
06.5.1414.260E	Broca Cilíndrica Ø 6,0x 70,0x C14
06.5.1514.260E	Broca Cilíndrica Ø 6,0x 75,0x C14
06.5.1614.260E	Broca Cilíndrica Ø 6,0x 80,0x C14
06.5.1714.260E	Broca Cilíndrica Ø 6,0x 85,0x C14
06.5.1814.260E	Broca Cilíndrica Ø 6,0x 90,0x C14
06.5.1914.260E	Broca Cilíndrica Ø 6,0x 95,0x C14
06.5.2014.260E	Broca Cilíndrica Ø 6,0x 100,0x C14
06.5.2114.260E	Broca Cilíndrica Ø 6,0x 105,0x C14
06.5.2214.260E	Broca Cilíndrica Ø 6,0x 110,0x C14
06.5.2314.260E	Broca Cilíndrica Ø 6,0x 115,0x C14
06.5.2414.260E	Broca Cilíndrica Ø 6,0x 120,0x C14
06.5.2514.260E	Broca Cilíndrica Ø 6,0x 125,0x C14

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

06.5.2614.260E	Broca Cilíndrica Ø 6,0x 130,0x C14
06.5.2714.260E	Broca Cilíndrica Ø 6,0x 135,0x C14
06.5.2814.260E	Broca Cilíndrica Ø 6,0x 140,0x C14
06.5.2914.260E	Broca Cilíndrica Ø 6,0x 145,0x C14
06.5.3014.260E	Broca Cilíndrica Ø 6,0x 150,0x C14
06.5.1016.270E	Broca Cilíndrica Ø 7,0x 50,0x C16
06.5.1116.270E	Broca Cilíndrica Ø 7,0x 55,0x C16
06.5.1216.270E	Broca Cilíndrica Ø 7,0x 60,0x C16
06.5.1316.270E	Broca Cilíndrica Ø 7,0x 65,0x C16
06.5.1416.270E	Broca Cilíndrica Ø 7,0x 70,0x C16
06.5.1516.270E	Broca Cilíndrica Ø 7,0x 75,0x C16
06.5.1616.270E	Broca Cilíndrica Ø 7,0x 80,0x C16
06.5.1716.270E	Broca Cilíndrica Ø 7,0x 85,0x C16
06.5.1816.270E	Broca Cilíndrica Ø 7,0x 90,0x C16
06.5.1916.270E	Broca Cilíndrica Ø 7,0x 95,0x C16
06.5.2016.270E	Broca Cilíndrica Ø 7,0x 100,0x C16
06.5.2116.270E	Broca Cilíndrica Ø 7,0x 105,0x C16
06.5.2216.270E	Broca Cilíndrica Ø 7,0x 110,0x C16
06.5.2316.270E	Broca Cilíndrica Ø 7,0x 115,0x C16
06.5.2416.270E	Broca Cilíndrica Ø 7,0x 120,0x C16
06.5.2516.270E	Broca Cilíndrica Ø 7,0x 125,0x C16
06.5.2616.270E	Broca Cilíndrica Ø 7,0x 130,0x C16
06.5.2716.270E	Broca Cilíndrica Ø 7,0x 135,0x C16
06.5.2816.270E	Broca Cilíndrica Ø 7,0x 140,0x C16
06.5.2916.270E	Broca Cilíndrica Ø 7,0x 145,0x C16
06.5.3016.270E	Broca Cilíndrica Ø 7,0x 150,0x C16
06.5.1016.280E	Broca Cilíndrica Ø 8,0x 50,0x C16
06.5.1116.280E	Broca Cilíndrica Ø 8,0x 55,0x C16
06.5.1216.280E	Broca Cilíndrica Ø 8,0x 60,0x C16
06.5.1316.280E	Broca Cilíndrica Ø 8,0x 65,0x C16
06.5.1416.280E	Broca Cilíndrica Ø 8,0x 70,0x C16
06.5.1516.280E	Broca Cilíndrica Ø 8,0x 75,0x C16
06.5.1616.280E	Broca Cilíndrica Ø 8,0x 80,0x C16
06.5.1716.280E	Broca Cilíndrica Ø 8,0x 85,0x C16
06.5.1816.280E	Broca Cilíndrica Ø 8,0x 90,0x C16
06.5.1916.280E	Broca Cilíndrica Ø 8,0x 95,0x C16
06.5.2016.280E	Broca Cilíndrica Ø 8,0x 100,0x C16
06.5.2116.280E	Broca Cilíndrica Ø 8,0x 105,0x C16
06.5.2216.280E	Broca Cilíndrica Ø 8,0x 110,0x C16
06.5.2316.280E	Broca Cilíndrica Ø 8,0x 115,0x C16
06.5.2416.280E	Broca Cilíndrica Ø 8,0x 120,0x C16
06.5.2516.280E	Broca Cilíndrica Ø 8,0x 125,0x C16
06.5.2616.280E	Broca Cilíndrica Ø 8,0x 130,0x C16
06.5.2716.280E	Broca Cilíndrica Ø 8,0x 135,0x C16
06.5.2816.280E	Broca Cilíndrica Ø 8,0x 140,0x C16
06.5.2916.280E	Broca Cilíndrica Ø 8,0x 145,0x C16
06.5.3016.280E	Broca Cilíndrica Ø 8,0x 150,0x C16

BROCA CORTANTE BUD ESTÉRIL	
	
CODIGO	DESCRIÇÃO
06.5.1010.340E	Broca Cortante BUD Ø 4,0x 50,0x C10
06.5.1110.340E	Broca Cortante BUD Ø 4,0x 55,0x C10
06.5.1210.340E	Broca Cortante BUD Ø 4,0x 60,0x C10
06.5.1310.340E	Broca Cortante BUD Ø 4,0x 65,0x C10
06.5.1410.340E	Broca Cortante BUD Ø 4,0x 70,0x C10
06.5.1510.340E	Broca Cortante BUD Ø 4,0x 75,0x C10
06.5.1610.340E	Broca Cortante BUD Ø 4,0x 80,0x C10
06.5.1710.340E	Broca Cortante BUD Ø 4,0x 85,0x C10
06.5.1810.340E	Broca Cortante BUD Ø 4,0x 90,0x C10
06.5.1910.340E	Broca Cortante BUD Ø 4,0x 95,0x C10
06.5.2010.340E	Broca Cortante BUD Ø 4,0x 100,0x C10
06.5.2110.340E	Broca Cortante BUD Ø 4,0x 105,0x C10
06.5.2210.340E	Broca Cortante BUD Ø 4,0x 110,0x C10
06.5.2310.340E	Broca Cortante BUD Ø 4,0x 115,0x C10
06.5.2410.340E	Broca Cortante BUD Ø 4,0x 120,0x C10
06.5.2510.340E	Broca Cortante BUD Ø 4,0x 125,0x C10
06.5.2610.340E	Broca Cortante BUD Ø 4,0x 130,0x C10
06.5.2710.340E	Broca Cortante BUD Ø 4,0x 135,0x C10
06.5.2810.340E	Broca Cortante BUD Ø 4,0x 140,0x C10
06.5.2910.340E	Broca Cortante BUD Ø 4,0x 145,0x C10
06.5.3010.340E	Broca Cortante BUD Ø 4,0x 150,0x C10
06.5.1012.350E	Broca Cortante BUD Ø 5,0x 50,0x C12
06.5.1112.350E	Broca Cortante BUD Ø 5,0x 55,0x C12
06.5.1212.350E	Broca Cortante BUD Ø 5,0x 60,0x C12
06.5.1312.350E	Broca Cortante BUD Ø 5,0x 65,0x C12
06.5.1412.350E	Broca Cortante BUD Ø 5,0x 70,0x C12
06.5.1512.350E	Broca Cortante BUD Ø 5,0x 75,0x C12
06.5.1612.350E	Broca Cortante BUD Ø 5,0x 80,0x C12
06.5.1712.350E	Broca Cortante BUD Ø 5,0x 85,0x C12
06.5.1812.350E	Broca Cortante BUD Ø 5,0x 90,0x C12
06.5.1912.350E	Broca Cortante BUD Ø 5,0x 95,0x C12
06.5.2012.350E	Broca Cortante BUD Ø 5,0x 100,0x C12
06.5.2112.350E	Broca Cortante BUD Ø 5,0x 105,0x C12
06.5.2212.350E	Broca Cortante BUD Ø 5,0x 110,0x C12
06.5.2312.350E	Broca Cortante BUD Ø 5,0x 115,0x C12
06.5.2412.350E	Broca Cortante BUD Ø 5,0x 120,0x C12
06.5.2512.350E	Broca Cortante BUD Ø 5,0x 125,0x C12
06.5.2612.350E	Broca Cortante BUD Ø 5,0x 130,0x C12
06.5.2712.350E	Broca Cortante BUD Ø 5,0x 135,0x C12
06.5.2812.350E	Broca Cortante BUD Ø 5,0x 140,0x C12
06.5.2912.350E	Broca Cortante BUD Ø 5,0x 145,0x C12
06.5.3012.350E	Broca Cortante BUD Ø 5,0x 150,0x C12
06.5.1014.360E	Broca Cortante BUD Ø 6,0x 50,0x C14

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003


06.5.1114.360E	Broca Cortante BUD Ø 6,0x 55,0x C14
06.5.1214.360E	Broca Cortante BUD Ø 6,0x 60,0x C14
06.5.1314.360E	Broca Cortante BUD Ø 6,0x 65,0x C14
06.5.1414.360E	Broca Cortante BUD Ø 6,0x 70,0x C14
06.5.1514.360E	Broca Cortante BUD Ø 6,0x 75,0x C14
06.5.1614.360E	Broca Cortante BUD Ø 6,0x 80,0x C14
06.5.1714.360E	Broca Cortante BUD Ø 6,0x 85,0x C14
06.5.1814.360E	Broca Cortante BUD Ø 6,0x 90,0x C14
06.5.1914.360E	Broca Cortante BUD Ø 6,0x 95,0x C14
06.5.2014.360E	Broca Cortante BUD Ø 6,0x 100,0x C14
06.5.2114.360E	Broca Cortante BUD Ø 6,0x 105,0x C14
06.5.2214.360E	Broca Cortante BUD Ø 6,0x 110,0x C14
06.5.2314.360E	Broca Cortante BUD Ø 6,0x 115,0x C14
06.5.2414.360E	Broca Cortante BUD Ø 6,0x 120,0x C14
06.5.2514.360E	Broca Cortante BUD Ø 6,0x 125,0x C14
06.5.2614.360E	Broca Cortante BUD Ø 6,0x 130,0x C14
06.5.2714.360E	Broca Cortante BUD Ø 6,0x 135,0x C14
06.5.2814.360E	Broca Cortante BUD Ø 6,0x 140,0x C14
06.5.2914.360E	Broca Cortante BUD Ø 6,0x 145,0x C14
06.5.3014.360E	Broca Cortante BUD Ø 6,0x 150,0x C14
06.5.1016.370E	Broca Cortante BUD Ø 7,0x 50,0x C16
06.5.1116.370E	Broca Cortante BUD Ø 7,0x 55,0x C16
06.5.1216.370E	Broca Cortante BUD Ø 7,0x 60,0x C16
06.5.1316.370E	Broca Cortante BUD Ø 7,0x 65,0x C16
06.5.1416.370E	Broca Cortante BUD Ø 7,0x 70,0x C16
06.5.1516.370E	Broca Cortante BUD Ø 7,0x 75,0x C16
06.5.1616.370E	Broca Cortante BUD Ø 7,0x 80,0x C16
06.5.1716.370E	Broca Cortante BUD Ø 7,0x 85,0x C16
06.5.1816.370E	Broca Cortante BUD Ø 7,0x 90,0x C16
06.5.1916.370E	Broca Cortante BUD Ø 7,0x 95,0x C16
06.5.2016.370E	Broca Cortante BUD Ø 7,0x 100,0x C16
06.5.2116.370E	Broca Cortante BUD Ø 7,0x 105,0x C16
06.5.2216.370E	Broca Cortante BUD Ø 7,0x 110,0x C16
06.5.2316.370E	Broca Cortante BUD Ø 7,0x 115,0x C16
06.5.2416.370E	Broca Cortante BUD Ø 7,0x 120,0x C16
06.5.2516.370E	Broca Cortante BUD Ø 7,0x 125,0x C16
06.5.2616.370E	Broca Cortante BUD Ø 7,0x 130,0x C16
06.5.2716.370E	Broca Cortante BUD Ø 7,0x 135,0x C16
06.5.2816.370E	Broca Cortante BUD Ø 7,0x 140,0x C16
06.5.2916.370E	Broca Cortante BUD Ø 7,0x 145,0x C16
06.5.3016.370E	Broca Cortante BUD Ø 7,0x 150,0x C16
06.5.1016.380E	Broca Cortante BUD Ø 8,0x 50,0x C16
06.5.1116.380E	Broca Cortante BUD Ø 8,0x 55,0x C16
06.5.1216.380E	Broca Cortante BUD Ø 8,0x 60,0x C16
06.5.1316.380E	Broca Cortante BUD Ø 8,0x 65,0x C16
06.5.1416.380E	Broca Cortante BUD Ø 8,0x 70,0x C16
06.5.1516.380E	Broca Cortante BUD Ø 8,0x 75,0x C16
06.5.1616.380E	Broca Cortante BUD Ø 8,0x 80,0x C16
06.5.1716.380E	Broca Cortante BUD Ø 8,0x 85,0x C16


06.5.1816.380E	Broca Cortante BUD Ø 8,0x 90,0x C16
06.5.1916.380E	Broca Cortante BUD Ø 8,0x 95,0x C16
06.5.2016.380E	Broca Cortante BUD Ø 8,0x 100,0x C16
06.5.2116.380E	Broca Cortante BUD Ø 8,0x 105,0x C16
06.5.2216.380E	Broca Cortante BUD Ø 8,0x 110,0x C16
06.5.2316.380E	Broca Cortante BUD Ø 8,0x 115,0x C16
06.5.2416.380E	Broca Cortante BUD Ø 8,0x 120,0x C16
06.5.2516.380E	Broca Cortante BUD Ø 8,0x 125,0x C16
06.5.2616.380E	Broca Cortante BUD Ø 8,0x 130,0x C16
06.5.2716.380E	Broca Cortante BUD Ø 8,0x 135,0x C16
06.5.2816.380E	Broca Cortante BUD Ø 8,0x 140,0x C16
06.5.2916.380E	Broca Cortante BUD Ø 8,0x 145,0x C16
06.5.3016.380E	Broca Cortante BUD Ø 8,0x 150,0x C16

BROCA CORTANTE BUD CROSS-CUT ESTÉRIL	
	
CODIGO	DESCRIÇÃO
06.5.1010.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 50,0x C10
06.5.1110.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 55,0x C10
06.5.1210.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 60,0x C10
06.5.1310.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 65,0x C10
06.5.1410.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 70,0x C10
06.5.1510.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 75,0x C10
06.5.1610.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 80,0x C10
06.5.1710.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 85,0x C10
06.5.1810.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 90,0x C10
06.5.1910.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 95,0x C10
06.5.2010.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 100,0x C10
06.5.2110.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 105,0x C10
06.5.2210.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 110,0x C10
06.5.2310.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 115,0x C10
06.5.2410.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 120,0x C10
06.5.2510.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 125,0x C10
06.5.2610.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 130,0x C10
06.5.2710.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 135,0x C10
06.5.2810.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 140,0x C10
06.5.2910.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 145,0x C10
06.5.3010.440E	Broca Cortante BUD Cross-Cut Ø 4,0x 150,0x C10
06.5.1012.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 50,0x C12
06.5.1112.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 55,0x C12
06.5.1212.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 60,0x C12
06.5.1312.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 65,0x C12
06.5.1412.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 70,0x C12
06.5.1512.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 75,0x C12
06.5.1612.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 80,0x C12
06.5.1712.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 85,0x C12

06.5.1812.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 90,0x C12
06.5.1912.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 95,0x C12
06.5.2012.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 100,0x C12
06.5.2112.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 105,0x C12
06.5.2212.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 110,0x C12
06.5.2312.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 115,0x C12
06.5.2412.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 120,0x C12
06.5.2512.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 125,0x C12
06.5.2612.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 130,0x C12
06.5.2712.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 135,0x C12
06.5.2812.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 140,0x C12
06.5.2912.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 145,0x C12
06.5.3012.450E	Broca Cortante BUD Cross-Cut Ø 5,0x 150,0x C12
06.5.1014.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 50,0x C14
06.5.1114.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 55,0x C14
06.5.1214.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 60,0x C14
06.5.1314.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 65,0x C14
06.5.1414.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 70,0x C14
06.5.1514.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 75,0x C14
06.5.1614.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 80,0x C14
06.5.1714.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 85,0x C14
06.5.1814.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 90,0x C14
06.5.1914.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 95,0x C14
06.5.2014.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 100,0x C14
06.5.2114.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 105,0x C14
06.5.2214.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 110,0x C14
06.5.2314.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 115,0x C14
06.5.2414.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 120,0x C14
06.5.2514.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 125,0x C14
06.5.2614.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 130,0x C14
06.5.2714.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 135,0x C14
06.5.2814.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 140,0x C14
06.5.2914.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 145,0x C14
06.5.3014.460E	Broca Cortante BUD Cross-Cut Ø 6,0x 150,0x C14
06.5.1016.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 50,0x C16
06.5.1116.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 55,0x C16
06.5.1216.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 60,0x C16
06.5.1316.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 65,0x C16
06.5.1416.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 70,0x C16
06.5.1516.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 75,0x C16
06.5.1616.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 80,0x C16
06.5.1716.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 85,0x C16
06.5.1816.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 90,0x C16
06.5.1916.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 95,0x C16
06.5.2016.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 100,0x C16
06.5.2116.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 105,0x C16
06.5.2216.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 110,0x C16
06.5.2316.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 115,0x C16
06.5.2416.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 120,0x C16

06.5.2516.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 125,0x C16
06.5.2616.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 130,0x C16
06.5.2716.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 135,0x C16
06.5.2816.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 140,0x C16
06.5.2916.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 145,0x C16
06.5.3016.470E	Broca Cortante BUD Cross-Cut Ø 7,0x 150,0x C16
06.5.1016.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 50,0x C16
06.5.1116.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 55,0x C16
06.5.1216.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 60,0x C16
06.5.1316.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 65,0x C16
06.5.1416.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 70,0x C16
06.5.1516.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 75,0x C16
06.5.1616.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 80,0x C16
06.5.1716.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 85,0x C16
06.5.1816.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 90,0x C16
06.5.1916.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 95,0x C16
06.5.2016.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 100,0x C16
06.5.2116.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 105,0x C16
06.5.2216.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 110,0x C16
06.5.2316.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 115,0x C16
06.5.2416.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 120,0x C16
06.5.2516.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 125,0x C16
06.5.2616.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 130,0x C16
06.5.2716.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 135,0x C16
06.5.2816.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 140,0x C16
06.5.2916.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 145,0x C16
06.5.3016.480E	Broca Cortante BUD Cross-Cut Ø 8,0x 150,0x C16

BROCA AUTOMÁTICA PARA CRANIOTOMIA ESTÉRIL	
	
12.5.0101.110E	Broca Automática para Craniotomia Ø 14.00/11.00
12.5.0101.080E	Broca Automática para Craniotomia Ø 11.00/8.00
12.5.0101.060E	Broca Automática para Craniotomia Ø 9.00/6.00

FRESA CRANIÓTOMO ESTÉRIL	
	
12.5.0701.020E	Fresa Craniotomo Ø2.00
12.5.0701.023E	Fresa Craniotomo Ø2.30
12.5.0701.025E	Fresa Craniotomo Ø2.50
12.5.0701.027E	Fresa Craniotomo Ø2.70

12.5.0701.030E	Fresa Craniotomo Ø3.00
12.5.0701.032E	Fresa Craniotomo Ø3.20
12.5.0701.037E	Fresa Craniotomo Ø3.70

LÂMINA REMOVEDORA DE ACETÁBULO ESTÉRIL	
	
03.5.2101.400E	Lâmina removedora de Acetábulo Inicial Ø 40.00
03.5.2102.400E	Lâmina removedora de Acetábulo Final Ø 40.00
03.5.2101.410E	Lâmina removedora de Acetábulo Inicial Ø 41.00
03.5.2102.410E	Lâmina removedora de Acetábulo Final Ø 41.00
03.5.2101.420E	Lâmina removedora de Acetábulo Inicial Ø 42.00
03.5.2102.420E	Lâmina removedora de Acetábulo Final Ø 42.00
03.5.2101.430E	Lâmina removedora de Acetábulo Inicial Ø 43.00
03.5.2102.430E	Lâmina removedora de Acetábulo Final Ø 43.00
03.5.2101.440E	Lâmina removedora de Acetábulo Inicial Ø 44.00
03.5.2102.440E	Lâmina removedora de Acetábulo Final Ø 44.00
03.5.2101.450E	Lâmina removedora de Acetábulo Inicial Ø 45.00
03.5.2102.450E	Lâmina removedora de Acetábulo Final Ø 45.00
03.5.2101.460E	Lâmina removedora de Acetábulo Inicial Ø 46.00
03.5.2102.460E	Lâmina removedora de Acetábulo Final Ø 46.00
03.5.2101.470E	Lâmina removedora de Acetábulo Inicial Ø 47.00
03.5.2102.470E	Lâmina removedora de Acetábulo Final Ø 47.00
03.5.2101.480E	Lâmina removedora de Acetábulo Inicial Ø 48.00
03.5.2102.480E	Lâmina removedora de Acetábulo Final Ø 48.00
03.5.2101.490E	Lâmina removedora de Acetábulo Inicial Ø 49.00
03.5.2102.490E	Lâmina removedora de Acetábulo Final Ø 49.00
03.5.2101.500E	Lâmina removedora de Acetábulo Inicial Ø 50.00
03.5.2102.500E	Lâmina removedora de Acetábulo Final Ø 50.00
03.5.2101.510E	Lâmina removedora de Acetábulo Inicial Ø 51.00
03.5.2102.510E	Lâmina removedora de Acetábulo Final Ø 51.00
03.5.2101.520E	Lâmina removedora de Acetábulo Inicial Ø 52.00
03.5.2102.520E	Lâmina removedora de Acetábulo Final Ø 52.00
03.5.2101.530E	Lâmina removedora de Acetábulo Inicial Ø 53.00
03.5.2102.530E	Lâmina removedora de Acetábulo Final Ø 53.00
03.5.2101.540E	Lâmina removedora de Acetábulo Inicial Ø 54.00
03.5.2102.540E	Lâmina removedora de Acetábulo Final Ø 54.00
03.5.2101.550E	Lâmina removedora de Acetábulo Inicial Ø 55.00
03.5.2102.550E	Lâmina removedora de Acetábulo Final Ø 55.00
03.5.2101.560E	Lâmina removedora de Acetábulo Inicial Ø 56.00
03.5.2102.560E	Lâmina removedora de Acetábulo Final Ø 56.00
03.5.2101.570E	Lâmina removedora de Acetábulo Inicial Ø 57.00
03.5.2102.570E	Lâmina removedora de Acetábulo Final Ø 57.00
03.5.2101.580E	Lâmina removedora de Acetábulo Inicial Ø 58.00
03.5.2102.580E	Lâmina removedora de Acetábulo Final Ø 58.00

03.5.2101.590E	Lâmina removedora de Acetábulo Inicial Ø 59.00
03.5.2102.590E	Lâmina removedora de Acetábulo Final Ø 59.00
03.5.2101.600E	Lâmina removedora de Acetábulo Inicial Ø 60.00
03.5.2102.600E	Lâmina removedora de Acetábulo Final Ø 60.00
03.5.2101.610E	Lâmina removedora de Acetábulo Inicial Ø 61.00
03.5.2102.610E	Lâmina removedora de Acetábulo Final Ø 61.00
03.5.2101.620E	Lâmina removedora de Acetábulo Inicial Ø 62.00
03.5.2102.620E	Lâmina removedora de Acetábulo Final Ø 62.00
03.5.2101.630E	Lâmina removedora de Acetábulo Inicial Ø 63.00
03.5.2102.630E	Lâmina removedora de Acetábulo Final Ø 63.00
03.5.2101.640E	Lâmina removedora de Acetábulo Inicial Ø 64.00
03.5.2102.640E	Lâmina removedora de Acetábulo Final Ø 64.00
03.5.2101.650E	Lâmina removedora de Acetábulo Inicial Ø 65.00
03.5.2102.650E	Lâmina removedora de Acetábulo Final Ø 65.00
03.5.2101.660E	Lâmina removedora de Acetábulo Inicial Ø 66.00
03.5.2102.660E	Lâmina removedora de Acetábulo Final Ø 66.00
03.5.2101.670E	Lâmina removedora de Acetábulo Inicial Ø 67.00
03.5.2102.670E	Lâmina removedora de Acetábulo Final Ø 67.00
03.5.2101.680E	Lâmina removedora de Acetábulo Inicial Ø 68.00
03.5.2102.680E	Lâmina removedora de Acetábulo Final Ø 68.00
03.5.2101.690E	Lâmina removedora de Acetábulo Inicial Ø 69.00
03.5.2102.690E	Lâmina removedora de Acetábulo Final Ø 69.00
03.5.2101.700E	Lâmina removedora de Acetábulo Inicial Ø 70.00
03.5.2102.700E	Lâmina removedora de Acetábulo Final Ø 70.00
03.5.2101.710E	Lâmina removedora de Acetábulo Inicial Ø 71.00
03.5.2102.710E	Lâmina removedora de Acetábulo Final Ø 71.00
03.5.2101.720E	Lâmina removedora de Acetábulo Inicial Ø 72.00
03.5.2102.720E	Lâmina removedora de Acetábulo Final Ø 72.00
03.5.2101.730E	Lâmina removedora de Acetábulo Inicial Ø 73.00
03.5.2102.730E	Lâmina removedora de Acetábulo Final Ø 73.00
03.5.2101.740E	Lâmina removedora de Acetábulo Inicial Ø 74.00
03.5.2102.740E	Lâmina removedora de Acetábulo Final Ø 74.00
03.5.2101.750E	Lâmina removedora de Acetábulo Inicial Ø 75.00
03.5.2102.750E	Lâmina removedora de Acetábulo Final Ø 75.00
03.5.2101.760E	Lâmina removedora de Acetábulo Inicial Ø 76.00
03.5.2102.760E	Lâmina removedora de Acetábulo Final Ø 76.00
03.5.2101.770E	Lâmina removedora de Acetábulo Inicial Ø 77.00
03.5.2102.770E	Lâmina removedora de Acetábulo Final Ø 77.00
03.5.2101.780E	Lâmina removedora de Acetábulo Inicial Ø 78.00
03.5.2102.780E	Lâmina removedora de Acetábulo Final Ø 78.00
03.5.2101.790E	Lâmina removedora de Acetábulo Inicial Ø 79.00
03.5.2102.790E	Lâmina removedora de Acetábulo Final Ø 79.00
03.5.2101.800E	Lâmina removedora de Acetábulo Inicial Ø 80.00
03.5.2102.800E	Lâmina removedora de Acetábulo Final Ø 80.00

BROCA BUCO MAXILO ESTÉRIL

	
CODIGO	DESCRIÇÃO
05.5.0001.109E	Broca Ø 0,9x Stop 4mmx Comprimento 50mm
05.5.0002.109E	Broca Ø 0,9x Stop 4mmx Comprimento 55mm
05.5.0003.109E	Broca Ø 0,9x Stop 4mmx Comprimento 60mm
05.5.0004.109E	Broca Ø 0,9x Stop 4mmx Comprimento 65mm
05.5.0005.109E	Broca Ø 0,9x Stop 4mmx Comprimento 70mm
05.5.0006.109E	Broca Ø 0,9x Stop 4mmx Comprimento 75mm
05.5.0007.109E	Broca Ø 0,9x Stop 4mmx Comprimento 80mm
05.5.0008.109E	Broca Ø 0,9x Stop 4mmx Comprimento 85mm
05.5.0009.109E	Broca Ø 0,9x Stop 4mmx Comprimento 90mm
05.5.0010.109E	Broca Ø 0,9x Stop 4mmx Comprimento 95mm
05.5.0011.109E	Broca Ø 0,9x Stop 4mmx Comprimento 100mm
05.5.0012.109E	Broca Ø 0,9x Stop 4mmx Comprimento 105mm
05.5.0013.109E	Broca Ø 0,9x Stop 4mmx Comprimento 110mm
05.5.0014.109E	Broca Ø 0,9x Stop 4mmx Comprimento 115mm
05.5.0015.109E	Broca Ø 0,9x Stop 4mmx Comprimento 120mm
05.5.0016.109E	Broca Ø 0,9x Stop 4mmx Comprimento 125mm
05.5.0017.109E	Broca Ø 0,9x Stop 4mmx Comprimento 130mm
05.5.0018.109E	Broca Ø 0,9x Stop 4mmx Comprimento 135mm
05.5.0019.109E	Broca Ø 0,9x Stop 4mmx Comprimento 140mm
05.5.0020.109E	Broca Ø 0,9x Stop 4mmx Comprimento 145mm
05.5.0021.109E	Broca Ø 0,9x Stop 4mmx Comprimento 150mm
05.5.0022.109E	Broca Ø 0,9x Stop 5mmx Comprimento 50mm
05.5.0023.109E	Broca Ø 0,9x Stop 5mmx Comprimento 55mm
05.5.0024.109E	Broca Ø 0,9x Stop 5mmx Comprimento 60mm
05.5.0025.109E	Broca Ø 0,9x Stop 5mmx Comprimento 65mm
05.5.0026.109E	Broca Ø 0,9x Stop 5mmx Comprimento 70mm
05.5.0027.109E	Broca Ø 0,9x Stop 5mmx Comprimento 75mm
05.5.0028.109E	Broca Ø 0,9x Stop 5mmx Comprimento 80mm
05.5.0029.109E	Broca Ø 0,9x Stop 5mmx Comprimento 85mm
05.5.0030.109E	Broca Ø 0,9x Stop 5mmx Comprimento 90mm
05.5.0031.109E	Broca Ø 0,9x Stop 5mmx Comprimento 95mm
05.5.0032.109E	Broca Ø 0,9x Stop 5mmx Comprimento 100mm
05.5.0033.109E	Broca Ø 0,9x Stop 5mmx Comprimento 105mm
05.5.0034.109E	Broca Ø 0,9x Stop 5mmx Comprimento 110mm
05.5.0035.109E	Broca Ø 0,9x Stop 5mmx Comprimento 115mm
05.5.0036.109E	Broca Ø 0,9x Stop 5mmx Comprimento 120mm
05.5.0037.109E	Broca Ø 0,9x Stop 5mmx Comprimento 125mm
05.5.0038.109E	Broca Ø 0,9x Stop 5mmx Comprimento 130mm
05.5.0039.109E	Broca Ø 0,9x Stop 5mmx Comprimento 135mm
05.5.0040.109E	Broca Ø 0,9x Stop 5mmx Comprimento 140mm
05.5.0041.109E	Broca Ø 0,9x Stop 5mmx Comprimento 145mm
05.5.0042.109E	Broca Ø 0,9x Stop 5mmx Comprimento 150mm
05.5.0043.109E	Broca Ø 0,9x Stop 6mmx Comprimento 50mm
05.5.0044.109E	Broca Ø 0,9x Stop 6mmx Comprimento 55mm
05.5.0045.109E	Broca Ø 0,9x Stop 6mmx Comprimento 60mm

05.5.0046.109E	Broca Ø 0,9x Stop 6mmx Comprimento 65mm
05.5.0047.109E	Broca Ø 0,9x Stop 6mmx Comprimento 70mm
05.5.0048.109E	Broca Ø 0,9x Stop 6mmx Comprimento 75mm
05.5.0049.109E	Broca Ø 0,9x Stop 6mmx Comprimento 80mm
05.5.0050.109E	Broca Ø 0,9x Stop 6mmx Comprimento 85mm
05.5.0051.109E	Broca Ø 0,9x Stop 6mmx Comprimento 90mm
05.5.0052.109E	Broca Ø 0,9x Stop 6mmx Comprimento 95mm
05.5.0053.109E	Broca Ø 0,9x Stop 6mmx Comprimento 100mm
05.5.0054.109E	Broca Ø 0,9x Stop 6mmx Comprimento 105mm
05.5.0055.109E	Broca Ø 0,9x Stop 6mmx Comprimento 110mm
05.5.0056.109E	Broca Ø 0,9x Stop 6mmx Comprimento 115mm
05.5.0057.109E	Broca Ø 0,9x Stop 6mmx Comprimento 120mm
05.5.0058.109E	Broca Ø 0,9x Stop 6mmx Comprimento 125mm
05.5.0059.109E	Broca Ø 0,9x Stop 6mmx Comprimento 130mm
05.5.0060.109E	Broca Ø 0,9x Stop 6mmx Comprimento 135mm
05.5.0061.109E	Broca Ø 0,9x Stop 6mmx Comprimento 140mm
05.5.0062.109E	Broca Ø 0,9x Stop 6mmx Comprimento 145mm
05.5.0063.109E	Broca Ø 0,9x Stop 6mmx Comprimento 150mm
05.5.0064.109E	Broca Ø 0,9x Stop 7mmx Comprimento 50mm
05.5.0065.109E	Broca Ø 0,9x Stop 7mmx Comprimento 55mm
05.5.0066.109E	Broca Ø 0,9x Stop 7mmx Comprimento 60mm
05.5.0067.109E	Broca Ø 0,9x Stop 7mmx Comprimento 65mm
05.5.0068.109E	Broca Ø 0,9x Stop 7mmx Comprimento 70mm
05.5.0069.109E	Broca Ø 0,9x Stop 7mmx Comprimento 75mm
05.5.0070.109E	Broca Ø 0,9x Stop 7mmx Comprimento 80mm
05.5.0071.109E	Broca Ø 0,9x Stop 7mmx Comprimento 85mm
05.5.0072.109E	Broca Ø 0,9x Stop 7mmx Comprimento 90mm
05.5.0073.109E	Broca Ø 0,9x Stop 7mmx Comprimento 95mm
05.5.0074.109E	Broca Ø 0,9x Stop 7mmx Comprimento 100mm
05.5.0075.109E	Broca Ø 0,9x Stop 7mmx Comprimento 105mm
05.5.0076.109E	Broca Ø 0,9x Stop 7mmx Comprimento 110mm
05.5.0077.109E	Broca Ø 0,9x Stop 7mmx Comprimento 115mm
05.5.0078.109E	Broca Ø 0,9x Stop 7mmx Comprimento 120mm
05.5.0079.109E	Broca Ø 0,9x Stop 7mmx Comprimento 125mm
05.5.0080.109E	Broca Ø 0,9x Stop 7mmx Comprimento 130mm
05.5.0081.109E	Broca Ø 0,9x Stop 7mmx Comprimento 135mm
05.5.0082.109E	Broca Ø 0,9x Stop 7mmx Comprimento 140mm
05.5.0083.109E	Broca Ø 0,9x Stop 7mmx Comprimento 145mm
05.5.0084.109E	Broca Ø 0,9x Stop 7mmx Comprimento 150mm
05.5.0085.109E	Broca Ø 0,9x Stop 8mmx Comprimento 50mm
05.5.0086.109E	Broca Ø 0,9x Stop 8mmx Comprimento 55mm
05.5.0087.109E	Broca Ø 0,9x Stop 8mmx Comprimento 60mm
05.5.0088.109E	Broca Ø 0,9x Stop 8mmx Comprimento 65mm
05.5.0089.109E	Broca Ø 0,9x Stop 8mmx Comprimento 70mm
05.5.0090.109E	Broca Ø 0,9x Stop 8mmx Comprimento 75mm
05.5.0091.109E	Broca Ø 0,9x Stop 8mmx Comprimento 80mm
05.5.0092.109E	Broca Ø 0,9x Stop 8mmx Comprimento 85mm
05.5.0093.109E	Broca Ø 0,9x Stop 8mmx Comprimento 90mm
05.5.0094.109E	Broca Ø 0,9x Stop 8mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0095.109E	Broca Ø 0,9x Stop 8mmx Comprimento 100mm
05.5.0096.109E	Broca Ø 0,9x Stop 8mmx Comprimento 105mm
05.5.0097.109E	Broca Ø 0,9x Stop 8mmx Comprimento 110mm
05.5.0098.109E	Broca Ø 0,9x Stop 8mmx Comprimento 115mm
05.5.0099.109E	Broca Ø 0,9x Stop 8mmx Comprimento 120mm
05.5.0100.109E	Broca Ø 0,9x Stop 8mmx Comprimento 125mm
05.5.0101.109E	Broca Ø 0,9x Stop 8mmx Comprimento 130mm
05.5.0102.109E	Broca Ø 0,9x Stop 8mmx Comprimento 135mm
05.5.0103.109E	Broca Ø 0,9x Stop 8mmx Comprimento 140mm
05.5.0104.109E	Broca Ø 0,9x Stop 8mmx Comprimento 145mm
05.5.0105.109E	Broca Ø 0,9x Stop 8mmx Comprimento 150mm
05.5.0106.109E	Broca Ø 0,9x Stop 9mmx Comprimento 50mm
05.5.0107.109E	Broca Ø 0,9x Stop 9mmx Comprimento 55mm
05.5.0108.109E	Broca Ø 0,9x Stop 9mmx Comprimento 60mm
05.5.0109.109E	Broca Ø 0,9x Stop 9mmx Comprimento 65mm
05.5.0110.109E	Broca Ø 0,9x Stop 9mmx Comprimento 70mm
05.5.0111.109E	Broca Ø 0,9x Stop 9mmx Comprimento 75mm
05.5.0112.109E	Broca Ø 0,9x Stop 9mmx Comprimento 80mm
05.5.0113.109E	Broca Ø 0,9x Stop 9mmx Comprimento 85mm
05.5.0114.109E	Broca Ø 0,9x Stop 9mmx Comprimento 90mm
05.5.0115.109E	Broca Ø 0,9x Stop 9mmx Comprimento 95mm
05.5.0116.109E	Broca Ø 0,9x Stop 9mmx Comprimento 100mm
05.5.0117.109E	Broca Ø 0,9x Stop 9mmx Comprimento 105mm
05.5.0118.109E	Broca Ø 0,9x Stop 9mmx Comprimento 110mm
05.5.0119.109E	Broca Ø 0,9x Stop 9mmx Comprimento 115mm
05.5.0120.109E	Broca Ø 0,9x Stop 9mmx Comprimento 120mm
05.5.0121.109E	Broca Ø 0,9x Stop 9mmx Comprimento 125mm
05.5.0122.109E	Broca Ø 0,9x Stop 9mmx Comprimento 130mm
05.5.0123.109E	Broca Ø 0,9x Stop 9mmx Comprimento 135mm
05.5.0124.109E	Broca Ø 0,9x Stop 9mmx Comprimento 140mm
05.5.0125.109E	Broca Ø 0,9x Stop 9mmx Comprimento 145mm
05.5.0126.109E	Broca Ø 0,9x Stop 9mmx Comprimento 150mm
05.5.0127.109E	Broca Ø 0,9x Stop 10mmx Comprimento 50mm
05.5.0128.109E	Broca Ø 0,9x Stop 10mmx Comprimento 55mm
05.5.0129.109E	Broca Ø 0,9x Stop 10mmx Comprimento 60mm
05.5.0130.109E	Broca Ø 0,9x Stop 10mmx Comprimento 65mm
05.5.0131.109E	Broca Ø 0,9x Stop 10mmx Comprimento 70mm
05.5.0132.109E	Broca Ø 0,9x Stop 10mmx Comprimento 75mm
05.5.0133.109E	Broca Ø 0,9x Stop 10mmx Comprimento 80mm
05.5.0134.109E	Broca Ø 0,9x Stop 10mmx Comprimento 85mm
05.5.0135.109E	Broca Ø 0,9x Stop 10mmx Comprimento 90mm
05.5.0136.109E	Broca Ø 0,9x Stop 10mmx Comprimento 95mm
05.5.0137.109E	Broca Ø 0,9x Stop 10mmx Comprimento 100mm
05.5.0138.109E	Broca Ø 0,9x Stop 10mmx Comprimento 105mm
05.5.0139.109E	Broca Ø 0,9x Stop 10mmx Comprimento 110mm
05.5.0140.109E	Broca Ø 0,9x Stop 10mmx Comprimento 115mm
05.5.0141.109E	Broca Ø 0,9x Stop 10mmx Comprimento 120mm
05.5.0142.109E	Broca Ø 0,9x Stop 10mmx Comprimento 125mm
05.5.0143.109E	Broca Ø 0,9x Stop 10mmx Comprimento 130mm

05.5.0144.109E	Broca Ø 0,9x Stop 10mmx Comprimento 135mm
05.5.0145.109E	Broca Ø 0,9x Stop 10mmx Comprimento 140mm
05.5.0146.109E	Broca Ø 0,9x Stop 10mmx Comprimento 145mm
05.5.0147.109E	Broca Ø 0,9x Stop 10mmx Comprimento 150mm
05.5.0148.109E	Broca Ø 0,9x Stop 11mmx Comprimento 50mm
05.5.0149.109E	Broca Ø 0,9x Stop 11mmx Comprimento 55mm
05.5.0150.109E	Broca Ø 0,9x Stop 11mmx Comprimento 60mm
05.5.0151.109E	Broca Ø 0,9x Stop 11mmx Comprimento 65mm
05.5.0152.109E	Broca Ø 0,9x Stop 11mmx Comprimento 70mm
05.5.0153.109E	Broca Ø 0,9x Stop 11mmx Comprimento 75mm
05.5.0154.109E	Broca Ø 0,9x Stop 11mmx Comprimento 80mm
05.5.0155.109E	Broca Ø 0,9x Stop 11mmx Comprimento 85mm
05.5.0156.109E	Broca Ø 0,9x Stop 11mmx Comprimento 90mm
05.5.0157.109E	Broca Ø 0,9x Stop 11mmx Comprimento 95mm
05.5.0158.109E	Broca Ø 0,9x Stop 11mmx Comprimento 100mm
05.5.0159.109E	Broca Ø 0,9x Stop 11mmx Comprimento 105mm
05.5.0160.109E	Broca Ø 0,9x Stop 11mmx Comprimento 110mm
05.5.0161.109E	Broca Ø 0,9x Stop 11mmx Comprimento 115mm
05.5.0162.109E	Broca Ø 0,9x Stop 11mmx Comprimento 120mm
05.5.0163.109E	Broca Ø 0,9x Stop 11mmx Comprimento 125mm
05.5.0164.109E	Broca Ø 0,9x Stop 11mmx Comprimento 130mm
05.5.0165.109E	Broca Ø 0,9x Stop 11mmx Comprimento 135mm
05.5.0166.109E	Broca Ø 0,9x Stop 11mmx Comprimento 140mm
05.5.0167.109E	Broca Ø 0,9x Stop 11mmx Comprimento 145mm
05.5.0168.109E	Broca Ø 0,9x Stop 11mmx Comprimento 150mm
05.5.0169.109E	Broca Ø 0,9x Stop 12mmx Comprimento 50mm
05.5.0170.109E	Broca Ø 0,9x Stop 12mmx Comprimento 55mm
05.5.0171.109E	Broca Ø 0,9x Stop 12mmx Comprimento 60mm
05.5.0172.109E	Broca Ø 0,9x Stop 12mmx Comprimento 65mm
05.5.0173.109E	Broca Ø 0,9x Stop 12mmx Comprimento 70mm
05.5.0174.109E	Broca Ø 0,9x Stop 12mmx Comprimento 75mm
05.5.0175.109E	Broca Ø 0,9x Stop 12mmx Comprimento 80mm
05.5.0176.109E	Broca Ø 0,9x Stop 12mmx Comprimento 85mm
05.5.0177.109E	Broca Ø 0,9x Stop 12mmx Comprimento 90mm
05.5.0178.109E	Broca Ø 0,9x Stop 12mmx Comprimento 95mm
05.5.0179.109E	Broca Ø 0,9x Stop 12mmx Comprimento 100mm
05.5.0180.109E	Broca Ø 0,9x Stop 12mmx Comprimento 105mm
05.5.0181.109E	Broca Ø 0,9x Stop 12mmx Comprimento 110mm
05.5.0182.109E	Broca Ø 0,9x Stop 12mmx Comprimento 115mm
05.5.0183.109E	Broca Ø 0,9x Stop 12mmx Comprimento 120mm
05.5.0184.109E	Broca Ø 0,9x Stop 12mmx Comprimento 125mm
05.5.0185.109E	Broca Ø 0,9x Stop 12mmx Comprimento 130mm
05.5.0186.109E	Broca Ø 0,9x Stop 12mmx Comprimento 135mm
05.5.0187.109E	Broca Ø 0,9x Stop 12mmx Comprimento 140mm
05.5.0188.109E	Broca Ø 0,9x Stop 12mmx Comprimento 145mm
05.5.0189.109E	Broca Ø 0,9x Stop 12mmx Comprimento 150mm
05.5.0190.109E	Broca Ø 0,9x Stop 13mmx Comprimento 50mm
05.5.0191.109E	Broca Ø 0,9x Stop 13mmx Comprimento 55mm
05.5.0192.109E	Broca Ø 0,9x Stop 13mmx Comprimento 60mm

05.5.0193.109E	Broca Ø 0,9x Stop 13mmx Comprimento 65mm
05.5.0194.109E	Broca Ø 0,9x Stop 13mmx Comprimento 70mm
05.5.0195.109E	Broca Ø 0,9x Stop 13mmx Comprimento 75mm
05.5.0196.109E	Broca Ø 0,9x Stop 13mmx Comprimento 80mm
05.5.0197.109E	Broca Ø 0,9x Stop 13mmx Comprimento 85mm
05.5.0198.109E	Broca Ø 0,9x Stop 13mmx Comprimento 90mm
05.5.0199.109E	Broca Ø 0,9x Stop 13mmx Comprimento 95mm
05.5.0200.109E	Broca Ø 0,9x Stop 13mmx Comprimento 100mm
05.5.0201.109E	Broca Ø 0,9x Stop 13mmx Comprimento 105mm
05.5.0202.109E	Broca Ø 0,9x Stop 13mmx Comprimento 110mm
05.5.0203.109E	Broca Ø 0,9x Stop 13mmx Comprimento 115mm
05.5.0204.109E	Broca Ø 0,9x Stop 13mmx Comprimento 120mm
05.5.0205.109E	Broca Ø 0,9x Stop 13mmx Comprimento 125mm
05.5.0206.109E	Broca Ø 0,9x Stop 13mmx Comprimento 130mm
05.5.0207.109E	Broca Ø 0,9x Stop 13mmx Comprimento 135mm
05.5.0208.109E	Broca Ø 0,9x Stop 13mmx Comprimento 140mm
05.5.0209.109E	Broca Ø 0,9x Stop 13mmx Comprimento 145mm
05.5.0210.109E	Broca Ø 0,9x Stop 13mmx Comprimento 150mm
05.5.0211.109E	Broca Ø 0,9x Stop 14mmx Comprimento 50mm
05.5.0212.109E	Broca Ø 0,9x Stop 14mmx Comprimento 55mm
05.5.0213.109E	Broca Ø 0,9x Stop 14mmx Comprimento 60mm
05.5.0214.109E	Broca Ø 0,9x Stop 14mmx Comprimento 65mm
05.5.0215.109E	Broca Ø 0,9x Stop 14mmx Comprimento 70mm
05.5.0216.109E	Broca Ø 0,9x Stop 14mmx Comprimento 75mm
05.5.0217.109E	Broca Ø 0,9x Stop 14mmx Comprimento 80mm
05.5.0218.109E	Broca Ø 0,9x Stop 14mmx Comprimento 85mm
05.5.0219.109E	Broca Ø 0,9x Stop 14mmx Comprimento 90mm
05.5.0220.109E	Broca Ø 0,9x Stop 14mmx Comprimento 95mm
05.5.0221.109E	Broca Ø 0,9x Stop 14mmx Comprimento 100mm
05.5.0222.109E	Broca Ø 0,9x Stop 14mmx Comprimento 105mm
05.5.0223.109E	Broca Ø 0,9x Stop 14mmx Comprimento 110mm
05.5.0224.109E	Broca Ø 0,9x Stop 14mmx Comprimento 115mm
05.5.0225.109E	Broca Ø 0,9x Stop 14mmx Comprimento 120mm
05.5.0226.109E	Broca Ø 0,9x Stop 14mmx Comprimento 125mm
05.5.0227.109E	Broca Ø 0,9x Stop 14mmx Comprimento 130mm
05.5.0228.109E	Broca Ø 0,9x Stop 14mmx Comprimento 135mm
05.5.0229.109E	Broca Ø 0,9x Stop 14mmx Comprimento 140mm
05.5.0230.109E	Broca Ø 0,9x Stop 14mmx Comprimento 145mm
05.5.0231.109E	Broca Ø 0,9x Stop 14mmx Comprimento 150mm
05.5.0232.109E	Broca Ø 0,9x Stop 15mmx Comprimento 50mm
05.5.0233.109E	Broca Ø 0,9x Stop 15mmx Comprimento 55mm
05.5.0234.109E	Broca Ø 0,9x Stop 15mmx Comprimento 60mm
05.5.0235.109E	Broca Ø 0,9x Stop 15mmx Comprimento 65mm
05.5.0236.109E	Broca Ø 0,9x Stop 15mmx Comprimento 70mm
05.5.0237.109E	Broca Ø 0,9x Stop 15mmx Comprimento 75mm
05.5.0238.109E	Broca Ø 0,9x Stop 15mmx Comprimento 80mm
05.5.0239.109E	Broca Ø 0,9x Stop 15mmx Comprimento 85mm
05.5.0240.109E	Broca Ø 0,9x Stop 15mmx Comprimento 90mm
05.5.0241.109E	Broca Ø 0,9x Stop 15mmx Comprimento 95mm

05.5.0242.109E	Broca Ø 0,9x Stop 15mmx Comprimento 100mm
05.5.0243.109E	Broca Ø 0,9x Stop 15mmx Comprimento 105mm
05.5.0244.109E	Broca Ø 0,9x Stop 15mmx Comprimento 110mm
05.5.0245.109E	Broca Ø 0,9x Stop 15mmx Comprimento 115mm
05.5.0246.109E	Broca Ø 0,9x Stop 15mmx Comprimento 120mm
05.5.0247.109E	Broca Ø 0,9x Stop 15mmx Comprimento 125mm
05.5.0248.109E	Broca Ø 0,9x Stop 15mmx Comprimento 130mm
05.5.0249.109E	Broca Ø 0,9x Stop 15mmx Comprimento 135mm
05.5.0250.109E	Broca Ø 0,9x Stop 15mmx Comprimento 140mm
05.5.0251.109E	Broca Ø 0,9x Stop 15mmx Comprimento 145mm
05.5.0252.109E	Broca Ø 0,9x Stop 15mmx Comprimento 150mm
05.5.0253.109E	Broca Ø 0,9x Stop 16mmx Comprimento 50mm
05.5.0254.109E	Broca Ø 0,9x Stop 16mmx Comprimento 55mm
05.5.0255.109E	Broca Ø 0,9x Stop 16mmx Comprimento 60mm
05.5.0256.109E	Broca Ø 0,9x Stop 16mmx Comprimento 65mm
05.5.0257.109E	Broca Ø 0,9x Stop 16mmx Comprimento 70mm
05.5.0258.109E	Broca Ø 0,9x Stop 16mmx Comprimento 75mm
05.5.0259.109E	Broca Ø 0,9x Stop 16mmx Comprimento 80mm
05.5.0260.109E	Broca Ø 0,9x Stop 16mmx Comprimento 85mm
05.5.0261.109E	Broca Ø 0,9x Stop 16mmx Comprimento 90mm
05.5.0262.109E	Broca Ø 0,9x Stop 16mmx Comprimento 95mm
05.5.0263.109E	Broca Ø 0,9x Stop 16mmx Comprimento 100mm
05.5.0264.109E	Broca Ø 0,9x Stop 16mmx Comprimento 105mm
05.5.0265.109E	Broca Ø 0,9x Stop 16mmx Comprimento 110mm
05.5.0266.109E	Broca Ø 0,9x Stop 16mmx Comprimento 115mm
05.5.0267.109E	Broca Ø 0,9x Stop 16mmx Comprimento 120mm
05.5.0268.109E	Broca Ø 0,9x Stop 16mmx Comprimento 125mm
05.5.0269.109E	Broca Ø 0,9x Stop 16mmx Comprimento 130mm
05.5.0270.109E	Broca Ø 0,9x Stop 16mmx Comprimento 135mm
05.5.0271.109E	Broca Ø 0,9x Stop 16mmx Comprimento 140mm
05.5.0272.109E	Broca Ø 0,9x Stop 16mmx Comprimento 145mm
05.5.0273.109E	Broca Ø 0,9x Stop 16mmx Comprimento 150mm
05.5.0274.109E	Broca Ø 0,9x Stop 17mmx Comprimento 50mm
05.5.0275.109E	Broca Ø 0,9x Stop 17mmx Comprimento 55mm
05.5.0276.109E	Broca Ø 0,9x Stop 17mmx Comprimento 60mm
05.5.0277.109E	Broca Ø 0,9x Stop 17mmx Comprimento 65mm
05.5.0278.109E	Broca Ø 0,9x Stop 17mmx Comprimento 70mm
05.5.0279.109E	Broca Ø 0,9x Stop 17mmx Comprimento 75mm
05.5.0280.109E	Broca Ø 0,9x Stop 17mmx Comprimento 80mm
05.5.0281.109E	Broca Ø 0,9x Stop 17mmx Comprimento 85mm
05.5.0282.109E	Broca Ø 0,9x Stop 17mmx Comprimento 90mm
05.5.0283.109E	Broca Ø 0,9x Stop 17mmx Comprimento 95mm
05.5.0284.109E	Broca Ø 0,9x Stop 17mmx Comprimento 100mm
05.5.0285.109E	Broca Ø 0,9x Stop 17mmx Comprimento 105mm
05.5.0286.109E	Broca Ø 0,9x Stop 17mmx Comprimento 110mm
05.5.0287.109E	Broca Ø 0,9x Stop 17mmx Comprimento 115mm
05.5.0288.109E	Broca Ø 0,9x Stop 17mmx Comprimento 120mm
05.5.0289.109E	Broca Ø 0,9x Stop 17mmx Comprimento 125mm
05.5.0290.109E	Broca Ø 0,9x Stop 17mmx Comprimento 130mm

05.5.0291.109E	Broca Ø 0,9x Stop 17mmx Comprimento 135mm
05.5.0292.109E	Broca Ø 0,9x Stop 17mmx Comprimento 140mm
05.5.0293.109E	Broca Ø 0,9x Stop 17mmx Comprimento 145mm
05.5.0294.109E	Broca Ø 0,9x Stop 17mmx Comprimento 150mm
05.5.0295.109E	Broca Ø 0,9x Stop 18mmx Comprimento 50mm
05.5.0296.109E	Broca Ø 0,9x Stop 18mmx Comprimento 55mm
05.5.0297.109E	Broca Ø 0,9x Stop 18mmx Comprimento 60mm
05.5.0298.109E	Broca Ø 0,9x Stop 18mmx Comprimento 65mm
05.5.0299.109E	Broca Ø 0,9x Stop 18mmx Comprimento 70mm
05.5.0300.109E	Broca Ø 0,9x Stop 18mmx Comprimento 75mm
05.5.0301.109E	Broca Ø 0,9x Stop 18mmx Comprimento 80mm
05.5.0302.109E	Broca Ø 0,9x Stop 18mmx Comprimento 85mm
05.5.0303.109E	Broca Ø 0,9x Stop 18mmx Comprimento 90mm
05.5.0304.109E	Broca Ø 0,9x Stop 18mmx Comprimento 95mm
05.5.0305.109E	Broca Ø 0,9x Stop 18mmx Comprimento 100mm
05.5.0306.109E	Broca Ø 0,9x Stop 18mmx Comprimento 105mm
05.5.0307.109E	Broca Ø 0,9x Stop 18mmx Comprimento 110mm
05.5.0308.109E	Broca Ø 0,9x Stop 18mmx Comprimento 115mm
05.5.0309.109E	Broca Ø 0,9x Stop 18mmx Comprimento 120mm
05.5.0310.109E	Broca Ø 0,9x Stop 18mmx Comprimento 125mm
05.5.0311.109E	Broca Ø 0,9x Stop 18mmx Comprimento 130mm
05.5.0312.109E	Broca Ø 0,9x Stop 18mmx Comprimento 135mm
05.5.0313.109E	Broca Ø 0,9x Stop 18mmx Comprimento 140mm
05.5.0314.109E	Broca Ø 0,9x Stop 18mmx Comprimento 145mm
05.5.0315.109E	Broca Ø 0,9x Stop 18mmx Comprimento 150mm
05.5.0316.109E	Broca Ø 0,9x Stop 19mmx Comprimento 50mm
05.5.0317.109E	Broca Ø 0,9x Stop 19mmx Comprimento 55mm
05.5.0318.109E	Broca Ø 0,9x Stop 19mmx Comprimento 60mm
05.5.0319.109E	Broca Ø 0,9x Stop 19mmx Comprimento 65mm
05.5.0320.109E	Broca Ø 0,9x Stop 19mmx Comprimento 70mm
05.5.0321.109E	Broca Ø 0,9x Stop 19mmx Comprimento 75mm
05.5.0322.109E	Broca Ø 0,9x Stop 19mmx Comprimento 80mm
05.5.0323.109E	Broca Ø 0,9x Stop 19mmx Comprimento 85mm
05.5.0324.109E	Broca Ø 0,9x Stop 19mmx Comprimento 90mm
05.5.0325.109E	Broca Ø 0,9x Stop 19mmx Comprimento 95mm
05.5.0326.109E	Broca Ø 0,9x Stop 19mmx Comprimento 100mm
05.5.0327.109E	Broca Ø 0,9x Stop 19mmx Comprimento 105mm
05.5.0328.109E	Broca Ø 0,9x Stop 19mmx Comprimento 110mm
05.5.0329.109E	Broca Ø 0,9x Stop 19mmx Comprimento 115mm
05.5.0330.109E	Broca Ø 0,9x Stop 19mmx Comprimento 120mm
05.5.0331.109E	Broca Ø 0,9x Stop 19mmx Comprimento 125mm
05.5.0332.109E	Broca Ø 0,9x Stop 19mmx Comprimento 130mm
05.5.0333.109E	Broca Ø 0,9x Stop 19mmx Comprimento 135mm
05.5.0334.109E	Broca Ø 0,9x Stop 19mmx Comprimento 140mm
05.5.0335.109E	Broca Ø 0,9x Stop 19mmx Comprimento 145mm
05.5.0336.109E	Broca Ø 0,9x Stop 19mmx Comprimento 150mm
05.5.0337.109E	Broca Ø 0,9x Stop 20mmx Comprimento 50mm
05.5.0338.109E	Broca Ø 0,9x Stop 20mmx Comprimento 55mm
05.5.0339.109E	Broca Ø 0,9x Stop 20mmx Comprimento 60mm

05.5.0340.109E	Broca Ø 0,9x Stop 20mmx Comprimento 65mm
05.5.0341.109E	Broca Ø 0,9x Stop 20mmx Comprimento 70mm
05.5.0342.109E	Broca Ø 0,9x Stop 20mmx Comprimento 75mm
05.5.0343.109E	Broca Ø 0,9x Stop 20mmx Comprimento 80mm
05.5.0344.109E	Broca Ø 0,9x Stop 20mmx Comprimento 85mm
05.5.0345.109E	Broca Ø 0,9x Stop 20mmx Comprimento 90mm
05.5.0346.109E	Broca Ø 0,9x Stop 20mmx Comprimento 95mm
05.5.0347.109E	Broca Ø 0,9x Stop 20mmx Comprimento 100mm
05.5.0348.109E	Broca Ø 0,9x Stop 20mmx Comprimento 105mm
05.5.0349.109E	Broca Ø 0,9x Stop 20mmx Comprimento 110mm
05.5.0350.109E	Broca Ø 0,9x Stop 20mmx Comprimento 115mm
05.5.0351.109E	Broca Ø 0,9x Stop 20mmx Comprimento 120mm
05.5.0352.109E	Broca Ø 0,9x Stop 20mmx Comprimento 125mm
05.5.0353.109E	Broca Ø 0,9x Stop 20mmx Comprimento 130mm
05.5.0354.109E	Broca Ø 0,9x Stop 20mmx Comprimento 135mm
05.5.0355.109E	Broca Ø 0,9x Stop 20mmx Comprimento 140mm
05.5.0356.109E	Broca Ø 0,9x Stop 20mmx Comprimento 145mm
05.5.0357.109E	Broca Ø 0,9x Stop 20mmx Comprimento 150mm
05.5.0358.109E	Broca Ø 0,9x Stop 21mmx Comprimento 50mm
05.5.0359.109E	Broca Ø 0,9x Stop 21mmx Comprimento 55mm
05.5.0360.109E	Broca Ø 0,9x Stop 21mmx Comprimento 60mm
05.5.0361.109E	Broca Ø 0,9x Stop 21mmx Comprimento 65mm
05.5.0362.109E	Broca Ø 0,9x Stop 21mmx Comprimento 70mm
05.5.0363.109E	Broca Ø 0,9x Stop 21mmx Comprimento 75mm
05.5.0364.109E	Broca Ø 0,9x Stop 21mmx Comprimento 80mm
05.5.0365.109E	Broca Ø 0,9x Stop 21mmx Comprimento 85mm
05.5.0366.109E	Broca Ø 0,9x Stop 21mmx Comprimento 90mm
05.5.0367.109E	Broca Ø 0,9x Stop 21mmx Comprimento 95mm
05.5.0368.109E	Broca Ø 0,9x Stop 21mmx Comprimento 100mm
05.5.0369.109E	Broca Ø 0,9x Stop 21mmx Comprimento 105mm
05.5.0370.109E	Broca Ø 0,9x Stop 21mmx Comprimento 110mm
05.5.0371.109E	Broca Ø 0,9x Stop 21mmx Comprimento 115mm
05.5.0372.109E	Broca Ø 0,9x Stop 21mmx Comprimento 120mm
05.5.0373.109E	Broca Ø 0,9x Stop 21mmx Comprimento 125mm
05.5.0374.109E	Broca Ø 0,9x Stop 21mmx Comprimento 130mm
05.5.0375.109E	Broca Ø 0,9x Stop 21mmx Comprimento 135mm
05.5.0376.109E	Broca Ø 0,9x Stop 21mmx Comprimento 140mm
05.5.0377.109E	Broca Ø 0,9x Stop 21mmx Comprimento 145mm
05.5.0378.109E	Broca Ø 0,9x Stop 21mmx Comprimento 150mm
05.5.0379.109E	Broca Ø 0,9x Stop 22mmx Comprimento 50mm
05.5.0380.109E	Broca Ø 0,9x Stop 22mmx Comprimento 55mm
05.5.0381.109E	Broca Ø 0,9x Stop 22mmx Comprimento 60mm
05.5.0382.109E	Broca Ø 0,9x Stop 22mmx Comprimento 65mm
05.5.0383.109E	Broca Ø 0,9x Stop 22mmx Comprimento 70mm
05.5.0384.109E	Broca Ø 0,9x Stop 22mmx Comprimento 75mm
05.5.0385.109E	Broca Ø 0,9x Stop 22mmx Comprimento 80mm
05.5.0386.109E	Broca Ø 0,9x Stop 22mmx Comprimento 85mm
05.5.0387.109E	Broca Ø 0,9x Stop 22mmx Comprimento 90mm
05.5.0388.109E	Broca Ø 0,9x Stop 22mmx Comprimento 95mm

05.5.0389.109E	Broca Ø 0,9x Stop 22mmx Comprimento 100mm
05.5.0390.109E	Broca Ø 0,9x Stop 22mmx Comprimento 105mm
05.5.0391.109E	Broca Ø 0,9x Stop 22mmx Comprimento 110mm
05.5.0392.109E	Broca Ø 0,9x Stop 22mmx Comprimento 115mm
05.5.0393.109E	Broca Ø 0,9x Stop 22mmx Comprimento 120mm
05.5.0394.109E	Broca Ø 0,9x Stop 22mmx Comprimento 125mm
05.5.0395.109E	Broca Ø 0,9x Stop 22mmx Comprimento 130mm
05.5.0396.109E	Broca Ø 0,9x Stop 22mmx Comprimento 135mm
05.5.0397.109E	Broca Ø 0,9x Stop 22mmx Comprimento 140mm
05.5.0398.109E	Broca Ø 0,9x Stop 22mmx Comprimento 145mm
05.5.0399.109E	Broca Ø 0,9x Stop 22mmx Comprimento 150mm
05.5.0400.109E	Broca Ø 0,9x Stop 23mmx Comprimento 50mm
05.5.0401.109E	Broca Ø 0,9x Stop 23mmx Comprimento 55mm
05.5.0402.109E	Broca Ø 0,9x Stop 23mmx Comprimento 60mm
05.5.0403.109E	Broca Ø 0,9x Stop 23mmx Comprimento 65mm
05.5.0404.109E	Broca Ø 0,9x Stop 23mmx Comprimento 70mm
05.5.0405.109E	Broca Ø 0,9x Stop 23mmx Comprimento 75mm
05.5.0406.109E	Broca Ø 0,9x Stop 23mmx Comprimento 80mm
05.5.0407.109E	Broca Ø 0,9x Stop 23mmx Comprimento 85mm
05.5.0408.109E	Broca Ø 0,9x Stop 23mmx Comprimento 90mm
05.5.0409.109E	Broca Ø 0,9x Stop 23mmx Comprimento 95mm
05.5.0410.109E	Broca Ø 0,9x Stop 23mmx Comprimento 100mm
05.5.0411.109E	Broca Ø 0,9x Stop 23mmx Comprimento 105mm
05.5.0412.109E	Broca Ø 0,9x Stop 23mmx Comprimento 110mm
05.5.0413.109E	Broca Ø 0,9x Stop 23mmx Comprimento 115mm
05.5.0414.109E	Broca Ø 0,9x Stop 23mmx Comprimento 120mm
05.5.0415.109E	Broca Ø 0,9x Stop 23mmx Comprimento 125mm
05.5.0416.109E	Broca Ø 0,9x Stop 23mmx Comprimento 130mm
05.5.0417.109E	Broca Ø 0,9x Stop 23mmx Comprimento 135mm
05.5.0418.109E	Broca Ø 0,9x Stop 23mmx Comprimento 140mm
05.5.0419.109E	Broca Ø 0,9x Stop 23mmx Comprimento 145mm
05.5.0420.109E	Broca Ø 0,9x Stop 23mmx Comprimento 150mm
05.5.0421.109E	Broca Ø 0,9x Stop 24mmx Comprimento 50mm
05.5.0422.109E	Broca Ø 0,9x Stop 24mmx Comprimento 55mm
05.5.0423.109E	Broca Ø 0,9x Stop 24mmx Comprimento 60mm
05.5.0424.109E	Broca Ø 0,9x Stop 24mmx Comprimento 65mm
05.5.0425.109E	Broca Ø 0,9x Stop 24mmx Comprimento 70mm
05.5.0426.109E	Broca Ø 0,9x Stop 24mmx Comprimento 75mm
05.5.0427.109E	Broca Ø 0,9x Stop 24mmx Comprimento 80mm
05.5.0428.109E	Broca Ø 0,9x Stop 24mmx Comprimento 85mm
05.5.0429.109E	Broca Ø 0,9x Stop 24mmx Comprimento 90mm
05.5.0430.109E	Broca Ø 0,9x Stop 24mmx Comprimento 95mm
05.5.0431.109E	Broca Ø 0,9x Stop 24mmx Comprimento 100mm
05.5.0432.109E	Broca Ø 0,9x Stop 24mmx Comprimento 105mm
05.5.0433.109E	Broca Ø 0,9x Stop 24mmx Comprimento 110mm
05.5.0434.109E	Broca Ø 0,9x Stop 24mmx Comprimento 115mm
05.5.0435.109E	Broca Ø 0,9x Stop 24mmx Comprimento 120mm
05.5.0436.109E	Broca Ø 0,9x Stop 24mmx Comprimento 125mm
05.5.0437.109E	Broca Ø 0,9x Stop 24mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0438.109E	Broca Ø 0,9x Stop 24mmx Comprimento 135mm
05.5.0439.109E	Broca Ø 0,9x Stop 24mmx Comprimento 140mm
05.5.0440.109E	Broca Ø 0,9x Stop 24mmx Comprimento 145mm
05.5.0441.109E	Broca Ø 0,9x Stop 24mmx Comprimento 150mm
05.5.0442.109E	Broca Ø 0,9x Stop 25mmx Comprimento 50mm
05.5.0443.109E	Broca Ø 0,9x Stop 25mmx Comprimento 55mm
05.5.0444.109E	Broca Ø 0,9x Stop 25mmx Comprimento 60mm
05.5.0445.109E	Broca Ø 0,9x Stop 25mmx Comprimento 65mm
05.5.0446.109E	Broca Ø 0,9x Stop 25mmx Comprimento 70mm
05.5.0447.109E	Broca Ø 0,9x Stop 25mmx Comprimento 75mm
05.5.0448.109E	Broca Ø 0,9x Stop 25mmx Comprimento 80mm
05.5.0449.109E	Broca Ø 0,9x Stop 25mmx Comprimento 85mm
05.5.0450.109E	Broca Ø 0,9x Stop 25mmx Comprimento 90mm
05.5.0451.109E	Broca Ø 0,9x Stop 25mmx Comprimento 95mm
05.5.0452.109E	Broca Ø 0,9x Stop 25mmx Comprimento 100mm
05.5.0453.109E	Broca Ø 0,9x Stop 25mmx Comprimento 105mm
05.5.0454.109E	Broca Ø 0,9x Stop 25mmx Comprimento 110mm
05.5.0455.109E	Broca Ø 0,9x Stop 25mmx Comprimento 115mm
05.5.0456.109E	Broca Ø 0,9x Stop 25mmx Comprimento 120mm
05.5.0457.109E	Broca Ø 0,9x Stop 25mmx Comprimento 125mm
05.5.0458.109E	Broca Ø 0,9x Stop 25mmx Comprimento 130mm
05.5.0459.109E	Broca Ø 0,9x Stop 25mmx Comprimento 135mm
05.5.0460.109E	Broca Ø 0,9x Stop 25mmx Comprimento 140mm
05.5.0461.109E	Broca Ø 0,9x Stop 25mmx Comprimento 145mm
05.5.0462.109E	Broca Ø 0,9x Stop 25mmx Comprimento 150mm
05.5.0463.109E	Broca Ø 0,9x Stop 30mmx Comprimento 50mm
05.5.0464.109E	Broca Ø 0,9x Stop 30mmx Comprimento 55mm
05.5.0465.109E	Broca Ø 0,9x Stop 30mmx Comprimento 60mm
05.5.0466.109E	Broca Ø 0,9x Stop 30mmx Comprimento 65mm
05.5.0467.109E	Broca Ø 0,9x Stop 30mmx Comprimento 70mm
05.5.0468.109E	Broca Ø 0,9x Stop 30mmx Comprimento 75mm
05.5.0469.109E	Broca Ø 0,9x Stop 30mmx Comprimento 80mm
05.5.0470.109E	Broca Ø 0,9x Stop 30mmx Comprimento 85mm
05.5.0471.109E	Broca Ø 0,9x Stop 30mmx Comprimento 90mm
05.5.0472.109E	Broca Ø 0,9x Stop 30mmx Comprimento 95mm
05.5.0473.109E	Broca Ø 0,9x Stop 30mmx Comprimento 100mm
05.5.0474.109E	Broca Ø 0,9x Stop 30mmx Comprimento 105mm
05.5.0475.109E	Broca Ø 0,9x Stop 30mmx Comprimento 110mm
05.5.0476.109E	Broca Ø 0,9x Stop 30mmx Comprimento 115mm
05.5.0477.109E	Broca Ø 0,9x Stop 30mmx Comprimento 120mm
05.5.0478.109E	Broca Ø 0,9x Stop 30mmx Comprimento 125mm
05.5.0479.109E	Broca Ø 0,9x Stop 30mmx Comprimento 130mm
05.5.0480.109E	Broca Ø 0,9x Stop 30mmx Comprimento 135mm
05.5.0481.109E	Broca Ø 0,9x Stop 30mmx Comprimento 140mm
05.5.0482.109E	Broca Ø 0,9x Stop 30mmx Comprimento 145mm
05.5.0483.109E	Broca Ø 0,9x Stop 30mmx Comprimento 150mm
05.5.0484.109E	Broca Ø 0,9x Stop 35mmx Comprimento 50mm
05.5.0485.109E	Broca Ø 0,9x Stop 35mmx Comprimento 55mm
05.5.0486.109E	Broca Ø 0,9x Stop 35mmx Comprimento 60mm

05.5.0487.109E	Broca Ø 0,9x Stop 35mmx Comprimento 65mm
05.5.0488.109E	Broca Ø 0,9x Stop 35mmx Comprimento 70mm
05.5.0489.109E	Broca Ø 0,9x Stop 35mmx Comprimento 75mm
05.5.0490.109E	Broca Ø 0,9x Stop 35mmx Comprimento 80mm
05.5.0491.109E	Broca Ø 0,9x Stop 35mmx Comprimento 85mm
05.5.0492.109E	Broca Ø 0,9x Stop 35mmx Comprimento 90mm
05.5.0493.109E	Broca Ø 0,9x Stop 35mmx Comprimento 95mm
05.5.0494.109E	Broca Ø 0,9x Stop 35mmx Comprimento 100mm
05.5.0495.109E	Broca Ø 0,9x Stop 35mmx Comprimento 105mm
05.5.0496.109E	Broca Ø 0,9x Stop 35mmx Comprimento 110mm
05.5.0497.109E	Broca Ø 0,9x Stop 35mmx Comprimento 115mm
05.5.0498.109E	Broca Ø 0,9x Stop 35mmx Comprimento 120mm
05.5.0499.109E	Broca Ø 0,9x Stop 35mmx Comprimento 125mm
05.5.0500.109E	Broca Ø 0,9x Stop 35mmx Comprimento 130mm
05.5.0501.109E	Broca Ø 0,9x Stop 35mmx Comprimento 135mm
05.5.0502.109E	Broca Ø 0,9x Stop 35mmx Comprimento 140mm
05.5.0503.109E	Broca Ø 0,9x Stop 35mmx Comprimento 145mm
05.5.0504.109E	Broca Ø 0,9x Stop 35mmx Comprimento 150mm
05.5.0001.110E	Broca Ø 1,0x Stop 4mmx Comprimento 50mm
05.5.0002.110E	Broca Ø 1,0x Stop 4mmx Comprimento 55mm
05.5.0003.110E	Broca Ø 1,0x Stop 4mmx Comprimento 60mm
05.5.0004.110E	Broca Ø 1,0x Stop 4mmx Comprimento 65mm
05.5.0005.110E	Broca Ø 1,0x Stop 4mmx Comprimento 70mm
05.5.0006.110E	Broca Ø 1,0x Stop 4mmx Comprimento 75mm
05.5.0007.110E	Broca Ø 1,0x Stop 4mmx Comprimento 80mm
05.5.0008.110E	Broca Ø 1,0x Stop 4mmx Comprimento 85mm
05.5.0009.110E	Broca Ø 1,0x Stop 4mmx Comprimento 90mm
05.5.0010.110E	Broca Ø 1,0x Stop 4mmx Comprimento 95mm
05.5.0011.110E	Broca Ø 1,0x Stop 4mmx Comprimento 100mm
05.5.0012.110E	Broca Ø 1,0x Stop 4mmx Comprimento 105mm
05.5.0013.110E	Broca Ø 1,0x Stop 4mmx Comprimento 110mm
05.5.0014.110E	Broca Ø 1,0x Stop 4mmx Comprimento 115mm
05.5.0015.110E	Broca Ø 1,0x Stop 4mmx Comprimento 120mm
05.5.0016.110E	Broca Ø 1,0x Stop 4mmx Comprimento 125mm
05.5.0017.110E	Broca Ø 1,0x Stop 4mmx Comprimento 130mm
05.5.0018.110E	Broca Ø 1,0x Stop 4mmx Comprimento 135mm
05.5.0019.110E	Broca Ø 1,0x Stop 4mmx Comprimento 140mm
05.5.0020.110E	Broca Ø 1,0x Stop 4mmx Comprimento 145mm
05.5.0021.110E	Broca Ø 1,0x Stop 4mmx Comprimento 150mm
05.5.0022.110E	Broca Ø 1,0x Stop 5mmx Comprimento 50mm
05.5.0023.110E	Broca Ø 1,0x Stop 5mmx Comprimento 55mm
05.5.0024.110E	Broca Ø 1,0x Stop 5mmx Comprimento 60mm
05.5.0025.110E	Broca Ø 1,0x Stop 5mmx Comprimento 65mm
05.5.0026.110E	Broca Ø 1,0x Stop 5mmx Comprimento 70mm
05.5.0027.110E	Broca Ø 1,0x Stop 5mmx Comprimento 75mm
05.5.0028.110E	Broca Ø 1,0x Stop 5mmx Comprimento 80mm
05.5.0029.110E	Broca Ø 1,0x Stop 5mmx Comprimento 85mm
05.5.0030.110E	Broca Ø 1,0x Stop 5mmx Comprimento 90mm
05.5.0031.110E	Broca Ø 1,0x Stop 5mmx Comprimento 95mm

05.5.0032.110E	Broca Ø 1,0x Stop 5mmx Comprimento 100mm
05.5.0033.110E	Broca Ø 1,0x Stop 5mmx Comprimento 105mm
05.5.0034.110E	Broca Ø 1,0x Stop 5mmx Comprimento 110mm
05.5.0035.110E	Broca Ø 1,0x Stop 5mmx Comprimento 115mm
05.5.0036.110E	Broca Ø 1,0x Stop 5mmx Comprimento 120mm
05.5.0037.110E	Broca Ø 1,0x Stop 5mmx Comprimento 125mm
05.5.0038.110E	Broca Ø 1,0x Stop 5mmx Comprimento 130mm
05.5.0039.110E	Broca Ø 1,0x Stop 5mmx Comprimento 135mm
05.5.0040.110E	Broca Ø 1,0x Stop 5mmx Comprimento 140mm
05.5.0041.110E	Broca Ø 1,0x Stop 5mmx Comprimento 145mm
05.5.0042.110E	Broca Ø 1,0x Stop 5mmx Comprimento 150mm
05.5.0043.110E	Broca Ø 1,0x Stop 6mmx Comprimento 50mm
05.5.0044.110E	Broca Ø 1,0x Stop 6mmx Comprimento 55mm
05.5.0045.110E	Broca Ø 1,0x Stop 6mmx Comprimento 60mm
05.5.0046.110E	Broca Ø 1,0x Stop 6mmx Comprimento 65mm
05.5.0047.110E	Broca Ø 1,0x Stop 6mmx Comprimento 70mm
05.5.0048.110E	Broca Ø 1,0x Stop 6mmx Comprimento 75mm
05.5.0049.110E	Broca Ø 1,0x Stop 6mmx Comprimento 80mm
05.5.0050.110E	Broca Ø 1,0x Stop 6mmx Comprimento 85mm
05.5.0051.110E	Broca Ø 1,0x Stop 6mmx Comprimento 90mm
05.5.0052.110E	Broca Ø 1,0x Stop 6mmx Comprimento 95mm
05.5.0053.110E	Broca Ø 1,0x Stop 6mmx Comprimento 100mm
05.5.0054.110E	Broca Ø 1,0x Stop 6mmx Comprimento 105mm
05.5.0055.110E	Broca Ø 1,0x Stop 6mmx Comprimento 110mm
05.5.0056.110E	Broca Ø 1,0x Stop 6mmx Comprimento 115mm
05.5.0057.110E	Broca Ø 1,0x Stop 6mmx Comprimento 120mm
05.5.0058.110E	Broca Ø 1,0x Stop 6mmx Comprimento 125mm
05.5.0059.110E	Broca Ø 1,0x Stop 6mmx Comprimento 130mm
05.5.0060.110E	Broca Ø 1,0x Stop 6mmx Comprimento 135mm
05.5.0061.110E	Broca Ø 1,0x Stop 6mmx Comprimento 140mm
05.5.0062.110E	Broca Ø 1,0x Stop 6mmx Comprimento 145mm
05.5.0063.110E	Broca Ø 1,0x Stop 6mmx Comprimento 150mm
05.5.0064.110E	Broca Ø 1,0x Stop 7mmx Comprimento 50mm
05.5.0065.110E	Broca Ø 1,0x Stop 7mmx Comprimento 55mm
05.5.0066.110E	Broca Ø 1,0x Stop 7mmx Comprimento 60mm
05.5.0067.110E	Broca Ø 1,0x Stop 7mmx Comprimento 65mm
05.5.0068.110E	Broca Ø 1,0x Stop 7mmx Comprimento 70mm
05.5.0069.110E	Broca Ø 1,0x Stop 7mmx Comprimento 75mm
05.5.0070.110E	Broca Ø 1,0x Stop 7mmx Comprimento 80mm
05.5.0071.110E	Broca Ø 1,0x Stop 7mmx Comprimento 85mm
05.5.0072.110E	Broca Ø 1,0x Stop 7mmx Comprimento 90mm
05.5.0073.110E	Broca Ø 1,0x Stop 7mmx Comprimento 95mm
05.5.0074.110E	Broca Ø 1,0x Stop 7mmx Comprimento 100mm
05.5.0075.110E	Broca Ø 1,0x Stop 7mmx Comprimento 105mm
05.5.0076.110E	Broca Ø 1,0x Stop 7mmx Comprimento 110mm
05.5.0077.110E	Broca Ø 1,0x Stop 7mmx Comprimento 115mm
05.5.0078.110E	Broca Ø 1,0x Stop 7mmx Comprimento 120mm
05.5.0079.110E	Broca Ø 1,0x Stop 7mmx Comprimento 125mm
05.5.0080.110E	Broca Ø 1,0x Stop 7mmx Comprimento 130mm

05.5.0081.110E	Broca Ø 1,0x Stop 7mmx Comprimento 135mm
05.5.0082.110E	Broca Ø 1,0x Stop 7mmx Comprimento 140mm
05.5.0083.110E	Broca Ø 1,0x Stop 7mmx Comprimento 145mm
05.5.0084.110E	Broca Ø 1,0x Stop 7mmx Comprimento 150mm
05.5.0085.110E	Broca Ø 1,0x Stop 8mmx Comprimento 50mm
05.5.0086.110E	Broca Ø 1,0x Stop 8mmx Comprimento 55mm
05.5.0087.110E	Broca Ø 1,0x Stop 8mmx Comprimento 60mm
05.5.0088.110E	Broca Ø 1,0x Stop 8mmx Comprimento 65mm
05.5.0089.110E	Broca Ø 1,0x Stop 8mmx Comprimento 70mm
05.5.0090.110E	Broca Ø 1,0x Stop 8mmx Comprimento 75mm
05.5.0091.110E	Broca Ø 1,0x Stop 8mmx Comprimento 80mm
05.5.0092.110E	Broca Ø 1,0x Stop 8mmx Comprimento 85mm
05.5.0093.110E	Broca Ø 1,0x Stop 8mmx Comprimento 90mm
05.5.0094.110E	Broca Ø 1,0x Stop 8mmx Comprimento 95mm
05.5.0095.110E	Broca Ø 1,0x Stop 8mmx Comprimento 100mm
05.5.0096.110E	Broca Ø 1,0x Stop 8mmx Comprimento 105mm
05.5.0097.110E	Broca Ø 1,0x Stop 8mmx Comprimento 110mm
05.5.0098.110E	Broca Ø 1,0x Stop 8mmx Comprimento 115mm
05.5.0099.110E	Broca Ø 1,0x Stop 8mmx Comprimento 120mm
05.5.0100.110E	Broca Ø 1,0x Stop 8mmx Comprimento 125mm
05.5.0101.110E	Broca Ø 1,0x Stop 8mmx Comprimento 130mm
05.5.0102.110E	Broca Ø 1,0x Stop 8mmx Comprimento 135mm
05.5.0103.110E	Broca Ø 1,0x Stop 8mmx Comprimento 140mm
05.5.0104.110E	Broca Ø 1,0x Stop 8mmx Comprimento 145mm
05.5.0105.110E	Broca Ø 1,0x Stop 8mmx Comprimento 150mm
05.5.0106.110E	Broca Ø 1,0x Stop 9mmx Comprimento 50mm
05.5.0107.110E	Broca Ø 1,0x Stop 9mmx Comprimento 55mm
05.5.0108.110E	Broca Ø 1,0x Stop 9mmx Comprimento 60mm
05.5.0109.110E	Broca Ø 1,0x Stop 9mmx Comprimento 65mm
05.5.0110.110E	Broca Ø 1,0x Stop 9mmx Comprimento 70mm
05.5.0111.110E	Broca Ø 1,0x Stop 9mmx Comprimento 75mm
05.5.0112.110E	Broca Ø 1,0x Stop 9mmx Comprimento 80mm
05.5.0113.110E	Broca Ø 1,0x Stop 9mmx Comprimento 85mm
05.5.0114.110E	Broca Ø 1,0x Stop 9mmx Comprimento 90mm
05.5.0115.110E	Broca Ø 1,0x Stop 9mmx Comprimento 95mm
05.5.0116.110E	Broca Ø 1,0x Stop 9mmx Comprimento 100mm
05.5.0117.110E	Broca Ø 1,0x Stop 9mmx Comprimento 105mm
05.5.0118.110E	Broca Ø 1,0x Stop 9mmx Comprimento 110mm
05.5.0119.110E	Broca Ø 1,0x Stop 9mmx Comprimento 115mm
05.5.0120.110E	Broca Ø 1,0x Stop 9mmx Comprimento 120mm
05.5.0121.110E	Broca Ø 1,0x Stop 9mmx Comprimento 125mm
05.5.0122.110E	Broca Ø 1,0x Stop 9mmx Comprimento 130mm
05.5.0123.110E	Broca Ø 1,0x Stop 9mmx Comprimento 135mm
05.5.0124.110E	Broca Ø 1,0x Stop 9mmx Comprimento 140mm
05.5.0125.110E	Broca Ø 1,0x Stop 9mmx Comprimento 145mm
05.5.0126.110E	Broca Ø 1,0x Stop 9mmx Comprimento 150mm
05.5.0127.110E	Broca Ø 1,0x Stop 10mmx Comprimento 50mm
05.5.0128.110E	Broca Ø 1,0x Stop 10mmx Comprimento 55mm
05.5.0129.110E	Broca Ø 1,0x Stop 10mmx Comprimento 60mm

05.5.0130.110E	Broca Ø 1,0x Stop 10mmx Comprimento 65mm
05.5.0131.110E	Broca Ø 1,0x Stop 10mmx Comprimento 70mm
05.5.0132.110E	Broca Ø 1,0x Stop 10mmx Comprimento 75mm
05.5.0133.110E	Broca Ø 1,0x Stop 10mmx Comprimento 80mm
05.5.0134.110E	Broca Ø 1,0x Stop 10mmx Comprimento 85mm
05.5.0135.110E	Broca Ø 1,0x Stop 10mmx Comprimento 90mm
05.5.0136.110E	Broca Ø 1,0x Stop 10mmx Comprimento 95mm
05.5.0137.110E	Broca Ø 1,0x Stop 10mmx Comprimento 100mm
05.5.0138.110E	Broca Ø 1,0x Stop 10mmx Comprimento 105mm
05.5.0139.110E	Broca Ø 1,0x Stop 10mmx Comprimento 110mm
05.5.0140.110E	Broca Ø 1,0x Stop 10mmx Comprimento 115mm
05.5.0141.110E	Broca Ø 1,0x Stop 10mmx Comprimento 120mm
05.5.0142.110E	Broca Ø 1,0x Stop 10mmx Comprimento 125mm
05.5.0143.110E	Broca Ø 1,0x Stop 10mmx Comprimento 130mm
05.5.0144.110E	Broca Ø 1,0x Stop 10mmx Comprimento 135mm
05.5.0145.110E	Broca Ø 1,0x Stop 10mmx Comprimento 140mm
05.5.0146.110E	Broca Ø 1,0x Stop 10mmx Comprimento 145mm
05.5.0147.110E	Broca Ø 1,0x Stop 10mmx Comprimento 150mm
05.5.0148.110E	Broca Ø 1,0x Stop 11mmx Comprimento 50mm
05.5.0149.110E	Broca Ø 1,0x Stop 11mmx Comprimento 55mm
05.5.0150.110E	Broca Ø 1,0x Stop 11mmx Comprimento 60mm
05.5.0151.110E	Broca Ø 1,0x Stop 11mmx Comprimento 65mm
05.5.0152.110E	Broca Ø 1,0x Stop 11mmx Comprimento 70mm
05.5.0153.110E	Broca Ø 1,0x Stop 11mmx Comprimento 75mm
05.5.0154.110E	Broca Ø 1,0x Stop 11mmx Comprimento 80mm
05.5.0155.110E	Broca Ø 1,0x Stop 11mmx Comprimento 85mm
05.5.0156.110E	Broca Ø 1,0x Stop 11mmx Comprimento 90mm
05.5.0157.110E	Broca Ø 1,0x Stop 11mmx Comprimento 95mm
05.5.0158.110E	Broca Ø 1,0x Stop 11mmx Comprimento 100mm
05.5.0159.110E	Broca Ø 1,0x Stop 11mmx Comprimento 105mm
05.5.0160.110E	Broca Ø 1,0x Stop 11mmx Comprimento 110mm
05.5.0161.110E	Broca Ø 1,0x Stop 11mmx Comprimento 115mm
05.5.0162.110E	Broca Ø 1,0x Stop 11mmx Comprimento 120mm
05.5.0163.110E	Broca Ø 1,0x Stop 11mmx Comprimento 125mm
05.5.0164.110E	Broca Ø 1,0x Stop 11mmx Comprimento 130mm
05.5.0165.110E	Broca Ø 1,0x Stop 11mmx Comprimento 135mm
05.5.0166.110E	Broca Ø 1,0x Stop 11mmx Comprimento 140mm
05.5.0167.110E	Broca Ø 1,0x Stop 11mmx Comprimento 145mm
05.5.0168.110E	Broca Ø 1,0x Stop 11mmx Comprimento 150mm
05.5.0169.110E	Broca Ø 1,0x Stop 12mmx Comprimento 50mm
05.5.0170.110E	Broca Ø 1,0x Stop 12mmx Comprimento 55mm
05.5.0171.110E	Broca Ø 1,0x Stop 12mmx Comprimento 60mm
05.5.0172.110E	Broca Ø 1,0x Stop 12mmx Comprimento 65mm
05.5.0173.110E	Broca Ø 1,0x Stop 12mmx Comprimento 70mm
05.5.0174.110E	Broca Ø 1,0x Stop 12mmx Comprimento 75mm
05.5.0175.110E	Broca Ø 1,0x Stop 12mmx Comprimento 80mm
05.5.0176.110E	Broca Ø 1,0x Stop 12mmx Comprimento 85mm
05.5.0177.110E	Broca Ø 1,0x Stop 12mmx Comprimento 90mm
05.5.0178.110E	Broca Ø 1,0x Stop 12mmx Comprimento 95mm

05.5.0179.110E	Broca Ø 1,0x Stop 12mmx Comprimento 100mm
05.5.0180.110E	Broca Ø 1,0x Stop 12mmx Comprimento 105mm
05.5.0181.110E	Broca Ø 1,0x Stop 12mmx Comprimento 110mm
05.5.0182.110E	Broca Ø 1,0x Stop 12mmx Comprimento 115mm
05.5.0183.110E	Broca Ø 1,0x Stop 12mmx Comprimento 120mm
05.5.0184.110E	Broca Ø 1,0x Stop 12mmx Comprimento 125mm
05.5.0185.110E	Broca Ø 1,0x Stop 12mmx Comprimento 130mm
05.5.0186.110E	Broca Ø 1,0x Stop 12mmx Comprimento 135mm
05.5.0187.110E	Broca Ø 1,0x Stop 12mmx Comprimento 140mm
05.5.0188.110E	Broca Ø 1,0x Stop 12mmx Comprimento 145mm
05.5.0189.110E	Broca Ø 1,0x Stop 12mmx Comprimento 150mm
05.5.0190.110E	Broca Ø 1,0x Stop 13mmx Comprimento 50mm
05.5.0191.110E	Broca Ø 1,0x Stop 13mmx Comprimento 55mm
05.5.0192.110E	Broca Ø 1,0x Stop 13mmx Comprimento 60mm
05.5.0193.110E	Broca Ø 1,0x Stop 13mmx Comprimento 65mm
05.5.0194.110E	Broca Ø 1,0x Stop 13mmx Comprimento 70mm
05.5.0195.110E	Broca Ø 1,0x Stop 13mmx Comprimento 75mm
05.5.0196.110E	Broca Ø 1,0x Stop 13mmx Comprimento 80mm
05.5.0197.110E	Broca Ø 1,0x Stop 13mmx Comprimento 85mm
05.5.0198.110E	Broca Ø 1,0x Stop 13mmx Comprimento 90mm
05.5.0199.110E	Broca Ø 1,0x Stop 13mmx Comprimento 95mm
05.5.0200.110E	Broca Ø 1,0x Stop 13mmx Comprimento 100mm
05.5.0201.110E	Broca Ø 1,0x Stop 13mmx Comprimento 105mm
05.5.0202.110E	Broca Ø 1,0x Stop 13mmx Comprimento 110mm
05.5.0203.110E	Broca Ø 1,0x Stop 13mmx Comprimento 115mm
05.5.0204.110E	Broca Ø 1,0x Stop 13mmx Comprimento 120mm
05.5.0205.110E	Broca Ø 1,0x Stop 13mmx Comprimento 125mm
05.5.0206.110E	Broca Ø 1,0x Stop 13mmx Comprimento 130mm
05.5.0207.110E	Broca Ø 1,0x Stop 13mmx Comprimento 135mm
05.5.0208.110E	Broca Ø 1,0x Stop 13mmx Comprimento 140mm
05.5.0209.110E	Broca Ø 1,0x Stop 13mmx Comprimento 145mm
05.5.0210.110E	Broca Ø 1,0x Stop 13mmx Comprimento 150mm
05.5.0211.110E	Broca Ø 1,0x Stop 14mmx Comprimento 50mm
05.5.0212.110E	Broca Ø 1,0x Stop 14mmx Comprimento 55mm
05.5.0213.110E	Broca Ø 1,0x Stop 14mmx Comprimento 60mm
05.5.0214.110E	Broca Ø 1,0x Stop 14mmx Comprimento 65mm
05.5.0215.110E	Broca Ø 1,0x Stop 14mmx Comprimento 70mm
05.5.0216.110E	Broca Ø 1,0x Stop 14mmx Comprimento 75mm
05.5.0217.110E	Broca Ø 1,0x Stop 14mmx Comprimento 80mm
05.5.0218.110E	Broca Ø 1,0x Stop 14mmx Comprimento 85mm
05.5.0219.110E	Broca Ø 1,0x Stop 14mmx Comprimento 90mm
05.5.0220.110E	Broca Ø 1,0x Stop 14mmx Comprimento 95mm
05.5.0221.110E	Broca Ø 1,0x Stop 14mmx Comprimento 100mm
05.5.0222.110E	Broca Ø 1,0x Stop 14mmx Comprimento 105mm
05.5.0223.110E	Broca Ø 1,0x Stop 14mmx Comprimento 110mm
05.5.0224.110E	Broca Ø 1,0x Stop 14mmx Comprimento 115mm
05.5.0225.110E	Broca Ø 1,0x Stop 14mmx Comprimento 120mm
05.5.0226.110E	Broca Ø 1,0x Stop 14mmx Comprimento 125mm
05.5.0227.110E	Broca Ø 1,0x Stop 14mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0228.110E	Broca Ø 1,0x Stop 14mmx Comprimento 135mm
05.5.0229.110E	Broca Ø 1,0x Stop 14mmx Comprimento 140mm
05.5.0230.110E	Broca Ø 1,0x Stop 14mmx Comprimento 145mm
05.5.0231.110E	Broca Ø 1,0x Stop 14mmx Comprimento 150mm
05.5.0232.110E	Broca Ø 1,0x Stop 15mmx Comprimento 50mm
05.5.0233.110E	Broca Ø 1,0x Stop 15mmx Comprimento 55mm
05.5.0234.110E	Broca Ø 1,0x Stop 15mmx Comprimento 60mm
05.5.0235.110E	Broca Ø 1,0x Stop 15mmx Comprimento 65mm
05.5.0236.110E	Broca Ø 1,0x Stop 15mmx Comprimento 70mm
05.5.0237.110E	Broca Ø 1,0x Stop 15mmx Comprimento 75mm
05.5.0238.110E	Broca Ø 1,0x Stop 15mmx Comprimento 80mm
05.5.0239.110E	Broca Ø 1,0x Stop 15mmx Comprimento 85mm
05.5.0240.110E	Broca Ø 1,0x Stop 15mmx Comprimento 90mm
05.5.0241.110E	Broca Ø 1,0x Stop 15mmx Comprimento 95mm
05.5.0242.110E	Broca Ø 1,0x Stop 15mmx Comprimento 100mm
05.5.0243.110E	Broca Ø 1,0x Stop 15mmx Comprimento 105mm
05.5.0244.110E	Broca Ø 1,0x Stop 15mmx Comprimento 110mm
05.5.0245.110E	Broca Ø 1,0x Stop 15mmx Comprimento 115mm
05.5.0246.110E	Broca Ø 1,0x Stop 15mmx Comprimento 120mm
05.5.0247.110E	Broca Ø 1,0x Stop 15mmx Comprimento 125mm
05.5.0248.110E	Broca Ø 1,0x Stop 15mmx Comprimento 130mm
05.5.0249.110E	Broca Ø 1,0x Stop 15mmx Comprimento 135mm
05.5.0250.110E	Broca Ø 1,0x Stop 15mmx Comprimento 140mm
05.5.0251.110E	Broca Ø 1,0x Stop 15mmx Comprimento 145mm
05.5.0252.110E	Broca Ø 1,0x Stop 15mmx Comprimento 150mm
05.5.0253.110E	Broca Ø 1,0x Stop 16mmx Comprimento 50mm
05.5.0254.110E	Broca Ø 1,0x Stop 16mmx Comprimento 55mm
05.5.0255.110E	Broca Ø 1,0x Stop 16mmx Comprimento 60mm
05.5.0256.110E	Broca Ø 1,0x Stop 16mmx Comprimento 65mm
05.5.0257.110E	Broca Ø 1,0x Stop 16mmx Comprimento 70mm
05.5.0258.110E	Broca Ø 1,0x Stop 16mmx Comprimento 75mm
05.5.0259.110E	Broca Ø 1,0x Stop 16mmx Comprimento 80mm
05.5.0260.110E	Broca Ø 1,0x Stop 16mmx Comprimento 85mm
05.5.0261.110E	Broca Ø 1,0x Stop 16mmx Comprimento 90mm
05.5.0262.110E	Broca Ø 1,0x Stop 16mmx Comprimento 95mm
05.5.0263.110E	Broca Ø 1,0x Stop 16mmx Comprimento 100mm
05.5.0264.110E	Broca Ø 1,0x Stop 16mmx Comprimento 105mm
05.5.0265.110E	Broca Ø 1,0x Stop 16mmx Comprimento 110mm
05.5.0266.110E	Broca Ø 1,0x Stop 16mmx Comprimento 115mm
05.5.0267.110E	Broca Ø 1,0x Stop 16mmx Comprimento 120mm
05.5.0268.110E	Broca Ø 1,0x Stop 16mmx Comprimento 125mm
05.5.0269.110E	Broca Ø 1,0x Stop 16mmx Comprimento 130mm
05.5.0270.110E	Broca Ø 1,0x Stop 16mmx Comprimento 135mm
05.5.0271.110E	Broca Ø 1,0x Stop 16mmx Comprimento 140mm
05.5.0272.110E	Broca Ø 1,0x Stop 16mmx Comprimento 145mm
05.5.0273.110E	Broca Ø 1,0x Stop 16mmx Comprimento 150mm
05.5.0274.110E	Broca Ø 1,0x Stop 17mmx Comprimento 50mm
05.5.0275.110E	Broca Ø 1,0x Stop 17mmx Comprimento 55mm
05.5.0276.110E	Broca Ø 1,0x Stop 17mmx Comprimento 60mm

05.5.0277.110E	Broca Ø 1,0x Stop 17mmx Comprimento 65mm
05.5.0278.110E	Broca Ø 1,0x Stop 17mmx Comprimento 70mm
05.5.0279.110E	Broca Ø 1,0x Stop 17mmx Comprimento 75mm
05.5.0280.110E	Broca Ø 1,0x Stop 17mmx Comprimento 80mm
05.5.0281.110E	Broca Ø 1,0x Stop 17mmx Comprimento 85mm
05.5.0282.110E	Broca Ø 1,0x Stop 17mmx Comprimento 90mm
05.5.0283.110E	Broca Ø 1,0x Stop 17mmx Comprimento 95mm
05.5.0284.110E	Broca Ø 1,0x Stop 17mmx Comprimento 100mm
05.5.0285.110E	Broca Ø 1,0x Stop 17mmx Comprimento 105mm
05.5.0286.110E	Broca Ø 1,0x Stop 17mmx Comprimento 110mm
05.5.0287.110E	Broca Ø 1,0x Stop 17mmx Comprimento 115mm
05.5.0288.110E	Broca Ø 1,0x Stop 17mmx Comprimento 120mm
05.5.0289.110E	Broca Ø 1,0x Stop 17mmx Comprimento 125mm
05.5.0290.110E	Broca Ø 1,0x Stop 17mmx Comprimento 130mm
05.5.0291.110E	Broca Ø 1,0x Stop 17mmx Comprimento 135mm
05.5.0292.110E	Broca Ø 1,0x Stop 17mmx Comprimento 140mm
05.5.0293.110E	Broca Ø 1,0x Stop 17mmx Comprimento 145mm
05.5.0294.110E	Broca Ø 1,0x Stop 17mmx Comprimento 150mm
05.5.0295.110E	Broca Ø 1,0x Stop 18mmx Comprimento 50mm
05.5.0296.110E	Broca Ø 1,0x Stop 18mmx Comprimento 55mm
05.5.0297.110E	Broca Ø 1,0x Stop 18mmx Comprimento 60mm
05.5.0298.110E	Broca Ø 1,0x Stop 18mmx Comprimento 65mm
05.5.0299.110E	Broca Ø 1,0x Stop 18mmx Comprimento 70mm
05.5.0300.110E	Broca Ø 1,0x Stop 18mmx Comprimento 75mm
05.5.0301.110E	Broca Ø 1,0x Stop 18mmx Comprimento 80mm
05.5.0302.110E	Broca Ø 1,0x Stop 18mmx Comprimento 85mm
05.5.0303.110E	Broca Ø 1,0x Stop 18mmx Comprimento 90mm
05.5.0304.110E	Broca Ø 1,0x Stop 18mmx Comprimento 95mm
05.5.0305.110E	Broca Ø 1,0x Stop 18mmx Comprimento 100mm
05.5.0306.110E	Broca Ø 1,0x Stop 18mmx Comprimento 105mm
05.5.0307.110E	Broca Ø 1,0x Stop 18mmx Comprimento 110mm
05.5.0308.110E	Broca Ø 1,0x Stop 18mmx Comprimento 115mm
05.5.0309.110E	Broca Ø 1,0x Stop 18mmx Comprimento 120mm
05.5.0310.110E	Broca Ø 1,0x Stop 18mmx Comprimento 125mm
05.5.0311.110E	Broca Ø 1,0x Stop 18mmx Comprimento 130mm
05.5.0312.110E	Broca Ø 1,0x Stop 18mmx Comprimento 135mm
05.5.0313.110E	Broca Ø 1,0x Stop 18mmx Comprimento 140mm
05.5.0314.110E	Broca Ø 1,0x Stop 18mmx Comprimento 145mm
05.5.0315.110E	Broca Ø 1,0x Stop 18mmx Comprimento 150mm
05.5.0316.110E	Broca Ø 1,0x Stop 19mmx Comprimento 50mm
05.5.0317.110E	Broca Ø 1,0x Stop 19mmx Comprimento 55mm
05.5.0318.110E	Broca Ø 1,0x Stop 19mmx Comprimento 60mm
05.5.0319.110E	Broca Ø 1,0x Stop 19mmx Comprimento 65mm
05.5.0320.110E	Broca Ø 1,0x Stop 19mmx Comprimento 70mm
05.5.0321.110E	Broca Ø 1,0x Stop 19mmx Comprimento 75mm
05.5.0322.110E	Broca Ø 1,0x Stop 19mmx Comprimento 80mm
05.5.0323.110E	Broca Ø 1,0x Stop 19mmx Comprimento 85mm
05.5.0324.110E	Broca Ø 1,0x Stop 19mmx Comprimento 90mm
05.5.0325.110E	Broca Ø 1,0x Stop 19mmx Comprimento 95mm

05.5.0326.110E	Broca Ø 1,0x Stop 19mmx Comprimento 100mm
05.5.0327.110E	Broca Ø 1,0x Stop 19mmx Comprimento 105mm
05.5.0328.110E	Broca Ø 1,0x Stop 19mmx Comprimento 110mm
05.5.0329.110E	Broca Ø 1,0x Stop 19mmx Comprimento 115mm
05.5.0330.110E	Broca Ø 1,0x Stop 19mmx Comprimento 120mm
05.5.0331.110E	Broca Ø 1,0x Stop 19mmx Comprimento 125mm
05.5.0332.110E	Broca Ø 1,0x Stop 19mmx Comprimento 130mm
05.5.0333.110E	Broca Ø 1,0x Stop 19mmx Comprimento 135mm
05.5.0334.110E	Broca Ø 1,0x Stop 19mmx Comprimento 140mm
05.5.0335.110E	Broca Ø 1,0x Stop 19mmx Comprimento 145mm
05.5.0336.110E	Broca Ø 1,0x Stop 19mmx Comprimento 150mm
05.5.0337.110E	Broca Ø 1,0x Stop 20mmx Comprimento 50mm
05.5.0338.110E	Broca Ø 1,0x Stop 20mmx Comprimento 55mm
05.5.0339.110E	Broca Ø 1,0x Stop 20mmx Comprimento 60mm
05.5.0340.110E	Broca Ø 1,0x Stop 20mmx Comprimento 65mm
05.5.0341.110E	Broca Ø 1,0x Stop 20mmx Comprimento 70mm
05.5.0342.110E	Broca Ø 1,0x Stop 20mmx Comprimento 75mm
05.5.0343.110E	Broca Ø 1,0x Stop 20mmx Comprimento 80mm
05.5.0344.110E	Broca Ø 1,0x Stop 20mmx Comprimento 85mm
05.5.0345.110E	Broca Ø 1,0x Stop 20mmx Comprimento 90mm
05.5.0346.110E	Broca Ø 1,0x Stop 20mmx Comprimento 95mm
05.5.0347.110E	Broca Ø 1,0x Stop 20mmx Comprimento 100mm
05.5.0348.110E	Broca Ø 1,0x Stop 20mmx Comprimento 105mm
05.5.0349.110E	Broca Ø 1,0x Stop 20mmx Comprimento 110mm
05.5.0350.110E	Broca Ø 1,0x Stop 20mmx Comprimento 115mm
05.5.0351.110E	Broca Ø 1,0x Stop 20mmx Comprimento 120mm
05.5.0352.110E	Broca Ø 1,0x Stop 20mmx Comprimento 125mm
05.5.0353.110E	Broca Ø 1,0x Stop 20mmx Comprimento 130mm
05.5.0354.110E	Broca Ø 1,0x Stop 20mmx Comprimento 135mm
05.5.0355.110E	Broca Ø 1,0x Stop 20mmx Comprimento 140mm
05.5.0356.110E	Broca Ø 1,0x Stop 20mmx Comprimento 145mm
05.5.0357.110E	Broca Ø 1,0x Stop 20mmx Comprimento 150mm
05.5.0358.110E	Broca Ø 1,0x Stop 21mmx Comprimento 50mm
05.5.0359.110E	Broca Ø 1,0x Stop 21mmx Comprimento 55mm
05.5.0360.110E	Broca Ø 1,0x Stop 21mmx Comprimento 60mm
05.5.0361.110E	Broca Ø 1,0x Stop 21mmx Comprimento 65mm
05.5.0362.110E	Broca Ø 1,0x Stop 21mmx Comprimento 70mm
05.5.0363.110E	Broca Ø 1,0x Stop 21mmx Comprimento 75mm
05.5.0364.110E	Broca Ø 1,0x Stop 21mmx Comprimento 80mm
05.5.0365.110E	Broca Ø 1,0x Stop 21mmx Comprimento 85mm
05.5.0366.110E	Broca Ø 1,0x Stop 21mmx Comprimento 90mm
05.5.0367.110E	Broca Ø 1,0x Stop 21mmx Comprimento 95mm
05.5.0368.110E	Broca Ø 1,0x Stop 21mmx Comprimento 100mm
05.5.0369.110E	Broca Ø 1,0x Stop 21mmx Comprimento 105mm
05.5.0370.110E	Broca Ø 1,0x Stop 21mmx Comprimento 110mm
05.5.0371.110E	Broca Ø 1,0x Stop 21mmx Comprimento 115mm
05.5.0372.110E	Broca Ø 1,0x Stop 21mmx Comprimento 120mm
05.5.0373.110E	Broca Ø 1,0x Stop 21mmx Comprimento 125mm
05.5.0374.110E	Broca Ø 1,0x Stop 21mmx Comprimento 130mm

05.5.0375.110E	Broca Ø 1,0x Stop 21mmx Comprimento 135mm
05.5.0376.110E	Broca Ø 1,0x Stop 21mmx Comprimento 140mm
05.5.0377.110E	Broca Ø 1,0x Stop 21mmx Comprimento 145mm
05.5.0378.110E	Broca Ø 1,0x Stop 21mmx Comprimento 150mm
05.5.0379.110E	Broca Ø 1,0x Stop 22mmx Comprimento 50mm
05.5.0380.110E	Broca Ø 1,0x Stop 22mmx Comprimento 55mm
05.5.0381.110E	Broca Ø 1,0x Stop 22mmx Comprimento 60mm
05.5.0382.110E	Broca Ø 1,0x Stop 22mmx Comprimento 65mm
05.5.0383.110E	Broca Ø 1,0x Stop 22mmx Comprimento 70mm
05.5.0384.110E	Broca Ø 1,0x Stop 22mmx Comprimento 75mm
05.5.0385.110E	Broca Ø 1,0x Stop 22mmx Comprimento 80mm
05.5.0386.110E	Broca Ø 1,0x Stop 22mmx Comprimento 85mm
05.5.0387.110E	Broca Ø 1,0x Stop 22mmx Comprimento 90mm
05.5.0388.110E	Broca Ø 1,0x Stop 22mmx Comprimento 95mm
05.5.0389.110E	Broca Ø 1,0x Stop 22mmx Comprimento 100mm
05.5.0390.110E	Broca Ø 1,0x Stop 22mmx Comprimento 105mm
05.5.0391.110E	Broca Ø 1,0x Stop 22mmx Comprimento 110mm
05.5.0392.110E	Broca Ø 1,0x Stop 22mmx Comprimento 115mm
05.5.0393.110E	Broca Ø 1,0x Stop 22mmx Comprimento 120mm
05.5.0394.110E	Broca Ø 1,0x Stop 22mmx Comprimento 125mm
05.5.0395.110E	Broca Ø 1,0x Stop 22mmx Comprimento 130mm
05.5.0396.110E	Broca Ø 1,0x Stop 22mmx Comprimento 135mm
05.5.0397.110E	Broca Ø 1,0x Stop 22mmx Comprimento 140mm
05.5.0398.110E	Broca Ø 1,0x Stop 22mmx Comprimento 145mm
05.5.0399.110E	Broca Ø 1,0x Stop 22mmx Comprimento 150mm
05.5.0400.110E	Broca Ø 1,0x Stop 23mmx Comprimento 50mm
05.5.0401.110E	Broca Ø 1,0x Stop 23mmx Comprimento 55mm
05.5.0402.110E	Broca Ø 1,0x Stop 23mmx Comprimento 60mm
05.5.0403.110E	Broca Ø 1,0x Stop 23mmx Comprimento 65mm
05.5.0404.110E	Broca Ø 1,0x Stop 23mmx Comprimento 70mm
05.5.0405.110E	Broca Ø 1,0x Stop 23mmx Comprimento 75mm
05.5.0406.110E	Broca Ø 1,0x Stop 23mmx Comprimento 80mm
05.5.0407.110E	Broca Ø 1,0x Stop 23mmx Comprimento 85mm
05.5.0408.110E	Broca Ø 1,0x Stop 23mmx Comprimento 90mm
05.5.0409.110E	Broca Ø 1,0x Stop 23mmx Comprimento 95mm
05.5.0410.110E	Broca Ø 1,0x Stop 23mmx Comprimento 100mm
05.5.0411.110E	Broca Ø 1,0x Stop 23mmx Comprimento 105mm
05.5.0412.110E	Broca Ø 1,0x Stop 23mmx Comprimento 110mm
05.5.0413.110E	Broca Ø 1,0x Stop 23mmx Comprimento 115mm
05.5.0414.110E	Broca Ø 1,0x Stop 23mmx Comprimento 120mm
05.5.0415.110E	Broca Ø 1,0x Stop 23mmx Comprimento 125mm
05.5.0416.110E	Broca Ø 1,0x Stop 23mmx Comprimento 130mm
05.5.0417.110E	Broca Ø 1,0x Stop 23mmx Comprimento 135mm
05.5.0418.110E	Broca Ø 1,0x Stop 23mmx Comprimento 140mm
05.5.0419.110E	Broca Ø 1,0x Stop 23mmx Comprimento 145mm
05.5.0420.110E	Broca Ø 1,0x Stop 23mmx Comprimento 150mm
05.5.0421.110E	Broca Ø 1,0x Stop 24mmx Comprimento 50mm
05.5.0422.110E	Broca Ø 1,0x Stop 24mmx Comprimento 55mm
05.5.0423.110E	Broca Ø 1,0x Stop 24mmx Comprimento 60mm

05.5.0424.110E	Broca Ø 1,0x Stop 24mmx Comprimento 65mm
05.5.0425.110E	Broca Ø 1,0x Stop 24mmx Comprimento 70mm
05.5.0426.110E	Broca Ø 1,0x Stop 24mmx Comprimento 75mm
05.5.0427.110E	Broca Ø 1,0x Stop 24mmx Comprimento 80mm
05.5.0428.110E	Broca Ø 1,0x Stop 24mmx Comprimento 85mm
05.5.0429.110E	Broca Ø 1,0x Stop 24mmx Comprimento 90mm
05.5.0430.110E	Broca Ø 1,0x Stop 24mmx Comprimento 95mm
05.5.0431.110E	Broca Ø 1,0x Stop 24mmx Comprimento 100mm
05.5.0432.110E	Broca Ø 1,0x Stop 24mmx Comprimento 105mm
05.5.0433.110E	Broca Ø 1,0x Stop 24mmx Comprimento 110mm
05.5.0434.110E	Broca Ø 1,0x Stop 24mmx Comprimento 115mm
05.5.0435.110E	Broca Ø 1,0x Stop 24mmx Comprimento 120mm
05.5.0436.110E	Broca Ø 1,0x Stop 24mmx Comprimento 125mm
05.5.0437.110E	Broca Ø 1,0x Stop 24mmx Comprimento 130mm
05.5.0438.110E	Broca Ø 1,0x Stop 24mmx Comprimento 135mm
05.5.0439.110E	Broca Ø 1,0x Stop 24mmx Comprimento 140mm
05.5.0440.110E	Broca Ø 1,0x Stop 24mmx Comprimento 145mm
05.5.0441.110E	Broca Ø 1,0x Stop 24mmx Comprimento 150mm
05.5.0442.110E	Broca Ø 1,0x Stop 25mmx Comprimento 50mm
05.5.0443.110E	Broca Ø 1,0x Stop 25mmx Comprimento 55mm
05.5.0444.110E	Broca Ø 1,0x Stop 25mmx Comprimento 60mm
05.5.0445.110E	Broca Ø 1,0x Stop 25mmx Comprimento 65mm
05.5.0446.110E	Broca Ø 1,0x Stop 25mmx Comprimento 70mm
05.5.0447.110E	Broca Ø 1,0x Stop 25mmx Comprimento 75mm
05.5.0448.110E	Broca Ø 1,0x Stop 25mmx Comprimento 80mm
05.5.0449.110E	Broca Ø 1,0x Stop 25mmx Comprimento 85mm
05.5.0450.110E	Broca Ø 1,0x Stop 25mmx Comprimento 90mm
05.5.0451.110E	Broca Ø 1,0x Stop 25mmx Comprimento 95mm
05.5.0452.110E	Broca Ø 1,0x Stop 25mmx Comprimento 100mm
05.5.0453.110E	Broca Ø 1,0x Stop 25mmx Comprimento 105mm
05.5.0454.110E	Broca Ø 1,0x Stop 25mmx Comprimento 110mm
05.5.0455.110E	Broca Ø 1,0x Stop 25mmx Comprimento 115mm
05.5.0456.110E	Broca Ø 1,0x Stop 25mmx Comprimento 120mm
05.5.0457.110E	Broca Ø 1,0x Stop 25mmx Comprimento 125mm
05.5.0458.110E	Broca Ø 1,0x Stop 25mmx Comprimento 130mm
05.5.0459.110E	Broca Ø 1,0x Stop 25mmx Comprimento 135mm
05.5.0460.110E	Broca Ø 1,0x Stop 25mmx Comprimento 140mm
05.5.0461.110E	Broca Ø 1,0x Stop 25mmx Comprimento 145mm
05.5.0462.110E	Broca Ø 1,0x Stop 25mmx Comprimento 150mm
05.5.0463.110E	Broca Ø 1,0x Stop 30mmx Comprimento 50mm
05.5.0464.110E	Broca Ø 1,0x Stop 30mmx Comprimento 55mm
05.5.0465.110E	Broca Ø 1,0x Stop 30mmx Comprimento 60mm
05.5.0466.110E	Broca Ø 1,0x Stop 30mmx Comprimento 65mm
05.5.0467.110E	Broca Ø 1,0x Stop 30mmx Comprimento 70mm
05.5.0468.110E	Broca Ø 1,0x Stop 30mmx Comprimento 75mm
05.5.0469.110E	Broca Ø 1,0x Stop 30mmx Comprimento 80mm
05.5.0470.110E	Broca Ø 1,0x Stop 30mmx Comprimento 85mm
05.5.0471.110E	Broca Ø 1,0x Stop 30mmx Comprimento 90mm
05.5.0472.110E	Broca Ø 1,0x Stop 30mmx Comprimento 95mm

05.5.0473.110E	Broca Ø 1,0x Stop 30mmx Comprimento 100mm
05.5.0474.110E	Broca Ø 1,0x Stop 30mmx Comprimento 105mm
05.5.0475.110E	Broca Ø 1,0x Stop 30mmx Comprimento 110mm
05.5.0476.110E	Broca Ø 1,0x Stop 30mmx Comprimento 115mm
05.5.0477.110E	Broca Ø 1,0x Stop 30mmx Comprimento 120mm
05.5.0478.110E	Broca Ø 1,0x Stop 30mmx Comprimento 125mm
05.5.0479.110E	Broca Ø 1,0x Stop 30mmx Comprimento 130mm
05.5.0480.110E	Broca Ø 1,0x Stop 30mmx Comprimento 135mm
05.5.0481.110E	Broca Ø 1,0x Stop 30mmx Comprimento 140mm
05.5.0482.110E	Broca Ø 1,0x Stop 30mmx Comprimento 145mm
05.5.0483.110E	Broca Ø 1,0x Stop 30mmx Comprimento 150mm
05.5.0484.110E	Broca Ø 1,0x Stop 35mmx Comprimento 50mm
05.5.0485.110E	Broca Ø 1,0x Stop 35mmx Comprimento 55mm
05.5.0486.110E	Broca Ø 1,0x Stop 35mmx Comprimento 60mm
05.5.0487.110E	Broca Ø 1,0x Stop 35mmx Comprimento 65mm
05.5.0488.110E	Broca Ø 1,0x Stop 35mmx Comprimento 70mm
05.5.0489.110E	Broca Ø 1,0x Stop 35mmx Comprimento 75mm
05.5.0490.110E	Broca Ø 1,0x Stop 35mmx Comprimento 80mm
05.5.0491.110E	Broca Ø 1,0x Stop 35mmx Comprimento 85mm
05.5.0492.110E	Broca Ø 1,0x Stop 35mmx Comprimento 90mm
05.5.0493.110E	Broca Ø 1,0x Stop 35mmx Comprimento 95mm
05.5.0494.110E	Broca Ø 1,0x Stop 35mmx Comprimento 100mm
05.5.0495.110E	Broca Ø 1,0x Stop 35mmx Comprimento 105mm
05.5.0496.110E	Broca Ø 1,0x Stop 35mmx Comprimento 110mm
05.5.0497.110E	Broca Ø 1,0x Stop 35mmx Comprimento 115mm
05.5.0498.110E	Broca Ø 1,0x Stop 35mmx Comprimento 120mm
05.5.0499.110E	Broca Ø 1,0x Stop 35mmx Comprimento 125mm
05.5.0500.110E	Broca Ø 1,0x Stop 35mmx Comprimento 130mm
05.5.0501.110E	Broca Ø 1,0x Stop 35mmx Comprimento 135mm
05.5.0502.110E	Broca Ø 1,0x Stop 35mmx Comprimento 140mm
05.5.0503.110E	Broca Ø 1,0x Stop 35mmx Comprimento 145mm
05.5.0504.110E	Broca Ø 1,0x Stop 35mmx Comprimento 150mm
05.5.0001.111E	Broca Ø 1,1x Stop 4mmx Comprimento 50mm
05.5.0002.111E	Broca Ø 1,1x Stop 4mmx Comprimento 55mm
05.5.0003.111E	Broca Ø 1,1x Stop 4mmx Comprimento 60mm
05.5.0004.111E	Broca Ø 1,1x Stop 4mmx Comprimento 65mm
05.5.0005.111E	Broca Ø 1,1x Stop 4mmx Comprimento 70mm
05.5.0006.111E	Broca Ø 1,1x Stop 4mmx Comprimento 75mm
05.5.0007.111E	Broca Ø 1,1x Stop 4mmx Comprimento 80mm
05.5.0008.111E	Broca Ø 1,1x Stop 4mmx Comprimento 85mm
05.5.0009.111E	Broca Ø 1,1x Stop 4mmx Comprimento 90mm
05.5.0010.111E	Broca Ø 1,1x Stop 4mmx Comprimento 95mm
05.5.0011.111E	Broca Ø 1,1x Stop 4mmx Comprimento 100mm
05.5.0012.111E	Broca Ø 1,1x Stop 4mmx Comprimento 105mm
05.5.0013.111E	Broca Ø 1,1x Stop 4mmx Comprimento 110mm
05.5.0014.111E	Broca Ø 1,1x Stop 4mmx Comprimento 115mm
05.5.0015.111E	Broca Ø 1,1x Stop 4mmx Comprimento 120mm
05.5.0016.111E	Broca Ø 1,1x Stop 4mmx Comprimento 125mm
05.5.0017.111E	Broca Ø 1,1x Stop 4mmx Comprimento 130mm

05.5.0018.111E	Broca Ø 1,1x Stop 4mmx Comprimento 135mm
05.5.0019.111E	Broca Ø 1,1x Stop 4mmx Comprimento 140mm
05.5.0020.111E	Broca Ø 1,1x Stop 4mmx Comprimento 145mm
05.5.0021.111E	Broca Ø 1,1x Stop 4mmx Comprimento 150mm
05.5.0022.111E	Broca Ø 1,1x Stop 5mmx Comprimento 50mm
05.5.0023.111E	Broca Ø 1,1x Stop 5mmx Comprimento 55mm
05.5.0024.111E	Broca Ø 1,1x Stop 5mmx Comprimento 60mm
05.5.0025.111E	Broca Ø 1,1x Stop 5mmx Comprimento 65mm
05.5.0026.111E	Broca Ø 1,1x Stop 5mmx Comprimento 70mm
05.5.0027.111E	Broca Ø 1,1x Stop 5mmx Comprimento 75mm
05.5.0028.111E	Broca Ø 1,1x Stop 5mmx Comprimento 80mm
05.5.0029.111E	Broca Ø 1,1x Stop 5mmx Comprimento 85mm
05.5.0030.111E	Broca Ø 1,1x Stop 5mmx Comprimento 90mm
05.5.0031.111E	Broca Ø 1,1x Stop 5mmx Comprimento 95mm
05.5.0032.111E	Broca Ø 1,1x Stop 5mmx Comprimento 100mm
05.5.0033.111E	Broca Ø 1,1x Stop 5mmx Comprimento 105mm
05.5.0034.111E	Broca Ø 1,1x Stop 5mmx Comprimento 110mm
05.5.0035.111E	Broca Ø 1,1x Stop 5mmx Comprimento 115mm
05.5.0036.111E	Broca Ø 1,1x Stop 5mmx Comprimento 120mm
05.5.0037.111E	Broca Ø 1,1x Stop 5mmx Comprimento 125mm
05.5.0038.111E	Broca Ø 1,1x Stop 5mmx Comprimento 130mm
05.5.0039.111E	Broca Ø 1,1x Stop 5mmx Comprimento 135mm
05.5.0040.111E	Broca Ø 1,1x Stop 5mmx Comprimento 140mm
05.5.0041.111E	Broca Ø 1,1x Stop 5mmx Comprimento 145mm
05.5.0042.111E	Broca Ø 1,1x Stop 5mmx Comprimento 150mm
05.5.0043.111E	Broca Ø 1,1x Stop 6mmx Comprimento 50mm
05.5.0044.111E	Broca Ø 1,1x Stop 6mmx Comprimento 55mm
05.5.0045.111E	Broca Ø 1,1x Stop 6mmx Comprimento 60mm
05.5.0046.111E	Broca Ø 1,1x Stop 6mmx Comprimento 65mm
05.5.0047.111E	Broca Ø 1,1x Stop 6mmx Comprimento 70mm
05.5.0048.111E	Broca Ø 1,1x Stop 6mmx Comprimento 75mm
05.5.0049.111E	Broca Ø 1,1x Stop 6mmx Comprimento 80mm
05.5.0050.111E	Broca Ø 1,1x Stop 6mmx Comprimento 85mm
05.5.0051.111E	Broca Ø 1,1x Stop 6mmx Comprimento 90mm
05.5.0052.111E	Broca Ø 1,1x Stop 6mmx Comprimento 95mm
05.5.0053.111E	Broca Ø 1,1x Stop 6mmx Comprimento 100mm
05.5.0054.111E	Broca Ø 1,1x Stop 6mmx Comprimento 105mm
05.5.0055.111E	Broca Ø 1,1x Stop 6mmx Comprimento 110mm
05.5.0056.111E	Broca Ø 1,1x Stop 6mmx Comprimento 115mm
05.5.0057.111E	Broca Ø 1,1x Stop 6mmx Comprimento 120mm
05.5.0058.111E	Broca Ø 1,1x Stop 6mmx Comprimento 125mm
05.5.0059.111E	Broca Ø 1,1x Stop 6mmx Comprimento 130mm
05.5.0060.111E	Broca Ø 1,1x Stop 6mmx Comprimento 135mm
05.5.0061.111E	Broca Ø 1,1x Stop 6mmx Comprimento 140mm
05.5.0062.111E	Broca Ø 1,1x Stop 6mmx Comprimento 145mm
05.5.0063.111E	Broca Ø 1,1x Stop 6mmx Comprimento 150mm
05.5.0064.111E	Broca Ø 1,1x Stop 7mmx Comprimento 50mm
05.5.0065.111E	Broca Ø 1,1x Stop 7mmx Comprimento 55mm
05.5.0066.111E	Broca Ø 1,1x Stop 7mmx Comprimento 60mm

05.5.0067.111E	Broca Ø 1,1x Stop 7mmx Comprimento 65mm
05.5.0068.111E	Broca Ø 1,1x Stop 7mmx Comprimento 70mm
05.5.0069.111E	Broca Ø 1,1x Stop 7mmx Comprimento 75mm
05.5.0070.111E	Broca Ø 1,1x Stop 7mmx Comprimento 80mm
05.5.0071.111E	Broca Ø 1,1x Stop 7mmx Comprimento 85mm
05.5.0072.111E	Broca Ø 1,1x Stop 7mmx Comprimento 90mm
05.5.0073.111E	Broca Ø 1,1x Stop 7mmx Comprimento 95mm
05.5.0074.111E	Broca Ø 1,1x Stop 7mmx Comprimento 100mm
05.5.0075.111E	Broca Ø 1,1x Stop 7mmx Comprimento 105mm
05.5.0076.111E	Broca Ø 1,1x Stop 7mmx Comprimento 110mm
05.5.0077.111E	Broca Ø 1,1x Stop 7mmx Comprimento 115mm
05.5.0078.111E	Broca Ø 1,1x Stop 7mmx Comprimento 120mm
05.5.0079.111E	Broca Ø 1,1x Stop 7mmx Comprimento 125mm
05.5.0080.111E	Broca Ø 1,1x Stop 7mmx Comprimento 130mm
05.5.0081.111E	Broca Ø 1,1x Stop 7mmx Comprimento 135mm
05.5.0082.111E	Broca Ø 1,1x Stop 7mmx Comprimento 140mm
05.5.0083.111E	Broca Ø 1,1x Stop 7mmx Comprimento 145mm
05.5.0084.111E	Broca Ø 1,1x Stop 7mmx Comprimento 150mm
05.5.0085.111E	Broca Ø 1,1x Stop 8mmx Comprimento 50mm
05.5.0086.111E	Broca Ø 1,1x Stop 8mmx Comprimento 55mm
05.5.0087.111E	Broca Ø 1,1x Stop 8mmx Comprimento 60mm
05.5.0088.111E	Broca Ø 1,1x Stop 8mmx Comprimento 65mm
05.5.0089.111E	Broca Ø 1,1x Stop 8mmx Comprimento 70mm
05.5.0090.111E	Broca Ø 1,1x Stop 8mmx Comprimento 75mm
05.5.0091.111E	Broca Ø 1,1x Stop 8mmx Comprimento 80mm
05.5.0092.111E	Broca Ø 1,1x Stop 8mmx Comprimento 85mm
05.5.0093.111E	Broca Ø 1,1x Stop 8mmx Comprimento 90mm
05.5.0094.111E	Broca Ø 1,1x Stop 8mmx Comprimento 95mm
05.5.0095.111E	Broca Ø 1,1x Stop 8mmx Comprimento 100mm
05.5.0096.111E	Broca Ø 1,1x Stop 8mmx Comprimento 105mm
05.5.0097.111E	Broca Ø 1,1x Stop 8mmx Comprimento 110mm
05.5.0098.111E	Broca Ø 1,1x Stop 8mmx Comprimento 115mm
05.5.0099.111E	Broca Ø 1,1x Stop 8mmx Comprimento 120mm
05.5.0100.111E	Broca Ø 1,1x Stop 8mmx Comprimento 125mm
05.5.0101.111E	Broca Ø 1,1x Stop 8mmx Comprimento 130mm
05.5.0102.111E	Broca Ø 1,1x Stop 8mmx Comprimento 135mm
05.5.0103.111E	Broca Ø 1,1x Stop 8mmx Comprimento 140mm
05.5.0104.111E	Broca Ø 1,1x Stop 8mmx Comprimento 145mm
05.5.0105.111E	Broca Ø 1,1x Stop 8mmx Comprimento 150mm
05.5.0106.111E	Broca Ø 1,1x Stop 9mmx Comprimento 50mm
05.5.0107.111E	Broca Ø 1,1x Stop 9mmx Comprimento 55mm
05.5.0108.111E	Broca Ø 1,1x Stop 9mmx Comprimento 60mm
05.5.0109.111E	Broca Ø 1,1x Stop 9mmx Comprimento 65mm
05.5.0110.111E	Broca Ø 1,1x Stop 9mmx Comprimento 70mm
05.5.0111.111E	Broca Ø 1,1x Stop 9mmx Comprimento 75mm
05.5.0112.111E	Broca Ø 1,1x Stop 9mmx Comprimento 80mm
05.5.0113.111E	Broca Ø 1,1x Stop 9mmx Comprimento 85mm
05.5.0114.111E	Broca Ø 1,1x Stop 9mmx Comprimento 90mm
05.5.0115.111E	Broca Ø 1,1x Stop 9mmx Comprimento 95mm

05.5.0116.111E	Broca Ø 1,1x Stop 9mmx Comprimento 100mm
05.5.0117.111E	Broca Ø 1,1x Stop 9mmx Comprimento 105mm
05.5.0118.111E	Broca Ø 1,1x Stop 9mmx Comprimento 110mm
05.5.0119.111E	Broca Ø 1,1x Stop 9mmx Comprimento 115mm
05.5.0120.111E	Broca Ø 1,1x Stop 9mmx Comprimento 120mm
05.5.0121.111E	Broca Ø 1,1x Stop 9mmx Comprimento 125mm
05.5.0122.111E	Broca Ø 1,1x Stop 9mmx Comprimento 130mm
05.5.0123.111E	Broca Ø 1,1x Stop 9mmx Comprimento 135mm
05.5.0124.111E	Broca Ø 1,1x Stop 9mmx Comprimento 140mm
05.5.0125.111E	Broca Ø 1,1x Stop 9mmx Comprimento 145mm
05.5.0126.111E	Broca Ø 1,1x Stop 9mmx Comprimento 150mm
05.5.0127.111E	Broca Ø 1,1x Stop 10mmx Comprimento 50mm
05.5.0128.111E	Broca Ø 1,1x Stop 10mmx Comprimento 55mm
05.5.0129.111E	Broca Ø 1,1x Stop 10mmx Comprimento 60mm
05.5.0130.111E	Broca Ø 1,1x Stop 10mmx Comprimento 65mm
05.5.0131.111E	Broca Ø 1,1x Stop 10mmx Comprimento 70mm
05.5.0132.111E	Broca Ø 1,1x Stop 10mmx Comprimento 75mm
05.5.0133.111E	Broca Ø 1,1x Stop 10mmx Comprimento 80mm
05.5.0134.111E	Broca Ø 1,1x Stop 10mmx Comprimento 85mm
05.5.0135.111E	Broca Ø 1,1x Stop 10mmx Comprimento 90mm
05.5.0136.111E	Broca Ø 1,1x Stop 10mmx Comprimento 95mm
05.5.0137.111E	Broca Ø 1,1x Stop 10mmx Comprimento 100mm
05.5.0138.111E	Broca Ø 1,1x Stop 10mmx Comprimento 105mm
05.5.0139.111E	Broca Ø 1,1x Stop 10mmx Comprimento 110mm
05.5.0140.111E	Broca Ø 1,1x Stop 10mmx Comprimento 115mm
05.5.0141.111E	Broca Ø 1,1x Stop 10mmx Comprimento 120mm
05.5.0142.111E	Broca Ø 1,1x Stop 10mmx Comprimento 125mm
05.5.0143.111E	Broca Ø 1,1x Stop 10mmx Comprimento 130mm
05.5.0144.111E	Broca Ø 1,1x Stop 10mmx Comprimento 135mm
05.5.0145.111E	Broca Ø 1,1x Stop 10mmx Comprimento 140mm
05.5.0146.111E	Broca Ø 1,1x Stop 10mmx Comprimento 145mm
05.5.0147.111E	Broca Ø 1,1x Stop 10mmx Comprimento 150mm
05.5.0148.111E	Broca Ø 1,1x Stop 11mmx Comprimento 50mm
05.5.0149.111E	Broca Ø 1,1x Stop 11mmx Comprimento 55mm
05.5.0150.111E	Broca Ø 1,1x Stop 11mmx Comprimento 60mm
05.5.0151.111E	Broca Ø 1,1x Stop 11mmx Comprimento 65mm
05.5.0152.111E	Broca Ø 1,1x Stop 11mmx Comprimento 70mm
05.5.0153.111E	Broca Ø 1,1x Stop 11mmx Comprimento 75mm
05.5.0154.111E	Broca Ø 1,1x Stop 11mmx Comprimento 80mm
05.5.0155.111E	Broca Ø 1,1x Stop 11mmx Comprimento 85mm
05.5.0156.111E	Broca Ø 1,1x Stop 11mmx Comprimento 90mm
05.5.0157.111E	Broca Ø 1,1x Stop 11mmx Comprimento 95mm
05.5.0158.111E	Broca Ø 1,1x Stop 11mmx Comprimento 100mm
05.5.0159.111E	Broca Ø 1,1x Stop 11mmx Comprimento 105mm
05.5.0160.111E	Broca Ø 1,1x Stop 11mmx Comprimento 110mm
05.5.0161.111E	Broca Ø 1,1x Stop 11mmx Comprimento 115mm
05.5.0162.111E	Broca Ø 1,1x Stop 11mmx Comprimento 120mm
05.5.0163.111E	Broca Ø 1,1x Stop 11mmx Comprimento 125mm
05.5.0164.111E	Broca Ø 1,1x Stop 11mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0165.111E	Broca Ø 1,1x Stop 11mmx Comprimento 135mm
05.5.0166.111E	Broca Ø 1,1x Stop 11mmx Comprimento 140mm
05.5.0167.111E	Broca Ø 1,1x Stop 11mmx Comprimento 145mm
05.5.0168.111E	Broca Ø 1,1x Stop 11mmx Comprimento 150mm
05.5.0169.111E	Broca Ø 1,1x Stop 12mmx Comprimento 50mm
05.5.0170.111E	Broca Ø 1,1x Stop 12mmx Comprimento 55mm
05.5.0171.111E	Broca Ø 1,1x Stop 12mmx Comprimento 60mm
05.5.0172.111E	Broca Ø 1,1x Stop 12mmx Comprimento 65mm
05.5.0173.111E	Broca Ø 1,1x Stop 12mmx Comprimento 70mm
05.5.0174.111E	Broca Ø 1,1x Stop 12mmx Comprimento 75mm
05.5.0175.111E	Broca Ø 1,1x Stop 12mmx Comprimento 80mm
05.5.0176.111E	Broca Ø 1,1x Stop 12mmx Comprimento 85mm
05.5.0177.111E	Broca Ø 1,1x Stop 12mmx Comprimento 90mm
05.5.0178.111E	Broca Ø 1,1x Stop 12mmx Comprimento 95mm
05.5.0179.111E	Broca Ø 1,1x Stop 12mmx Comprimento 100mm
05.5.0180.111E	Broca Ø 1,1x Stop 12mmx Comprimento 105mm
05.5.0181.111E	Broca Ø 1,1x Stop 12mmx Comprimento 110mm
05.5.0182.111E	Broca Ø 1,1x Stop 12mmx Comprimento 115mm
05.5.0183.111E	Broca Ø 1,1x Stop 12mmx Comprimento 120mm
05.5.0184.111E	Broca Ø 1,1x Stop 12mmx Comprimento 125mm
05.5.0185.111E	Broca Ø 1,1x Stop 12mmx Comprimento 130mm
05.5.0186.111E	Broca Ø 1,1x Stop 12mmx Comprimento 135mm
05.5.0187.111E	Broca Ø 1,1x Stop 12mmx Comprimento 140mm
05.5.0188.111E	Broca Ø 1,1x Stop 12mmx Comprimento 145mm
05.5.0189.111E	Broca Ø 1,1x Stop 12mmx Comprimento 150mm
05.5.0190.111E	Broca Ø 1,1x Stop 13mmx Comprimento 50mm
05.5.0191.111E	Broca Ø 1,1x Stop 13mmx Comprimento 55mm
05.5.0192.111E	Broca Ø 1,1x Stop 13mmx Comprimento 60mm
05.5.0193.111E	Broca Ø 1,1x Stop 13mmx Comprimento 65mm
05.5.0194.111E	Broca Ø 1,1x Stop 13mmx Comprimento 70mm
05.5.0195.111E	Broca Ø 1,1x Stop 13mmx Comprimento 75mm
05.5.0196.111E	Broca Ø 1,1x Stop 13mmx Comprimento 80mm
05.5.0197.111E	Broca Ø 1,1x Stop 13mmx Comprimento 85mm
05.5.0198.111E	Broca Ø 1,1x Stop 13mmx Comprimento 90mm
05.5.0199.111E	Broca Ø 1,1x Stop 13mmx Comprimento 95mm
05.5.0200.111E	Broca Ø 1,1x Stop 13mmx Comprimento 100mm
05.5.0201.111E	Broca Ø 1,1x Stop 13mmx Comprimento 105mm
05.5.0202.111E	Broca Ø 1,1x Stop 13mmx Comprimento 110mm
05.5.0203.111E	Broca Ø 1,1x Stop 13mmx Comprimento 115mm
05.5.0204.111E	Broca Ø 1,1x Stop 13mmx Comprimento 120mm
05.5.0205.111E	Broca Ø 1,1x Stop 13mmx Comprimento 125mm
05.5.0206.111E	Broca Ø 1,1x Stop 13mmx Comprimento 130mm
05.5.0207.111E	Broca Ø 1,1x Stop 13mmx Comprimento 135mm
05.5.0208.111E	Broca Ø 1,1x Stop 13mmx Comprimento 140mm
05.5.0209.111E	Broca Ø 1,1x Stop 13mmx Comprimento 145mm
05.5.0210.111E	Broca Ø 1,1x Stop 13mmx Comprimento 150mm
05.5.0211.111E	Broca Ø 1,1x Stop 14mmx Comprimento 50mm
05.5.0212.111E	Broca Ø 1,1x Stop 14mmx Comprimento 55mm
05.5.0213.111E	Broca Ø 1,1x Stop 14mmx Comprimento 60mm

05.5.0214.111E	Broca Ø 1,1x Stop 14mmx Comprimento 65mm
05.5.0215.111E	Broca Ø 1,1x Stop 14mmx Comprimento 70mm
05.5.0216.111E	Broca Ø 1,1x Stop 14mmx Comprimento 75mm
05.5.0217.111E	Broca Ø 1,1x Stop 14mmx Comprimento 80mm
05.5.0218.111E	Broca Ø 1,1x Stop 14mmx Comprimento 85mm
05.5.0219.111E	Broca Ø 1,1x Stop 14mmx Comprimento 90mm
05.5.0220.111E	Broca Ø 1,1x Stop 14mmx Comprimento 95mm
05.5.0221.111E	Broca Ø 1,1x Stop 14mmx Comprimento 100mm
05.5.0222.111E	Broca Ø 1,1x Stop 14mmx Comprimento 105mm
05.5.0223.111E	Broca Ø 1,1x Stop 14mmx Comprimento 110mm
05.5.0224.111E	Broca Ø 1,1x Stop 14mmx Comprimento 115mm
05.5.0225.111E	Broca Ø 1,1x Stop 14mmx Comprimento 120mm
05.5.0226.111E	Broca Ø 1,1x Stop 14mmx Comprimento 125mm
05.5.0227.111E	Broca Ø 1,1x Stop 14mmx Comprimento 130mm
05.5.0228.111E	Broca Ø 1,1x Stop 14mmx Comprimento 135mm
05.5.0229.111E	Broca Ø 1,1x Stop 14mmx Comprimento 140mm
05.5.0230.111E	Broca Ø 1,1x Stop 14mmx Comprimento 145mm
05.5.0231.111E	Broca Ø 1,1x Stop 14mmx Comprimento 150mm
05.5.0232.111E	Broca Ø 1,1x Stop 15mmx Comprimento 50mm
05.5.0233.111E	Broca Ø 1,1x Stop 15mmx Comprimento 55mm
05.5.0234.111E	Broca Ø 1,1x Stop 15mmx Comprimento 60mm
05.5.0235.111E	Broca Ø 1,1x Stop 15mmx Comprimento 65mm
05.5.0236.111E	Broca Ø 1,1x Stop 15mmx Comprimento 70mm
05.5.0237.111E	Broca Ø 1,1x Stop 15mmx Comprimento 75mm
05.5.0238.111E	Broca Ø 1,1x Stop 15mmx Comprimento 80mm
05.5.0239.111E	Broca Ø 1,1x Stop 15mmx Comprimento 85mm
05.5.0240.111E	Broca Ø 1,1x Stop 15mmx Comprimento 90mm
05.5.0241.111E	Broca Ø 1,1x Stop 15mmx Comprimento 95mm
05.5.0242.111E	Broca Ø 1,1x Stop 15mmx Comprimento 100mm
05.5.0243.111E	Broca Ø 1,1x Stop 15mmx Comprimento 105mm
05.5.0244.111E	Broca Ø 1,1x Stop 15mmx Comprimento 110mm
05.5.0245.111E	Broca Ø 1,1x Stop 15mmx Comprimento 115mm
05.5.0246.111E	Broca Ø 1,1x Stop 15mmx Comprimento 120mm
05.5.0247.111E	Broca Ø 1,1x Stop 15mmx Comprimento 125mm
05.5.0248.111E	Broca Ø 1,1x Stop 15mmx Comprimento 130mm
05.5.0249.111E	Broca Ø 1,1x Stop 15mmx Comprimento 135mm
05.5.0250.111E	Broca Ø 1,1x Stop 15mmx Comprimento 140mm
05.5.0251.111E	Broca Ø 1,1x Stop 15mmx Comprimento 145mm
05.5.0252.111E	Broca Ø 1,1x Stop 15mmx Comprimento 150mm
05.5.0253.111E	Broca Ø 1,1x Stop 16mmx Comprimento 50mm
05.5.0254.111E	Broca Ø 1,1x Stop 16mmx Comprimento 55mm
05.5.0255.111E	Broca Ø 1,1x Stop 16mmx Comprimento 60mm
05.5.0256.111E	Broca Ø 1,1x Stop 16mmx Comprimento 65mm
05.5.0257.111E	Broca Ø 1,1x Stop 16mmx Comprimento 70mm
05.5.0258.111E	Broca Ø 1,1x Stop 16mmx Comprimento 75mm
05.5.0259.111E	Broca Ø 1,1x Stop 16mmx Comprimento 80mm
05.5.0260.111E	Broca Ø 1,1x Stop 16mmx Comprimento 85mm
05.5.0261.111E	Broca Ø 1,1x Stop 16mmx Comprimento 90mm
05.5.0262.111E	Broca Ø 1,1x Stop 16mmx Comprimento 95mm

05.5.0263.111E	Broca Ø 1,1x Stop 16mmx Comprimento 100mm
05.5.0264.111E	Broca Ø 1,1x Stop 16mmx Comprimento 105mm
05.5.0265.111E	Broca Ø 1,1x Stop 16mmx Comprimento 110mm
05.5.0266.111E	Broca Ø 1,1x Stop 16mmx Comprimento 115mm
05.5.0267.111E	Broca Ø 1,1x Stop 16mmx Comprimento 120mm
05.5.0268.111E	Broca Ø 1,1x Stop 16mmx Comprimento 125mm
05.5.0269.111E	Broca Ø 1,1x Stop 16mmx Comprimento 130mm
05.5.0270.111E	Broca Ø 1,1x Stop 16mmx Comprimento 135mm
05.5.0271.111E	Broca Ø 1,1x Stop 16mmx Comprimento 140mm
05.5.0272.111E	Broca Ø 1,1x Stop 16mmx Comprimento 145mm
05.5.0273.111E	Broca Ø 1,1x Stop 16mmx Comprimento 150mm
05.5.0274.111E	Broca Ø 1,1x Stop 17mmx Comprimento 50mm
05.5.0275.111E	Broca Ø 1,1x Stop 17mmx Comprimento 55mm
05.5.0276.111E	Broca Ø 1,1x Stop 17mmx Comprimento 60mm
05.5.0277.111E	Broca Ø 1,1x Stop 17mmx Comprimento 65mm
05.5.0278.111E	Broca Ø 1,1x Stop 17mmx Comprimento 70mm
05.5.0279.111E	Broca Ø 1,1x Stop 17mmx Comprimento 75mm
05.5.0280.111E	Broca Ø 1,1x Stop 17mmx Comprimento 80mm
05.5.0281.111E	Broca Ø 1,1x Stop 17mmx Comprimento 85mm
05.5.0282.111E	Broca Ø 1,1x Stop 17mmx Comprimento 90mm
05.5.0283.111E	Broca Ø 1,1x Stop 17mmx Comprimento 95mm
05.5.0284.111E	Broca Ø 1,1x Stop 17mmx Comprimento 100mm
05.5.0285.111E	Broca Ø 1,1x Stop 17mmx Comprimento 105mm
05.5.0286.111E	Broca Ø 1,1x Stop 17mmx Comprimento 110mm
05.5.0287.111E	Broca Ø 1,1x Stop 17mmx Comprimento 115mm
05.5.0288.111E	Broca Ø 1,1x Stop 17mmx Comprimento 120mm
05.5.0289.111E	Broca Ø 1,1x Stop 17mmx Comprimento 125mm
05.5.0290.111E	Broca Ø 1,1x Stop 17mmx Comprimento 130mm
05.5.0291.111E	Broca Ø 1,1x Stop 17mmx Comprimento 135mm
05.5.0292.111E	Broca Ø 1,1x Stop 17mmx Comprimento 140mm
05.5.0293.111E	Broca Ø 1,1x Stop 17mmx Comprimento 145mm
05.5.0294.111E	Broca Ø 1,1x Stop 17mmx Comprimento 150mm
05.5.0295.111E	Broca Ø 1,1x Stop 18mmx Comprimento 50mm
05.5.0296.111E	Broca Ø 1,1x Stop 18mmx Comprimento 55mm
05.5.0297.111E	Broca Ø 1,1x Stop 18mmx Comprimento 60mm
05.5.0298.111E	Broca Ø 1,1x Stop 18mmx Comprimento 65mm
05.5.0299.111E	Broca Ø 1,1x Stop 18mmx Comprimento 70mm
05.5.0300.111E	Broca Ø 1,1x Stop 18mmx Comprimento 75mm
05.5.0301.111E	Broca Ø 1,1x Stop 18mmx Comprimento 80mm
05.5.0302.111E	Broca Ø 1,1x Stop 18mmx Comprimento 85mm
05.5.0303.111E	Broca Ø 1,1x Stop 18mmx Comprimento 90mm
05.5.0304.111E	Broca Ø 1,1x Stop 18mmx Comprimento 95mm
05.5.0305.111E	Broca Ø 1,1x Stop 18mmx Comprimento 100mm
05.5.0306.111E	Broca Ø 1,1x Stop 18mmx Comprimento 105mm
05.5.0307.111E	Broca Ø 1,1x Stop 18mmx Comprimento 110mm
05.5.0308.111E	Broca Ø 1,1x Stop 18mmx Comprimento 115mm
05.5.0309.111E	Broca Ø 1,1x Stop 18mmx Comprimento 120mm
05.5.0310.111E	Broca Ø 1,1x Stop 18mmx Comprimento 125mm
05.5.0311.111E	Broca Ø 1,1x Stop 18mmx Comprimento 130mm

05.5.0312.111E	Broca Ø 1,1x Stop 18mmx Comprimento 135mm
05.5.0313.111E	Broca Ø 1,1x Stop 18mmx Comprimento 140mm
05.5.0314.111E	Broca Ø 1,1x Stop 18mmx Comprimento 145mm
05.5.0315.111E	Broca Ø 1,1x Stop 18mmx Comprimento 150mm
05.5.0316.111E	Broca Ø 1,1x Stop 19mmx Comprimento 50mm
05.5.0317.111E	Broca Ø 1,1x Stop 19mmx Comprimento 55mm
05.5.0318.111E	Broca Ø 1,1x Stop 19mmx Comprimento 60mm
05.5.0319.111E	Broca Ø 1,1x Stop 19mmx Comprimento 65mm
05.5.0320.111E	Broca Ø 1,1x Stop 19mmx Comprimento 70mm
05.5.0321.111E	Broca Ø 1,1x Stop 19mmx Comprimento 75mm
05.5.0322.111E	Broca Ø 1,1x Stop 19mmx Comprimento 80mm
05.5.0323.111E	Broca Ø 1,1x Stop 19mmx Comprimento 85mm
05.5.0324.111E	Broca Ø 1,1x Stop 19mmx Comprimento 90mm
05.5.0325.111E	Broca Ø 1,1x Stop 19mmx Comprimento 95mm
05.5.0326.111E	Broca Ø 1,1x Stop 19mmx Comprimento 100mm
05.5.0327.111E	Broca Ø 1,1x Stop 19mmx Comprimento 105mm
05.5.0328.111E	Broca Ø 1,1x Stop 19mmx Comprimento 110mm
05.5.0329.111E	Broca Ø 1,1x Stop 19mmx Comprimento 115mm
05.5.0330.111E	Broca Ø 1,1x Stop 19mmx Comprimento 120mm
05.5.0331.111E	Broca Ø 1,1x Stop 19mmx Comprimento 125mm
05.5.0332.111E	Broca Ø 1,1x Stop 19mmx Comprimento 130mm
05.5.0333.111E	Broca Ø 1,1x Stop 19mmx Comprimento 135mm
05.5.0334.111E	Broca Ø 1,1x Stop 19mmx Comprimento 140mm
05.5.0335.111E	Broca Ø 1,1x Stop 19mmx Comprimento 145mm
05.5.0336.111E	Broca Ø 1,1x Stop 19mmx Comprimento 150mm
05.5.0337.111E	Broca Ø 1,1x Stop 20mmx Comprimento 50mm
05.5.0338.111E	Broca Ø 1,1x Stop 20mmx Comprimento 55mm
05.5.0339.111E	Broca Ø 1,1x Stop 20mmx Comprimento 60mm
05.5.0340.111E	Broca Ø 1,1x Stop 20mmx Comprimento 65mm
05.5.0341.111E	Broca Ø 1,1x Stop 20mmx Comprimento 70mm
05.5.0342.111E	Broca Ø 1,1x Stop 20mmx Comprimento 75mm
05.5.0343.111E	Broca Ø 1,1x Stop 20mmx Comprimento 80mm
05.5.0344.111E	Broca Ø 1,1x Stop 20mmx Comprimento 85mm
05.5.0345.111E	Broca Ø 1,1x Stop 20mmx Comprimento 90mm
05.5.0346.111E	Broca Ø 1,1x Stop 20mmx Comprimento 95mm
05.5.0347.111E	Broca Ø 1,1x Stop 20mmx Comprimento 100mm
05.5.0348.111E	Broca Ø 1,1x Stop 20mmx Comprimento 105mm
05.5.0349.111E	Broca Ø 1,1x Stop 20mmx Comprimento 110mm
05.5.0350.111E	Broca Ø 1,1x Stop 20mmx Comprimento 115mm
05.5.0351.111E	Broca Ø 1,1x Stop 20mmx Comprimento 120mm
05.5.0352.111E	Broca Ø 1,1x Stop 20mmx Comprimento 125mm
05.5.0353.111E	Broca Ø 1,1x Stop 20mmx Comprimento 130mm
05.5.0354.111E	Broca Ø 1,1x Stop 20mmx Comprimento 135mm
05.5.0355.111E	Broca Ø 1,1x Stop 20mmx Comprimento 140mm
05.5.0356.111E	Broca Ø 1,1x Stop 20mmx Comprimento 145mm
05.5.0357.111E	Broca Ø 1,1x Stop 20mmx Comprimento 150mm
05.5.0358.111E	Broca Ø 1,1x Stop 21mmx Comprimento 50mm
05.5.0359.111E	Broca Ø 1,1x Stop 21mmx Comprimento 55mm
05.5.0360.111E	Broca Ø 1,1x Stop 21mmx Comprimento 60mm

05.5.0361.111E	Broca Ø 1,1x Stop 21mmx Comprimento 65mm
05.5.0362.111E	Broca Ø 1,1x Stop 21mmx Comprimento 70mm
05.5.0363.111E	Broca Ø 1,1x Stop 21mmx Comprimento 75mm
05.5.0364.111E	Broca Ø 1,1x Stop 21mmx Comprimento 80mm
05.5.0365.111E	Broca Ø 1,1x Stop 21mmx Comprimento 85mm
05.5.0366.111E	Broca Ø 1,1x Stop 21mmx Comprimento 90mm
05.5.0367.111E	Broca Ø 1,1x Stop 21mmx Comprimento 95mm
05.5.0368.111E	Broca Ø 1,1x Stop 21mmx Comprimento 100mm
05.5.0369.111E	Broca Ø 1,1x Stop 21mmx Comprimento 105mm
05.5.0370.111E	Broca Ø 1,1x Stop 21mmx Comprimento 110mm
05.5.0371.111E	Broca Ø 1,1x Stop 21mmx Comprimento 115mm
05.5.0372.111E	Broca Ø 1,1x Stop 21mmx Comprimento 120mm
05.5.0373.111E	Broca Ø 1,1x Stop 21mmx Comprimento 125mm
05.5.0374.111E	Broca Ø 1,1x Stop 21mmx Comprimento 130mm
05.5.0375.111E	Broca Ø 1,1x Stop 21mmx Comprimento 135mm
05.5.0376.111E	Broca Ø 1,1x Stop 21mmx Comprimento 140mm
05.5.0377.111E	Broca Ø 1,1x Stop 21mmx Comprimento 145mm
05.5.0378.111E	Broca Ø 1,1x Stop 21mmx Comprimento 150mm
05.5.0379.111E	Broca Ø 1,1x Stop 22mmx Comprimento 50mm
05.5.0380.111E	Broca Ø 1,1x Stop 22mmx Comprimento 55mm
05.5.0381.111E	Broca Ø 1,1x Stop 22mmx Comprimento 60mm
05.5.0382.111E	Broca Ø 1,1x Stop 22mmx Comprimento 65mm
05.5.0383.111E	Broca Ø 1,1x Stop 22mmx Comprimento 70mm
05.5.0384.111E	Broca Ø 1,1x Stop 22mmx Comprimento 75mm
05.5.0385.111E	Broca Ø 1,1x Stop 22mmx Comprimento 80mm
05.5.0386.111E	Broca Ø 1,1x Stop 22mmx Comprimento 85mm
05.5.0387.111E	Broca Ø 1,1x Stop 22mmx Comprimento 90mm
05.5.0388.111E	Broca Ø 1,1x Stop 22mmx Comprimento 95mm
05.5.0389.111E	Broca Ø 1,1x Stop 22mmx Comprimento 100mm
05.5.0390.111E	Broca Ø 1,1x Stop 22mmx Comprimento 105mm
05.5.0391.111E	Broca Ø 1,1x Stop 22mmx Comprimento 110mm
05.5.0392.111E	Broca Ø 1,1x Stop 22mmx Comprimento 115mm
05.5.0393.111E	Broca Ø 1,1x Stop 22mmx Comprimento 120mm
05.5.0394.111E	Broca Ø 1,1x Stop 22mmx Comprimento 125mm
05.5.0395.111E	Broca Ø 1,1x Stop 22mmx Comprimento 130mm
05.5.0396.111E	Broca Ø 1,1x Stop 22mmx Comprimento 135mm
05.5.0397.111E	Broca Ø 1,1x Stop 22mmx Comprimento 140mm
05.5.0398.111E	Broca Ø 1,1x Stop 22mmx Comprimento 145mm
05.5.0399.111E	Broca Ø 1,1x Stop 22mmx Comprimento 150mm
05.5.0400.111E	Broca Ø 1,1x Stop 23mmx Comprimento 50mm
05.5.0401.111E	Broca Ø 1,1x Stop 23mmx Comprimento 55mm
05.5.0402.111E	Broca Ø 1,1x Stop 23mmx Comprimento 60mm
05.5.0403.111E	Broca Ø 1,1x Stop 23mmx Comprimento 65mm
05.5.0404.111E	Broca Ø 1,1x Stop 23mmx Comprimento 70mm
05.5.0405.111E	Broca Ø 1,1x Stop 23mmx Comprimento 75mm
05.5.0406.111E	Broca Ø 1,1x Stop 23mmx Comprimento 80mm
05.5.0407.111E	Broca Ø 1,1x Stop 23mmx Comprimento 85mm
05.5.0408.111E	Broca Ø 1,1x Stop 23mmx Comprimento 90mm
05.5.0409.111E	Broca Ø 1,1x Stop 23mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0410.111E	Broca Ø 1,1x Stop 23mmx Comprimento 100mm
05.5.0411.111E	Broca Ø 1,1x Stop 23mmx Comprimento 105mm
05.5.0412.111E	Broca Ø 1,1x Stop 23mmx Comprimento 110mm
05.5.0413.111E	Broca Ø 1,1x Stop 23mmx Comprimento 115mm
05.5.0414.111E	Broca Ø 1,1x Stop 23mmx Comprimento 120mm
05.5.0415.111E	Broca Ø 1,1x Stop 23mmx Comprimento 125mm
05.5.0416.111E	Broca Ø 1,1x Stop 23mmx Comprimento 130mm
05.5.0417.111E	Broca Ø 1,1x Stop 23mmx Comprimento 135mm
05.5.0418.111E	Broca Ø 1,1x Stop 23mmx Comprimento 140mm
05.5.0419.111E	Broca Ø 1,1x Stop 23mmx Comprimento 145mm
05.5.0420.111E	Broca Ø 1,1x Stop 23mmx Comprimento 150mm
05.5.0421.111E	Broca Ø 1,1x Stop 24mmx Comprimento 50mm
05.5.0422.111E	Broca Ø 1,1x Stop 24mmx Comprimento 55mm
05.5.0423.111E	Broca Ø 1,1x Stop 24mmx Comprimento 60mm
05.5.0424.111E	Broca Ø 1,1x Stop 24mmx Comprimento 65mm
05.5.0425.111E	Broca Ø 1,1x Stop 24mmx Comprimento 70mm
05.5.0426.111E	Broca Ø 1,1x Stop 24mmx Comprimento 75mm
05.5.0427.111E	Broca Ø 1,1x Stop 24mmx Comprimento 80mm
05.5.0428.111E	Broca Ø 1,1x Stop 24mmx Comprimento 85mm
05.5.0429.111E	Broca Ø 1,1x Stop 24mmx Comprimento 90mm
05.5.0430.111E	Broca Ø 1,1x Stop 24mmx Comprimento 95mm
05.5.0431.111E	Broca Ø 1,1x Stop 24mmx Comprimento 100mm
05.5.0432.111E	Broca Ø 1,1x Stop 24mmx Comprimento 105mm
05.5.0433.111E	Broca Ø 1,1x Stop 24mmx Comprimento 110mm
05.5.0434.111E	Broca Ø 1,1x Stop 24mmx Comprimento 115mm
05.5.0435.111E	Broca Ø 1,1x Stop 24mmx Comprimento 120mm
05.5.0436.111E	Broca Ø 1,1x Stop 24mmx Comprimento 125mm
05.5.0437.111E	Broca Ø 1,1x Stop 24mmx Comprimento 130mm
05.5.0438.111E	Broca Ø 1,1x Stop 24mmx Comprimento 135mm
05.5.0439.111E	Broca Ø 1,1x Stop 24mmx Comprimento 140mm
05.5.0440.111E	Broca Ø 1,1x Stop 24mmx Comprimento 145mm
05.5.0441.111E	Broca Ø 1,1x Stop 24mmx Comprimento 150mm
05.5.0442.111E	Broca Ø 1,1x Stop 25mmx Comprimento 50mm
05.5.0443.111E	Broca Ø 1,1x Stop 25mmx Comprimento 55mm
05.5.0444.111E	Broca Ø 1,1x Stop 25mmx Comprimento 60mm
05.5.0445.111E	Broca Ø 1,1x Stop 25mmx Comprimento 65mm
05.5.0446.111E	Broca Ø 1,1x Stop 25mmx Comprimento 70mm
05.5.0447.111E	Broca Ø 1,1x Stop 25mmx Comprimento 75mm
05.5.0448.111E	Broca Ø 1,1x Stop 25mmx Comprimento 80mm
05.5.0449.111E	Broca Ø 1,1x Stop 25mmx Comprimento 85mm
05.5.0450.111E	Broca Ø 1,1x Stop 25mmx Comprimento 90mm
05.5.0451.111E	Broca Ø 1,1x Stop 25mmx Comprimento 95mm
05.5.0452.111E	Broca Ø 1,1x Stop 25mmx Comprimento 100mm
05.5.0453.111E	Broca Ø 1,1x Stop 25mmx Comprimento 105mm
05.5.0454.111E	Broca Ø 1,1x Stop 25mmx Comprimento 110mm
05.5.0455.111E	Broca Ø 1,1x Stop 25mmx Comprimento 115mm
05.5.0456.111E	Broca Ø 1,1x Stop 25mmx Comprimento 120mm
05.5.0457.111E	Broca Ø 1,1x Stop 25mmx Comprimento 125mm
05.5.0458.111E	Broca Ø 1,1x Stop 25mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0459.111E	Broca Ø 1,1x Stop 25mmx Comprimento 135mm
05.5.0460.111E	Broca Ø 1,1x Stop 25mmx Comprimento 140mm
05.5.0461.111E	Broca Ø 1,1x Stop 25mmx Comprimento 145mm
05.5.0462.111E	Broca Ø 1,1x Stop 25mmx Comprimento 150mm
05.5.0463.111E	Broca Ø 1,1x Stop 30mmx Comprimento 50mm
05.5.0464.111E	Broca Ø 1,1x Stop 30mmx Comprimento 55mm
05.5.0465.111E	Broca Ø 1,1x Stop 30mmx Comprimento 60mm
05.5.0466.111E	Broca Ø 1,1x Stop 30mmx Comprimento 65mm
05.5.0467.111E	Broca Ø 1,1x Stop 30mmx Comprimento 70mm
05.5.0468.111E	Broca Ø 1,1x Stop 30mmx Comprimento 75mm
05.5.0469.111E	Broca Ø 1,1x Stop 30mmx Comprimento 80mm
05.5.0470.111E	Broca Ø 1,1x Stop 30mmx Comprimento 85mm
05.5.0471.111E	Broca Ø 1,1x Stop 30mmx Comprimento 90mm
05.5.0472.111E	Broca Ø 1,1x Stop 30mmx Comprimento 95mm
05.5.0473.111E	Broca Ø 1,1x Stop 30mmx Comprimento 100mm
05.5.0474.111E	Broca Ø 1,1x Stop 30mmx Comprimento 105mm
05.5.0475.111E	Broca Ø 1,1x Stop 30mmx Comprimento 110mm
05.5.0476.111E	Broca Ø 1,1x Stop 30mmx Comprimento 115mm
05.5.0477.111E	Broca Ø 1,1x Stop 30mmx Comprimento 120mm
05.5.0478.111E	Broca Ø 1,1x Stop 30mmx Comprimento 125mm
05.5.0479.111E	Broca Ø 1,1x Stop 30mmx Comprimento 130mm
05.5.0480.111E	Broca Ø 1,1x Stop 30mmx Comprimento 135mm
05.5.0481.111E	Broca Ø 1,1x Stop 30mmx Comprimento 140mm
05.5.0482.111E	Broca Ø 1,1x Stop 30mmx Comprimento 145mm
05.5.0483.111E	Broca Ø 1,1x Stop 30mmx Comprimento 150mm
05.5.0484.111E	Broca Ø 1,1x Stop 35mmx Comprimento 50mm
05.5.0485.111E	Broca Ø 1,1x Stop 35mmx Comprimento 55mm
05.5.0486.111E	Broca Ø 1,1x Stop 35mmx Comprimento 60mm
05.5.0487.111E	Broca Ø 1,1x Stop 35mmx Comprimento 65mm
05.5.0488.111E	Broca Ø 1,1x Stop 35mmx Comprimento 70mm
05.5.0489.111E	Broca Ø 1,1x Stop 35mmx Comprimento 75mm
05.5.0490.111E	Broca Ø 1,1x Stop 35mmx Comprimento 80mm
05.5.0491.111E	Broca Ø 1,1x Stop 35mmx Comprimento 85mm
05.5.0492.111E	Broca Ø 1,1x Stop 35mmx Comprimento 90mm
05.5.0493.111E	Broca Ø 1,1x Stop 35mmx Comprimento 95mm
05.5.0494.111E	Broca Ø 1,1x Stop 35mmx Comprimento 100mm
05.5.0495.111E	Broca Ø 1,1x Stop 35mmx Comprimento 105mm
05.5.0496.111E	Broca Ø 1,1x Stop 35mmx Comprimento 110mm
05.5.0497.111E	Broca Ø 1,1x Stop 35mmx Comprimento 115mm
05.5.0498.111E	Broca Ø 1,1x Stop 35mmx Comprimento 120mm
05.5.0499.111E	Broca Ø 1,1x Stop 35mmx Comprimento 125mm
05.5.0500.111E	Broca Ø 1,1x Stop 35mmx Comprimento 130mm
05.5.0501.111E	Broca Ø 1,1x Stop 35mmx Comprimento 135mm
05.5.0502.111E	Broca Ø 1,1x Stop 35mmx Comprimento 140mm
05.5.0503.111E	Broca Ø 1,1x Stop 35mmx Comprimento 145mm
05.5.0504.111E	Broca Ø 1,1x Stop 35mmx Comprimento 150mm
05.5.0001.112E	Broca Ø 1,2x Stop 4mmx Comprimento 50mm
05.5.0002.112E	Broca Ø 1,2x Stop 4mmx Comprimento 55mm
05.5.0003.112E	Broca Ø 1,2x Stop 4mmx Comprimento 60mm

05.5.0004.112E	Broca Ø 1,2x Stop 4mmx Comprimento 65mm
05.5.0005.112E	Broca Ø 1,2x Stop 4mmx Comprimento 70mm
05.5.0006.112E	Broca Ø 1,2x Stop 4mmx Comprimento 75mm
05.5.0007.112E	Broca Ø 1,2x Stop 4mmx Comprimento 80mm
05.5.0008.112E	Broca Ø 1,2x Stop 4mmx Comprimento 85mm
05.5.0009.112E	Broca Ø 1,2x Stop 4mmx Comprimento 90mm
05.5.0010.112E	Broca Ø 1,2x Stop 4mmx Comprimento 95mm
05.5.0011.112E	Broca Ø 1,2x Stop 4mmx Comprimento 100mm
05.5.0012.112E	Broca Ø 1,2x Stop 4mmx Comprimento 105mm
05.5.0013.112E	Broca Ø 1,2x Stop 4mmx Comprimento 110mm
05.5.0014.112E	Broca Ø 1,2x Stop 4mmx Comprimento 115mm
05.5.0015.112E	Broca Ø 1,2x Stop 4mmx Comprimento 120mm
05.5.0016.112E	Broca Ø 1,2x Stop 4mmx Comprimento 125mm
05.5.0017.112E	Broca Ø 1,2x Stop 4mmx Comprimento 130mm
05.5.0018.112E	Broca Ø 1,2x Stop 4mmx Comprimento 135mm
05.5.0019.112E	Broca Ø 1,2x Stop 4mmx Comprimento 140mm
05.5.0020.112E	Broca Ø 1,2x Stop 4mmx Comprimento 145mm
05.5.0021.112E	Broca Ø 1,2x Stop 4mmx Comprimento 150mm
05.5.0022.112E	Broca Ø 1,2x Stop 5mmx Comprimento 50mm
05.5.0023.112E	Broca Ø 1,2x Stop 5mmx Comprimento 55mm
05.5.0024.112E	Broca Ø 1,2x Stop 5mmx Comprimento 60mm
05.5.0025.112E	Broca Ø 1,2x Stop 5mmx Comprimento 65mm
05.5.0026.112E	Broca Ø 1,2x Stop 5mmx Comprimento 70mm
05.5.0027.112E	Broca Ø 1,2x Stop 5mmx Comprimento 75mm
05.5.0028.112E	Broca Ø 1,2x Stop 5mmx Comprimento 80mm
05.5.0029.112E	Broca Ø 1,2x Stop 5mmx Comprimento 85mm
05.5.0030.112E	Broca Ø 1,2x Stop 5mmx Comprimento 90mm
05.5.0031.112E	Broca Ø 1,2x Stop 5mmx Comprimento 95mm
05.5.0032.112E	Broca Ø 1,2x Stop 5mmx Comprimento 100mm
05.5.0033.112E	Broca Ø 1,2x Stop 5mmx Comprimento 105mm
05.5.0034.112E	Broca Ø 1,2x Stop 5mmx Comprimento 110mm
05.5.0035.112E	Broca Ø 1,2x Stop 5mmx Comprimento 115mm
05.5.0036.112E	Broca Ø 1,2x Stop 5mmx Comprimento 120mm
05.5.0037.112E	Broca Ø 1,2x Stop 5mmx Comprimento 125mm
05.5.0038.112E	Broca Ø 1,2x Stop 5mmx Comprimento 130mm
05.5.0039.112E	Broca Ø 1,2x Stop 5mmx Comprimento 135mm
05.5.0040.112E	Broca Ø 1,2x Stop 5mmx Comprimento 140mm
05.5.0041.112E	Broca Ø 1,2x Stop 5mmx Comprimento 145mm
05.5.0042.112E	Broca Ø 1,2x Stop 5mmx Comprimento 150mm
05.5.0043.112E	Broca Ø 1,2x Stop 6mmx Comprimento 50mm
05.5.0044.112E	Broca Ø 1,2x Stop 6mmx Comprimento 55mm
05.5.0045.112E	Broca Ø 1,2x Stop 6mmx Comprimento 60mm
05.5.0046.112E	Broca Ø 1,2x Stop 6mmx Comprimento 65mm
05.5.0047.112E	Broca Ø 1,2x Stop 6mmx Comprimento 70mm
05.5.0048.112E	Broca Ø 1,2x Stop 6mmx Comprimento 75mm
05.5.0049.112E	Broca Ø 1,2x Stop 6mmx Comprimento 80mm
05.5.0050.112E	Broca Ø 1,2x Stop 6mmx Comprimento 85mm
05.5.0051.112E	Broca Ø 1,2x Stop 6mmx Comprimento 90mm
05.5.0052.112E	Broca Ø 1,2x Stop 6mmx Comprimento 95mm

05.5.0053.112E	Broca Ø 1,2x Stop 6mmx Comprimento 100mm
05.5.0054.112E	Broca Ø 1,2x Stop 6mmx Comprimento 105mm
05.5.0055.112E	Broca Ø 1,2x Stop 6mmx Comprimento 110mm
05.5.0056.112E	Broca Ø 1,2x Stop 6mmx Comprimento 115mm
05.5.0057.112E	Broca Ø 1,2x Stop 6mmx Comprimento 120mm
05.5.0058.112E	Broca Ø 1,2x Stop 6mmx Comprimento 125mm
05.5.0059.112E	Broca Ø 1,2x Stop 6mmx Comprimento 130mm
05.5.0060.112E	Broca Ø 1,2x Stop 6mmx Comprimento 135mm
05.5.0061.112E	Broca Ø 1,2x Stop 6mmx Comprimento 140mm
05.5.0062.112E	Broca Ø 1,2x Stop 6mmx Comprimento 145mm
05.5.0063.112E	Broca Ø 1,2x Stop 6mmx Comprimento 150mm
05.5.0064.112E	Broca Ø 1,2x Stop 7mmx Comprimento 50mm
05.5.0065.112E	Broca Ø 1,2x Stop 7mmx Comprimento 55mm
05.5.0066.112E	Broca Ø 1,2x Stop 7mmx Comprimento 60mm
05.5.0067.112E	Broca Ø 1,2x Stop 7mmx Comprimento 65mm
05.5.0068.112E	Broca Ø 1,2x Stop 7mmx Comprimento 70mm
05.5.0069.112E	Broca Ø 1,2x Stop 7mmx Comprimento 75mm
05.5.0070.112E	Broca Ø 1,2x Stop 7mmx Comprimento 80mm
05.5.0071.112E	Broca Ø 1,2x Stop 7mmx Comprimento 85mm
05.5.0072.112E	Broca Ø 1,2x Stop 7mmx Comprimento 90mm
05.5.0073.112E	Broca Ø 1,2x Stop 7mmx Comprimento 95mm
05.5.0074.112E	Broca Ø 1,2x Stop 7mmx Comprimento 100mm
05.5.0075.112E	Broca Ø 1,2x Stop 7mmx Comprimento 105mm
05.5.0076.112E	Broca Ø 1,2x Stop 7mmx Comprimento 110mm
05.5.0077.112E	Broca Ø 1,2x Stop 7mmx Comprimento 115mm
05.5.0078.112E	Broca Ø 1,2x Stop 7mmx Comprimento 120mm
05.5.0079.112E	Broca Ø 1,2x Stop 7mmx Comprimento 125mm
05.5.0080.112E	Broca Ø 1,2x Stop 7mmx Comprimento 130mm
05.5.0081.112E	Broca Ø 1,2x Stop 7mmx Comprimento 135mm
05.5.0082.112E	Broca Ø 1,2x Stop 7mmx Comprimento 140mm
05.5.0083.112E	Broca Ø 1,2x Stop 7mmx Comprimento 145mm
05.5.0084.112E	Broca Ø 1,2x Stop 7mmx Comprimento 150mm
05.5.0085.112E	Broca Ø 1,2x Stop 8mmx Comprimento 50mm
05.5.0086.112E	Broca Ø 1,2x Stop 8mmx Comprimento 55mm
05.5.0087.112E	Broca Ø 1,2x Stop 8mmx Comprimento 60mm
05.5.0088.112E	Broca Ø 1,2x Stop 8mmx Comprimento 65mm
05.5.0089.112E	Broca Ø 1,2x Stop 8mmx Comprimento 70mm
05.5.0090.112E	Broca Ø 1,2x Stop 8mmx Comprimento 75mm
05.5.0091.112E	Broca Ø 1,2x Stop 8mmx Comprimento 80mm
05.5.0092.112E	Broca Ø 1,2x Stop 8mmx Comprimento 85mm
05.5.0093.112E	Broca Ø 1,2x Stop 8mmx Comprimento 90mm
05.5.0094.112E	Broca Ø 1,2x Stop 8mmx Comprimento 95mm
05.5.0095.112E	Broca Ø 1,2x Stop 8mmx Comprimento 100mm
05.5.0096.112E	Broca Ø 1,2x Stop 8mmx Comprimento 105mm
05.5.0097.112E	Broca Ø 1,2x Stop 8mmx Comprimento 110mm
05.5.0098.112E	Broca Ø 1,2x Stop 8mmx Comprimento 115mm
05.5.0099.112E	Broca Ø 1,2x Stop 8mmx Comprimento 120mm
05.5.0100.112E	Broca Ø 1,2x Stop 8mmx Comprimento 125mm
05.5.0101.112E	Broca Ø 1,2x Stop 8mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0102.112E	Broca Ø 1,2x Stop 8mmx Comprimento 135mm
05.5.0103.112E	Broca Ø 1,2x Stop 8mmx Comprimento 140mm
05.5.0104.112E	Broca Ø 1,2x Stop 8mmx Comprimento 145mm
05.5.0105.112E	Broca Ø 1,2x Stop 8mmx Comprimento 150mm
05.5.0106.112E	Broca Ø 1,2x Stop 9mmx Comprimento 50mm
05.5.0107.112E	Broca Ø 1,2x Stop 9mmx Comprimento 55mm
05.5.0108.112E	Broca Ø 1,2x Stop 9mmx Comprimento 60mm
05.5.0109.112E	Broca Ø 1,2x Stop 9mmx Comprimento 65mm
05.5.0110.112E	Broca Ø 1,2x Stop 9mmx Comprimento 70mm
05.5.0111.112E	Broca Ø 1,2x Stop 9mmx Comprimento 75mm
05.5.0112.112E	Broca Ø 1,2x Stop 9mmx Comprimento 80mm
05.5.0113.112E	Broca Ø 1,2x Stop 9mmx Comprimento 85mm
05.5.0114.112E	Broca Ø 1,2x Stop 9mmx Comprimento 90mm
05.5.0115.112E	Broca Ø 1,2x Stop 9mmx Comprimento 95mm
05.5.0116.112E	Broca Ø 1,2x Stop 9mmx Comprimento 100mm
05.5.0117.112E	Broca Ø 1,2x Stop 9mmx Comprimento 105mm
05.5.0118.112E	Broca Ø 1,2x Stop 9mmx Comprimento 110mm
05.5.0119.112E	Broca Ø 1,2x Stop 9mmx Comprimento 115mm
05.5.0120.112E	Broca Ø 1,2x Stop 9mmx Comprimento 120mm
05.5.0121.112E	Broca Ø 1,2x Stop 9mmx Comprimento 125mm
05.5.0122.112E	Broca Ø 1,2x Stop 9mmx Comprimento 130mm
05.5.0123.112E	Broca Ø 1,2x Stop 9mmx Comprimento 135mm
05.5.0124.112E	Broca Ø 1,2x Stop 9mmx Comprimento 140mm
05.5.0125.112E	Broca Ø 1,2x Stop 9mmx Comprimento 145mm
05.5.0126.112E	Broca Ø 1,2x Stop 9mmx Comprimento 150mm
05.5.0127.112E	Broca Ø 1,2x Stop 10mmx Comprimento 50mm
05.5.0128.112E	Broca Ø 1,2x Stop 10mmx Comprimento 55mm
05.5.0129.112E	Broca Ø 1,2x Stop 10mmx Comprimento 60mm
05.5.0130.112E	Broca Ø 1,2x Stop 10mmx Comprimento 65mm
05.5.0131.112E	Broca Ø 1,2x Stop 10mmx Comprimento 70mm
05.5.0132.112E	Broca Ø 1,2x Stop 10mmx Comprimento 75mm
05.5.0133.112E	Broca Ø 1,2x Stop 10mmx Comprimento 80mm
05.5.0134.112E	Broca Ø 1,2x Stop 10mmx Comprimento 85mm
05.5.0135.112E	Broca Ø 1,2x Stop 10mmx Comprimento 90mm
05.5.0136.112E	Broca Ø 1,2x Stop 10mmx Comprimento 95mm
05.5.0137.112E	Broca Ø 1,2x Stop 10mmx Comprimento 100mm
05.5.0138.112E	Broca Ø 1,2x Stop 10mmx Comprimento 105mm
05.5.0139.112E	Broca Ø 1,2x Stop 10mmx Comprimento 110mm
05.5.0140.112E	Broca Ø 1,2x Stop 10mmx Comprimento 115mm
05.5.0141.112E	Broca Ø 1,2x Stop 10mmx Comprimento 120mm
05.5.0142.112E	Broca Ø 1,2x Stop 10mmx Comprimento 125mm
05.5.0143.112E	Broca Ø 1,2x Stop 10mmx Comprimento 130mm
05.5.0144.112E	Broca Ø 1,2x Stop 10mmx Comprimento 135mm
05.5.0145.112E	Broca Ø 1,2x Stop 10mmx Comprimento 140mm
05.5.0146.112E	Broca Ø 1,2x Stop 10mmx Comprimento 145mm
05.5.0147.112E	Broca Ø 1,2x Stop 10mmx Comprimento 150mm
05.5.0148.112E	Broca Ø 1,2x Stop 11mmx Comprimento 50mm
05.5.0149.112E	Broca Ø 1,2x Stop 11mmx Comprimento 55mm
05.5.0150.112E	Broca Ø 1,2x Stop 11mmx Comprimento 60mm

05.5.0151.112E	Broca Ø 1,2x Stop 11mmx Comprimento 65mm
05.5.0152.112E	Broca Ø 1,2x Stop 11mmx Comprimento 70mm
05.5.0153.112E	Broca Ø 1,2x Stop 11mmx Comprimento 75mm
05.5.0154.112E	Broca Ø 1,2x Stop 11mmx Comprimento 80mm
05.5.0155.112E	Broca Ø 1,2x Stop 11mmx Comprimento 85mm
05.5.0156.112E	Broca Ø 1,2x Stop 11mmx Comprimento 90mm
05.5.0157.112E	Broca Ø 1,2x Stop 11mmx Comprimento 95mm
05.5.0158.112E	Broca Ø 1,2x Stop 11mmx Comprimento 100mm
05.5.0159.112E	Broca Ø 1,2x Stop 11mmx Comprimento 105mm
05.5.0160.112E	Broca Ø 1,2x Stop 11mmx Comprimento 110mm
05.5.0161.112E	Broca Ø 1,2x Stop 11mmx Comprimento 115mm
05.5.0162.112E	Broca Ø 1,2x Stop 11mmx Comprimento 120mm
05.5.0163.112E	Broca Ø 1,2x Stop 11mmx Comprimento 125mm
05.5.0164.112E	Broca Ø 1,2x Stop 11mmx Comprimento 130mm
05.5.0165.112E	Broca Ø 1,2x Stop 11mmx Comprimento 135mm
05.5.0166.112E	Broca Ø 1,2x Stop 11mmx Comprimento 140mm
05.5.0167.112E	Broca Ø 1,2x Stop 11mmx Comprimento 145mm
05.5.0168.112E	Broca Ø 1,2x Stop 11mmx Comprimento 150mm
05.5.0169.112E	Broca Ø 1,2x Stop 12mmx Comprimento 50mm
05.5.0170.112E	Broca Ø 1,2x Stop 12mmx Comprimento 55mm
05.5.0171.112E	Broca Ø 1,2x Stop 12mmx Comprimento 60mm
05.5.0172.112E	Broca Ø 1,2x Stop 12mmx Comprimento 65mm
05.5.0173.112E	Broca Ø 1,2x Stop 12mmx Comprimento 70mm
05.5.0174.112E	Broca Ø 1,2x Stop 12mmx Comprimento 75mm
05.5.0175.112E	Broca Ø 1,2x Stop 12mmx Comprimento 80mm
05.5.0176.112E	Broca Ø 1,2x Stop 12mmx Comprimento 85mm
05.5.0177.112E	Broca Ø 1,2x Stop 12mmx Comprimento 90mm
05.5.0178.112E	Broca Ø 1,2x Stop 12mmx Comprimento 95mm
05.5.0179.112E	Broca Ø 1,2x Stop 12mmx Comprimento 100mm
05.5.0180.112E	Broca Ø 1,2x Stop 12mmx Comprimento 105mm
05.5.0181.112E	Broca Ø 1,2x Stop 12mmx Comprimento 110mm
05.5.0182.112E	Broca Ø 1,2x Stop 12mmx Comprimento 115mm
05.5.0183.112E	Broca Ø 1,2x Stop 12mmx Comprimento 120mm
05.5.0184.112E	Broca Ø 1,2x Stop 12mmx Comprimento 125mm
05.5.0185.112E	Broca Ø 1,2x Stop 12mmx Comprimento 130mm
05.5.0186.112E	Broca Ø 1,2x Stop 12mmx Comprimento 135mm
05.5.0187.112E	Broca Ø 1,2x Stop 12mmx Comprimento 140mm
05.5.0188.112E	Broca Ø 1,2x Stop 12mmx Comprimento 145mm
05.5.0189.112E	Broca Ø 1,2x Stop 12mmx Comprimento 150mm
05.5.0190.112E	Broca Ø 1,2x Stop 13mmx Comprimento 50mm
05.5.0191.112E	Broca Ø 1,2x Stop 13mmx Comprimento 55mm
05.5.0192.112E	Broca Ø 1,2x Stop 13mmx Comprimento 60mm
05.5.0193.112E	Broca Ø 1,2x Stop 13mmx Comprimento 65mm
05.5.0194.112E	Broca Ø 1,2x Stop 13mmx Comprimento 70mm
05.5.0195.112E	Broca Ø 1,2x Stop 13mmx Comprimento 75mm
05.5.0196.112E	Broca Ø 1,2x Stop 13mmx Comprimento 80mm
05.5.0197.112E	Broca Ø 1,2x Stop 13mmx Comprimento 85mm
05.5.0198.112E	Broca Ø 1,2x Stop 13mmx Comprimento 90mm
05.5.0199.112E	Broca Ø 1,2x Stop 13mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0200.112E	Broca Ø 1,2x Stop 13mmx Comprimento 100mm
05.5.0201.112E	Broca Ø 1,2x Stop 13mmx Comprimento 105mm
05.5.0202.112E	Broca Ø 1,2x Stop 13mmx Comprimento 110mm
05.5.0203.112E	Broca Ø 1,2x Stop 13mmx Comprimento 115mm
05.5.0204.112E	Broca Ø 1,2x Stop 13mmx Comprimento 120mm
05.5.0205.112E	Broca Ø 1,2x Stop 13mmx Comprimento 125mm
05.5.0206.112E	Broca Ø 1,2x Stop 13mmx Comprimento 130mm
05.5.0207.112E	Broca Ø 1,2x Stop 13mmx Comprimento 135mm
05.5.0208.112E	Broca Ø 1,2x Stop 13mmx Comprimento 140mm
05.5.0209.112E	Broca Ø 1,2x Stop 13mmx Comprimento 145mm
05.5.0210.112E	Broca Ø 1,2x Stop 13mmx Comprimento 150mm
05.5.0211.112E	Broca Ø 1,2x Stop 14mmx Comprimento 50mm
05.5.0212.112E	Broca Ø 1,2x Stop 14mmx Comprimento 55mm
05.5.0213.112E	Broca Ø 1,2x Stop 14mmx Comprimento 60mm
05.5.0214.112E	Broca Ø 1,2x Stop 14mmx Comprimento 65mm
05.5.0215.112E	Broca Ø 1,2x Stop 14mmx Comprimento 70mm
05.5.0216.112E	Broca Ø 1,2x Stop 14mmx Comprimento 75mm
05.5.0217.112E	Broca Ø 1,2x Stop 14mmx Comprimento 80mm
05.5.0218.112E	Broca Ø 1,2x Stop 14mmx Comprimento 85mm
05.5.0219.112E	Broca Ø 1,2x Stop 14mmx Comprimento 90mm
05.5.0220.112E	Broca Ø 1,2x Stop 14mmx Comprimento 95mm
05.5.0221.112E	Broca Ø 1,2x Stop 14mmx Comprimento 100mm
05.5.0222.112E	Broca Ø 1,2x Stop 14mmx Comprimento 105mm
05.5.0223.112E	Broca Ø 1,2x Stop 14mmx Comprimento 110mm
05.5.0224.112E	Broca Ø 1,2x Stop 14mmx Comprimento 115mm
05.5.0225.112E	Broca Ø 1,2x Stop 14mmx Comprimento 120mm
05.5.0226.112E	Broca Ø 1,2x Stop 14mmx Comprimento 125mm
05.5.0227.112E	Broca Ø 1,2x Stop 14mmx Comprimento 130mm
05.5.0228.112E	Broca Ø 1,2x Stop 14mmx Comprimento 135mm
05.5.0229.112E	Broca Ø 1,2x Stop 14mmx Comprimento 140mm
05.5.0230.112E	Broca Ø 1,2x Stop 14mmx Comprimento 145mm
05.5.0231.112E	Broca Ø 1,2x Stop 14mmx Comprimento 150mm
05.5.0232.112E	Broca Ø 1,2x Stop 15mmx Comprimento 50mm
05.5.0233.112E	Broca Ø 1,2x Stop 15mmx Comprimento 55mm
05.5.0234.112E	Broca Ø 1,2x Stop 15mmx Comprimento 60mm
05.5.0235.112E	Broca Ø 1,2x Stop 15mmx Comprimento 65mm
05.5.0236.112E	Broca Ø 1,2x Stop 15mmx Comprimento 70mm
05.5.0237.112E	Broca Ø 1,2x Stop 15mmx Comprimento 75mm
05.5.0238.112E	Broca Ø 1,2x Stop 15mmx Comprimento 80mm
05.5.0239.112E	Broca Ø 1,2x Stop 15mmx Comprimento 85mm
05.5.0240.112E	Broca Ø 1,2x Stop 15mmx Comprimento 90mm
05.5.0241.112E	Broca Ø 1,2x Stop 15mmx Comprimento 95mm
05.5.0242.112E	Broca Ø 1,2x Stop 15mmx Comprimento 100mm
05.5.0243.112E	Broca Ø 1,2x Stop 15mmx Comprimento 105mm
05.5.0244.112E	Broca Ø 1,2x Stop 15mmx Comprimento 110mm
05.5.0245.112E	Broca Ø 1,2x Stop 15mmx Comprimento 115mm
05.5.0246.112E	Broca Ø 1,2x Stop 15mmx Comprimento 120mm
05.5.0247.112E	Broca Ø 1,2x Stop 15mmx Comprimento 125mm
05.5.0248.112E	Broca Ø 1,2x Stop 15mmx Comprimento 130mm

05.5.0249.112E	Broca Ø 1,2x Stop 15mmx Comprimento 135mm
05.5.0250.112E	Broca Ø 1,2x Stop 15mmx Comprimento 140mm
05.5.0251.112E	Broca Ø 1,2x Stop 15mmx Comprimento 145mm
05.5.0252.112E	Broca Ø 1,2x Stop 15mmx Comprimento 150mm
05.5.0253.112E	Broca Ø 1,2x Stop 16mmx Comprimento 50mm
05.5.0254.112E	Broca Ø 1,2x Stop 16mmx Comprimento 55mm
05.5.0255.112E	Broca Ø 1,2x Stop 16mmx Comprimento 60mm
05.5.0256.112E	Broca Ø 1,2x Stop 16mmx Comprimento 65mm
05.5.0257.112E	Broca Ø 1,2x Stop 16mmx Comprimento 70mm
05.5.0258.112E	Broca Ø 1,2x Stop 16mmx Comprimento 75mm
05.5.0259.112E	Broca Ø 1,2x Stop 16mmx Comprimento 80mm
05.5.0260.112E	Broca Ø 1,2x Stop 16mmx Comprimento 85mm
05.5.0261.112E	Broca Ø 1,2x Stop 16mmx Comprimento 90mm
05.5.0262.112E	Broca Ø 1,2x Stop 16mmx Comprimento 95mm
05.5.0263.112E	Broca Ø 1,2x Stop 16mmx Comprimento 100mm
05.5.0264.112E	Broca Ø 1,2x Stop 16mmx Comprimento 105mm
05.5.0265.112E	Broca Ø 1,2x Stop 16mmx Comprimento 110mm
05.5.0266.112E	Broca Ø 1,2x Stop 16mmx Comprimento 115mm
05.5.0267.112E	Broca Ø 1,2x Stop 16mmx Comprimento 120mm
05.5.0268.112E	Broca Ø 1,2x Stop 16mmx Comprimento 125mm
05.5.0269.112E	Broca Ø 1,2x Stop 16mmx Comprimento 130mm
05.5.0270.112E	Broca Ø 1,2x Stop 16mmx Comprimento 135mm
05.5.0271.112E	Broca Ø 1,2x Stop 16mmx Comprimento 140mm
05.5.0272.112E	Broca Ø 1,2x Stop 16mmx Comprimento 145mm
05.5.0273.112E	Broca Ø 1,2x Stop 16mmx Comprimento 150mm
05.5.0274.112E	Broca Ø 1,2x Stop 17mmx Comprimento 50mm
05.5.0275.112E	Broca Ø 1,2x Stop 17mmx Comprimento 55mm
05.5.0276.112E	Broca Ø 1,2x Stop 17mmx Comprimento 60mm
05.5.0277.112E	Broca Ø 1,2x Stop 17mmx Comprimento 65mm
05.5.0278.112E	Broca Ø 1,2x Stop 17mmx Comprimento 70mm
05.5.0279.112E	Broca Ø 1,2x Stop 17mmx Comprimento 75mm
05.5.0280.112E	Broca Ø 1,2x Stop 17mmx Comprimento 80mm
05.5.0281.112E	Broca Ø 1,2x Stop 17mmx Comprimento 85mm
05.5.0282.112E	Broca Ø 1,2x Stop 17mmx Comprimento 90mm
05.5.0283.112E	Broca Ø 1,2x Stop 17mmx Comprimento 95mm
05.5.0284.112E	Broca Ø 1,2x Stop 17mmx Comprimento 100mm
05.5.0285.112E	Broca Ø 1,2x Stop 17mmx Comprimento 105mm
05.5.0286.112E	Broca Ø 1,2x Stop 17mmx Comprimento 110mm
05.5.0287.112E	Broca Ø 1,2x Stop 17mmx Comprimento 115mm
05.5.0288.112E	Broca Ø 1,2x Stop 17mmx Comprimento 120mm
05.5.0289.112E	Broca Ø 1,2x Stop 17mmx Comprimento 125mm
05.5.0290.112E	Broca Ø 1,2x Stop 17mmx Comprimento 130mm
05.5.0291.112E	Broca Ø 1,2x Stop 17mmx Comprimento 135mm
05.5.0292.112E	Broca Ø 1,2x Stop 17mmx Comprimento 140mm
05.5.0293.112E	Broca Ø 1,2x Stop 17mmx Comprimento 145mm
05.5.0294.112E	Broca Ø 1,2x Stop 17mmx Comprimento 150mm
05.5.0295.112E	Broca Ø 1,2x Stop 18mmx Comprimento 50mm
05.5.0296.112E	Broca Ø 1,2x Stop 18mmx Comprimento 55mm
05.5.0297.112E	Broca Ø 1,2x Stop 18mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0298.112E	Broca Ø 1,2x Stop 18mmx Comprimento 65mm
05.5.0299.112E	Broca Ø 1,2x Stop 18mmx Comprimento 70mm
05.5.0300.112E	Broca Ø 1,2x Stop 18mmx Comprimento 75mm
05.5.0301.112E	Broca Ø 1,2x Stop 18mmx Comprimento 80mm
05.5.0302.112E	Broca Ø 1,2x Stop 18mmx Comprimento 85mm
05.5.0303.112E	Broca Ø 1,2x Stop 18mmx Comprimento 90mm
05.5.0304.112E	Broca Ø 1,2x Stop 18mmx Comprimento 95mm
05.5.0305.112E	Broca Ø 1,2x Stop 18mmx Comprimento 100mm
05.5.0306.112E	Broca Ø 1,2x Stop 18mmx Comprimento 105mm
05.5.0307.112E	Broca Ø 1,2x Stop 18mmx Comprimento 110mm
05.5.0308.112E	Broca Ø 1,2x Stop 18mmx Comprimento 115mm
05.5.0309.112E	Broca Ø 1,2x Stop 18mmx Comprimento 120mm
05.5.0310.112E	Broca Ø 1,2x Stop 18mmx Comprimento 125mm
05.5.0311.112E	Broca Ø 1,2x Stop 18mmx Comprimento 130mm
05.5.0312.112E	Broca Ø 1,2x Stop 18mmx Comprimento 135mm
05.5.0313.112E	Broca Ø 1,2x Stop 18mmx Comprimento 140mm
05.5.0314.112E	Broca Ø 1,2x Stop 18mmx Comprimento 145mm
05.5.0315.112E	Broca Ø 1,2x Stop 18mmx Comprimento 150mm
05.5.0316.112E	Broca Ø 1,2x Stop 19mmx Comprimento 50mm
05.5.0317.112E	Broca Ø 1,2x Stop 19mmx Comprimento 55mm
05.5.0318.112E	Broca Ø 1,2x Stop 19mmx Comprimento 60mm
05.5.0319.112E	Broca Ø 1,2x Stop 19mmx Comprimento 65mm
05.5.0320.112E	Broca Ø 1,2x Stop 19mmx Comprimento 70mm
05.5.0321.112E	Broca Ø 1,2x Stop 19mmx Comprimento 75mm
05.5.0322.112E	Broca Ø 1,2x Stop 19mmx Comprimento 80mm
05.5.0323.112E	Broca Ø 1,2x Stop 19mmx Comprimento 85mm
05.5.0324.112E	Broca Ø 1,2x Stop 19mmx Comprimento 90mm
05.5.0325.112E	Broca Ø 1,2x Stop 19mmx Comprimento 95mm
05.5.0326.112E	Broca Ø 1,2x Stop 19mmx Comprimento 100mm
05.5.0327.112E	Broca Ø 1,2x Stop 19mmx Comprimento 105mm
05.5.0328.112E	Broca Ø 1,2x Stop 19mmx Comprimento 110mm
05.5.0329.112E	Broca Ø 1,2x Stop 19mmx Comprimento 115mm
05.5.0330.112E	Broca Ø 1,2x Stop 19mmx Comprimento 120mm
05.5.0331.112E	Broca Ø 1,2x Stop 19mmx Comprimento 125mm
05.5.0332.112E	Broca Ø 1,2x Stop 19mmx Comprimento 130mm
05.5.0333.112E	Broca Ø 1,2x Stop 19mmx Comprimento 135mm
05.5.0334.112E	Broca Ø 1,2x Stop 19mmx Comprimento 140mm
05.5.0335.112E	Broca Ø 1,2x Stop 19mmx Comprimento 145mm
05.5.0336.112E	Broca Ø 1,2x Stop 19mmx Comprimento 150mm
05.5.0337.112E	Broca Ø 1,2x Stop 20mmx Comprimento 50mm
05.5.0338.112E	Broca Ø 1,2x Stop 20mmx Comprimento 55mm
05.5.0339.112E	Broca Ø 1,2x Stop 20mmx Comprimento 60mm
05.5.0340.112E	Broca Ø 1,2x Stop 20mmx Comprimento 65mm
05.5.0341.112E	Broca Ø 1,2x Stop 20mmx Comprimento 70mm
05.5.0342.112E	Broca Ø 1,2x Stop 20mmx Comprimento 75mm
05.5.0343.112E	Broca Ø 1,2x Stop 20mmx Comprimento 80mm
05.5.0344.112E	Broca Ø 1,2x Stop 20mmx Comprimento 85mm
05.5.0345.112E	Broca Ø 1,2x Stop 20mmx Comprimento 90mm
05.5.0346.112E	Broca Ø 1,2x Stop 20mmx Comprimento 95mm

05.5.0347.112E	Broca Ø 1,2x Stop 20mmx Comprimento 100mm
05.5.0348.112E	Broca Ø 1,2x Stop 20mmx Comprimento 105mm
05.5.0349.112E	Broca Ø 1,2x Stop 20mmx Comprimento 110mm
05.5.0350.112E	Broca Ø 1,2x Stop 20mmx Comprimento 115mm
05.5.0351.112E	Broca Ø 1,2x Stop 20mmx Comprimento 120mm
05.5.0352.112E	Broca Ø 1,2x Stop 20mmx Comprimento 125mm
05.5.0353.112E	Broca Ø 1,2x Stop 20mmx Comprimento 130mm
05.5.0354.112E	Broca Ø 1,2x Stop 20mmx Comprimento 135mm
05.5.0355.112E	Broca Ø 1,2x Stop 20mmx Comprimento 140mm
05.5.0356.112E	Broca Ø 1,2x Stop 20mmx Comprimento 145mm
05.5.0357.112E	Broca Ø 1,2x Stop 20mmx Comprimento 150mm
05.5.0358.112E	Broca Ø 1,2x Stop 21mmx Comprimento 50mm
05.5.0359.112E	Broca Ø 1,2x Stop 21mmx Comprimento 55mm
05.5.0360.112E	Broca Ø 1,2x Stop 21mmx Comprimento 60mm
05.5.0361.112E	Broca Ø 1,2x Stop 21mmx Comprimento 65mm
05.5.0362.112E	Broca Ø 1,2x Stop 21mmx Comprimento 70mm
05.5.0363.112E	Broca Ø 1,2x Stop 21mmx Comprimento 75mm
05.5.0364.112E	Broca Ø 1,2x Stop 21mmx Comprimento 80mm
05.5.0365.112E	Broca Ø 1,2x Stop 21mmx Comprimento 85mm
05.5.0366.112E	Broca Ø 1,2x Stop 21mmx Comprimento 90mm
05.5.0367.112E	Broca Ø 1,2x Stop 21mmx Comprimento 95mm
05.5.0368.112E	Broca Ø 1,2x Stop 21mmx Comprimento 100mm
05.5.0369.112E	Broca Ø 1,2x Stop 21mmx Comprimento 105mm
05.5.0370.112E	Broca Ø 1,2x Stop 21mmx Comprimento 110mm
05.5.0371.112E	Broca Ø 1,2x Stop 21mmx Comprimento 115mm
05.5.0372.112E	Broca Ø 1,2x Stop 21mmx Comprimento 120mm
05.5.0373.112E	Broca Ø 1,2x Stop 21mmx Comprimento 125mm
05.5.0374.112E	Broca Ø 1,2x Stop 21mmx Comprimento 130mm
05.5.0375.112E	Broca Ø 1,2x Stop 21mmx Comprimento 135mm
05.5.0376.112E	Broca Ø 1,2x Stop 21mmx Comprimento 140mm
05.5.0377.112E	Broca Ø 1,2x Stop 21mmx Comprimento 145mm
05.5.0378.112E	Broca Ø 1,2x Stop 21mmx Comprimento 150mm
05.5.0379.112E	Broca Ø 1,2x Stop 22mmx Comprimento 50mm
05.5.0380.112E	Broca Ø 1,2x Stop 22mmx Comprimento 55mm
05.5.0381.112E	Broca Ø 1,2x Stop 22mmx Comprimento 60mm
05.5.0382.112E	Broca Ø 1,2x Stop 22mmx Comprimento 65mm
05.5.0383.112E	Broca Ø 1,2x Stop 22mmx Comprimento 70mm
05.5.0384.112E	Broca Ø 1,2x Stop 22mmx Comprimento 75mm
05.5.0385.112E	Broca Ø 1,2x Stop 22mmx Comprimento 80mm
05.5.0386.112E	Broca Ø 1,2x Stop 22mmx Comprimento 85mm
05.5.0387.112E	Broca Ø 1,2x Stop 22mmx Comprimento 90mm
05.5.0388.112E	Broca Ø 1,2x Stop 22mmx Comprimento 95mm
05.5.0389.112E	Broca Ø 1,2x Stop 22mmx Comprimento 100mm
05.5.0390.112E	Broca Ø 1,2x Stop 22mmx Comprimento 105mm
05.5.0391.112E	Broca Ø 1,2x Stop 22mmx Comprimento 110mm
05.5.0392.112E	Broca Ø 1,2x Stop 22mmx Comprimento 115mm
05.5.0393.112E	Broca Ø 1,2x Stop 22mmx Comprimento 120mm
05.5.0394.112E	Broca Ø 1,2x Stop 22mmx Comprimento 125mm
05.5.0395.112E	Broca Ø 1,2x Stop 22mmx Comprimento 130mm

05.5.0396.112E	Broca Ø 1,2x Stop 22mmx Comprimento 135mm
05.5.0397.112E	Broca Ø 1,2x Stop 22mmx Comprimento 140mm
05.5.0398.112E	Broca Ø 1,2x Stop 22mmx Comprimento 145mm
05.5.0399.112E	Broca Ø 1,2x Stop 22mmx Comprimento 150mm
05.5.0400.112E	Broca Ø 1,2x Stop 23mmx Comprimento 50mm
05.5.0401.112E	Broca Ø 1,2x Stop 23mmx Comprimento 55mm
05.5.0402.112E	Broca Ø 1,2x Stop 23mmx Comprimento 60mm
05.5.0403.112E	Broca Ø 1,2x Stop 23mmx Comprimento 65mm
05.5.0404.112E	Broca Ø 1,2x Stop 23mmx Comprimento 70mm
05.5.0405.112E	Broca Ø 1,2x Stop 23mmx Comprimento 75mm
05.5.0406.112E	Broca Ø 1,2x Stop 23mmx Comprimento 80mm
05.5.0407.112E	Broca Ø 1,2x Stop 23mmx Comprimento 85mm
05.5.0408.112E	Broca Ø 1,2x Stop 23mmx Comprimento 90mm
05.5.0409.112E	Broca Ø 1,2x Stop 23mmx Comprimento 95mm
05.5.0410.112E	Broca Ø 1,2x Stop 23mmx Comprimento 100mm
05.5.0411.112E	Broca Ø 1,2x Stop 23mmx Comprimento 105mm
05.5.0412.112E	Broca Ø 1,2x Stop 23mmx Comprimento 110mm
05.5.0413.112E	Broca Ø 1,2x Stop 23mmx Comprimento 115mm
05.5.0414.112E	Broca Ø 1,2x Stop 23mmx Comprimento 120mm
05.5.0415.112E	Broca Ø 1,2x Stop 23mmx Comprimento 125mm
05.5.0416.112E	Broca Ø 1,2x Stop 23mmx Comprimento 130mm
05.5.0417.112E	Broca Ø 1,2x Stop 23mmx Comprimento 135mm
05.5.0418.112E	Broca Ø 1,2x Stop 23mmx Comprimento 140mm
05.5.0419.112E	Broca Ø 1,2x Stop 23mmx Comprimento 145mm
05.5.0420.112E	Broca Ø 1,2x Stop 23mmx Comprimento 150mm
05.5.0421.112E	Broca Ø 1,2x Stop 24mmx Comprimento 50mm
05.5.0422.112E	Broca Ø 1,2x Stop 24mmx Comprimento 55mm
05.5.0423.112E	Broca Ø 1,2x Stop 24mmx Comprimento 60mm
05.5.0424.112E	Broca Ø 1,2x Stop 24mmx Comprimento 65mm
05.5.0425.112E	Broca Ø 1,2x Stop 24mmx Comprimento 70mm
05.5.0426.112E	Broca Ø 1,2x Stop 24mmx Comprimento 75mm
05.5.0427.112E	Broca Ø 1,2x Stop 24mmx Comprimento 80mm
05.5.0428.112E	Broca Ø 1,2x Stop 24mmx Comprimento 85mm
05.5.0429.112E	Broca Ø 1,2x Stop 24mmx Comprimento 90mm
05.5.0430.112E	Broca Ø 1,2x Stop 24mmx Comprimento 95mm
05.5.0431.112E	Broca Ø 1,2x Stop 24mmx Comprimento 100mm
05.5.0432.112E	Broca Ø 1,2x Stop 24mmx Comprimento 105mm
05.5.0433.112E	Broca Ø 1,2x Stop 24mmx Comprimento 110mm
05.5.0434.112E	Broca Ø 1,2x Stop 24mmx Comprimento 115mm
05.5.0435.112E	Broca Ø 1,2x Stop 24mmx Comprimento 120mm
05.5.0436.112E	Broca Ø 1,2x Stop 24mmx Comprimento 125mm
05.5.0437.112E	Broca Ø 1,2x Stop 24mmx Comprimento 130mm
05.5.0438.112E	Broca Ø 1,2x Stop 24mmx Comprimento 135mm
05.5.0439.112E	Broca Ø 1,2x Stop 24mmx Comprimento 140mm
05.5.0440.112E	Broca Ø 1,2x Stop 24mmx Comprimento 145mm
05.5.0441.112E	Broca Ø 1,2x Stop 24mmx Comprimento 150mm
05.5.0442.112E	Broca Ø 1,2x Stop 25mmx Comprimento 50mm
05.5.0443.112E	Broca Ø 1,2x Stop 25mmx Comprimento 55mm
05.5.0444.112E	Broca Ø 1,2x Stop 25mmx Comprimento 60mm

05.5.0445.112E	Broca Ø 1,2x Stop 25mmx Comprimento 65mm
05.5.0446.112E	Broca Ø 1,2x Stop 25mmx Comprimento 70mm
05.5.0447.112E	Broca Ø 1,2x Stop 25mmx Comprimento 75mm
05.5.0448.112E	Broca Ø 1,2x Stop 25mmx Comprimento 80mm
05.5.0449.112E	Broca Ø 1,2x Stop 25mmx Comprimento 85mm
05.5.0450.112E	Broca Ø 1,2x Stop 25mmx Comprimento 90mm
05.5.0451.112E	Broca Ø 1,2x Stop 25mmx Comprimento 95mm
05.5.0452.112E	Broca Ø 1,2x Stop 25mmx Comprimento 100mm
05.5.0453.112E	Broca Ø 1,2x Stop 25mmx Comprimento 105mm
05.5.0454.112E	Broca Ø 1,2x Stop 25mmx Comprimento 110mm
05.5.0455.112E	Broca Ø 1,2x Stop 25mmx Comprimento 115mm
05.5.0456.112E	Broca Ø 1,2x Stop 25mmx Comprimento 120mm
05.5.0457.112E	Broca Ø 1,2x Stop 25mmx Comprimento 125mm
05.5.0458.112E	Broca Ø 1,2x Stop 25mmx Comprimento 130mm
05.5.0459.112E	Broca Ø 1,2x Stop 25mmx Comprimento 135mm
05.5.0460.112E	Broca Ø 1,2x Stop 25mmx Comprimento 140mm
05.5.0461.112E	Broca Ø 1,2x Stop 25mmx Comprimento 145mm
05.5.0462.112E	Broca Ø 1,2x Stop 25mmx Comprimento 150mm
05.5.0463.112E	Broca Ø 1,2x Stop 30mmx Comprimento 50mm
05.5.0464.112E	Broca Ø 1,2x Stop 30mmx Comprimento 55mm
05.5.0465.112E	Broca Ø 1,2x Stop 30mmx Comprimento 60mm
05.5.0466.112E	Broca Ø 1,2x Stop 30mmx Comprimento 65mm
05.5.0467.112E	Broca Ø 1,2x Stop 30mmx Comprimento 70mm
05.5.0468.112E	Broca Ø 1,2x Stop 30mmx Comprimento 75mm
05.5.0469.112E	Broca Ø 1,2x Stop 30mmx Comprimento 80mm
05.5.0470.112E	Broca Ø 1,2x Stop 30mmx Comprimento 85mm
05.5.0471.112E	Broca Ø 1,2x Stop 30mmx Comprimento 90mm
05.5.0472.112E	Broca Ø 1,2x Stop 30mmx Comprimento 95mm
05.5.0473.112E	Broca Ø 1,2x Stop 30mmx Comprimento 100mm
05.5.0474.112E	Broca Ø 1,2x Stop 30mmx Comprimento 105mm
05.5.0475.112E	Broca Ø 1,2x Stop 30mmx Comprimento 110mm
05.5.0476.112E	Broca Ø 1,2x Stop 30mmx Comprimento 115mm
05.5.0477.112E	Broca Ø 1,2x Stop 30mmx Comprimento 120mm
05.5.0478.112E	Broca Ø 1,2x Stop 30mmx Comprimento 125mm
05.5.0479.112E	Broca Ø 1,2x Stop 30mmx Comprimento 130mm
05.5.0480.112E	Broca Ø 1,2x Stop 30mmx Comprimento 135mm
05.5.0481.112E	Broca Ø 1,2x Stop 30mmx Comprimento 140mm
05.5.0482.112E	Broca Ø 1,2x Stop 30mmx Comprimento 145mm
05.5.0483.112E	Broca Ø 1,2x Stop 30mmx Comprimento 150mm
05.5.0484.112E	Broca Ø 1,2x Stop 35mmx Comprimento 50mm
05.5.0485.112E	Broca Ø 1,2x Stop 35mmx Comprimento 55mm
05.5.0486.112E	Broca Ø 1,2x Stop 35mmx Comprimento 60mm
05.5.0487.112E	Broca Ø 1,2x Stop 35mmx Comprimento 65mm
05.5.0488.112E	Broca Ø 1,2x Stop 35mmx Comprimento 70mm
05.5.0489.112E	Broca Ø 1,2x Stop 35mmx Comprimento 75mm
05.5.0490.112E	Broca Ø 1,2x Stop 35mmx Comprimento 80mm
05.5.0491.112E	Broca Ø 1,2x Stop 35mmx Comprimento 85mm
05.5.0492.112E	Broca Ø 1,2x Stop 35mmx Comprimento 90mm
05.5.0493.112E	Broca Ø 1,2x Stop 35mmx Comprimento 95mm

05.5.0494.112E	Broca Ø 1,2x Stop 35mmx Comprimento 100mm
05.5.0495.112E	Broca Ø 1,2x Stop 35mmx Comprimento 105mm
05.5.0496.112E	Broca Ø 1,2x Stop 35mmx Comprimento 110mm
05.5.0497.112E	Broca Ø 1,2x Stop 35mmx Comprimento 115mm
05.5.0498.112E	Broca Ø 1,2x Stop 35mmx Comprimento 120mm
05.5.0499.112E	Broca Ø 1,2x Stop 35mmx Comprimento 125mm
05.5.0500.112E	Broca Ø 1,2x Stop 35mmx Comprimento 130mm
05.5.0501.112E	Broca Ø 1,2x Stop 35mmx Comprimento 135mm
05.5.0502.112E	Broca Ø 1,2x Stop 35mmx Comprimento 140mm
05.5.0503.112E	Broca Ø 1,2x Stop 35mmx Comprimento 145mm
05.5.0504.112E	Broca Ø 1,2x Stop 35mmx Comprimento 150mm
05.5.0001.113E	Broca Ø 1,3x Stop 4mmx Comprimento 50mm
05.5.0002.113E	Broca Ø 1,3x Stop 4mmx Comprimento 55mm
05.5.0003.113E	Broca Ø 1,3x Stop 4mmx Comprimento 60mm
05.5.0004.113E	Broca Ø 1,3x Stop 4mmx Comprimento 65mm
05.5.0005.113E	Broca Ø 1,3x Stop 4mmx Comprimento 70mm
05.5.0006.113E	Broca Ø 1,3x Stop 4mmx Comprimento 75mm
05.5.0007.113E	Broca Ø 1,3x Stop 4mmx Comprimento 80mm
05.5.0008.113E	Broca Ø 1,3x Stop 4mmx Comprimento 85mm
05.5.0009.113E	Broca Ø 1,3x Stop 4mmx Comprimento 90mm
05.5.0010.113E	Broca Ø 1,3x Stop 4mmx Comprimento 95mm
05.5.0011.113E	Broca Ø 1,3x Stop 4mmx Comprimento 100mm
05.5.0012.113E	Broca Ø 1,3x Stop 4mmx Comprimento 105mm
05.5.0013.113E	Broca Ø 1,3x Stop 4mmx Comprimento 110mm
05.5.0014.113E	Broca Ø 1,3x Stop 4mmx Comprimento 115mm
05.5.0015.113E	Broca Ø 1,3x Stop 4mmx Comprimento 120mm
05.5.0016.113E	Broca Ø 1,3x Stop 4mmx Comprimento 125mm
05.5.0017.113E	Broca Ø 1,3x Stop 4mmx Comprimento 130mm
05.5.0018.113E	Broca Ø 1,3x Stop 4mmx Comprimento 135mm
05.5.0019.113E	Broca Ø 1,3x Stop 4mmx Comprimento 140mm
05.5.0020.113E	Broca Ø 1,3x Stop 4mmx Comprimento 145mm
05.5.0021.113E	Broca Ø 1,3x Stop 4mmx Comprimento 150mm
05.5.0022.113E	Broca Ø 1,3x Stop 5mmx Comprimento 50mm
05.5.0023.113E	Broca Ø 1,3x Stop 5mmx Comprimento 55mm
05.5.0024.113E	Broca Ø 1,3x Stop 5mmx Comprimento 60mm
05.5.0025.113E	Broca Ø 1,3x Stop 5mmx Comprimento 65mm
05.5.0026.113E	Broca Ø 1,3x Stop 5mmx Comprimento 70mm
05.5.0027.113E	Broca Ø 1,3x Stop 5mmx Comprimento 75mm
05.5.0028.113E	Broca Ø 1,3x Stop 5mmx Comprimento 80mm
05.5.0029.113E	Broca Ø 1,3x Stop 5mmx Comprimento 85mm
05.5.0030.113E	Broca Ø 1,3x Stop 5mmx Comprimento 90mm
05.5.0031.113E	Broca Ø 1,3x Stop 5mmx Comprimento 95mm
05.5.0032.113E	Broca Ø 1,3x Stop 5mmx Comprimento 100mm
05.5.0033.113E	Broca Ø 1,3x Stop 5mmx Comprimento 105mm
05.5.0034.113E	Broca Ø 1,3x Stop 5mmx Comprimento 110mm
05.5.0035.113E	Broca Ø 1,3x Stop 5mmx Comprimento 115mm
05.5.0036.113E	Broca Ø 1,3x Stop 5mmx Comprimento 120mm
05.5.0037.113E	Broca Ø 1,3x Stop 5mmx Comprimento 125mm
05.5.0038.113E	Broca Ø 1,3x Stop 5mmx Comprimento 130mm

05.5.0039.113E	Broca Ø 1,3x Stop 5mmx Comprimento 135mm
05.5.0040.113E	Broca Ø 1,3x Stop 5mmx Comprimento 140mm
05.5.0041.113E	Broca Ø 1,3x Stop 5mmx Comprimento 145mm
05.5.0042.113E	Broca Ø 1,3x Stop 5mmx Comprimento 150mm
05.5.0043.113E	Broca Ø 1,3x Stop 6mmx Comprimento 50mm
05.5.0044.113E	Broca Ø 1,3x Stop 6mmx Comprimento 55mm
05.5.0045.113E	Broca Ø 1,3x Stop 6mmx Comprimento 60mm
05.5.0046.113E	Broca Ø 1,3x Stop 6mmx Comprimento 65mm
05.5.0047.113E	Broca Ø 1,3x Stop 6mmx Comprimento 70mm
05.5.0048.113E	Broca Ø 1,3x Stop 6mmx Comprimento 75mm
05.5.0049.113E	Broca Ø 1,3x Stop 6mmx Comprimento 80mm
05.5.0050.113E	Broca Ø 1,3x Stop 6mmx Comprimento 85mm
05.5.0051.113E	Broca Ø 1,3x Stop 6mmx Comprimento 90mm
05.5.0052.113E	Broca Ø 1,3x Stop 6mmx Comprimento 95mm
05.5.0053.113E	Broca Ø 1,3x Stop 6mmx Comprimento 100mm
05.5.0054.113E	Broca Ø 1,3x Stop 6mmx Comprimento 105mm
05.5.0055.113E	Broca Ø 1,3x Stop 6mmx Comprimento 110mm
05.5.0056.113E	Broca Ø 1,3x Stop 6mmx Comprimento 115mm
05.5.0057.113E	Broca Ø 1,3x Stop 6mmx Comprimento 120mm
05.5.0058.113E	Broca Ø 1,3x Stop 6mmx Comprimento 125mm
05.5.0059.113E	Broca Ø 1,3x Stop 6mmx Comprimento 130mm
05.5.0060.113E	Broca Ø 1,3x Stop 6mmx Comprimento 135mm
05.5.0061.113E	Broca Ø 1,3x Stop 6mmx Comprimento 140mm
05.5.0062.113E	Broca Ø 1,3x Stop 6mmx Comprimento 145mm
05.5.0063.113E	Broca Ø 1,3x Stop 6mmx Comprimento 150mm
05.5.0064.113E	Broca Ø 1,3x Stop 7mmx Comprimento 50mm
05.5.0065.113E	Broca Ø 1,3x Stop 7mmx Comprimento 55mm
05.5.0066.113E	Broca Ø 1,3x Stop 7mmx Comprimento 60mm
05.5.0067.113E	Broca Ø 1,3x Stop 7mmx Comprimento 65mm
05.5.0068.113E	Broca Ø 1,3x Stop 7mmx Comprimento 70mm
05.5.0069.113E	Broca Ø 1,3x Stop 7mmx Comprimento 75mm
05.5.0070.113E	Broca Ø 1,3x Stop 7mmx Comprimento 80mm
05.5.0071.113E	Broca Ø 1,3x Stop 7mmx Comprimento 85mm
05.5.0072.113E	Broca Ø 1,3x Stop 7mmx Comprimento 90mm
05.5.0073.113E	Broca Ø 1,3x Stop 7mmx Comprimento 95mm
05.5.0074.113E	Broca Ø 1,3x Stop 7mmx Comprimento 100mm
05.5.0075.113E	Broca Ø 1,3x Stop 7mmx Comprimento 105mm
05.5.0076.113E	Broca Ø 1,3x Stop 7mmx Comprimento 110mm
05.5.0077.113E	Broca Ø 1,3x Stop 7mmx Comprimento 115mm
05.5.0078.113E	Broca Ø 1,3x Stop 7mmx Comprimento 120mm
05.5.0079.113E	Broca Ø 1,3x Stop 7mmx Comprimento 125mm
05.5.0080.113E	Broca Ø 1,3x Stop 7mmx Comprimento 130mm
05.5.0081.113E	Broca Ø 1,3x Stop 7mmx Comprimento 135mm
05.5.0082.113E	Broca Ø 1,3x Stop 7mmx Comprimento 140mm
05.5.0083.113E	Broca Ø 1,3x Stop 7mmx Comprimento 145mm
05.5.0084.113E	Broca Ø 1,3x Stop 7mmx Comprimento 150mm
05.5.0085.113E	Broca Ø 1,3x Stop 8mmx Comprimento 50mm
05.5.0086.113E	Broca Ø 1,3x Stop 8mmx Comprimento 55mm
05.5.0087.113E	Broca Ø 1,3x Stop 8mmx Comprimento 60mm

05.5.0088.113E	Broca Ø 1,3x Stop 8mmx Comprimento 65mm
05.5.0089.113E	Broca Ø 1,3x Stop 8mmx Comprimento 70mm
05.5.0090.113E	Broca Ø 1,3x Stop 8mmx Comprimento 75mm
05.5.0091.113E	Broca Ø 1,3x Stop 8mmx Comprimento 80mm
05.5.0092.113E	Broca Ø 1,3x Stop 8mmx Comprimento 85mm
05.5.0093.113E	Broca Ø 1,3x Stop 8mmx Comprimento 90mm
05.5.0094.113E	Broca Ø 1,3x Stop 8mmx Comprimento 95mm
05.5.0095.113E	Broca Ø 1,3x Stop 8mmx Comprimento 100mm
05.5.0096.113E	Broca Ø 1,3x Stop 8mmx Comprimento 105mm
05.5.0097.113E	Broca Ø 1,3x Stop 8mmx Comprimento 110mm
05.5.0098.113E	Broca Ø 1,3x Stop 8mmx Comprimento 115mm
05.5.0099.113E	Broca Ø 1,3x Stop 8mmx Comprimento 120mm
05.5.0100.113E	Broca Ø 1,3x Stop 8mmx Comprimento 125mm
05.5.0101.113E	Broca Ø 1,3x Stop 8mmx Comprimento 130mm
05.5.0102.113E	Broca Ø 1,3x Stop 8mmx Comprimento 135mm
05.5.0103.113E	Broca Ø 1,3x Stop 8mmx Comprimento 140mm
05.5.0104.113E	Broca Ø 1,3x Stop 8mmx Comprimento 145mm
05.5.0105.113E	Broca Ø 1,3x Stop 8mmx Comprimento 150mm
05.5.0106.113E	Broca Ø 1,3x Stop 9mmx Comprimento 50mm
05.5.0107.113E	Broca Ø 1,3x Stop 9mmx Comprimento 55mm
05.5.0108.113E	Broca Ø 1,3x Stop 9mmx Comprimento 60mm
05.5.0109.113E	Broca Ø 1,3x Stop 9mmx Comprimento 65mm
05.5.0110.113E	Broca Ø 1,3x Stop 9mmx Comprimento 70mm
05.5.0111.113E	Broca Ø 1,3x Stop 9mmx Comprimento 75mm
05.5.0112.113E	Broca Ø 1,3x Stop 9mmx Comprimento 80mm
05.5.0113.113E	Broca Ø 1,3x Stop 9mmx Comprimento 85mm
05.5.0114.113E	Broca Ø 1,3x Stop 9mmx Comprimento 90mm
05.5.0115.113E	Broca Ø 1,3x Stop 9mmx Comprimento 95mm
05.5.0116.113E	Broca Ø 1,3x Stop 9mmx Comprimento 100mm
05.5.0117.113E	Broca Ø 1,3x Stop 9mmx Comprimento 105mm
05.5.0118.113E	Broca Ø 1,3x Stop 9mmx Comprimento 110mm
05.5.0119.113E	Broca Ø 1,3x Stop 9mmx Comprimento 115mm
05.5.0120.113E	Broca Ø 1,3x Stop 9mmx Comprimento 120mm
05.5.0121.113E	Broca Ø 1,3x Stop 9mmx Comprimento 125mm
05.5.0122.113E	Broca Ø 1,3x Stop 9mmx Comprimento 130mm
05.5.0123.113E	Broca Ø 1,3x Stop 9mmx Comprimento 135mm
05.5.0124.113E	Broca Ø 1,3x Stop 9mmx Comprimento 140mm
05.5.0125.113E	Broca Ø 1,3x Stop 9mmx Comprimento 145mm
05.5.0126.113E	Broca Ø 1,3x Stop 9mmx Comprimento 150mm
05.5.0127.113E	Broca Ø 1,3x Stop 10mmx Comprimento 50mm
05.5.0128.113E	Broca Ø 1,3x Stop 10mmx Comprimento 55mm
05.5.0129.113E	Broca Ø 1,3x Stop 10mmx Comprimento 60mm
05.5.0130.113E	Broca Ø 1,3x Stop 10mmx Comprimento 65mm
05.5.0131.113E	Broca Ø 1,3x Stop 10mmx Comprimento 70mm
05.5.0132.113E	Broca Ø 1,3x Stop 10mmx Comprimento 75mm
05.5.0133.113E	Broca Ø 1,3x Stop 10mmx Comprimento 80mm
05.5.0134.113E	Broca Ø 1,3x Stop 10mmx Comprimento 85mm
05.5.0135.113E	Broca Ø 1,3x Stop 10mmx Comprimento 90mm
05.5.0136.113E	Broca Ø 1,3x Stop 10mmx Comprimento 95mm

05.5.0137.113E	Broca Ø 1,3x Stop 10mmx Comprimento 100mm
05.5.0138.113E	Broca Ø 1,3x Stop 10mmx Comprimento 105mm
05.5.0139.113E	Broca Ø 1,3x Stop 10mmx Comprimento 110mm
05.5.0140.113E	Broca Ø 1,3x Stop 10mmx Comprimento 115mm
05.5.0141.113E	Broca Ø 1,3x Stop 10mmx Comprimento 120mm
05.5.0142.113E	Broca Ø 1,3x Stop 10mmx Comprimento 125mm
05.5.0143.113E	Broca Ø 1,3x Stop 10mmx Comprimento 130mm
05.5.0144.113E	Broca Ø 1,3x Stop 10mmx Comprimento 135mm
05.5.0145.113E	Broca Ø 1,3x Stop 10mmx Comprimento 140mm
05.5.0146.113E	Broca Ø 1,3x Stop 10mmx Comprimento 145mm
05.5.0147.113E	Broca Ø 1,3x Stop 10mmx Comprimento 150mm
05.5.0148.113E	Broca Ø 1,3x Stop 11mmx Comprimento 50mm
05.5.0149.113E	Broca Ø 1,3x Stop 11mmx Comprimento 55mm
05.5.0150.113E	Broca Ø 1,3x Stop 11mmx Comprimento 60mm
05.5.0151.113E	Broca Ø 1,3x Stop 11mmx Comprimento 65mm
05.5.0152.113E	Broca Ø 1,3x Stop 11mmx Comprimento 70mm
05.5.0153.113E	Broca Ø 1,3x Stop 11mmx Comprimento 75mm
05.5.0154.113E	Broca Ø 1,3x Stop 11mmx Comprimento 80mm
05.5.0155.113E	Broca Ø 1,3x Stop 11mmx Comprimento 85mm
05.5.0156.113E	Broca Ø 1,3x Stop 11mmx Comprimento 90mm
05.5.0157.113E	Broca Ø 1,3x Stop 11mmx Comprimento 95mm
05.5.0158.113E	Broca Ø 1,3x Stop 11mmx Comprimento 100mm
05.5.0159.113E	Broca Ø 1,3x Stop 11mmx Comprimento 105mm
05.5.0160.113E	Broca Ø 1,3x Stop 11mmx Comprimento 110mm
05.5.0161.113E	Broca Ø 1,3x Stop 11mmx Comprimento 115mm
05.5.0162.113E	Broca Ø 1,3x Stop 11mmx Comprimento 120mm
05.5.0163.113E	Broca Ø 1,3x Stop 11mmx Comprimento 125mm
05.5.0164.113E	Broca Ø 1,3x Stop 11mmx Comprimento 130mm
05.5.0165.113E	Broca Ø 1,3x Stop 11mmx Comprimento 135mm
05.5.0166.113E	Broca Ø 1,3x Stop 11mmx Comprimento 140mm
05.5.0167.113E	Broca Ø 1,3x Stop 11mmx Comprimento 145mm
05.5.0168.113E	Broca Ø 1,3x Stop 11mmx Comprimento 150mm
05.5.0169.113E	Broca Ø 1,3x Stop 12mmx Comprimento 50mm
05.5.0170.113E	Broca Ø 1,3x Stop 12mmx Comprimento 55mm
05.5.0171.113E	Broca Ø 1,3x Stop 12mmx Comprimento 60mm
05.5.0172.113E	Broca Ø 1,3x Stop 12mmx Comprimento 65mm
05.5.0173.113E	Broca Ø 1,3x Stop 12mmx Comprimento 70mm
05.5.0174.113E	Broca Ø 1,3x Stop 12mmx Comprimento 75mm
05.5.0175.113E	Broca Ø 1,3x Stop 12mmx Comprimento 80mm
05.5.0176.113E	Broca Ø 1,3x Stop 12mmx Comprimento 85mm
05.5.0177.113E	Broca Ø 1,3x Stop 12mmx Comprimento 90mm
05.5.0178.113E	Broca Ø 1,3x Stop 12mmx Comprimento 95mm
05.5.0179.113E	Broca Ø 1,3x Stop 12mmx Comprimento 100mm
05.5.0180.113E	Broca Ø 1,3x Stop 12mmx Comprimento 105mm
05.5.0181.113E	Broca Ø 1,3x Stop 12mmx Comprimento 110mm
05.5.0182.113E	Broca Ø 1,3x Stop 12mmx Comprimento 115mm
05.5.0183.113E	Broca Ø 1,3x Stop 12mmx Comprimento 120mm
05.5.0184.113E	Broca Ø 1,3x Stop 12mmx Comprimento 125mm
05.5.0185.113E	Broca Ø 1,3x Stop 12mmx Comprimento 130mm

05.5.0186.113E	Broca Ø 1,3x Stop 12mmx Comprimento 135mm
05.5.0187.113E	Broca Ø 1,3x Stop 12mmx Comprimento 140mm
05.5.0188.113E	Broca Ø 1,3x Stop 12mmx Comprimento 145mm
05.5.0189.113E	Broca Ø 1,3x Stop 12mmx Comprimento 150mm
05.5.0190.113E	Broca Ø 1,3x Stop 13mmx Comprimento 50mm
05.5.0191.113E	Broca Ø 1,3x Stop 13mmx Comprimento 55mm
05.5.0192.113E	Broca Ø 1,3x Stop 13mmx Comprimento 60mm
05.5.0193.113E	Broca Ø 1,3x Stop 13mmx Comprimento 65mm
05.5.0194.113E	Broca Ø 1,3x Stop 13mmx Comprimento 70mm
05.5.0195.113E	Broca Ø 1,3x Stop 13mmx Comprimento 75mm
05.5.0196.113E	Broca Ø 1,3x Stop 13mmx Comprimento 80mm
05.5.0197.113E	Broca Ø 1,3x Stop 13mmx Comprimento 85mm
05.5.0198.113E	Broca Ø 1,3x Stop 13mmx Comprimento 90mm
05.5.0199.113E	Broca Ø 1,3x Stop 13mmx Comprimento 95mm
05.5.0200.113E	Broca Ø 1,3x Stop 13mmx Comprimento 100mm
05.5.0201.113E	Broca Ø 1,3x Stop 13mmx Comprimento 105mm
05.5.0202.113E	Broca Ø 1,3x Stop 13mmx Comprimento 110mm
05.5.0203.113E	Broca Ø 1,3x Stop 13mmx Comprimento 115mm
05.5.0204.113E	Broca Ø 1,3x Stop 13mmx Comprimento 120mm
05.5.0205.113E	Broca Ø 1,3x Stop 13mmx Comprimento 125mm
05.5.0206.113E	Broca Ø 1,3x Stop 13mmx Comprimento 130mm
05.5.0207.113E	Broca Ø 1,3x Stop 13mmx Comprimento 135mm
05.5.0208.113E	Broca Ø 1,3x Stop 13mmx Comprimento 140mm
05.5.0209.113E	Broca Ø 1,3x Stop 13mmx Comprimento 145mm
05.5.0210.113E	Broca Ø 1,3x Stop 13mmx Comprimento 150mm
05.5.0211.113E	Broca Ø 1,3x Stop 14mmx Comprimento 50mm
05.5.0212.113E	Broca Ø 1,3x Stop 14mmx Comprimento 55mm
05.5.0213.113E	Broca Ø 1,3x Stop 14mmx Comprimento 60mm
05.5.0214.113E	Broca Ø 1,3x Stop 14mmx Comprimento 65mm
05.5.0215.113E	Broca Ø 1,3x Stop 14mmx Comprimento 70mm
05.5.0216.113E	Broca Ø 1,3x Stop 14mmx Comprimento 75mm
05.5.0217.113E	Broca Ø 1,3x Stop 14mmx Comprimento 80mm
05.5.0218.113E	Broca Ø 1,3x Stop 14mmx Comprimento 85mm
05.5.0219.113E	Broca Ø 1,3x Stop 14mmx Comprimento 90mm
05.5.0220.113E	Broca Ø 1,3x Stop 14mmx Comprimento 95mm
05.5.0221.113E	Broca Ø 1,3x Stop 14mmx Comprimento 100mm
05.5.0222.113E	Broca Ø 1,3x Stop 14mmx Comprimento 105mm
05.5.0223.113E	Broca Ø 1,3x Stop 14mmx Comprimento 110mm
05.5.0224.113E	Broca Ø 1,3x Stop 14mmx Comprimento 115mm
05.5.0225.113E	Broca Ø 1,3x Stop 14mmx Comprimento 120mm
05.5.0226.113E	Broca Ø 1,3x Stop 14mmx Comprimento 125mm
05.5.0227.113E	Broca Ø 1,3x Stop 14mmx Comprimento 130mm
05.5.0228.113E	Broca Ø 1,3x Stop 14mmx Comprimento 135mm
05.5.0229.113E	Broca Ø 1,3x Stop 14mmx Comprimento 140mm
05.5.0230.113E	Broca Ø 1,3x Stop 14mmx Comprimento 145mm
05.5.0231.113E	Broca Ø 1,3x Stop 14mmx Comprimento 150mm
05.5.0232.113E	Broca Ø 1,3x Stop 15mmx Comprimento 50mm
05.5.0233.113E	Broca Ø 1,3x Stop 15mmx Comprimento 55mm
05.5.0234.113E	Broca Ø 1,3x Stop 15mmx Comprimento 60mm

05.5.0235.113E	Broca Ø 1,3x Stop 15mmx Comprimento 65mm
05.5.0236.113E	Broca Ø 1,3x Stop 15mmx Comprimento 70mm
05.5.0237.113E	Broca Ø 1,3x Stop 15mmx Comprimento 75mm
05.5.0238.113E	Broca Ø 1,3x Stop 15mmx Comprimento 80mm
05.5.0239.113E	Broca Ø 1,3x Stop 15mmx Comprimento 85mm
05.5.0240.113E	Broca Ø 1,3x Stop 15mmx Comprimento 90mm
05.5.0241.113E	Broca Ø 1,3x Stop 15mmx Comprimento 95mm
05.5.0242.113E	Broca Ø 1,3x Stop 15mmx Comprimento 100mm
05.5.0243.113E	Broca Ø 1,3x Stop 15mmx Comprimento 105mm
05.5.0244.113E	Broca Ø 1,3x Stop 15mmx Comprimento 110mm
05.5.0245.113E	Broca Ø 1,3x Stop 15mmx Comprimento 115mm
05.5.0246.113E	Broca Ø 1,3x Stop 15mmx Comprimento 120mm
05.5.0247.113E	Broca Ø 1,3x Stop 15mmx Comprimento 125mm
05.5.0248.113E	Broca Ø 1,3x Stop 15mmx Comprimento 130mm
05.5.0249.113E	Broca Ø 1,3x Stop 15mmx Comprimento 135mm
05.5.0250.113E	Broca Ø 1,3x Stop 15mmx Comprimento 140mm
05.5.0251.113E	Broca Ø 1,3x Stop 15mmx Comprimento 145mm
05.5.0252.113E	Broca Ø 1,3x Stop 15mmx Comprimento 150mm
05.5.0253.113E	Broca Ø 1,3x Stop 16mmx Comprimento 50mm
05.5.0254.113E	Broca Ø 1,3x Stop 16mmx Comprimento 55mm
05.5.0255.113E	Broca Ø 1,3x Stop 16mmx Comprimento 60mm
05.5.0256.113E	Broca Ø 1,3x Stop 16mmx Comprimento 65mm
05.5.0257.113E	Broca Ø 1,3x Stop 16mmx Comprimento 70mm
05.5.0258.113E	Broca Ø 1,3x Stop 16mmx Comprimento 75mm
05.5.0259.113E	Broca Ø 1,3x Stop 16mmx Comprimento 80mm
05.5.0260.113E	Broca Ø 1,3x Stop 16mmx Comprimento 85mm
05.5.0261.113E	Broca Ø 1,3x Stop 16mmx Comprimento 90mm
05.5.0262.113E	Broca Ø 1,3x Stop 16mmx Comprimento 95mm
05.5.0263.113E	Broca Ø 1,3x Stop 16mmx Comprimento 100mm
05.5.0264.113E	Broca Ø 1,3x Stop 16mmx Comprimento 105mm
05.5.0265.113E	Broca Ø 1,3x Stop 16mmx Comprimento 110mm
05.5.0266.113E	Broca Ø 1,3x Stop 16mmx Comprimento 115mm
05.5.0267.113E	Broca Ø 1,3x Stop 16mmx Comprimento 120mm
05.5.0268.113E	Broca Ø 1,3x Stop 16mmx Comprimento 125mm
05.5.0269.113E	Broca Ø 1,3x Stop 16mmx Comprimento 130mm
05.5.0270.113E	Broca Ø 1,3x Stop 16mmx Comprimento 135mm
05.5.0271.113E	Broca Ø 1,3x Stop 16mmx Comprimento 140mm
05.5.0272.113E	Broca Ø 1,3x Stop 16mmx Comprimento 145mm
05.5.0273.113E	Broca Ø 1,3x Stop 16mmx Comprimento 150mm
05.5.0274.113E	Broca Ø 1,3x Stop 17mmx Comprimento 50mm
05.5.0275.113E	Broca Ø 1,3x Stop 17mmx Comprimento 55mm
05.5.0276.113E	Broca Ø 1,3x Stop 17mmx Comprimento 60mm
05.5.0277.113E	Broca Ø 1,3x Stop 17mmx Comprimento 65mm
05.5.0278.113E	Broca Ø 1,3x Stop 17mmx Comprimento 70mm
05.5.0279.113E	Broca Ø 1,3x Stop 17mmx Comprimento 75mm
05.5.0280.113E	Broca Ø 1,3x Stop 17mmx Comprimento 80mm
05.5.0281.113E	Broca Ø 1,3x Stop 17mmx Comprimento 85mm
05.5.0282.113E	Broca Ø 1,3x Stop 17mmx Comprimento 90mm
05.5.0283.113E	Broca Ø 1,3x Stop 17mmx Comprimento 95mm

05.5.0284.113E	Broca Ø 1,3x Stop 17mmx Comprimento 100mm
05.5.0285.113E	Broca Ø 1,3x Stop 17mmx Comprimento 105mm
05.5.0286.113E	Broca Ø 1,3x Stop 17mmx Comprimento 110mm
05.5.0287.113E	Broca Ø 1,3x Stop 17mmx Comprimento 115mm
05.5.0288.113E	Broca Ø 1,3x Stop 17mmx Comprimento 120mm
05.5.0289.113E	Broca Ø 1,3x Stop 17mmx Comprimento 125mm
05.5.0290.113E	Broca Ø 1,3x Stop 17mmx Comprimento 130mm
05.5.0291.113E	Broca Ø 1,3x Stop 17mmx Comprimento 135mm
05.5.0292.113E	Broca Ø 1,3x Stop 17mmx Comprimento 140mm
05.5.0293.113E	Broca Ø 1,3x Stop 17mmx Comprimento 145mm
05.5.0294.113E	Broca Ø 1,3x Stop 17mmx Comprimento 150mm
05.5.0295.113E	Broca Ø 1,3x Stop 18mmx Comprimento 50mm
05.5.0296.113E	Broca Ø 1,3x Stop 18mmx Comprimento 55mm
05.5.0297.113E	Broca Ø 1,3x Stop 18mmx Comprimento 60mm
05.5.0298.113E	Broca Ø 1,3x Stop 18mmx Comprimento 65mm
05.5.0299.113E	Broca Ø 1,3x Stop 18mmx Comprimento 70mm
05.5.0300.113E	Broca Ø 1,3x Stop 18mmx Comprimento 75mm
05.5.0301.113E	Broca Ø 1,3x Stop 18mmx Comprimento 80mm
05.5.0302.113E	Broca Ø 1,3x Stop 18mmx Comprimento 85mm
05.5.0303.113E	Broca Ø 1,3x Stop 18mmx Comprimento 90mm
05.5.0304.113E	Broca Ø 1,3x Stop 18mmx Comprimento 95mm
05.5.0305.113E	Broca Ø 1,3x Stop 18mmx Comprimento 100mm
05.5.0306.113E	Broca Ø 1,3x Stop 18mmx Comprimento 105mm
05.5.0307.113E	Broca Ø 1,3x Stop 18mmx Comprimento 110mm
05.5.0308.113E	Broca Ø 1,3x Stop 18mmx Comprimento 115mm
05.5.0309.113E	Broca Ø 1,3x Stop 18mmx Comprimento 120mm
05.5.0310.113E	Broca Ø 1,3x Stop 18mmx Comprimento 125mm
05.5.0311.113E	Broca Ø 1,3x Stop 18mmx Comprimento 130mm
05.5.0312.113E	Broca Ø 1,3x Stop 18mmx Comprimento 135mm
05.5.0313.113E	Broca Ø 1,3x Stop 18mmx Comprimento 140mm
05.5.0314.113E	Broca Ø 1,3x Stop 18mmx Comprimento 145mm
05.5.0315.113E	Broca Ø 1,3x Stop 18mmx Comprimento 150mm
05.5.0316.113E	Broca Ø 1,3x Stop 19mmx Comprimento 50mm
05.5.0317.113E	Broca Ø 1,3x Stop 19mmx Comprimento 55mm
05.5.0318.113E	Broca Ø 1,3x Stop 19mmx Comprimento 60mm
05.5.0319.113E	Broca Ø 1,3x Stop 19mmx Comprimento 65mm
05.5.0320.113E	Broca Ø 1,3x Stop 19mmx Comprimento 70mm
05.5.0321.113E	Broca Ø 1,3x Stop 19mmx Comprimento 75mm
05.5.0322.113E	Broca Ø 1,3x Stop 19mmx Comprimento 80mm
05.5.0323.113E	Broca Ø 1,3x Stop 19mmx Comprimento 85mm
05.5.0324.113E	Broca Ø 1,3x Stop 19mmx Comprimento 90mm
05.5.0325.113E	Broca Ø 1,3x Stop 19mmx Comprimento 95mm
05.5.0326.113E	Broca Ø 1,3x Stop 19mmx Comprimento 100mm
05.5.0327.113E	Broca Ø 1,3x Stop 19mmx Comprimento 105mm
05.5.0328.113E	Broca Ø 1,3x Stop 19mmx Comprimento 110mm
05.5.0329.113E	Broca Ø 1,3x Stop 19mmx Comprimento 115mm
05.5.0330.113E	Broca Ø 1,3x Stop 19mmx Comprimento 120mm
05.5.0331.113E	Broca Ø 1,3x Stop 19mmx Comprimento 125mm
05.5.0332.113E	Broca Ø 1,3x Stop 19mmx Comprimento 130mm

05.5.0333.113E	Broca Ø 1,3x Stop 19mmx Comprimento 135mm
05.5.0334.113E	Broca Ø 1,3x Stop 19mmx Comprimento 140mm
05.5.0335.113E	Broca Ø 1,3x Stop 19mmx Comprimento 145mm
05.5.0336.113E	Broca Ø 1,3x Stop 19mmx Comprimento 150mm
05.5.0337.113E	Broca Ø 1,3x Stop 20mmx Comprimento 50mm
05.5.0338.113E	Broca Ø 1,3x Stop 20mmx Comprimento 55mm
05.5.0339.113E	Broca Ø 1,3x Stop 20mmx Comprimento 60mm
05.5.0340.113E	Broca Ø 1,3x Stop 20mmx Comprimento 65mm
05.5.0341.113E	Broca Ø 1,3x Stop 20mmx Comprimento 70mm
05.5.0342.113E	Broca Ø 1,3x Stop 20mmx Comprimento 75mm
05.5.0343.113E	Broca Ø 1,3x Stop 20mmx Comprimento 80mm
05.5.0344.113E	Broca Ø 1,3x Stop 20mmx Comprimento 85mm
05.5.0345.113E	Broca Ø 1,3x Stop 20mmx Comprimento 90mm
05.5.0346.113E	Broca Ø 1,3x Stop 20mmx Comprimento 95mm
05.5.0347.113E	Broca Ø 1,3x Stop 20mmx Comprimento 100mm
05.5.0348.113E	Broca Ø 1,3x Stop 20mmx Comprimento 105mm
05.5.0349.113E	Broca Ø 1,3x Stop 20mmx Comprimento 110mm
05.5.0350.113E	Broca Ø 1,3x Stop 20mmx Comprimento 115mm
05.5.0351.113E	Broca Ø 1,3x Stop 20mmx Comprimento 120mm
05.5.0352.113E	Broca Ø 1,3x Stop 20mmx Comprimento 125mm
05.5.0353.113E	Broca Ø 1,3x Stop 20mmx Comprimento 130mm
05.5.0354.113E	Broca Ø 1,3x Stop 20mmx Comprimento 135mm
05.5.0355.113E	Broca Ø 1,3x Stop 20mmx Comprimento 140mm
05.5.0356.113E	Broca Ø 1,3x Stop 20mmx Comprimento 145mm
05.5.0357.113E	Broca Ø 1,3x Stop 20mmx Comprimento 150mm
05.5.0358.113E	Broca Ø 1,3x Stop 21mmx Comprimento 50mm
05.5.0359.113E	Broca Ø 1,3x Stop 21mmx Comprimento 55mm
05.5.0360.113E	Broca Ø 1,3x Stop 21mmx Comprimento 60mm
05.5.0361.113E	Broca Ø 1,3x Stop 21mmx Comprimento 65mm
05.5.0362.113E	Broca Ø 1,3x Stop 21mmx Comprimento 70mm
05.5.0363.113E	Broca Ø 1,3x Stop 21mmx Comprimento 75mm
05.5.0364.113E	Broca Ø 1,3x Stop 21mmx Comprimento 80mm
05.5.0365.113E	Broca Ø 1,3x Stop 21mmx Comprimento 85mm
05.5.0366.113E	Broca Ø 1,3x Stop 21mmx Comprimento 90mm
05.5.0367.113E	Broca Ø 1,3x Stop 21mmx Comprimento 95mm
05.5.0368.113E	Broca Ø 1,3x Stop 21mmx Comprimento 100mm
05.5.0369.113E	Broca Ø 1,3x Stop 21mmx Comprimento 105mm
05.5.0370.113E	Broca Ø 1,3x Stop 21mmx Comprimento 110mm
05.5.0371.113E	Broca Ø 1,3x Stop 21mmx Comprimento 115mm
05.5.0372.113E	Broca Ø 1,3x Stop 21mmx Comprimento 120mm
05.5.0373.113E	Broca Ø 1,3x Stop 21mmx Comprimento 125mm
05.5.0374.113E	Broca Ø 1,3x Stop 21mmx Comprimento 130mm
05.5.0375.113E	Broca Ø 1,3x Stop 21mmx Comprimento 135mm
05.5.0376.113E	Broca Ø 1,3x Stop 21mmx Comprimento 140mm
05.5.0377.113E	Broca Ø 1,3x Stop 21mmx Comprimento 145mm
05.5.0378.113E	Broca Ø 1,3x Stop 21mmx Comprimento 150mm
05.5.0379.113E	Broca Ø 1,3x Stop 22mmx Comprimento 50mm
05.5.0380.113E	Broca Ø 1,3x Stop 22mmx Comprimento 55mm
05.5.0381.113E	Broca Ø 1,3x Stop 22mmx Comprimento 60mm

05.5.0382.113E	Broca Ø 1,3x Stop 22mmx Comprimento 65mm
05.5.0383.113E	Broca Ø 1,3x Stop 22mmx Comprimento 70mm
05.5.0384.113E	Broca Ø 1,3x Stop 22mmx Comprimento 75mm
05.5.0385.113E	Broca Ø 1,3x Stop 22mmx Comprimento 80mm
05.5.0386.113E	Broca Ø 1,3x Stop 22mmx Comprimento 85mm
05.5.0387.113E	Broca Ø 1,3x Stop 22mmx Comprimento 90mm
05.5.0388.113E	Broca Ø 1,3x Stop 22mmx Comprimento 95mm
05.5.0389.113E	Broca Ø 1,3x Stop 22mmx Comprimento 100mm
05.5.0390.113E	Broca Ø 1,3x Stop 22mmx Comprimento 105mm
05.5.0391.113E	Broca Ø 1,3x Stop 22mmx Comprimento 110mm
05.5.0392.113E	Broca Ø 1,3x Stop 22mmx Comprimento 115mm
05.5.0393.113E	Broca Ø 1,3x Stop 22mmx Comprimento 120mm
05.5.0394.113E	Broca Ø 1,3x Stop 22mmx Comprimento 125mm
05.5.0395.113E	Broca Ø 1,3x Stop 22mmx Comprimento 130mm
05.5.0396.113E	Broca Ø 1,3x Stop 22mmx Comprimento 135mm
05.5.0397.113E	Broca Ø 1,3x Stop 22mmx Comprimento 140mm
05.5.0398.113E	Broca Ø 1,3x Stop 22mmx Comprimento 145mm
05.5.0399.113E	Broca Ø 1,3x Stop 22mmx Comprimento 150mm
05.5.0400.113E	Broca Ø 1,3x Stop 23mmx Comprimento 50mm
05.5.0401.113E	Broca Ø 1,3x Stop 23mmx Comprimento 55mm
05.5.0402.113E	Broca Ø 1,3x Stop 23mmx Comprimento 60mm
05.5.0403.113E	Broca Ø 1,3x Stop 23mmx Comprimento 65mm
05.5.0404.113E	Broca Ø 1,3x Stop 23mmx Comprimento 70mm
05.5.0405.113E	Broca Ø 1,3x Stop 23mmx Comprimento 75mm
05.5.0406.113E	Broca Ø 1,3x Stop 23mmx Comprimento 80mm
05.5.0407.113E	Broca Ø 1,3x Stop 23mmx Comprimento 85mm
05.5.0408.113E	Broca Ø 1,3x Stop 23mmx Comprimento 90mm
05.5.0409.113E	Broca Ø 1,3x Stop 23mmx Comprimento 95mm
05.5.0410.113E	Broca Ø 1,3x Stop 23mmx Comprimento 100mm
05.5.0411.113E	Broca Ø 1,3x Stop 23mmx Comprimento 105mm
05.5.0412.113E	Broca Ø 1,3x Stop 23mmx Comprimento 110mm
05.5.0413.113E	Broca Ø 1,3x Stop 23mmx Comprimento 115mm
05.5.0414.113E	Broca Ø 1,3x Stop 23mmx Comprimento 120mm
05.5.0415.113E	Broca Ø 1,3x Stop 23mmx Comprimento 125mm
05.5.0416.113E	Broca Ø 1,3x Stop 23mmx Comprimento 130mm
05.5.0417.113E	Broca Ø 1,3x Stop 23mmx Comprimento 135mm
05.5.0418.113E	Broca Ø 1,3x Stop 23mmx Comprimento 140mm
05.5.0419.113E	Broca Ø 1,3x Stop 23mmx Comprimento 145mm
05.5.0420.113E	Broca Ø 1,3x Stop 23mmx Comprimento 150mm
05.5.0421.113E	Broca Ø 1,3x Stop 24mmx Comprimento 50mm
05.5.0422.113E	Broca Ø 1,3x Stop 24mmx Comprimento 55mm
05.5.0423.113E	Broca Ø 1,3x Stop 24mmx Comprimento 60mm
05.5.0424.113E	Broca Ø 1,3x Stop 24mmx Comprimento 65mm
05.5.0425.113E	Broca Ø 1,3x Stop 24mmx Comprimento 70mm
05.5.0426.113E	Broca Ø 1,3x Stop 24mmx Comprimento 75mm
05.5.0427.113E	Broca Ø 1,3x Stop 24mmx Comprimento 80mm
05.5.0428.113E	Broca Ø 1,3x Stop 24mmx Comprimento 85mm
05.5.0429.113E	Broca Ø 1,3x Stop 24mmx Comprimento 90mm
05.5.0430.113E	Broca Ø 1,3x Stop 24mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0431.113E	Broca Ø 1,3x Stop 24mmx Comprimento 100mm
05.5.0432.113E	Broca Ø 1,3x Stop 24mmx Comprimento 105mm
05.5.0433.113E	Broca Ø 1,3x Stop 24mmx Comprimento 110mm
05.5.0434.113E	Broca Ø 1,3x Stop 24mmx Comprimento 115mm
05.5.0435.113E	Broca Ø 1,3x Stop 24mmx Comprimento 120mm
05.5.0436.113E	Broca Ø 1,3x Stop 24mmx Comprimento 125mm
05.5.0437.113E	Broca Ø 1,3x Stop 24mmx Comprimento 130mm
05.5.0438.113E	Broca Ø 1,3x Stop 24mmx Comprimento 135mm
05.5.0439.113E	Broca Ø 1,3x Stop 24mmx Comprimento 140mm
05.5.0440.113E	Broca Ø 1,3x Stop 24mmx Comprimento 145mm
05.5.0441.113E	Broca Ø 1,3x Stop 24mmx Comprimento 150mm
05.5.0442.113E	Broca Ø 1,3x Stop 25mmx Comprimento 50mm
05.5.0443.113E	Broca Ø 1,3x Stop 25mmx Comprimento 55mm
05.5.0444.113E	Broca Ø 1,3x Stop 25mmx Comprimento 60mm
05.5.0445.113E	Broca Ø 1,3x Stop 25mmx Comprimento 65mm
05.5.0446.113E	Broca Ø 1,3x Stop 25mmx Comprimento 70mm
05.5.0447.113E	Broca Ø 1,3x Stop 25mmx Comprimento 75mm
05.5.0448.113E	Broca Ø 1,3x Stop 25mmx Comprimento 80mm
05.5.0449.113E	Broca Ø 1,3x Stop 25mmx Comprimento 85mm
05.5.0450.113E	Broca Ø 1,3x Stop 25mmx Comprimento 90mm
05.5.0451.113E	Broca Ø 1,3x Stop 25mmx Comprimento 95mm
05.5.0452.113E	Broca Ø 1,3x Stop 25mmx Comprimento 100mm
05.5.0453.113E	Broca Ø 1,3x Stop 25mmx Comprimento 105mm
05.5.0454.113E	Broca Ø 1,3x Stop 25mmx Comprimento 110mm
05.5.0455.113E	Broca Ø 1,3x Stop 25mmx Comprimento 115mm
05.5.0456.113E	Broca Ø 1,3x Stop 25mmx Comprimento 120mm
05.5.0457.113E	Broca Ø 1,3x Stop 25mmx Comprimento 125mm
05.5.0458.113E	Broca Ø 1,3x Stop 25mmx Comprimento 130mm
05.5.0459.113E	Broca Ø 1,3x Stop 25mmx Comprimento 135mm
05.5.0460.113E	Broca Ø 1,3x Stop 25mmx Comprimento 140mm
05.5.0461.113E	Broca Ø 1,3x Stop 25mmx Comprimento 145mm
05.5.0462.113E	Broca Ø 1,3x Stop 25mmx Comprimento 150mm
05.5.0463.113E	Broca Ø 1,3x Stop 30mmx Comprimento 50mm
05.5.0464.113E	Broca Ø 1,3x Stop 30mmx Comprimento 55mm
05.5.0465.113E	Broca Ø 1,3x Stop 30mmx Comprimento 60mm
05.5.0466.113E	Broca Ø 1,3x Stop 30mmx Comprimento 65mm
05.5.0467.113E	Broca Ø 1,3x Stop 30mmx Comprimento 70mm
05.5.0468.113E	Broca Ø 1,3x Stop 30mmx Comprimento 75mm
05.5.0469.113E	Broca Ø 1,3x Stop 30mmx Comprimento 80mm
05.5.0470.113E	Broca Ø 1,3x Stop 30mmx Comprimento 85mm
05.5.0471.113E	Broca Ø 1,3x Stop 30mmx Comprimento 90mm
05.5.0472.113E	Broca Ø 1,3x Stop 30mmx Comprimento 95mm
05.5.0473.113E	Broca Ø 1,3x Stop 30mmx Comprimento 100mm
05.5.0474.113E	Broca Ø 1,3x Stop 30mmx Comprimento 105mm
05.5.0475.113E	Broca Ø 1,3x Stop 30mmx Comprimento 110mm
05.5.0476.113E	Broca Ø 1,3x Stop 30mmx Comprimento 115mm
05.5.0477.113E	Broca Ø 1,3x Stop 30mmx Comprimento 120mm
05.5.0478.113E	Broca Ø 1,3x Stop 30mmx Comprimento 125mm
05.5.0479.113E	Broca Ø 1,3x Stop 30mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0480.113E	Broca Ø 1,3x Stop 30mmx Comprimento 135mm
05.5.0481.113E	Broca Ø 1,3x Stop 30mmx Comprimento 140mm
05.5.0482.113E	Broca Ø 1,3x Stop 30mmx Comprimento 145mm
05.5.0483.113E	Broca Ø 1,3x Stop 30mmx Comprimento 150mm
05.5.0484.113E	Broca Ø 1,3x Stop 35mmx Comprimento 50mm
05.5.0485.113E	Broca Ø 1,3x Stop 35mmx Comprimento 55mm
05.5.0486.113E	Broca Ø 1,3x Stop 35mmx Comprimento 60mm
05.5.0487.113E	Broca Ø 1,3x Stop 35mmx Comprimento 65mm
05.5.0488.113E	Broca Ø 1,3x Stop 35mmx Comprimento 70mm
05.5.0489.113E	Broca Ø 1,3x Stop 35mmx Comprimento 75mm
05.5.0490.113E	Broca Ø 1,3x Stop 35mmx Comprimento 80mm
05.5.0491.113E	Broca Ø 1,3x Stop 35mmx Comprimento 85mm
05.5.0492.113E	Broca Ø 1,3x Stop 35mmx Comprimento 90mm
05.5.0493.113E	Broca Ø 1,3x Stop 35mmx Comprimento 95mm
05.5.0494.113E	Broca Ø 1,3x Stop 35mmx Comprimento 100mm
05.5.0495.113E	Broca Ø 1,3x Stop 35mmx Comprimento 105mm
05.5.0496.113E	Broca Ø 1,3x Stop 35mmx Comprimento 110mm
05.5.0497.113E	Broca Ø 1,3x Stop 35mmx Comprimento 115mm
05.5.0498.113E	Broca Ø 1,3x Stop 35mmx Comprimento 120mm
05.5.0499.113E	Broca Ø 1,3x Stop 35mmx Comprimento 125mm
05.5.0500.113E	Broca Ø 1,3x Stop 35mmx Comprimento 130mm
05.5.0501.113E	Broca Ø 1,3x Stop 35mmx Comprimento 135mm
05.5.0502.113E	Broca Ø 1,3x Stop 35mmx Comprimento 140mm
05.5.0503.113E	Broca Ø 1,3x Stop 35mmx Comprimento 145mm
05.5.0504.113E	Broca Ø 1,3x Stop 35mmx Comprimento 150mm
05.5.0001.114E	Broca Ø 1,4x Stop 4mmx Comprimento 50mm
05.5.0002.114E	Broca Ø 1,4x Stop 4mmx Comprimento 55mm
05.5.0003.114E	Broca Ø 1,4x Stop 4mmx Comprimento 60mm
05.5.0004.114E	Broca Ø 1,4x Stop 4mmx Comprimento 65mm
05.5.0005.114E	Broca Ø 1,4x Stop 4mmx Comprimento 70mm
05.5.0006.114E	Broca Ø 1,4x Stop 4mmx Comprimento 75mm
05.5.0007.114E	Broca Ø 1,4x Stop 4mmx Comprimento 80mm
05.5.0008.114E	Broca Ø 1,4x Stop 4mmx Comprimento 85mm
05.5.0009.114E	Broca Ø 1,4x Stop 4mmx Comprimento 90mm
05.5.0010.114E	Broca Ø 1,4x Stop 4mmx Comprimento 95mm
05.5.0011.114E	Broca Ø 1,4x Stop 4mmx Comprimento 100mm
05.5.0012.114E	Broca Ø 1,4x Stop 4mmx Comprimento 105mm
05.5.0013.114E	Broca Ø 1,4x Stop 4mmx Comprimento 110mm
05.5.0014.114E	Broca Ø 1,4x Stop 4mmx Comprimento 115mm
05.5.0015.114E	Broca Ø 1,4x Stop 4mmx Comprimento 120mm
05.5.0016.114E	Broca Ø 1,4x Stop 4mmx Comprimento 125mm
05.5.0017.114E	Broca Ø 1,4x Stop 4mmx Comprimento 130mm
05.5.0018.114E	Broca Ø 1,4x Stop 4mmx Comprimento 135mm
05.5.0019.114E	Broca Ø 1,4x Stop 4mmx Comprimento 140mm
05.5.0020.114E	Broca Ø 1,4x Stop 4mmx Comprimento 145mm
05.5.0021.114E	Broca Ø 1,4x Stop 4mmx Comprimento 150mm
05.5.0022.114E	Broca Ø 1,4x Stop 5mmx Comprimento 50mm
05.5.0023.114E	Broca Ø 1,4x Stop 5mmx Comprimento 55mm
05.5.0024.114E	Broca Ø 1,4x Stop 5mmx Comprimento 60mm

05.5.0025.114E	Broca Ø 1,4x Stop 5mmx Comprimento 65mm
05.5.0026.114E	Broca Ø 1,4x Stop 5mmx Comprimento 70mm
05.5.0027.114E	Broca Ø 1,4x Stop 5mmx Comprimento 75mm
05.5.0028.114E	Broca Ø 1,4x Stop 5mmx Comprimento 80mm
05.5.0029.114E	Broca Ø 1,4x Stop 5mmx Comprimento 85mm
05.5.0030.114E	Broca Ø 1,4x Stop 5mmx Comprimento 90mm
05.5.0031.114E	Broca Ø 1,4x Stop 5mmx Comprimento 95mm
05.5.0032.114E	Broca Ø 1,4x Stop 5mmx Comprimento 100mm
05.5.0033.114E	Broca Ø 1,4x Stop 5mmx Comprimento 105mm
05.5.0034.114E	Broca Ø 1,4x Stop 5mmx Comprimento 110mm
05.5.0035.114E	Broca Ø 1,4x Stop 5mmx Comprimento 115mm
05.5.0036.114E	Broca Ø 1,4x Stop 5mmx Comprimento 120mm
05.5.0037.114E	Broca Ø 1,4x Stop 5mmx Comprimento 125mm
05.5.0038.114E	Broca Ø 1,4x Stop 5mmx Comprimento 130mm
05.5.0039.114E	Broca Ø 1,4x Stop 5mmx Comprimento 135mm
05.5.0040.114E	Broca Ø 1,4x Stop 5mmx Comprimento 140mm
05.5.0041.114E	Broca Ø 1,4x Stop 5mmx Comprimento 145mm
05.5.0042.114E	Broca Ø 1,4x Stop 5mmx Comprimento 150mm
05.5.0043.114E	Broca Ø 1,4x Stop 6mmx Comprimento 50mm
05.5.0044.114E	Broca Ø 1,4x Stop 6mmx Comprimento 55mm
05.5.0045.114E	Broca Ø 1,4x Stop 6mmx Comprimento 60mm
05.5.0046.114E	Broca Ø 1,4x Stop 6mmx Comprimento 65mm
05.5.0047.114E	Broca Ø 1,4x Stop 6mmx Comprimento 70mm
05.5.0048.114E	Broca Ø 1,4x Stop 6mmx Comprimento 75mm
05.5.0049.114E	Broca Ø 1,4x Stop 6mmx Comprimento 80mm
05.5.0050.114E	Broca Ø 1,4x Stop 6mmx Comprimento 85mm
05.5.0051.114E	Broca Ø 1,4x Stop 6mmx Comprimento 90mm
05.5.0052.114E	Broca Ø 1,4x Stop 6mmx Comprimento 95mm
05.5.0053.114E	Broca Ø 1,4x Stop 6mmx Comprimento 100mm
05.5.0054.114E	Broca Ø 1,4x Stop 6mmx Comprimento 105mm
05.5.0055.114E	Broca Ø 1,4x Stop 6mmx Comprimento 110mm
05.5.0056.114E	Broca Ø 1,4x Stop 6mmx Comprimento 115mm
05.5.0057.114E	Broca Ø 1,4x Stop 6mmx Comprimento 120mm
05.5.0058.114E	Broca Ø 1,4x Stop 6mmx Comprimento 125mm
05.5.0059.114E	Broca Ø 1,4x Stop 6mmx Comprimento 130mm
05.5.0060.114E	Broca Ø 1,4x Stop 6mmx Comprimento 135mm
05.5.0061.114E	Broca Ø 1,4x Stop 6mmx Comprimento 140mm
05.5.0062.114E	Broca Ø 1,4x Stop 6mmx Comprimento 145mm
05.5.0063.114E	Broca Ø 1,4x Stop 6mmx Comprimento 150mm
05.5.0064.114E	Broca Ø 1,4x Stop 7mmx Comprimento 50mm
05.5.0065.114E	Broca Ø 1,4x Stop 7mmx Comprimento 55mm
05.5.0066.114E	Broca Ø 1,4x Stop 7mmx Comprimento 60mm
05.5.0067.114E	Broca Ø 1,4x Stop 7mmx Comprimento 65mm
05.5.0068.114E	Broca Ø 1,4x Stop 7mmx Comprimento 70mm
05.5.0069.114E	Broca Ø 1,4x Stop 7mmx Comprimento 75mm
05.5.0070.114E	Broca Ø 1,4x Stop 7mmx Comprimento 80mm
05.5.0071.114E	Broca Ø 1,4x Stop 7mmx Comprimento 85mm
05.5.0072.114E	Broca Ø 1,4x Stop 7mmx Comprimento 90mm
05.5.0073.114E	Broca Ø 1,4x Stop 7mmx Comprimento 95mm

05.5.0074.114E	Broca Ø 1,4x Stop 7mmx Comprimento 100mm
05.5.0075.114E	Broca Ø 1,4x Stop 7mmx Comprimento 105mm
05.5.0076.114E	Broca Ø 1,4x Stop 7mmx Comprimento 110mm
05.5.0077.114E	Broca Ø 1,4x Stop 7mmx Comprimento 115mm
05.5.0078.114E	Broca Ø 1,4x Stop 7mmx Comprimento 120mm
05.5.0079.114E	Broca Ø 1,4x Stop 7mmx Comprimento 125mm
05.5.0080.114E	Broca Ø 1,4x Stop 7mmx Comprimento 130mm
05.5.0081.114E	Broca Ø 1,4x Stop 7mmx Comprimento 135mm
05.5.0082.114E	Broca Ø 1,4x Stop 7mmx Comprimento 140mm
05.5.0083.114E	Broca Ø 1,4x Stop 7mmx Comprimento 145mm
05.5.0084.114E	Broca Ø 1,4x Stop 7mmx Comprimento 150mm
05.5.0085.114E	Broca Ø 1,4x Stop 8mmx Comprimento 50mm
05.5.0086.114E	Broca Ø 1,4x Stop 8mmx Comprimento 55mm
05.5.0087.114E	Broca Ø 1,4x Stop 8mmx Comprimento 60mm
05.5.0088.114E	Broca Ø 1,4x Stop 8mmx Comprimento 65mm
05.5.0089.114E	Broca Ø 1,4x Stop 8mmx Comprimento 70mm
05.5.0090.114E	Broca Ø 1,4x Stop 8mmx Comprimento 75mm
05.5.0091.114E	Broca Ø 1,4x Stop 8mmx Comprimento 80mm
05.5.0092.114E	Broca Ø 1,4x Stop 8mmx Comprimento 85mm
05.5.0093.114E	Broca Ø 1,4x Stop 8mmx Comprimento 90mm
05.5.0094.114E	Broca Ø 1,4x Stop 8mmx Comprimento 95mm
05.5.0095.114E	Broca Ø 1,4x Stop 8mmx Comprimento 100mm
05.5.0096.114E	Broca Ø 1,4x Stop 8mmx Comprimento 105mm
05.5.0097.114E	Broca Ø 1,4x Stop 8mmx Comprimento 110mm
05.5.0098.114E	Broca Ø 1,4x Stop 8mmx Comprimento 115mm
05.5.0099.114E	Broca Ø 1,4x Stop 8mmx Comprimento 120mm
05.5.0100.114E	Broca Ø 1,4x Stop 8mmx Comprimento 125mm
05.5.0101.114E	Broca Ø 1,4x Stop 8mmx Comprimento 130mm
05.5.0102.114E	Broca Ø 1,4x Stop 8mmx Comprimento 135mm
05.5.0103.114E	Broca Ø 1,4x Stop 8mmx Comprimento 140mm
05.5.0104.114E	Broca Ø 1,4x Stop 8mmx Comprimento 145mm
05.5.0105.114E	Broca Ø 1,4x Stop 8mmx Comprimento 150mm
05.5.0106.114E	Broca Ø 1,4x Stop 9mmx Comprimento 50mm
05.5.0107.114E	Broca Ø 1,4x Stop 9mmx Comprimento 55mm
05.5.0108.114E	Broca Ø 1,4x Stop 9mmx Comprimento 60mm
05.5.0109.114E	Broca Ø 1,4x Stop 9mmx Comprimento 65mm
05.5.0110.114E	Broca Ø 1,4x Stop 9mmx Comprimento 70mm
05.5.0111.114E	Broca Ø 1,4x Stop 9mmx Comprimento 75mm
05.5.0112.114E	Broca Ø 1,4x Stop 9mmx Comprimento 80mm
05.5.0113.114E	Broca Ø 1,4x Stop 9mmx Comprimento 85mm
05.5.0114.114E	Broca Ø 1,4x Stop 9mmx Comprimento 90mm
05.5.0115.114E	Broca Ø 1,4x Stop 9mmx Comprimento 95mm
05.5.0116.114E	Broca Ø 1,4x Stop 9mmx Comprimento 100mm
05.5.0117.114E	Broca Ø 1,4x Stop 9mmx Comprimento 105mm
05.5.0118.114E	Broca Ø 1,4x Stop 9mmx Comprimento 110mm
05.5.0119.114E	Broca Ø 1,4x Stop 9mmx Comprimento 115mm
05.5.0120.114E	Broca Ø 1,4x Stop 9mmx Comprimento 120mm
05.5.0121.114E	Broca Ø 1,4x Stop 9mmx Comprimento 125mm
05.5.0122.114E	Broca Ø 1,4x Stop 9mmx Comprimento 130mm

05.5.0123.114E	Broca Ø 1,4x Stop 9mmx Comprimento 135mm
05.5.0124.114E	Broca Ø 1,4x Stop 9mmx Comprimento 140mm
05.5.0125.114E	Broca Ø 1,4x Stop 9mmx Comprimento 145mm
05.5.0126.114E	Broca Ø 1,4x Stop 9mmx Comprimento 150mm
05.5.0127.114E	Broca Ø 1,4x Stop 10mmx Comprimento 50mm
05.5.0128.114E	Broca Ø 1,4x Stop 10mmx Comprimento 55mm
05.5.0129.114E	Broca Ø 1,4x Stop 10mmx Comprimento 60mm
05.5.0130.114E	Broca Ø 1,4x Stop 10mmx Comprimento 65mm
05.5.0131.114E	Broca Ø 1,4x Stop 10mmx Comprimento 70mm
05.5.0132.114E	Broca Ø 1,4x Stop 10mmx Comprimento 75mm
05.5.0133.114E	Broca Ø 1,4x Stop 10mmx Comprimento 80mm
05.5.0134.114E	Broca Ø 1,4x Stop 10mmx Comprimento 85mm
05.5.0135.114E	Broca Ø 1,4x Stop 10mmx Comprimento 90mm
05.5.0136.114E	Broca Ø 1,4x Stop 10mmx Comprimento 95mm
05.5.0137.114E	Broca Ø 1,4x Stop 10mmx Comprimento 100mm
05.5.0138.114E	Broca Ø 1,4x Stop 10mmx Comprimento 105mm
05.5.0139.114E	Broca Ø 1,4x Stop 10mmx Comprimento 110mm
05.5.0140.114E	Broca Ø 1,4x Stop 10mmx Comprimento 115mm
05.5.0141.114E	Broca Ø 1,4x Stop 10mmx Comprimento 120mm
05.5.0142.114E	Broca Ø 1,4x Stop 10mmx Comprimento 125mm
05.5.0143.114E	Broca Ø 1,4x Stop 10mmx Comprimento 130mm
05.5.0144.114E	Broca Ø 1,4x Stop 10mmx Comprimento 135mm
05.5.0145.114E	Broca Ø 1,4x Stop 10mmx Comprimento 140mm
05.5.0146.114E	Broca Ø 1,4x Stop 10mmx Comprimento 145mm
05.5.0147.114E	Broca Ø 1,4x Stop 10mmx Comprimento 150mm
05.5.0148.114E	Broca Ø 1,4x Stop 11mmx Comprimento 50mm
05.5.0149.114E	Broca Ø 1,4x Stop 11mmx Comprimento 55mm
05.5.0150.114E	Broca Ø 1,4x Stop 11mmx Comprimento 60mm
05.5.0151.114E	Broca Ø 1,4x Stop 11mmx Comprimento 65mm
05.5.0152.114E	Broca Ø 1,4x Stop 11mmx Comprimento 70mm
05.5.0153.114E	Broca Ø 1,4x Stop 11mmx Comprimento 75mm
05.5.0154.114E	Broca Ø 1,4x Stop 11mmx Comprimento 80mm
05.5.0155.114E	Broca Ø 1,4x Stop 11mmx Comprimento 85mm
05.5.0156.114E	Broca Ø 1,4x Stop 11mmx Comprimento 90mm
05.5.0157.114E	Broca Ø 1,4x Stop 11mmx Comprimento 95mm
05.5.0158.114E	Broca Ø 1,4x Stop 11mmx Comprimento 100mm
05.5.0159.114E	Broca Ø 1,4x Stop 11mmx Comprimento 105mm
05.5.0160.114E	Broca Ø 1,4x Stop 11mmx Comprimento 110mm
05.5.0161.114E	Broca Ø 1,4x Stop 11mmx Comprimento 115mm
05.5.0162.114E	Broca Ø 1,4x Stop 11mmx Comprimento 120mm
05.5.0163.114E	Broca Ø 1,4x Stop 11mmx Comprimento 125mm
05.5.0164.114E	Broca Ø 1,4x Stop 11mmx Comprimento 130mm
05.5.0165.114E	Broca Ø 1,4x Stop 11mmx Comprimento 135mm
05.5.0166.114E	Broca Ø 1,4x Stop 11mmx Comprimento 140mm
05.5.0167.114E	Broca Ø 1,4x Stop 11mmx Comprimento 145mm
05.5.0168.114E	Broca Ø 1,4x Stop 11mmx Comprimento 150mm
05.5.0169.114E	Broca Ø 1,4x Stop 12mmx Comprimento 50mm
05.5.0170.114E	Broca Ø 1,4x Stop 12mmx Comprimento 55mm
05.5.0171.114E	Broca Ø 1,4x Stop 12mmx Comprimento 60mm

05.5.0172.114E	Broca Ø 1,4x Stop 12mmx Comprimento 65mm
05.5.0173.114E	Broca Ø 1,4x Stop 12mmx Comprimento 70mm
05.5.0174.114E	Broca Ø 1,4x Stop 12mmx Comprimento 75mm
05.5.0175.114E	Broca Ø 1,4x Stop 12mmx Comprimento 80mm
05.5.0176.114E	Broca Ø 1,4x Stop 12mmx Comprimento 85mm
05.5.0177.114E	Broca Ø 1,4x Stop 12mmx Comprimento 90mm
05.5.0178.114E	Broca Ø 1,4x Stop 12mmx Comprimento 95mm
05.5.0179.114E	Broca Ø 1,4x Stop 12mmx Comprimento 100mm
05.5.0180.114E	Broca Ø 1,4x Stop 12mmx Comprimento 105mm
05.5.0181.114E	Broca Ø 1,4x Stop 12mmx Comprimento 110mm
05.5.0182.114E	Broca Ø 1,4x Stop 12mmx Comprimento 115mm
05.5.0183.114E	Broca Ø 1,4x Stop 12mmx Comprimento 120mm
05.5.0184.114E	Broca Ø 1,4x Stop 12mmx Comprimento 125mm
05.5.0185.114E	Broca Ø 1,4x Stop 12mmx Comprimento 130mm
05.5.0186.114E	Broca Ø 1,4x Stop 12mmx Comprimento 135mm
05.5.0187.114E	Broca Ø 1,4x Stop 12mmx Comprimento 140mm
05.5.0188.114E	Broca Ø 1,4x Stop 12mmx Comprimento 145mm
05.5.0189.114E	Broca Ø 1,4x Stop 12mmx Comprimento 150mm
05.5.0190.114E	Broca Ø 1,4x Stop 13mmx Comprimento 50mm
05.5.0191.114E	Broca Ø 1,4x Stop 13mmx Comprimento 55mm
05.5.0192.114E	Broca Ø 1,4x Stop 13mmx Comprimento 60mm
05.5.0193.114E	Broca Ø 1,4x Stop 13mmx Comprimento 65mm
05.5.0194.114E	Broca Ø 1,4x Stop 13mmx Comprimento 70mm
05.5.0195.114E	Broca Ø 1,4x Stop 13mmx Comprimento 75mm
05.5.0196.114E	Broca Ø 1,4x Stop 13mmx Comprimento 80mm
05.5.0197.114E	Broca Ø 1,4x Stop 13mmx Comprimento 85mm
05.5.0198.114E	Broca Ø 1,4x Stop 13mmx Comprimento 90mm
05.5.0199.114E	Broca Ø 1,4x Stop 13mmx Comprimento 95mm
05.5.0200.114E	Broca Ø 1,4x Stop 13mmx Comprimento 100mm
05.5.0201.114E	Broca Ø 1,4x Stop 13mmx Comprimento 105mm
05.5.0202.114E	Broca Ø 1,4x Stop 13mmx Comprimento 110mm
05.5.0203.114E	Broca Ø 1,4x Stop 13mmx Comprimento 115mm
05.5.0204.114E	Broca Ø 1,4x Stop 13mmx Comprimento 120mm
05.5.0205.114E	Broca Ø 1,4x Stop 13mmx Comprimento 125mm
05.5.0206.114E	Broca Ø 1,4x Stop 13mmx Comprimento 130mm
05.5.0207.114E	Broca Ø 1,4x Stop 13mmx Comprimento 135mm
05.5.0208.114E	Broca Ø 1,4x Stop 13mmx Comprimento 140mm
05.5.0209.114E	Broca Ø 1,4x Stop 13mmx Comprimento 145mm
05.5.0210.114E	Broca Ø 1,4x Stop 13mmx Comprimento 150mm
05.5.0211.114E	Broca Ø 1,4x Stop 14mmx Comprimento 50mm
05.5.0212.114E	Broca Ø 1,4x Stop 14mmx Comprimento 55mm
05.5.0213.114E	Broca Ø 1,4x Stop 14mmx Comprimento 60mm
05.5.0214.114E	Broca Ø 1,4x Stop 14mmx Comprimento 65mm
05.5.0215.114E	Broca Ø 1,4x Stop 14mmx Comprimento 70mm
05.5.0216.114E	Broca Ø 1,4x Stop 14mmx Comprimento 75mm
05.5.0217.114E	Broca Ø 1,4x Stop 14mmx Comprimento 80mm
05.5.0218.114E	Broca Ø 1,4x Stop 14mmx Comprimento 85mm
05.5.0219.114E	Broca Ø 1,4x Stop 14mmx Comprimento 90mm
05.5.0220.114E	Broca Ø 1,4x Stop 14mmx Comprimento 95mm

05.5.0221.114E	Broca Ø 1,4x Stop 14mmx Comprimento 100mm
05.5.0222.114E	Broca Ø 1,4x Stop 14mmx Comprimento 105mm
05.5.0223.114E	Broca Ø 1,4x Stop 14mmx Comprimento 110mm
05.5.0224.114E	Broca Ø 1,4x Stop 14mmx Comprimento 115mm
05.5.0225.114E	Broca Ø 1,4x Stop 14mmx Comprimento 120mm
05.5.0226.114E	Broca Ø 1,4x Stop 14mmx Comprimento 125mm
05.5.0227.114E	Broca Ø 1,4x Stop 14mmx Comprimento 130mm
05.5.0228.114E	Broca Ø 1,4x Stop 14mmx Comprimento 135mm
05.5.0229.114E	Broca Ø 1,4x Stop 14mmx Comprimento 140mm
05.5.0230.114E	Broca Ø 1,4x Stop 14mmx Comprimento 145mm
05.5.0231.114E	Broca Ø 1,4x Stop 14mmx Comprimento 150mm
05.5.0232.114E	Broca Ø 1,4x Stop 15mmx Comprimento 50mm
05.5.0233.114E	Broca Ø 1,4x Stop 15mmx Comprimento 55mm
05.5.0234.114E	Broca Ø 1,4x Stop 15mmx Comprimento 60mm
05.5.0235.114E	Broca Ø 1,4x Stop 15mmx Comprimento 65mm
05.5.0236.114E	Broca Ø 1,4x Stop 15mmx Comprimento 70mm
05.5.0237.114E	Broca Ø 1,4x Stop 15mmx Comprimento 75mm
05.5.0238.114E	Broca Ø 1,4x Stop 15mmx Comprimento 80mm
05.5.0239.114E	Broca Ø 1,4x Stop 15mmx Comprimento 85mm
05.5.0240.114E	Broca Ø 1,4x Stop 15mmx Comprimento 90mm
05.5.0241.114E	Broca Ø 1,4x Stop 15mmx Comprimento 95mm
05.5.0242.114E	Broca Ø 1,4x Stop 15mmx Comprimento 100mm
05.5.0243.114E	Broca Ø 1,4x Stop 15mmx Comprimento 105mm
05.5.0244.114E	Broca Ø 1,4x Stop 15mmx Comprimento 110mm
05.5.0245.114E	Broca Ø 1,4x Stop 15mmx Comprimento 115mm
05.5.0246.114E	Broca Ø 1,4x Stop 15mmx Comprimento 120mm
05.5.0247.114E	Broca Ø 1,4x Stop 15mmx Comprimento 125mm
05.5.0248.114E	Broca Ø 1,4x Stop 15mmx Comprimento 130mm
05.5.0249.114E	Broca Ø 1,4x Stop 15mmx Comprimento 135mm
05.5.0250.114E	Broca Ø 1,4x Stop 15mmx Comprimento 140mm
05.5.0251.114E	Broca Ø 1,4x Stop 15mmx Comprimento 145mm
05.5.0252.114E	Broca Ø 1,4x Stop 15mmx Comprimento 150mm
05.5.0253.114E	Broca Ø 1,4x Stop 16mmx Comprimento 50mm
05.5.0254.114E	Broca Ø 1,4x Stop 16mmx Comprimento 55mm
05.5.0255.114E	Broca Ø 1,4x Stop 16mmx Comprimento 60mm
05.5.0256.114E	Broca Ø 1,4x Stop 16mmx Comprimento 65mm
05.5.0257.114E	Broca Ø 1,4x Stop 16mmx Comprimento 70mm
05.5.0258.114E	Broca Ø 1,4x Stop 16mmx Comprimento 75mm
05.5.0259.114E	Broca Ø 1,4x Stop 16mmx Comprimento 80mm
05.5.0260.114E	Broca Ø 1,4x Stop 16mmx Comprimento 85mm
05.5.0261.114E	Broca Ø 1,4x Stop 16mmx Comprimento 90mm
05.5.0262.114E	Broca Ø 1,4x Stop 16mmx Comprimento 95mm
05.5.0263.114E	Broca Ø 1,4x Stop 16mmx Comprimento 100mm
05.5.0264.114E	Broca Ø 1,4x Stop 16mmx Comprimento 105mm
05.5.0265.114E	Broca Ø 1,4x Stop 16mmx Comprimento 110mm
05.5.0266.114E	Broca Ø 1,4x Stop 16mmx Comprimento 115mm
05.5.0267.114E	Broca Ø 1,4x Stop 16mmx Comprimento 120mm
05.5.0268.114E	Broca Ø 1,4x Stop 16mmx Comprimento 125mm
05.5.0269.114E	Broca Ø 1,4x Stop 16mmx Comprimento 130mm

05.5.0270.114E	Broca Ø 1,4x Stop 16mmx Comprimento 135mm
05.5.0271.114E	Broca Ø 1,4x Stop 16mmx Comprimento 140mm
05.5.0272.114E	Broca Ø 1,4x Stop 16mmx Comprimento 145mm
05.5.0273.114E	Broca Ø 1,4x Stop 16mmx Comprimento 150mm
05.5.0274.114E	Broca Ø 1,4x Stop 17mmx Comprimento 50mm
05.5.0275.114E	Broca Ø 1,4x Stop 17mmx Comprimento 55mm
05.5.0276.114E	Broca Ø 1,4x Stop 17mmx Comprimento 60mm
05.5.0277.114E	Broca Ø 1,4x Stop 17mmx Comprimento 65mm
05.5.0278.114E	Broca Ø 1,4x Stop 17mmx Comprimento 70mm
05.5.0279.114E	Broca Ø 1,4x Stop 17mmx Comprimento 75mm
05.5.0280.114E	Broca Ø 1,4x Stop 17mmx Comprimento 80mm
05.5.0281.114E	Broca Ø 1,4x Stop 17mmx Comprimento 85mm
05.5.0282.114E	Broca Ø 1,4x Stop 17mmx Comprimento 90mm
05.5.0283.114E	Broca Ø 1,4x Stop 17mmx Comprimento 95mm
05.5.0284.114E	Broca Ø 1,4x Stop 17mmx Comprimento 100mm
05.5.0285.114E	Broca Ø 1,4x Stop 17mmx Comprimento 105mm
05.5.0286.114E	Broca Ø 1,4x Stop 17mmx Comprimento 110mm
05.5.0287.114E	Broca Ø 1,4x Stop 17mmx Comprimento 115mm
05.5.0288.114E	Broca Ø 1,4x Stop 17mmx Comprimento 120mm
05.5.0289.114E	Broca Ø 1,4x Stop 17mmx Comprimento 125mm
05.5.0290.114E	Broca Ø 1,4x Stop 17mmx Comprimento 130mm
05.5.0291.114E	Broca Ø 1,4x Stop 17mmx Comprimento 135mm
05.5.0292.114E	Broca Ø 1,4x Stop 17mmx Comprimento 140mm
05.5.0293.114E	Broca Ø 1,4x Stop 17mmx Comprimento 145mm
05.5.0294.114E	Broca Ø 1,4x Stop 17mmx Comprimento 150mm
05.5.0295.114E	Broca Ø 1,4x Stop 18mmx Comprimento 50mm
05.5.0296.114E	Broca Ø 1,4x Stop 18mmx Comprimento 55mm
05.5.0297.114E	Broca Ø 1,4x Stop 18mmx Comprimento 60mm
05.5.0298.114E	Broca Ø 1,4x Stop 18mmx Comprimento 65mm
05.5.0299.114E	Broca Ø 1,4x Stop 18mmx Comprimento 70mm
05.5.0300.114E	Broca Ø 1,4x Stop 18mmx Comprimento 75mm
05.5.0301.114E	Broca Ø 1,4x Stop 18mmx Comprimento 80mm
05.5.0302.114E	Broca Ø 1,4x Stop 18mmx Comprimento 85mm
05.5.0303.114E	Broca Ø 1,4x Stop 18mmx Comprimento 90mm
05.5.0304.114E	Broca Ø 1,4x Stop 18mmx Comprimento 95mm
05.5.0305.114E	Broca Ø 1,4x Stop 18mmx Comprimento 100mm
05.5.0306.114E	Broca Ø 1,4x Stop 18mmx Comprimento 105mm
05.5.0307.114E	Broca Ø 1,4x Stop 18mmx Comprimento 110mm
05.5.0308.114E	Broca Ø 1,4x Stop 18mmx Comprimento 115mm
05.5.0309.114E	Broca Ø 1,4x Stop 18mmx Comprimento 120mm
05.5.0310.114E	Broca Ø 1,4x Stop 18mmx Comprimento 125mm
05.5.0311.114E	Broca Ø 1,4x Stop 18mmx Comprimento 130mm
05.5.0312.114E	Broca Ø 1,4x Stop 18mmx Comprimento 135mm
05.5.0313.114E	Broca Ø 1,4x Stop 18mmx Comprimento 140mm
05.5.0314.114E	Broca Ø 1,4x Stop 18mmx Comprimento 145mm
05.5.0315.114E	Broca Ø 1,4x Stop 18mmx Comprimento 150mm
05.5.0316.114E	Broca Ø 1,4x Stop 19mmx Comprimento 50mm
05.5.0317.114E	Broca Ø 1,4x Stop 19mmx Comprimento 55mm
05.5.0318.114E	Broca Ø 1,4x Stop 19mmx Comprimento 60mm

05.5.0319.114E	Broca Ø 1,4x Stop 19mmx Comprimento 65mm
05.5.0320.114E	Broca Ø 1,4x Stop 19mmx Comprimento 70mm
05.5.0321.114E	Broca Ø 1,4x Stop 19mmx Comprimento 75mm
05.5.0322.114E	Broca Ø 1,4x Stop 19mmx Comprimento 80mm
05.5.0323.114E	Broca Ø 1,4x Stop 19mmx Comprimento 85mm
05.5.0324.114E	Broca Ø 1,4x Stop 19mmx Comprimento 90mm
05.5.0325.114E	Broca Ø 1,4x Stop 19mmx Comprimento 95mm
05.5.0326.114E	Broca Ø 1,4x Stop 19mmx Comprimento 100mm
05.5.0327.114E	Broca Ø 1,4x Stop 19mmx Comprimento 105mm
05.5.0328.114E	Broca Ø 1,4x Stop 19mmx Comprimento 110mm
05.5.0329.114E	Broca Ø 1,4x Stop 19mmx Comprimento 115mm
05.5.0330.114E	Broca Ø 1,4x Stop 19mmx Comprimento 120mm
05.5.0331.114E	Broca Ø 1,4x Stop 19mmx Comprimento 125mm
05.5.0332.114E	Broca Ø 1,4x Stop 19mmx Comprimento 130mm
05.5.0333.114E	Broca Ø 1,4x Stop 19mmx Comprimento 135mm
05.5.0334.114E	Broca Ø 1,4x Stop 19mmx Comprimento 140mm
05.5.0335.114E	Broca Ø 1,4x Stop 19mmx Comprimento 145mm
05.5.0336.114E	Broca Ø 1,4x Stop 19mmx Comprimento 150mm
05.5.0337.114E	Broca Ø 1,4x Stop 20mmx Comprimento 50mm
05.5.0338.114E	Broca Ø 1,4x Stop 20mmx Comprimento 55mm
05.5.0339.114E	Broca Ø 1,4x Stop 20mmx Comprimento 60mm
05.5.0340.114E	Broca Ø 1,4x Stop 20mmx Comprimento 65mm
05.5.0341.114E	Broca Ø 1,4x Stop 20mmx Comprimento 70mm
05.5.0342.114E	Broca Ø 1,4x Stop 20mmx Comprimento 75mm
05.5.0343.114E	Broca Ø 1,4x Stop 20mmx Comprimento 80mm
05.5.0344.114E	Broca Ø 1,4x Stop 20mmx Comprimento 85mm
05.5.0345.114E	Broca Ø 1,4x Stop 20mmx Comprimento 90mm
05.5.0346.114E	Broca Ø 1,4x Stop 20mmx Comprimento 95mm
05.5.0347.114E	Broca Ø 1,4x Stop 20mmx Comprimento 100mm
05.5.0348.114E	Broca Ø 1,4x Stop 20mmx Comprimento 105mm
05.5.0349.114E	Broca Ø 1,4x Stop 20mmx Comprimento 110mm
05.5.0350.114E	Broca Ø 1,4x Stop 20mmx Comprimento 115mm
05.5.0351.114E	Broca Ø 1,4x Stop 20mmx Comprimento 120mm
05.5.0352.114E	Broca Ø 1,4x Stop 20mmx Comprimento 125mm
05.5.0353.114E	Broca Ø 1,4x Stop 20mmx Comprimento 130mm
05.5.0354.114E	Broca Ø 1,4x Stop 20mmx Comprimento 135mm
05.5.0355.114E	Broca Ø 1,4x Stop 20mmx Comprimento 140mm
05.5.0356.114E	Broca Ø 1,4x Stop 20mmx Comprimento 145mm
05.5.0357.114E	Broca Ø 1,4x Stop 20mmx Comprimento 150mm
05.5.0358.114E	Broca Ø 1,4x Stop 21mmx Comprimento 50mm
05.5.0359.114E	Broca Ø 1,4x Stop 21mmx Comprimento 55mm
05.5.0360.114E	Broca Ø 1,4x Stop 21mmx Comprimento 60mm
05.5.0361.114E	Broca Ø 1,4x Stop 21mmx Comprimento 65mm
05.5.0362.114E	Broca Ø 1,4x Stop 21mmx Comprimento 70mm
05.5.0363.114E	Broca Ø 1,4x Stop 21mmx Comprimento 75mm
05.5.0364.114E	Broca Ø 1,4x Stop 21mmx Comprimento 80mm
05.5.0365.114E	Broca Ø 1,4x Stop 21mmx Comprimento 85mm
05.5.0366.114E	Broca Ø 1,4x Stop 21mmx Comprimento 90mm
05.5.0367.114E	Broca Ø 1,4x Stop 21mmx Comprimento 95mm

05.5.0368.114E	Broca Ø 1,4x Stop 21mmx Comprimento 100mm
05.5.0369.114E	Broca Ø 1,4x Stop 21mmx Comprimento 105mm
05.5.0370.114E	Broca Ø 1,4x Stop 21mmx Comprimento 110mm
05.5.0371.114E	Broca Ø 1,4x Stop 21mmx Comprimento 115mm
05.5.0372.114E	Broca Ø 1,4x Stop 21mmx Comprimento 120mm
05.5.0373.114E	Broca Ø 1,4x Stop 21mmx Comprimento 125mm
05.5.0374.114E	Broca Ø 1,4x Stop 21mmx Comprimento 130mm
05.5.0375.114E	Broca Ø 1,4x Stop 21mmx Comprimento 135mm
05.5.0376.114E	Broca Ø 1,4x Stop 21mmx Comprimento 140mm
05.5.0377.114E	Broca Ø 1,4x Stop 21mmx Comprimento 145mm
05.5.0378.114E	Broca Ø 1,4x Stop 21mmx Comprimento 150mm
05.5.0379.114E	Broca Ø 1,4x Stop 22mmx Comprimento 50mm
05.5.0380.114E	Broca Ø 1,4x Stop 22mmx Comprimento 55mm
05.5.0381.114E	Broca Ø 1,4x Stop 22mmx Comprimento 60mm
05.5.0382.114E	Broca Ø 1,4x Stop 22mmx Comprimento 65mm
05.5.0383.114E	Broca Ø 1,4x Stop 22mmx Comprimento 70mm
05.5.0384.114E	Broca Ø 1,4x Stop 22mmx Comprimento 75mm
05.5.0385.114E	Broca Ø 1,4x Stop 22mmx Comprimento 80mm
05.5.0386.114E	Broca Ø 1,4x Stop 22mmx Comprimento 85mm
05.5.0387.114E	Broca Ø 1,4x Stop 22mmx Comprimento 90mm
05.5.0388.114E	Broca Ø 1,4x Stop 22mmx Comprimento 95mm
05.5.0389.114E	Broca Ø 1,4x Stop 22mmx Comprimento 100mm
05.5.0390.114E	Broca Ø 1,4x Stop 22mmx Comprimento 105mm
05.5.0391.114E	Broca Ø 1,4x Stop 22mmx Comprimento 110mm
05.5.0392.114E	Broca Ø 1,4x Stop 22mmx Comprimento 115mm
05.5.0393.114E	Broca Ø 1,4x Stop 22mmx Comprimento 120mm
05.5.0394.114E	Broca Ø 1,4x Stop 22mmx Comprimento 125mm
05.5.0395.114E	Broca Ø 1,4x Stop 22mmx Comprimento 130mm
05.5.0396.114E	Broca Ø 1,4x Stop 22mmx Comprimento 135mm
05.5.0397.114E	Broca Ø 1,4x Stop 22mmx Comprimento 140mm
05.5.0398.114E	Broca Ø 1,4x Stop 22mmx Comprimento 145mm
05.5.0399.114E	Broca Ø 1,4x Stop 22mmx Comprimento 150mm
05.5.0400.114E	Broca Ø 1,4x Stop 23mmx Comprimento 50mm
05.5.0401.114E	Broca Ø 1,4x Stop 23mmx Comprimento 55mm
05.5.0402.114E	Broca Ø 1,4x Stop 23mmx Comprimento 60mm
05.5.0403.114E	Broca Ø 1,4x Stop 23mmx Comprimento 65mm
05.5.0404.114E	Broca Ø 1,4x Stop 23mmx Comprimento 70mm
05.5.0405.114E	Broca Ø 1,4x Stop 23mmx Comprimento 75mm
05.5.0406.114E	Broca Ø 1,4x Stop 23mmx Comprimento 80mm
05.5.0407.114E	Broca Ø 1,4x Stop 23mmx Comprimento 85mm
05.5.0408.114E	Broca Ø 1,4x Stop 23mmx Comprimento 90mm
05.5.0409.114E	Broca Ø 1,4x Stop 23mmx Comprimento 95mm
05.5.0410.114E	Broca Ø 1,4x Stop 23mmx Comprimento 100mm
05.5.0411.114E	Broca Ø 1,4x Stop 23mmx Comprimento 105mm
05.5.0412.114E	Broca Ø 1,4x Stop 23mmx Comprimento 110mm
05.5.0413.114E	Broca Ø 1,4x Stop 23mmx Comprimento 115mm
05.5.0414.114E	Broca Ø 1,4x Stop 23mmx Comprimento 120mm
05.5.0415.114E	Broca Ø 1,4x Stop 23mmx Comprimento 125mm
05.5.0416.114E	Broca Ø 1,4x Stop 23mmx Comprimento 130mm

05.5.0417.114E	Broca Ø 1,4x Stop 23mmx Comprimento 135mm
05.5.0418.114E	Broca Ø 1,4x Stop 23mmx Comprimento 140mm
05.5.0419.114E	Broca Ø 1,4x Stop 23mmx Comprimento 145mm
05.5.0420.114E	Broca Ø 1,4x Stop 23mmx Comprimento 150mm
05.5.0421.114E	Broca Ø 1,4x Stop 24mmx Comprimento 50mm
05.5.0422.114E	Broca Ø 1,4x Stop 24mmx Comprimento 55mm
05.5.0423.114E	Broca Ø 1,4x Stop 24mmx Comprimento 60mm
05.5.0424.114E	Broca Ø 1,4x Stop 24mmx Comprimento 65mm
05.5.0425.114E	Broca Ø 1,4x Stop 24mmx Comprimento 70mm
05.5.0426.114E	Broca Ø 1,4x Stop 24mmx Comprimento 75mm
05.5.0427.114E	Broca Ø 1,4x Stop 24mmx Comprimento 80mm
05.5.0428.114E	Broca Ø 1,4x Stop 24mmx Comprimento 85mm
05.5.0429.114E	Broca Ø 1,4x Stop 24mmx Comprimento 90mm
05.5.0430.114E	Broca Ø 1,4x Stop 24mmx Comprimento 95mm
05.5.0431.114E	Broca Ø 1,4x Stop 24mmx Comprimento 100mm
05.5.0432.114E	Broca Ø 1,4x Stop 24mmx Comprimento 105mm
05.5.0433.114E	Broca Ø 1,4x Stop 24mmx Comprimento 110mm
05.5.0434.114E	Broca Ø 1,4x Stop 24mmx Comprimento 115mm
05.5.0435.114E	Broca Ø 1,4x Stop 24mmx Comprimento 120mm
05.5.0436.114E	Broca Ø 1,4x Stop 24mmx Comprimento 125mm
05.5.0437.114E	Broca Ø 1,4x Stop 24mmx Comprimento 130mm
05.5.0438.114E	Broca Ø 1,4x Stop 24mmx Comprimento 135mm
05.5.0439.114E	Broca Ø 1,4x Stop 24mmx Comprimento 140mm
05.5.0440.114E	Broca Ø 1,4x Stop 24mmx Comprimento 145mm
05.5.0441.114E	Broca Ø 1,4x Stop 24mmx Comprimento 150mm
05.5.0442.114E	Broca Ø 1,4x Stop 25mmx Comprimento 50mm
05.5.0443.114E	Broca Ø 1,4x Stop 25mmx Comprimento 55mm
05.5.0444.114E	Broca Ø 1,4x Stop 25mmx Comprimento 60mm
05.5.0445.114E	Broca Ø 1,4x Stop 25mmx Comprimento 65mm
05.5.0446.114E	Broca Ø 1,4x Stop 25mmx Comprimento 70mm
05.5.0447.114E	Broca Ø 1,4x Stop 25mmx Comprimento 75mm
05.5.0448.114E	Broca Ø 1,4x Stop 25mmx Comprimento 80mm
05.5.0449.114E	Broca Ø 1,4x Stop 25mmx Comprimento 85mm
05.5.0450.114E	Broca Ø 1,4x Stop 25mmx Comprimento 90mm
05.5.0451.114E	Broca Ø 1,4x Stop 25mmx Comprimento 95mm
05.5.0452.114E	Broca Ø 1,4x Stop 25mmx Comprimento 100mm
05.5.0453.114E	Broca Ø 1,4x Stop 25mmx Comprimento 105mm
05.5.0454.114E	Broca Ø 1,4x Stop 25mmx Comprimento 110mm
05.5.0455.114E	Broca Ø 1,4x Stop 25mmx Comprimento 115mm
05.5.0456.114E	Broca Ø 1,4x Stop 25mmx Comprimento 120mm
05.5.0457.114E	Broca Ø 1,4x Stop 25mmx Comprimento 125mm
05.5.0458.114E	Broca Ø 1,4x Stop 25mmx Comprimento 130mm
05.5.0459.114E	Broca Ø 1,4x Stop 25mmx Comprimento 135mm
05.5.0460.114E	Broca Ø 1,4x Stop 25mmx Comprimento 140mm
05.5.0461.114E	Broca Ø 1,4x Stop 25mmx Comprimento 145mm
05.5.0462.114E	Broca Ø 1,4x Stop 25mmx Comprimento 150mm
05.5.0463.114E	Broca Ø 1,4x Stop 30mmx Comprimento 50mm
05.5.0464.114E	Broca Ø 1,4x Stop 30mmx Comprimento 55mm
05.5.0465.114E	Broca Ø 1,4x Stop 30mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0466.114E	Broca Ø 1,4x Stop 30mmx Comprimento 65mm
05.5.0467.114E	Broca Ø 1,4x Stop 30mmx Comprimento 70mm
05.5.0468.114E	Broca Ø 1,4x Stop 30mmx Comprimento 75mm
05.5.0469.114E	Broca Ø 1,4x Stop 30mmx Comprimento 80mm
05.5.0470.114E	Broca Ø 1,4x Stop 30mmx Comprimento 85mm
05.5.0471.114E	Broca Ø 1,4x Stop 30mmx Comprimento 90mm
05.5.0472.114E	Broca Ø 1,4x Stop 30mmx Comprimento 95mm
05.5.0473.114E	Broca Ø 1,4x Stop 30mmx Comprimento 100mm
05.5.0474.114E	Broca Ø 1,4x Stop 30mmx Comprimento 105mm
05.5.0475.114E	Broca Ø 1,4x Stop 30mmx Comprimento 110mm
05.5.0476.114E	Broca Ø 1,4x Stop 30mmx Comprimento 115mm
05.5.0477.114E	Broca Ø 1,4x Stop 30mmx Comprimento 120mm
05.5.0478.114E	Broca Ø 1,4x Stop 30mmx Comprimento 125mm
05.5.0479.114E	Broca Ø 1,4x Stop 30mmx Comprimento 130mm
05.5.0480.114E	Broca Ø 1,4x Stop 30mmx Comprimento 135mm
05.5.0481.114E	Broca Ø 1,4x Stop 30mmx Comprimento 140mm
05.5.0482.114E	Broca Ø 1,4x Stop 30mmx Comprimento 145mm
05.5.0483.114E	Broca Ø 1,4x Stop 30mmx Comprimento 150mm
05.5.0484.114E	Broca Ø 1,4x Stop 35mmx Comprimento 50mm
05.5.0485.114E	Broca Ø 1,4x Stop 35mmx Comprimento 55mm
05.5.0486.114E	Broca Ø 1,4x Stop 35mmx Comprimento 60mm
05.5.0487.114E	Broca Ø 1,4x Stop 35mmx Comprimento 65mm
05.5.0488.114E	Broca Ø 1,4x Stop 35mmx Comprimento 70mm
05.5.0489.114E	Broca Ø 1,4x Stop 35mmx Comprimento 75mm
05.5.0490.114E	Broca Ø 1,4x Stop 35mmx Comprimento 80mm
05.5.0491.114E	Broca Ø 1,4x Stop 35mmx Comprimento 85mm
05.5.0492.114E	Broca Ø 1,4x Stop 35mmx Comprimento 90mm
05.5.0493.114E	Broca Ø 1,4x Stop 35mmx Comprimento 95mm
05.5.0494.114E	Broca Ø 1,4x Stop 35mmx Comprimento 100mm
05.5.0495.114E	Broca Ø 1,4x Stop 35mmx Comprimento 105mm
05.5.0496.114E	Broca Ø 1,4x Stop 35mmx Comprimento 110mm
05.5.0497.114E	Broca Ø 1,4x Stop 35mmx Comprimento 115mm
05.5.0498.114E	Broca Ø 1,4x Stop 35mmx Comprimento 120mm
05.5.0499.114E	Broca Ø 1,4x Stop 35mmx Comprimento 125mm
05.5.0500.114E	Broca Ø 1,4x Stop 35mmx Comprimento 130mm
05.5.0501.114E	Broca Ø 1,4x Stop 35mmx Comprimento 135mm
05.5.0502.114E	Broca Ø 1,4x Stop 35mmx Comprimento 140mm
05.5.0503.114E	Broca Ø 1,4x Stop 35mmx Comprimento 145mm
05.5.0504.114E	Broca Ø 1,4x Stop 35mmx Comprimento 150mm
05.5.0001.115E	Broca Ø 1,5x Stop 4mmx Comprimento 50mm
05.5.0002.115E	Broca Ø 1,5x Stop 4mmx Comprimento 55mm
05.5.0003.115E	Broca Ø 1,5x Stop 4mmx Comprimento 60mm
05.5.0004.115E	Broca Ø 1,5x Stop 4mmx Comprimento 65mm
05.5.0005.115E	Broca Ø 1,5x Stop 4mmx Comprimento 70mm
05.5.0006.115E	Broca Ø 1,5x Stop 4mmx Comprimento 75mm
05.5.0007.115E	Broca Ø 1,5x Stop 4mmx Comprimento 80mm
05.5.0008.115E	Broca Ø 1,5x Stop 4mmx Comprimento 85mm
05.5.0009.115E	Broca Ø 1,5x Stop 4mmx Comprimento 90mm
05.5.0010.115E	Broca Ø 1,5x Stop 4mmx Comprimento 95mm

05.5.0011.115E	Broca Ø 1,5x Stop 4mmx Comprimento 100mm
05.5.0012.115E	Broca Ø 1,5x Stop 4mmx Comprimento 105mm
05.5.0013.115E	Broca Ø 1,5x Stop 4mmx Comprimento 110mm
05.5.0014.115E	Broca Ø 1,5x Stop 4mmx Comprimento 115mm
05.5.0015.115E	Broca Ø 1,5x Stop 4mmx Comprimento 120mm
05.5.0016.115E	Broca Ø 1,5x Stop 4mmx Comprimento 125mm
05.5.0017.115E	Broca Ø 1,5x Stop 4mmx Comprimento 130mm
05.5.0018.115E	Broca Ø 1,5x Stop 4mmx Comprimento 135mm
05.5.0019.115E	Broca Ø 1,5x Stop 4mmx Comprimento 140mm
05.5.0020.115E	Broca Ø 1,5x Stop 4mmx Comprimento 145mm
05.5.0021.115E	Broca Ø 1,5x Stop 4mmx Comprimento 150mm
05.5.0022.115E	Broca Ø 1,5x Stop 5mmx Comprimento 50mm
05.5.0023.115E	Broca Ø 1,5x Stop 5mmx Comprimento 55mm
05.5.0024.115E	Broca Ø 1,5x Stop 5mmx Comprimento 60mm
05.5.0025.115E	Broca Ø 1,5x Stop 5mmx Comprimento 65mm
05.5.0026.115E	Broca Ø 1,5x Stop 5mmx Comprimento 70mm
05.5.0027.115E	Broca Ø 1,5x Stop 5mmx Comprimento 75mm
05.5.0028.115E	Broca Ø 1,5x Stop 5mmx Comprimento 80mm
05.5.0029.115E	Broca Ø 1,5x Stop 5mmx Comprimento 85mm
05.5.0030.115E	Broca Ø 1,5x Stop 5mmx Comprimento 90mm
05.5.0031.115E	Broca Ø 1,5x Stop 5mmx Comprimento 95mm
05.5.0032.115E	Broca Ø 1,5x Stop 5mmx Comprimento 100mm
05.5.0033.115E	Broca Ø 1,5x Stop 5mmx Comprimento 105mm
05.5.0034.115E	Broca Ø 1,5x Stop 5mmx Comprimento 110mm
05.5.0035.115E	Broca Ø 1,5x Stop 5mmx Comprimento 115mm
05.5.0036.115E	Broca Ø 1,5x Stop 5mmx Comprimento 120mm
05.5.0037.115E	Broca Ø 1,5x Stop 5mmx Comprimento 125mm
05.5.0038.115E	Broca Ø 1,5x Stop 5mmx Comprimento 130mm
05.5.0039.115E	Broca Ø 1,5x Stop 5mmx Comprimento 135mm
05.5.0040.115E	Broca Ø 1,5x Stop 5mmx Comprimento 140mm
05.5.0041.115E	Broca Ø 1,5x Stop 5mmx Comprimento 145mm
05.5.0042.115E	Broca Ø 1,5x Stop 5mmx Comprimento 150mm
05.5.0043.115E	Broca Ø 1,5x Stop 6mmx Comprimento 50mm
05.5.0044.115E	Broca Ø 1,5x Stop 6mmx Comprimento 55mm
05.5.0045.115E	Broca Ø 1,5x Stop 6mmx Comprimento 60mm
05.5.0046.115E	Broca Ø 1,5x Stop 6mmx Comprimento 65mm
05.5.0047.115E	Broca Ø 1,5x Stop 6mmx Comprimento 70mm
05.5.0048.115E	Broca Ø 1,5x Stop 6mmx Comprimento 75mm
05.5.0049.115E	Broca Ø 1,5x Stop 6mmx Comprimento 80mm
05.5.0050.115E	Broca Ø 1,5x Stop 6mmx Comprimento 85mm
05.5.0051.115E	Broca Ø 1,5x Stop 6mmx Comprimento 90mm
05.5.0052.115E	Broca Ø 1,5x Stop 6mmx Comprimento 95mm
05.5.0053.115E	Broca Ø 1,5x Stop 6mmx Comprimento 100mm
05.5.0054.115E	Broca Ø 1,5x Stop 6mmx Comprimento 105mm
05.5.0055.115E	Broca Ø 1,5x Stop 6mmx Comprimento 110mm
05.5.0056.115E	Broca Ø 1,5x Stop 6mmx Comprimento 115mm
05.5.0057.115E	Broca Ø 1,5x Stop 6mmx Comprimento 120mm
05.5.0058.115E	Broca Ø 1,5x Stop 6mmx Comprimento 125mm
05.5.0059.115E	Broca Ø 1,5x Stop 6mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0060.115E	Broca Ø 1,5x Stop 6mmx Comprimento 135mm
05.5.0061.115E	Broca Ø 1,5x Stop 6mmx Comprimento 140mm
05.5.0062.115E	Broca Ø 1,5x Stop 6mmx Comprimento 145mm
05.5.0063.115E	Broca Ø 1,5x Stop 6mmx Comprimento 150mm
05.5.0064.115E	Broca Ø 1,5x Stop 7mmx Comprimento 50mm
05.5.0065.115E	Broca Ø 1,5x Stop 7mmx Comprimento 55mm
05.5.0066.115E	Broca Ø 1,5x Stop 7mmx Comprimento 60mm
05.5.0067.115E	Broca Ø 1,5x Stop 7mmx Comprimento 65mm
05.5.0068.115E	Broca Ø 1,5x Stop 7mmx Comprimento 70mm
05.5.0069.115E	Broca Ø 1,5x Stop 7mmx Comprimento 75mm
05.5.0070.115E	Broca Ø 1,5x Stop 7mmx Comprimento 80mm
05.5.0071.115E	Broca Ø 1,5x Stop 7mmx Comprimento 85mm
05.5.0072.115E	Broca Ø 1,5x Stop 7mmx Comprimento 90mm
05.5.0073.115E	Broca Ø 1,5x Stop 7mmx Comprimento 95mm
05.5.0074.115E	Broca Ø 1,5x Stop 7mmx Comprimento 100mm
05.5.0075.115E	Broca Ø 1,5x Stop 7mmx Comprimento 105mm
05.5.0076.115E	Broca Ø 1,5x Stop 7mmx Comprimento 110mm
05.5.0077.115E	Broca Ø 1,5x Stop 7mmx Comprimento 115mm
05.5.0078.115E	Broca Ø 1,5x Stop 7mmx Comprimento 120mm
05.5.0079.115E	Broca Ø 1,5x Stop 7mmx Comprimento 125mm
05.5.0080.115E	Broca Ø 1,5x Stop 7mmx Comprimento 130mm
05.5.0081.115E	Broca Ø 1,5x Stop 7mmx Comprimento 135mm
05.5.0082.115E	Broca Ø 1,5x Stop 7mmx Comprimento 140mm
05.5.0083.115E	Broca Ø 1,5x Stop 7mmx Comprimento 145mm
05.5.0084.115E	Broca Ø 1,5x Stop 7mmx Comprimento 150mm
05.5.0085.115E	Broca Ø 1,5x Stop 8mmx Comprimento 50mm
05.5.0086.115E	Broca Ø 1,5x Stop 8mmx Comprimento 55mm
05.5.0087.115E	Broca Ø 1,5x Stop 8mmx Comprimento 60mm
05.5.0088.115E	Broca Ø 1,5x Stop 8mmx Comprimento 65mm
05.5.0089.115E	Broca Ø 1,5x Stop 8mmx Comprimento 70mm
05.5.0090.115E	Broca Ø 1,5x Stop 8mmx Comprimento 75mm
05.5.0091.115E	Broca Ø 1,5x Stop 8mmx Comprimento 80mm
05.5.0092.115E	Broca Ø 1,5x Stop 8mmx Comprimento 85mm
05.5.0093.115E	Broca Ø 1,5x Stop 8mmx Comprimento 90mm
05.5.0094.115E	Broca Ø 1,5x Stop 8mmx Comprimento 95mm
05.5.0095.115E	Broca Ø 1,5x Stop 8mmx Comprimento 100mm
05.5.0096.115E	Broca Ø 1,5x Stop 8mmx Comprimento 105mm
05.5.0097.115E	Broca Ø 1,5x Stop 8mmx Comprimento 110mm
05.5.0098.115E	Broca Ø 1,5x Stop 8mmx Comprimento 115mm
05.5.0099.115E	Broca Ø 1,5x Stop 8mmx Comprimento 120mm
05.5.0100.115E	Broca Ø 1,5x Stop 8mmx Comprimento 125mm
05.5.0101.115E	Broca Ø 1,5x Stop 8mmx Comprimento 130mm
05.5.0102.115E	Broca Ø 1,5x Stop 8mmx Comprimento 135mm
05.5.0103.115E	Broca Ø 1,5x Stop 8mmx Comprimento 140mm
05.5.0104.115E	Broca Ø 1,5x Stop 8mmx Comprimento 145mm
05.5.0105.115E	Broca Ø 1,5x Stop 8mmx Comprimento 150mm
05.5.0106.115E	Broca Ø 1,5x Stop 9mmx Comprimento 50mm
05.5.0107.115E	Broca Ø 1,5x Stop 9mmx Comprimento 55mm
05.5.0108.115E	Broca Ø 1,5x Stop 9mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0109.115E	Broca Ø 1,5x Stop 9mmx Comprimento 65mm
05.5.0110.115E	Broca Ø 1,5x Stop 9mmx Comprimento 70mm
05.5.0111.115E	Broca Ø 1,5x Stop 9mmx Comprimento 75mm
05.5.0112.115E	Broca Ø 1,5x Stop 9mmx Comprimento 80mm
05.5.0113.115E	Broca Ø 1,5x Stop 9mmx Comprimento 85mm
05.5.0114.115E	Broca Ø 1,5x Stop 9mmx Comprimento 90mm
05.5.0115.115E	Broca Ø 1,5x Stop 9mmx Comprimento 95mm
05.5.0116.115E	Broca Ø 1,5x Stop 9mmx Comprimento 100mm
05.5.0117.115E	Broca Ø 1,5x Stop 9mmx Comprimento 105mm
05.5.0118.115E	Broca Ø 1,5x Stop 9mmx Comprimento 110mm
05.5.0119.115E	Broca Ø 1,5x Stop 9mmx Comprimento 115mm
05.5.0120.115E	Broca Ø 1,5x Stop 9mmx Comprimento 120mm
05.5.0121.115E	Broca Ø 1,5x Stop 9mmx Comprimento 125mm
05.5.0122.115E	Broca Ø 1,5x Stop 9mmx Comprimento 130mm
05.5.0123.115E	Broca Ø 1,5x Stop 9mmx Comprimento 135mm
05.5.0124.115E	Broca Ø 1,5x Stop 9mmx Comprimento 140mm
05.5.0125.115E	Broca Ø 1,5x Stop 9mmx Comprimento 145mm
05.5.0126.115E	Broca Ø 1,5x Stop 9mmx Comprimento 150mm
05.5.0127.115E	Broca Ø 1,5x Stop 10mmx Comprimento 50mm
05.5.0128.115E	Broca Ø 1,5x Stop 10mmx Comprimento 55mm
05.5.0129.115E	Broca Ø 1,5x Stop 10mmx Comprimento 60mm
05.5.0130.115E	Broca Ø 1,5x Stop 10mmx Comprimento 65mm
05.5.0131.115E	Broca Ø 1,5x Stop 10mmx Comprimento 70mm
05.5.0132.115E	Broca Ø 1,5x Stop 10mmx Comprimento 75mm
05.5.0133.115E	Broca Ø 1,5x Stop 10mmx Comprimento 80mm
05.5.0134.115E	Broca Ø 1,5x Stop 10mmx Comprimento 85mm
05.5.0135.115E	Broca Ø 1,5x Stop 10mmx Comprimento 90mm
05.5.0136.115E	Broca Ø 1,5x Stop 10mmx Comprimento 95mm
05.5.0137.115E	Broca Ø 1,5x Stop 10mmx Comprimento 100mm
05.5.0138.115E	Broca Ø 1,5x Stop 10mmx Comprimento 105mm
05.5.0139.115E	Broca Ø 1,5x Stop 10mmx Comprimento 110mm
05.5.0140.115E	Broca Ø 1,5x Stop 10mmx Comprimento 115mm
05.5.0141.115E	Broca Ø 1,5x Stop 10mmx Comprimento 120mm
05.5.0142.115E	Broca Ø 1,5x Stop 10mmx Comprimento 125mm
05.5.0143.115E	Broca Ø 1,5x Stop 10mmx Comprimento 130mm
05.5.0144.115E	Broca Ø 1,5x Stop 10mmx Comprimento 135mm
05.5.0145.115E	Broca Ø 1,5x Stop 10mmx Comprimento 140mm
05.5.0146.115E	Broca Ø 1,5x Stop 10mmx Comprimento 145mm
05.5.0147.115E	Broca Ø 1,5x Stop 10mmx Comprimento 150mm
05.5.0148.115E	Broca Ø 1,5x Stop 11mmx Comprimento 50mm
05.5.0149.115E	Broca Ø 1,5x Stop 11mmx Comprimento 55mm
05.5.0150.115E	Broca Ø 1,5x Stop 11mmx Comprimento 60mm
05.5.0151.115E	Broca Ø 1,5x Stop 11mmx Comprimento 65mm
05.5.0152.115E	Broca Ø 1,5x Stop 11mmx Comprimento 70mm
05.5.0153.115E	Broca Ø 1,5x Stop 11mmx Comprimento 75mm
05.5.0154.115E	Broca Ø 1,5x Stop 11mmx Comprimento 80mm
05.5.0155.115E	Broca Ø 1,5x Stop 11mmx Comprimento 85mm
05.5.0156.115E	Broca Ø 1,5x Stop 11mmx Comprimento 90mm
05.5.0157.115E	Broca Ø 1,5x Stop 11mmx Comprimento 95mm

05.5.0158.115E	Broca Ø 1,5x Stop 11mmx Comprimento 100mm
05.5.0159.115E	Broca Ø 1,5x Stop 11mmx Comprimento 105mm
05.5.0160.115E	Broca Ø 1,5x Stop 11mmx Comprimento 110mm
05.5.0161.115E	Broca Ø 1,5x Stop 11mmx Comprimento 115mm
05.5.0162.115E	Broca Ø 1,5x Stop 11mmx Comprimento 120mm
05.5.0163.115E	Broca Ø 1,5x Stop 11mmx Comprimento 125mm
05.5.0164.115E	Broca Ø 1,5x Stop 11mmx Comprimento 130mm
05.5.0165.115E	Broca Ø 1,5x Stop 11mmx Comprimento 135mm
05.5.0166.115E	Broca Ø 1,5x Stop 11mmx Comprimento 140mm
05.5.0167.115E	Broca Ø 1,5x Stop 11mmx Comprimento 145mm
05.5.0168.115E	Broca Ø 1,5x Stop 11mmx Comprimento 150mm
05.5.0169.115E	Broca Ø 1,5x Stop 12mmx Comprimento 50mm
05.5.0170.115E	Broca Ø 1,5x Stop 12mmx Comprimento 55mm
05.5.0171.115E	Broca Ø 1,5x Stop 12mmx Comprimento 60mm
05.5.0172.115E	Broca Ø 1,5x Stop 12mmx Comprimento 65mm
05.5.0173.115E	Broca Ø 1,5x Stop 12mmx Comprimento 70mm
05.5.0174.115E	Broca Ø 1,5x Stop 12mmx Comprimento 75mm
05.5.0175.115E	Broca Ø 1,5x Stop 12mmx Comprimento 80mm
05.5.0176.115E	Broca Ø 1,5x Stop 12mmx Comprimento 85mm
05.5.0177.115E	Broca Ø 1,5x Stop 12mmx Comprimento 90mm
05.5.0178.115E	Broca Ø 1,5x Stop 12mmx Comprimento 95mm
05.5.0179.115E	Broca Ø 1,5x Stop 12mmx Comprimento 100mm
05.5.0180.115E	Broca Ø 1,5x Stop 12mmx Comprimento 105mm
05.5.0181.115E	Broca Ø 1,5x Stop 12mmx Comprimento 110mm
05.5.0182.115E	Broca Ø 1,5x Stop 12mmx Comprimento 115mm
05.5.0183.115E	Broca Ø 1,5x Stop 12mmx Comprimento 120mm
05.5.0184.115E	Broca Ø 1,5x Stop 12mmx Comprimento 125mm
05.5.0185.115E	Broca Ø 1,5x Stop 12mmx Comprimento 130mm
05.5.0186.115E	Broca Ø 1,5x Stop 12mmx Comprimento 135mm
05.5.0187.115E	Broca Ø 1,5x Stop 12mmx Comprimento 140mm
05.5.0188.115E	Broca Ø 1,5x Stop 12mmx Comprimento 145mm
05.5.0189.115E	Broca Ø 1,5x Stop 12mmx Comprimento 150mm
05.5.0190.115E	Broca Ø 1,5x Stop 13mmx Comprimento 50mm
05.5.0191.115E	Broca Ø 1,5x Stop 13mmx Comprimento 55mm
05.5.0192.115E	Broca Ø 1,5x Stop 13mmx Comprimento 60mm
05.5.0193.115E	Broca Ø 1,5x Stop 13mmx Comprimento 65mm
05.5.0194.115E	Broca Ø 1,5x Stop 13mmx Comprimento 70mm
05.5.0195.115E	Broca Ø 1,5x Stop 13mmx Comprimento 75mm
05.5.0196.115E	Broca Ø 1,5x Stop 13mmx Comprimento 80mm
05.5.0197.115E	Broca Ø 1,5x Stop 13mmx Comprimento 85mm
05.5.0198.115E	Broca Ø 1,5x Stop 13mmx Comprimento 90mm
05.5.0199.115E	Broca Ø 1,5x Stop 13mmx Comprimento 95mm
05.5.0200.115E	Broca Ø 1,5x Stop 13mmx Comprimento 100mm
05.5.0201.115E	Broca Ø 1,5x Stop 13mmx Comprimento 105mm
05.5.0202.115E	Broca Ø 1,5x Stop 13mmx Comprimento 110mm
05.5.0203.115E	Broca Ø 1,5x Stop 13mmx Comprimento 115mm
05.5.0204.115E	Broca Ø 1,5x Stop 13mmx Comprimento 120mm
05.5.0205.115E	Broca Ø 1,5x Stop 13mmx Comprimento 125mm
05.5.0206.115E	Broca Ø 1,5x Stop 13mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0207.115E	Broca Ø 1,5x Stop 13mmx Comprimento 135mm
05.5.0208.115E	Broca Ø 1,5x Stop 13mmx Comprimento 140mm
05.5.0209.115E	Broca Ø 1,5x Stop 13mmx Comprimento 145mm
05.5.0210.115E	Broca Ø 1,5x Stop 13mmx Comprimento 150mm
05.5.0211.115E	Broca Ø 1,5x Stop 14mmx Comprimento 50mm
05.5.0212.115E	Broca Ø 1,5x Stop 14mmx Comprimento 55mm
05.5.0213.115E	Broca Ø 1,5x Stop 14mmx Comprimento 60mm
05.5.0214.115E	Broca Ø 1,5x Stop 14mmx Comprimento 65mm
05.5.0215.115E	Broca Ø 1,5x Stop 14mmx Comprimento 70mm
05.5.0216.115E	Broca Ø 1,5x Stop 14mmx Comprimento 75mm
05.5.0217.115E	Broca Ø 1,5x Stop 14mmx Comprimento 80mm
05.5.0218.115E	Broca Ø 1,5x Stop 14mmx Comprimento 85mm
05.5.0219.115E	Broca Ø 1,5x Stop 14mmx Comprimento 90mm
05.5.0220.115E	Broca Ø 1,5x Stop 14mmx Comprimento 95mm
05.5.0221.115E	Broca Ø 1,5x Stop 14mmx Comprimento 100mm
05.5.0222.115E	Broca Ø 1,5x Stop 14mmx Comprimento 105mm
05.5.0223.115E	Broca Ø 1,5x Stop 14mmx Comprimento 110mm
05.5.0224.115E	Broca Ø 1,5x Stop 14mmx Comprimento 115mm
05.5.0225.115E	Broca Ø 1,5x Stop 14mmx Comprimento 120mm
05.5.0226.115E	Broca Ø 1,5x Stop 14mmx Comprimento 125mm
05.5.0227.115E	Broca Ø 1,5x Stop 14mmx Comprimento 130mm
05.5.0228.115E	Broca Ø 1,5x Stop 14mmx Comprimento 135mm
05.5.0229.115E	Broca Ø 1,5x Stop 14mmx Comprimento 140mm
05.5.0230.115E	Broca Ø 1,5x Stop 14mmx Comprimento 145mm
05.5.0231.115E	Broca Ø 1,5x Stop 14mmx Comprimento 150mm
05.5.0232.115E	Broca Ø 1,5x Stop 15mmx Comprimento 50mm
05.5.0233.115E	Broca Ø 1,5x Stop 15mmx Comprimento 55mm
05.5.0234.115E	Broca Ø 1,5x Stop 15mmx Comprimento 60mm
05.5.0235.115E	Broca Ø 1,5x Stop 15mmx Comprimento 65mm
05.5.0236.115E	Broca Ø 1,5x Stop 15mmx Comprimento 70mm
05.5.0237.115E	Broca Ø 1,5x Stop 15mmx Comprimento 75mm
05.5.0238.115E	Broca Ø 1,5x Stop 15mmx Comprimento 80mm
05.5.0239.115E	Broca Ø 1,5x Stop 15mmx Comprimento 85mm
05.5.0240.115E	Broca Ø 1,5x Stop 15mmx Comprimento 90mm
05.5.0241.115E	Broca Ø 1,5x Stop 15mmx Comprimento 95mm
05.5.0242.115E	Broca Ø 1,5x Stop 15mmx Comprimento 100mm
05.5.0243.115E	Broca Ø 1,5x Stop 15mmx Comprimento 105mm
05.5.0244.115E	Broca Ø 1,5x Stop 15mmx Comprimento 110mm
05.5.0245.115E	Broca Ø 1,5x Stop 15mmx Comprimento 115mm
05.5.0246.115E	Broca Ø 1,5x Stop 15mmx Comprimento 120mm
05.5.0247.115E	Broca Ø 1,5x Stop 15mmx Comprimento 125mm
05.5.0248.115E	Broca Ø 1,5x Stop 15mmx Comprimento 130mm
05.5.0249.115E	Broca Ø 1,5x Stop 15mmx Comprimento 135mm
05.5.0250.115E	Broca Ø 1,5x Stop 15mmx Comprimento 140mm
05.5.0251.115E	Broca Ø 1,5x Stop 15mmx Comprimento 145mm
05.5.0252.115E	Broca Ø 1,5x Stop 15mmx Comprimento 150mm
05.5.0253.115E	Broca Ø 1,5x Stop 16mmx Comprimento 50mm
05.5.0254.115E	Broca Ø 1,5x Stop 16mmx Comprimento 55mm
05.5.0255.115E	Broca Ø 1,5x Stop 16mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0256.115E	Broca Ø 1,5x Stop 16mmx Comprimento 65mm
05.5.0257.115E	Broca Ø 1,5x Stop 16mmx Comprimento 70mm
05.5.0258.115E	Broca Ø 1,5x Stop 16mmx Comprimento 75mm
05.5.0259.115E	Broca Ø 1,5x Stop 16mmx Comprimento 80mm
05.5.0260.115E	Broca Ø 1,5x Stop 16mmx Comprimento 85mm
05.5.0261.115E	Broca Ø 1,5x Stop 16mmx Comprimento 90mm
05.5.0262.115E	Broca Ø 1,5x Stop 16mmx Comprimento 95mm
05.5.0263.115E	Broca Ø 1,5x Stop 16mmx Comprimento 100mm
05.5.0264.115E	Broca Ø 1,5x Stop 16mmx Comprimento 105mm
05.5.0265.115E	Broca Ø 1,5x Stop 16mmx Comprimento 110mm
05.5.0266.115E	Broca Ø 1,5x Stop 16mmx Comprimento 115mm
05.5.0267.115E	Broca Ø 1,5x Stop 16mmx Comprimento 120mm
05.5.0268.115E	Broca Ø 1,5x Stop 16mmx Comprimento 125mm
05.5.0269.115E	Broca Ø 1,5x Stop 16mmx Comprimento 130mm
05.5.0270.115E	Broca Ø 1,5x Stop 16mmx Comprimento 135mm
05.5.0271.115E	Broca Ø 1,5x Stop 16mmx Comprimento 140mm
05.5.0272.115E	Broca Ø 1,5x Stop 16mmx Comprimento 145mm
05.5.0273.115E	Broca Ø 1,5x Stop 16mmx Comprimento 150mm
05.5.0274.115E	Broca Ø 1,5x Stop 17mmx Comprimento 50mm
05.5.0275.115E	Broca Ø 1,5x Stop 17mmx Comprimento 55mm
05.5.0276.115E	Broca Ø 1,5x Stop 17mmx Comprimento 60mm
05.5.0277.115E	Broca Ø 1,5x Stop 17mmx Comprimento 65mm
05.5.0278.115E	Broca Ø 1,5x Stop 17mmx Comprimento 70mm
05.5.0279.115E	Broca Ø 1,5x Stop 17mmx Comprimento 75mm
05.5.0280.115E	Broca Ø 1,5x Stop 17mmx Comprimento 80mm
05.5.0281.115E	Broca Ø 1,5x Stop 17mmx Comprimento 85mm
05.5.0282.115E	Broca Ø 1,5x Stop 17mmx Comprimento 90mm
05.5.0283.115E	Broca Ø 1,5x Stop 17mmx Comprimento 95mm
05.5.0284.115E	Broca Ø 1,5x Stop 17mmx Comprimento 100mm
05.5.0285.115E	Broca Ø 1,5x Stop 17mmx Comprimento 105mm
05.5.0286.115E	Broca Ø 1,5x Stop 17mmx Comprimento 110mm
05.5.0287.115E	Broca Ø 1,5x Stop 17mmx Comprimento 115mm
05.5.0288.115E	Broca Ø 1,5x Stop 17mmx Comprimento 120mm
05.5.0289.115E	Broca Ø 1,5x Stop 17mmx Comprimento 125mm
05.5.0290.115E	Broca Ø 1,5x Stop 17mmx Comprimento 130mm
05.5.0291.115E	Broca Ø 1,5x Stop 17mmx Comprimento 135mm
05.5.0292.115E	Broca Ø 1,5x Stop 17mmx Comprimento 140mm
05.5.0293.115E	Broca Ø 1,5x Stop 17mmx Comprimento 145mm
05.5.0294.115E	Broca Ø 1,5x Stop 17mmx Comprimento 150mm
05.5.0295.115E	Broca Ø 1,5x Stop 18mmx Comprimento 50mm
05.5.0296.115E	Broca Ø 1,5x Stop 18mmx Comprimento 55mm
05.5.0297.115E	Broca Ø 1,5x Stop 18mmx Comprimento 60mm
05.5.0298.115E	Broca Ø 1,5x Stop 18mmx Comprimento 65mm
05.5.0299.115E	Broca Ø 1,5x Stop 18mmx Comprimento 70mm
05.5.0300.115E	Broca Ø 1,5x Stop 18mmx Comprimento 75mm
05.5.0301.115E	Broca Ø 1,5x Stop 18mmx Comprimento 80mm
05.5.0302.115E	Broca Ø 1,5x Stop 18mmx Comprimento 85mm
05.5.0303.115E	Broca Ø 1,5x Stop 18mmx Comprimento 90mm
05.5.0304.115E	Broca Ø 1,5x Stop 18mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0305.115E	Broca Ø 1,5x Stop 18mmx Comprimento 100mm
05.5.0306.115E	Broca Ø 1,5x Stop 18mmx Comprimento 105mm
05.5.0307.115E	Broca Ø 1,5x Stop 18mmx Comprimento 110mm
05.5.0308.115E	Broca Ø 1,5x Stop 18mmx Comprimento 115mm
05.5.0309.115E	Broca Ø 1,5x Stop 18mmx Comprimento 120mm
05.5.0310.115E	Broca Ø 1,5x Stop 18mmx Comprimento 125mm
05.5.0311.115E	Broca Ø 1,5x Stop 18mmx Comprimento 130mm
05.5.0312.115E	Broca Ø 1,5x Stop 18mmx Comprimento 135mm
05.5.0313.115E	Broca Ø 1,5x Stop 18mmx Comprimento 140mm
05.5.0314.115E	Broca Ø 1,5x Stop 18mmx Comprimento 145mm
05.5.0315.115E	Broca Ø 1,5x Stop 18mmx Comprimento 150mm
05.5.0316.115E	Broca Ø 1,5x Stop 19mmx Comprimento 50mm
05.5.0317.115E	Broca Ø 1,5x Stop 19mmx Comprimento 55mm
05.5.0318.115E	Broca Ø 1,5x Stop 19mmx Comprimento 60mm
05.5.0319.115E	Broca Ø 1,5x Stop 19mmx Comprimento 65mm
05.5.0320.115E	Broca Ø 1,5x Stop 19mmx Comprimento 70mm
05.5.0321.115E	Broca Ø 1,5x Stop 19mmx Comprimento 75mm
05.5.0322.115E	Broca Ø 1,5x Stop 19mmx Comprimento 80mm
05.5.0323.115E	Broca Ø 1,5x Stop 19mmx Comprimento 85mm
05.5.0324.115E	Broca Ø 1,5x Stop 19mmx Comprimento 90mm
05.5.0325.115E	Broca Ø 1,5x Stop 19mmx Comprimento 95mm
05.5.0326.115E	Broca Ø 1,5x Stop 19mmx Comprimento 100mm
05.5.0327.115E	Broca Ø 1,5x Stop 19mmx Comprimento 105mm
05.5.0328.115E	Broca Ø 1,5x Stop 19mmx Comprimento 110mm
05.5.0329.115E	Broca Ø 1,5x Stop 19mmx Comprimento 115mm
05.5.0330.115E	Broca Ø 1,5x Stop 19mmx Comprimento 120mm
05.5.0331.115E	Broca Ø 1,5x Stop 19mmx Comprimento 125mm
05.5.0332.115E	Broca Ø 1,5x Stop 19mmx Comprimento 130mm
05.5.0333.115E	Broca Ø 1,5x Stop 19mmx Comprimento 135mm
05.5.0334.115E	Broca Ø 1,5x Stop 19mmx Comprimento 140mm
05.5.0335.115E	Broca Ø 1,5x Stop 19mmx Comprimento 145mm
05.5.0336.115E	Broca Ø 1,5x Stop 19mmx Comprimento 150mm
05.5.0337.115E	Broca Ø 1,5x Stop 20mmx Comprimento 50mm
05.5.0338.115E	Broca Ø 1,5x Stop 20mmx Comprimento 55mm
05.5.0339.115E	Broca Ø 1,5x Stop 20mmx Comprimento 60mm
05.5.0340.115E	Broca Ø 1,5x Stop 20mmx Comprimento 65mm
05.5.0341.115E	Broca Ø 1,5x Stop 20mmx Comprimento 70mm
05.5.0342.115E	Broca Ø 1,5x Stop 20mmx Comprimento 75mm
05.5.0343.115E	Broca Ø 1,5x Stop 20mmx Comprimento 80mm
05.5.0344.115E	Broca Ø 1,5x Stop 20mmx Comprimento 85mm
05.5.0345.115E	Broca Ø 1,5x Stop 20mmx Comprimento 90mm
05.5.0346.115E	Broca Ø 1,5x Stop 20mmx Comprimento 95mm
05.5.0347.115E	Broca Ø 1,5x Stop 20mmx Comprimento 100mm
05.5.0348.115E	Broca Ø 1,5x Stop 20mmx Comprimento 105mm
05.5.0349.115E	Broca Ø 1,5x Stop 20mmx Comprimento 110mm
05.5.0350.115E	Broca Ø 1,5x Stop 20mmx Comprimento 115mm
05.5.0351.115E	Broca Ø 1,5x Stop 20mmx Comprimento 120mm
05.5.0352.115E	Broca Ø 1,5x Stop 20mmx Comprimento 125mm
05.5.0353.115E	Broca Ø 1,5x Stop 20mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0354.115E	Broca Ø 1,5x Stop 20mmx Comprimento 135mm
05.5.0355.115E	Broca Ø 1,5x Stop 20mmx Comprimento 140mm
05.5.0356.115E	Broca Ø 1,5x Stop 20mmx Comprimento 145mm
05.5.0357.115E	Broca Ø 1,5x Stop 20mmx Comprimento 150mm
05.5.0358.115E	Broca Ø 1,5x Stop 21mmx Comprimento 50mm
05.5.0359.115E	Broca Ø 1,5x Stop 21mmx Comprimento 55mm
05.5.0360.115E	Broca Ø 1,5x Stop 21mmx Comprimento 60mm
05.5.0361.115E	Broca Ø 1,5x Stop 21mmx Comprimento 65mm
05.5.0362.115E	Broca Ø 1,5x Stop 21mmx Comprimento 70mm
05.5.0363.115E	Broca Ø 1,5x Stop 21mmx Comprimento 75mm
05.5.0364.115E	Broca Ø 1,5x Stop 21mmx Comprimento 80mm
05.5.0365.115E	Broca Ø 1,5x Stop 21mmx Comprimento 85mm
05.5.0366.115E	Broca Ø 1,5x Stop 21mmx Comprimento 90mm
05.5.0367.115E	Broca Ø 1,5x Stop 21mmx Comprimento 95mm
05.5.0368.115E	Broca Ø 1,5x Stop 21mmx Comprimento 100mm
05.5.0369.115E	Broca Ø 1,5x Stop 21mmx Comprimento 105mm
05.5.0370.115E	Broca Ø 1,5x Stop 21mmx Comprimento 110mm
05.5.0371.115E	Broca Ø 1,5x Stop 21mmx Comprimento 115mm
05.5.0372.115E	Broca Ø 1,5x Stop 21mmx Comprimento 120mm
05.5.0373.115E	Broca Ø 1,5x Stop 21mmx Comprimento 125mm
05.5.0374.115E	Broca Ø 1,5x Stop 21mmx Comprimento 130mm
05.5.0375.115E	Broca Ø 1,5x Stop 21mmx Comprimento 135mm
05.5.0376.115E	Broca Ø 1,5x Stop 21mmx Comprimento 140mm
05.5.0377.115E	Broca Ø 1,5x Stop 21mmx Comprimento 145mm
05.5.0378.115E	Broca Ø 1,5x Stop 21mmx Comprimento 150mm
05.5.0379.115E	Broca Ø 1,5x Stop 22mmx Comprimento 50mm
05.5.0380.115E	Broca Ø 1,5x Stop 22mmx Comprimento 55mm
05.5.0381.115E	Broca Ø 1,5x Stop 22mmx Comprimento 60mm
05.5.0382.115E	Broca Ø 1,5x Stop 22mmx Comprimento 65mm
05.5.0383.115E	Broca Ø 1,5x Stop 22mmx Comprimento 70mm
05.5.0384.115E	Broca Ø 1,5x Stop 22mmx Comprimento 75mm
05.5.0385.115E	Broca Ø 1,5x Stop 22mmx Comprimento 80mm
05.5.0386.115E	Broca Ø 1,5x Stop 22mmx Comprimento 85mm
05.5.0387.115E	Broca Ø 1,5x Stop 22mmx Comprimento 90mm
05.5.0388.115E	Broca Ø 1,5x Stop 22mmx Comprimento 95mm
05.5.0389.115E	Broca Ø 1,5x Stop 22mmx Comprimento 100mm
05.5.0390.115E	Broca Ø 1,5x Stop 22mmx Comprimento 105mm
05.5.0391.115E	Broca Ø 1,5x Stop 22mmx Comprimento 110mm
05.5.0392.115E	Broca Ø 1,5x Stop 22mmx Comprimento 115mm
05.5.0393.115E	Broca Ø 1,5x Stop 22mmx Comprimento 120mm
05.5.0394.115E	Broca Ø 1,5x Stop 22mmx Comprimento 125mm
05.5.0395.115E	Broca Ø 1,5x Stop 22mmx Comprimento 130mm
05.5.0396.115E	Broca Ø 1,5x Stop 22mmx Comprimento 135mm
05.5.0397.115E	Broca Ø 1,5x Stop 22mmx Comprimento 140mm
05.5.0398.115E	Broca Ø 1,5x Stop 22mmx Comprimento 145mm
05.5.0399.115E	Broca Ø 1,5x Stop 22mmx Comprimento 150mm
05.5.0400.115E	Broca Ø 1,5x Stop 23mmx Comprimento 50mm
05.5.0401.115E	Broca Ø 1,5x Stop 23mmx Comprimento 55mm
05.5.0402.115E	Broca Ø 1,5x Stop 23mmx Comprimento 60mm

05.5.0403.115E	Broca Ø 1,5x Stop 23mmx Comprimento 65mm
05.5.0404.115E	Broca Ø 1,5x Stop 23mmx Comprimento 70mm
05.5.0405.115E	Broca Ø 1,5x Stop 23mmx Comprimento 75mm
05.5.0406.115E	Broca Ø 1,5x Stop 23mmx Comprimento 80mm
05.5.0407.115E	Broca Ø 1,5x Stop 23mmx Comprimento 85mm
05.5.0408.115E	Broca Ø 1,5x Stop 23mmx Comprimento 90mm
05.5.0409.115E	Broca Ø 1,5x Stop 23mmx Comprimento 95mm
05.5.0410.115E	Broca Ø 1,5x Stop 23mmx Comprimento 100mm
05.5.0411.115E	Broca Ø 1,5x Stop 23mmx Comprimento 105mm
05.5.0412.115E	Broca Ø 1,5x Stop 23mmx Comprimento 110mm
05.5.0413.115E	Broca Ø 1,5x Stop 23mmx Comprimento 115mm
05.5.0414.115E	Broca Ø 1,5x Stop 23mmx Comprimento 120mm
05.5.0415.115E	Broca Ø 1,5x Stop 23mmx Comprimento 125mm
05.5.0416.115E	Broca Ø 1,5x Stop 23mmx Comprimento 130mm
05.5.0417.115E	Broca Ø 1,5x Stop 23mmx Comprimento 135mm
05.5.0418.115E	Broca Ø 1,5x Stop 23mmx Comprimento 140mm
05.5.0419.115E	Broca Ø 1,5x Stop 23mmx Comprimento 145mm
05.5.0420.115E	Broca Ø 1,5x Stop 23mmx Comprimento 150mm
05.5.0421.115E	Broca Ø 1,5x Stop 24mmx Comprimento 50mm
05.5.0422.115E	Broca Ø 1,5x Stop 24mmx Comprimento 55mm
05.5.0423.115E	Broca Ø 1,5x Stop 24mmx Comprimento 60mm
05.5.0424.115E	Broca Ø 1,5x Stop 24mmx Comprimento 65mm
05.5.0425.115E	Broca Ø 1,5x Stop 24mmx Comprimento 70mm
05.5.0426.115E	Broca Ø 1,5x Stop 24mmx Comprimento 75mm
05.5.0427.115E	Broca Ø 1,5x Stop 24mmx Comprimento 80mm
05.5.0428.115E	Broca Ø 1,5x Stop 24mmx Comprimento 85mm
05.5.0429.115E	Broca Ø 1,5x Stop 24mmx Comprimento 90mm
05.5.0430.115E	Broca Ø 1,5x Stop 24mmx Comprimento 95mm
05.5.0431.115E	Broca Ø 1,5x Stop 24mmx Comprimento 100mm
05.5.0432.115E	Broca Ø 1,5x Stop 24mmx Comprimento 105mm
05.5.0433.115E	Broca Ø 1,5x Stop 24mmx Comprimento 110mm
05.5.0434.115E	Broca Ø 1,5x Stop 24mmx Comprimento 115mm
05.5.0435.115E	Broca Ø 1,5x Stop 24mmx Comprimento 120mm
05.5.0436.115E	Broca Ø 1,5x Stop 24mmx Comprimento 125mm
05.5.0437.115E	Broca Ø 1,5x Stop 24mmx Comprimento 130mm
05.5.0438.115E	Broca Ø 1,5x Stop 24mmx Comprimento 135mm
05.5.0439.115E	Broca Ø 1,5x Stop 24mmx Comprimento 140mm
05.5.0440.115E	Broca Ø 1,5x Stop 24mmx Comprimento 145mm
05.5.0441.115E	Broca Ø 1,5x Stop 24mmx Comprimento 150mm
05.5.0442.115E	Broca Ø 1,5x Stop 25mmx Comprimento 50mm
05.5.0443.115E	Broca Ø 1,5x Stop 25mmx Comprimento 55mm
05.5.0444.115E	Broca Ø 1,5x Stop 25mmx Comprimento 60mm
05.5.0445.115E	Broca Ø 1,5x Stop 25mmx Comprimento 65mm
05.5.0446.115E	Broca Ø 1,5x Stop 25mmx Comprimento 70mm
05.5.0447.115E	Broca Ø 1,5x Stop 25mmx Comprimento 75mm
05.5.0448.115E	Broca Ø 1,5x Stop 25mmx Comprimento 80mm
05.5.0449.115E	Broca Ø 1,5x Stop 25mmx Comprimento 85mm
05.5.0450.115E	Broca Ø 1,5x Stop 25mmx Comprimento 90mm
05.5.0451.115E	Broca Ø 1,5x Stop 25mmx Comprimento 95mm

05.5.0452.115E	Broca Ø 1,5x Stop 25mmx Comprimento 100mm
05.5.0453.115E	Broca Ø 1,5x Stop 25mmx Comprimento 105mm
05.5.0454.115E	Broca Ø 1,5x Stop 25mmx Comprimento 110mm
05.5.0455.115E	Broca Ø 1,5x Stop 25mmx Comprimento 115mm
05.5.0456.115E	Broca Ø 1,5x Stop 25mmx Comprimento 120mm
05.5.0457.115E	Broca Ø 1,5x Stop 25mmx Comprimento 125mm
05.5.0458.115E	Broca Ø 1,5x Stop 25mmx Comprimento 130mm
05.5.0459.115E	Broca Ø 1,5x Stop 25mmx Comprimento 135mm
05.5.0460.115E	Broca Ø 1,5x Stop 25mmx Comprimento 140mm
05.5.0461.115E	Broca Ø 1,5x Stop 25mmx Comprimento 145mm
05.5.0462.115E	Broca Ø 1,5x Stop 25mmx Comprimento 150mm
05.5.0463.115E	Broca Ø 1,5x Stop 30mmx Comprimento 50mm
05.5.0464.115E	Broca Ø 1,5x Stop 30mmx Comprimento 55mm
05.5.0465.115E	Broca Ø 1,5x Stop 30mmx Comprimento 60mm
05.5.0466.115E	Broca Ø 1,5x Stop 30mmx Comprimento 65mm
05.5.0467.115E	Broca Ø 1,5x Stop 30mmx Comprimento 70mm
05.5.0468.115E	Broca Ø 1,5x Stop 30mmx Comprimento 75mm
05.5.0469.115E	Broca Ø 1,5x Stop 30mmx Comprimento 80mm
05.5.0470.115E	Broca Ø 1,5x Stop 30mmx Comprimento 85mm
05.5.0471.115E	Broca Ø 1,5x Stop 30mmx Comprimento 90mm
05.5.0472.115E	Broca Ø 1,5x Stop 30mmx Comprimento 95mm
05.5.0473.115E	Broca Ø 1,5x Stop 30mmx Comprimento 100mm
05.5.0474.115E	Broca Ø 1,5x Stop 30mmx Comprimento 105mm
05.5.0475.115E	Broca Ø 1,5x Stop 30mmx Comprimento 110mm
05.5.0476.115E	Broca Ø 1,5x Stop 30mmx Comprimento 115mm
05.5.0477.115E	Broca Ø 1,5x Stop 30mmx Comprimento 120mm
05.5.0478.115E	Broca Ø 1,5x Stop 30mmx Comprimento 125mm
05.5.0479.115E	Broca Ø 1,5x Stop 30mmx Comprimento 130mm
05.5.0480.115E	Broca Ø 1,5x Stop 30mmx Comprimento 135mm
05.5.0481.115E	Broca Ø 1,5x Stop 30mmx Comprimento 140mm
05.5.0482.115E	Broca Ø 1,5x Stop 30mmx Comprimento 145mm
05.5.0483.115E	Broca Ø 1,5x Stop 30mmx Comprimento 150mm
05.5.0484.115E	Broca Ø 1,5x Stop 35mmx Comprimento 50mm
05.5.0485.115E	Broca Ø 1,5x Stop 35mmx Comprimento 55mm
05.5.0486.115E	Broca Ø 1,5x Stop 35mmx Comprimento 60mm
05.5.0487.115E	Broca Ø 1,5x Stop 35mmx Comprimento 65mm
05.5.0488.115E	Broca Ø 1,5x Stop 35mmx Comprimento 70mm
05.5.0489.115E	Broca Ø 1,5x Stop 35mmx Comprimento 75mm
05.5.0490.115E	Broca Ø 1,5x Stop 35mmx Comprimento 80mm
05.5.0491.115E	Broca Ø 1,5x Stop 35mmx Comprimento 85mm
05.5.0492.115E	Broca Ø 1,5x Stop 35mmx Comprimento 90mm
05.5.0493.115E	Broca Ø 1,5x Stop 35mmx Comprimento 95mm
05.5.0494.115E	Broca Ø 1,5x Stop 35mmx Comprimento 100mm
05.5.0495.115E	Broca Ø 1,5x Stop 35mmx Comprimento 105mm
05.5.0496.115E	Broca Ø 1,5x Stop 35mmx Comprimento 110mm
05.5.0497.115E	Broca Ø 1,5x Stop 35mmx Comprimento 115mm
05.5.0498.115E	Broca Ø 1,5x Stop 35mmx Comprimento 120mm
05.5.0499.115E	Broca Ø 1,5x Stop 35mmx Comprimento 125mm
05.5.0500.115E	Broca Ø 1,5x Stop 35mmx Comprimento 130mm

05.5.0501.115E	Broca Ø 1,5x Stop 35mmx Comprimento 135mm
05.5.0502.115E	Broca Ø 1,5x Stop 35mmx Comprimento 140mm
05.5.0503.115E	Broca Ø 1,5x Stop 35mmx Comprimento 145mm
05.5.0504.115E	Broca Ø 1,5x Stop 35mmx Comprimento 150mm
05.5.0001.116E	Broca Ø 1,6x Stop 4mmx Comprimento 50mm
05.5.0002.116E	Broca Ø 1,6x Stop 4mmx Comprimento 55mm
05.5.0003.116E	Broca Ø 1,6x Stop 4mmx Comprimento 60mm
05.5.0004.116E	Broca Ø 1,6x Stop 4mmx Comprimento 65mm
05.5.0005.116E	Broca Ø 1,6x Stop 4mmx Comprimento 70mm
05.5.0006.116E	Broca Ø 1,6x Stop 4mmx Comprimento 75mm
05.5.0007.116E	Broca Ø 1,6x Stop 4mmx Comprimento 80mm
05.5.0008.116E	Broca Ø 1,6x Stop 4mmx Comprimento 85mm
05.5.0009.116E	Broca Ø 1,6x Stop 4mmx Comprimento 90mm
05.5.0010.116E	Broca Ø 1,6x Stop 4mmx Comprimento 95mm
05.5.0011.116E	Broca Ø 1,6x Stop 4mmx Comprimento 100mm
05.5.0012.116E	Broca Ø 1,6x Stop 4mmx Comprimento 105mm
05.5.0013.116E	Broca Ø 1,6x Stop 4mmx Comprimento 110mm
05.5.0014.116E	Broca Ø 1,6x Stop 4mmx Comprimento 115mm
05.5.0015.116E	Broca Ø 1,6x Stop 4mmx Comprimento 120mm
05.5.0016.116E	Broca Ø 1,6x Stop 4mmx Comprimento 125mm
05.5.0017.116E	Broca Ø 1,6x Stop 4mmx Comprimento 130mm
05.5.0018.116E	Broca Ø 1,6x Stop 4mmx Comprimento 135mm
05.5.0019.116E	Broca Ø 1,6x Stop 4mmx Comprimento 140mm
05.5.0020.116E	Broca Ø 1,6x Stop 4mmx Comprimento 145mm
05.5.0021.116E	Broca Ø 1,6x Stop 4mmx Comprimento 150mm
05.5.0022.116E	Broca Ø 1,6x Stop 5mmx Comprimento 50mm
05.5.0023.116E	Broca Ø 1,6x Stop 5mmx Comprimento 55mm
05.5.0024.116E	Broca Ø 1,6x Stop 5mmx Comprimento 60mm
05.5.0025.116E	Broca Ø 1,6x Stop 5mmx Comprimento 65mm
05.5.0026.116E	Broca Ø 1,6x Stop 5mmx Comprimento 70mm
05.5.0027.116E	Broca Ø 1,6x Stop 5mmx Comprimento 75mm
05.5.0028.116E	Broca Ø 1,6x Stop 5mmx Comprimento 80mm
05.5.0029.116E	Broca Ø 1,6x Stop 5mmx Comprimento 85mm
05.5.0030.116E	Broca Ø 1,6x Stop 5mmx Comprimento 90mm
05.5.0031.116E	Broca Ø 1,6x Stop 5mmx Comprimento 95mm
05.5.0032.116E	Broca Ø 1,6x Stop 5mmx Comprimento 100mm
05.5.0033.116E	Broca Ø 1,6x Stop 5mmx Comprimento 105mm
05.5.0034.116E	Broca Ø 1,6x Stop 5mmx Comprimento 110mm
05.5.0035.116E	Broca Ø 1,6x Stop 5mmx Comprimento 115mm
05.5.0036.116E	Broca Ø 1,6x Stop 5mmx Comprimento 120mm
05.5.0037.116E	Broca Ø 1,6x Stop 5mmx Comprimento 125mm
05.5.0038.116E	Broca Ø 1,6x Stop 5mmx Comprimento 130mm
05.5.0039.116E	Broca Ø 1,6x Stop 5mmx Comprimento 135mm
05.5.0040.116E	Broca Ø 1,6x Stop 5mmx Comprimento 140mm
05.5.0041.116E	Broca Ø 1,6x Stop 5mmx Comprimento 145mm
05.5.0042.116E	Broca Ø 1,6x Stop 5mmx Comprimento 150mm
05.5.0043.116E	Broca Ø 1,6x Stop 6mmx Comprimento 50mm
05.5.0044.116E	Broca Ø 1,6x Stop 6mmx Comprimento 55mm
05.5.0045.116E	Broca Ø 1,6x Stop 6mmx Comprimento 60mm

05.5.0046.116E	Broca Ø 1,6x Stop 6mmx Comprimento 65mm
05.5.0047.116E	Broca Ø 1,6x Stop 6mmx Comprimento 70mm
05.5.0048.116E	Broca Ø 1,6x Stop 6mmx Comprimento 75mm
05.5.0049.116E	Broca Ø 1,6x Stop 6mmx Comprimento 80mm
05.5.0050.116E	Broca Ø 1,6x Stop 6mmx Comprimento 85mm
05.5.0051.116E	Broca Ø 1,6x Stop 6mmx Comprimento 90mm
05.5.0052.116E	Broca Ø 1,6x Stop 6mmx Comprimento 95mm
05.5.0053.116E	Broca Ø 1,6x Stop 6mmx Comprimento 100mm
05.5.0054.116E	Broca Ø 1,6x Stop 6mmx Comprimento 105mm
05.5.0055.116E	Broca Ø 1,6x Stop 6mmx Comprimento 110mm
05.5.0056.116E	Broca Ø 1,6x Stop 6mmx Comprimento 115mm
05.5.0057.116E	Broca Ø 1,6x Stop 6mmx Comprimento 120mm
05.5.0058.116E	Broca Ø 1,6x Stop 6mmx Comprimento 125mm
05.5.0059.116E	Broca Ø 1,6x Stop 6mmx Comprimento 130mm
05.5.0060.116E	Broca Ø 1,6x Stop 6mmx Comprimento 135mm
05.5.0061.116E	Broca Ø 1,6x Stop 6mmx Comprimento 140mm
05.5.0062.116E	Broca Ø 1,6x Stop 6mmx Comprimento 145mm
05.5.0063.116E	Broca Ø 1,6x Stop 6mmx Comprimento 150mm
05.5.0064.116E	Broca Ø 1,6x Stop 7mmx Comprimento 50mm
05.5.0065.116E	Broca Ø 1,6x Stop 7mmx Comprimento 55mm
05.5.0066.116E	Broca Ø 1,6x Stop 7mmx Comprimento 60mm
05.5.0067.116E	Broca Ø 1,6x Stop 7mmx Comprimento 65mm
05.5.0068.116E	Broca Ø 1,6x Stop 7mmx Comprimento 70mm
05.5.0069.116E	Broca Ø 1,6x Stop 7mmx Comprimento 75mm
05.5.0070.116E	Broca Ø 1,6x Stop 7mmx Comprimento 80mm
05.5.0071.116E	Broca Ø 1,6x Stop 7mmx Comprimento 85mm
05.5.0072.116E	Broca Ø 1,6x Stop 7mmx Comprimento 90mm
05.5.0073.116E	Broca Ø 1,6x Stop 7mmx Comprimento 95mm
05.5.0074.116E	Broca Ø 1,6x Stop 7mmx Comprimento 100mm
05.5.0075.116E	Broca Ø 1,6x Stop 7mmx Comprimento 105mm
05.5.0076.116E	Broca Ø 1,6x Stop 7mmx Comprimento 110mm
05.5.0077.116E	Broca Ø 1,6x Stop 7mmx Comprimento 115mm
05.5.0078.116E	Broca Ø 1,6x Stop 7mmx Comprimento 120mm
05.5.0079.116E	Broca Ø 1,6x Stop 7mmx Comprimento 125mm
05.5.0080.116E	Broca Ø 1,6x Stop 7mmx Comprimento 130mm
05.5.0081.116E	Broca Ø 1,6x Stop 7mmx Comprimento 135mm
05.5.0082.116E	Broca Ø 1,6x Stop 7mmx Comprimento 140mm
05.5.0083.116E	Broca Ø 1,6x Stop 7mmx Comprimento 145mm
05.5.0084.116E	Broca Ø 1,6x Stop 7mmx Comprimento 150mm
05.5.0085.116E	Broca Ø 1,6x Stop 8mmx Comprimento 50mm
05.5.0086.116E	Broca Ø 1,6x Stop 8mmx Comprimento 55mm
05.5.0087.116E	Broca Ø 1,6x Stop 8mmx Comprimento 60mm
05.5.0088.116E	Broca Ø 1,6x Stop 8mmx Comprimento 65mm
05.5.0089.116E	Broca Ø 1,6x Stop 8mmx Comprimento 70mm
05.5.0090.116E	Broca Ø 1,6x Stop 8mmx Comprimento 75mm
05.5.0091.116E	Broca Ø 1,6x Stop 8mmx Comprimento 80mm
05.5.0092.116E	Broca Ø 1,6x Stop 8mmx Comprimento 85mm
05.5.0093.116E	Broca Ø 1,6x Stop 8mmx Comprimento 90mm
05.5.0094.116E	Broca Ø 1,6x Stop 8mmx Comprimento 95mm

05.5.0095.116E	Broca Ø 1,6x Stop 8mmx Comprimento 100mm
05.5.0096.116E	Broca Ø 1,6x Stop 8mmx Comprimento 105mm
05.5.0097.116E	Broca Ø 1,6x Stop 8mmx Comprimento 110mm
05.5.0098.116E	Broca Ø 1,6x Stop 8mmx Comprimento 115mm
05.5.0099.116E	Broca Ø 1,6x Stop 8mmx Comprimento 120mm
05.5.0100.116E	Broca Ø 1,6x Stop 8mmx Comprimento 125mm
05.5.0101.116E	Broca Ø 1,6x Stop 8mmx Comprimento 130mm
05.5.0102.116E	Broca Ø 1,6x Stop 8mmx Comprimento 135mm
05.5.0103.116E	Broca Ø 1,6x Stop 8mmx Comprimento 140mm
05.5.0104.116E	Broca Ø 1,6x Stop 8mmx Comprimento 145mm
05.5.0105.116E	Broca Ø 1,6x Stop 8mmx Comprimento 150mm
05.5.0106.116E	Broca Ø 1,6x Stop 9mmx Comprimento 50mm
05.5.0107.116E	Broca Ø 1,6x Stop 9mmx Comprimento 55mm
05.5.0108.116E	Broca Ø 1,6x Stop 9mmx Comprimento 60mm
05.5.0109.116E	Broca Ø 1,6x Stop 9mmx Comprimento 65mm
05.5.0110.116E	Broca Ø 1,6x Stop 9mmx Comprimento 70mm
05.5.0111.116E	Broca Ø 1,6x Stop 9mmx Comprimento 75mm
05.5.0112.116E	Broca Ø 1,6x Stop 9mmx Comprimento 80mm
05.5.0113.116E	Broca Ø 1,6x Stop 9mmx Comprimento 85mm
05.5.0114.116E	Broca Ø 1,6x Stop 9mmx Comprimento 90mm
05.5.0115.116E	Broca Ø 1,6x Stop 9mmx Comprimento 95mm
05.5.0116.116E	Broca Ø 1,6x Stop 9mmx Comprimento 100mm
05.5.0117.116E	Broca Ø 1,6x Stop 9mmx Comprimento 105mm
05.5.0118.116E	Broca Ø 1,6x Stop 9mmx Comprimento 110mm
05.5.0119.116E	Broca Ø 1,6x Stop 9mmx Comprimento 115mm
05.5.0120.116E	Broca Ø 1,6x Stop 9mmx Comprimento 120mm
05.5.0121.116E	Broca Ø 1,6x Stop 9mmx Comprimento 125mm
05.5.0122.116E	Broca Ø 1,6x Stop 9mmx Comprimento 130mm
05.5.0123.116E	Broca Ø 1,6x Stop 9mmx Comprimento 135mm
05.5.0124.116E	Broca Ø 1,6x Stop 9mmx Comprimento 140mm
05.5.0125.116E	Broca Ø 1,6x Stop 9mmx Comprimento 145mm
05.5.0126.116E	Broca Ø 1,6x Stop 9mmx Comprimento 150mm
05.5.0127.116E	Broca Ø 1,6x Stop 10mmx Comprimento 50mm
05.5.0128.116E	Broca Ø 1,6x Stop 10mmx Comprimento 55mm
05.5.0129.116E	Broca Ø 1,6x Stop 10mmx Comprimento 60mm
05.5.0130.116E	Broca Ø 1,6x Stop 10mmx Comprimento 65mm
05.5.0131.116E	Broca Ø 1,6x Stop 10mmx Comprimento 70mm
05.5.0132.116E	Broca Ø 1,6x Stop 10mmx Comprimento 75mm
05.5.0133.116E	Broca Ø 1,6x Stop 10mmx Comprimento 80mm
05.5.0134.116E	Broca Ø 1,6x Stop 10mmx Comprimento 85mm
05.5.0135.116E	Broca Ø 1,6x Stop 10mmx Comprimento 90mm
05.5.0136.116E	Broca Ø 1,6x Stop 10mmx Comprimento 95mm
05.5.0137.116E	Broca Ø 1,6x Stop 10mmx Comprimento 100mm
05.5.0138.116E	Broca Ø 1,6x Stop 10mmx Comprimento 105mm
05.5.0139.116E	Broca Ø 1,6x Stop 10mmx Comprimento 110mm
05.5.0140.116E	Broca Ø 1,6x Stop 10mmx Comprimento 115mm
05.5.0141.116E	Broca Ø 1,6x Stop 10mmx Comprimento 120mm
05.5.0142.116E	Broca Ø 1,6x Stop 10mmx Comprimento 125mm
05.5.0143.116E	Broca Ø 1,6x Stop 10mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0144.116E	Broca Ø 1,6x Stop 10mmx Comprimento 135mm
05.5.0145.116E	Broca Ø 1,6x Stop 10mmx Comprimento 140mm
05.5.0146.116E	Broca Ø 1,6x Stop 10mmx Comprimento 145mm
05.5.0147.116E	Broca Ø 1,6x Stop 10mmx Comprimento 150mm
05.5.0148.116E	Broca Ø 1,6x Stop 11mmx Comprimento 50mm
05.5.0149.116E	Broca Ø 1,6x Stop 11mmx Comprimento 55mm
05.5.0150.116E	Broca Ø 1,6x Stop 11mmx Comprimento 60mm
05.5.0151.116E	Broca Ø 1,6x Stop 11mmx Comprimento 65mm
05.5.0152.116E	Broca Ø 1,6x Stop 11mmx Comprimento 70mm
05.5.0153.116E	Broca Ø 1,6x Stop 11mmx Comprimento 75mm
05.5.0154.116E	Broca Ø 1,6x Stop 11mmx Comprimento 80mm
05.5.0155.116E	Broca Ø 1,6x Stop 11mmx Comprimento 85mm
05.5.0156.116E	Broca Ø 1,6x Stop 11mmx Comprimento 90mm
05.5.0157.116E	Broca Ø 1,6x Stop 11mmx Comprimento 95mm
05.5.0158.116E	Broca Ø 1,6x Stop 11mmx Comprimento 100mm
05.5.0159.116E	Broca Ø 1,6x Stop 11mmx Comprimento 105mm
05.5.0160.116E	Broca Ø 1,6x Stop 11mmx Comprimento 110mm
05.5.0161.116E	Broca Ø 1,6x Stop 11mmx Comprimento 115mm
05.5.0162.116E	Broca Ø 1,6x Stop 11mmx Comprimento 120mm
05.5.0163.116E	Broca Ø 1,6x Stop 11mmx Comprimento 125mm
05.5.0164.116E	Broca Ø 1,6x Stop 11mmx Comprimento 130mm
05.5.0165.116E	Broca Ø 1,6x Stop 11mmx Comprimento 135mm
05.5.0166.116E	Broca Ø 1,6x Stop 11mmx Comprimento 140mm
05.5.0167.116E	Broca Ø 1,6x Stop 11mmx Comprimento 145mm
05.5.0168.116E	Broca Ø 1,6x Stop 11mmx Comprimento 150mm
05.5.0169.116E	Broca Ø 1,6x Stop 12mmx Comprimento 50mm
05.5.0170.116E	Broca Ø 1,6x Stop 12mmx Comprimento 55mm
05.5.0171.116E	Broca Ø 1,6x Stop 12mmx Comprimento 60mm
05.5.0172.116E	Broca Ø 1,6x Stop 12mmx Comprimento 65mm
05.5.0173.116E	Broca Ø 1,6x Stop 12mmx Comprimento 70mm
05.5.0174.116E	Broca Ø 1,6x Stop 12mmx Comprimento 75mm
05.5.0175.116E	Broca Ø 1,6x Stop 12mmx Comprimento 80mm
05.5.0176.116E	Broca Ø 1,6x Stop 12mmx Comprimento 85mm
05.5.0177.116E	Broca Ø 1,6x Stop 12mmx Comprimento 90mm
05.5.0178.116E	Broca Ø 1,6x Stop 12mmx Comprimento 95mm
05.5.0179.116E	Broca Ø 1,6x Stop 12mmx Comprimento 100mm
05.5.0180.116E	Broca Ø 1,6x Stop 12mmx Comprimento 105mm
05.5.0181.116E	Broca Ø 1,6x Stop 12mmx Comprimento 110mm
05.5.0182.116E	Broca Ø 1,6x Stop 12mmx Comprimento 115mm
05.5.0183.116E	Broca Ø 1,6x Stop 12mmx Comprimento 120mm
05.5.0184.116E	Broca Ø 1,6x Stop 12mmx Comprimento 125mm
05.5.0185.116E	Broca Ø 1,6x Stop 12mmx Comprimento 130mm
05.5.0186.116E	Broca Ø 1,6x Stop 12mmx Comprimento 135mm
05.5.0187.116E	Broca Ø 1,6x Stop 12mmx Comprimento 140mm
05.5.0188.116E	Broca Ø 1,6x Stop 12mmx Comprimento 145mm
05.5.0189.116E	Broca Ø 1,6x Stop 12mmx Comprimento 150mm
05.5.0190.116E	Broca Ø 1,6x Stop 13mmx Comprimento 50mm
05.5.0191.116E	Broca Ø 1,6x Stop 13mmx Comprimento 55mm
05.5.0192.116E	Broca Ø 1,6x Stop 13mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0193.116E	Broca Ø 1,6x Stop 13mmx Comprimento 65mm
05.5.0194.116E	Broca Ø 1,6x Stop 13mmx Comprimento 70mm
05.5.0195.116E	Broca Ø 1,6x Stop 13mmx Comprimento 75mm
05.5.0196.116E	Broca Ø 1,6x Stop 13mmx Comprimento 80mm
05.5.0197.116E	Broca Ø 1,6x Stop 13mmx Comprimento 85mm
05.5.0198.116E	Broca Ø 1,6x Stop 13mmx Comprimento 90mm
05.5.0199.116E	Broca Ø 1,6x Stop 13mmx Comprimento 95mm
05.5.0200.116E	Broca Ø 1,6x Stop 13mmx Comprimento 100mm
05.5.0201.116E	Broca Ø 1,6x Stop 13mmx Comprimento 105mm
05.5.0202.116E	Broca Ø 1,6x Stop 13mmx Comprimento 110mm
05.5.0203.116E	Broca Ø 1,6x Stop 13mmx Comprimento 115mm
05.5.0204.116E	Broca Ø 1,6x Stop 13mmx Comprimento 120mm
05.5.0205.116E	Broca Ø 1,6x Stop 13mmx Comprimento 125mm
05.5.0206.116E	Broca Ø 1,6x Stop 13mmx Comprimento 130mm
05.5.0207.116E	Broca Ø 1,6x Stop 13mmx Comprimento 135mm
05.5.0208.116E	Broca Ø 1,6x Stop 13mmx Comprimento 140mm
05.5.0209.116E	Broca Ø 1,6x Stop 13mmx Comprimento 145mm
05.5.0210.116E	Broca Ø 1,6x Stop 13mmx Comprimento 150mm
05.5.0211.116E	Broca Ø 1,6x Stop 14mmx Comprimento 50mm
05.5.0212.116E	Broca Ø 1,6x Stop 14mmx Comprimento 55mm
05.5.0213.116E	Broca Ø 1,6x Stop 14mmx Comprimento 60mm
05.5.0214.116E	Broca Ø 1,6x Stop 14mmx Comprimento 65mm
05.5.0215.116E	Broca Ø 1,6x Stop 14mmx Comprimento 70mm
05.5.0216.116E	Broca Ø 1,6x Stop 14mmx Comprimento 75mm
05.5.0217.116E	Broca Ø 1,6x Stop 14mmx Comprimento 80mm
05.5.0218.116E	Broca Ø 1,6x Stop 14mmx Comprimento 85mm
05.5.0219.116E	Broca Ø 1,6x Stop 14mmx Comprimento 90mm
05.5.0220.116E	Broca Ø 1,6x Stop 14mmx Comprimento 95mm
05.5.0221.116E	Broca Ø 1,6x Stop 14mmx Comprimento 100mm
05.5.0222.116E	Broca Ø 1,6x Stop 14mmx Comprimento 105mm
05.5.0223.116E	Broca Ø 1,6x Stop 14mmx Comprimento 110mm
05.5.0224.116E	Broca Ø 1,6x Stop 14mmx Comprimento 115mm
05.5.0225.116E	Broca Ø 1,6x Stop 14mmx Comprimento 120mm
05.5.0226.116E	Broca Ø 1,6x Stop 14mmx Comprimento 125mm
05.5.0227.116E	Broca Ø 1,6x Stop 14mmx Comprimento 130mm
05.5.0228.116E	Broca Ø 1,6x Stop 14mmx Comprimento 135mm
05.5.0229.116E	Broca Ø 1,6x Stop 14mmx Comprimento 140mm
05.5.0230.116E	Broca Ø 1,6x Stop 14mmx Comprimento 145mm
05.5.0231.116E	Broca Ø 1,6x Stop 14mmx Comprimento 150mm
05.5.0232.116E	Broca Ø 1,6x Stop 15mmx Comprimento 50mm
05.5.0233.116E	Broca Ø 1,6x Stop 15mmx Comprimento 55mm
05.5.0234.116E	Broca Ø 1,6x Stop 15mmx Comprimento 60mm
05.5.0235.116E	Broca Ø 1,6x Stop 15mmx Comprimento 65mm
05.5.0236.116E	Broca Ø 1,6x Stop 15mmx Comprimento 70mm
05.5.0237.116E	Broca Ø 1,6x Stop 15mmx Comprimento 75mm
05.5.0238.116E	Broca Ø 1,6x Stop 15mmx Comprimento 80mm
05.5.0239.116E	Broca Ø 1,6x Stop 15mmx Comprimento 85mm
05.5.0240.116E	Broca Ø 1,6x Stop 15mmx Comprimento 90mm
05.5.0241.116E	Broca Ø 1,6x Stop 15mmx Comprimento 95mm

05.5.0242.116E	Broca Ø 1,6x Stop 15mmx Comprimento 100mm
05.5.0243.116E	Broca Ø 1,6x Stop 15mmx Comprimento 105mm
05.5.0244.116E	Broca Ø 1,6x Stop 15mmx Comprimento 110mm
05.5.0245.116E	Broca Ø 1,6x Stop 15mmx Comprimento 115mm
05.5.0246.116E	Broca Ø 1,6x Stop 15mmx Comprimento 120mm
05.5.0247.116E	Broca Ø 1,6x Stop 15mmx Comprimento 125mm
05.5.0248.116E	Broca Ø 1,6x Stop 15mmx Comprimento 130mm
05.5.0249.116E	Broca Ø 1,6x Stop 15mmx Comprimento 135mm
05.5.0250.116E	Broca Ø 1,6x Stop 15mmx Comprimento 140mm
05.5.0251.116E	Broca Ø 1,6x Stop 15mmx Comprimento 145mm
05.5.0252.116E	Broca Ø 1,6x Stop 15mmx Comprimento 150mm
05.5.0253.116E	Broca Ø 1,6x Stop 16mmx Comprimento 50mm
05.5.0254.116E	Broca Ø 1,6x Stop 16mmx Comprimento 55mm
05.5.0255.116E	Broca Ø 1,6x Stop 16mmx Comprimento 60mm
05.5.0256.116E	Broca Ø 1,6x Stop 16mmx Comprimento 65mm
05.5.0257.116E	Broca Ø 1,6x Stop 16mmx Comprimento 70mm
05.5.0258.116E	Broca Ø 1,6x Stop 16mmx Comprimento 75mm
05.5.0259.116E	Broca Ø 1,6x Stop 16mmx Comprimento 80mm
05.5.0260.116E	Broca Ø 1,6x Stop 16mmx Comprimento 85mm
05.5.0261.116E	Broca Ø 1,6x Stop 16mmx Comprimento 90mm
05.5.0262.116E	Broca Ø 1,6x Stop 16mmx Comprimento 95mm
05.5.0263.116E	Broca Ø 1,6x Stop 16mmx Comprimento 100mm
05.5.0264.116E	Broca Ø 1,6x Stop 16mmx Comprimento 105mm
05.5.0265.116E	Broca Ø 1,6x Stop 16mmx Comprimento 110mm
05.5.0266.116E	Broca Ø 1,6x Stop 16mmx Comprimento 115mm
05.5.0267.116E	Broca Ø 1,6x Stop 16mmx Comprimento 120mm
05.5.0268.116E	Broca Ø 1,6x Stop 16mmx Comprimento 125mm
05.5.0269.116E	Broca Ø 1,6x Stop 16mmx Comprimento 130mm
05.5.0270.116E	Broca Ø 1,6x Stop 16mmx Comprimento 135mm
05.5.0271.116E	Broca Ø 1,6x Stop 16mmx Comprimento 140mm
05.5.0272.116E	Broca Ø 1,6x Stop 16mmx Comprimento 145mm
05.5.0273.116E	Broca Ø 1,6x Stop 16mmx Comprimento 150mm
05.5.0274.116E	Broca Ø 1,6x Stop 17mmx Comprimento 50mm
05.5.0275.116E	Broca Ø 1,6x Stop 17mmx Comprimento 55mm
05.5.0276.116E	Broca Ø 1,6x Stop 17mmx Comprimento 60mm
05.5.0277.116E	Broca Ø 1,6x Stop 17mmx Comprimento 65mm
05.5.0278.116E	Broca Ø 1,6x Stop 17mmx Comprimento 70mm
05.5.0279.116E	Broca Ø 1,6x Stop 17mmx Comprimento 75mm
05.5.0280.116E	Broca Ø 1,6x Stop 17mmx Comprimento 80mm
05.5.0281.116E	Broca Ø 1,6x Stop 17mmx Comprimento 85mm
05.5.0282.116E	Broca Ø 1,6x Stop 17mmx Comprimento 90mm
05.5.0283.116E	Broca Ø 1,6x Stop 17mmx Comprimento 95mm
05.5.0284.116E	Broca Ø 1,6x Stop 17mmx Comprimento 100mm
05.5.0285.116E	Broca Ø 1,6x Stop 17mmx Comprimento 105mm
05.5.0286.116E	Broca Ø 1,6x Stop 17mmx Comprimento 110mm
05.5.0287.116E	Broca Ø 1,6x Stop 17mmx Comprimento 115mm
05.5.0288.116E	Broca Ø 1,6x Stop 17mmx Comprimento 120mm
05.5.0289.116E	Broca Ø 1,6x Stop 17mmx Comprimento 125mm
05.5.0290.116E	Broca Ø 1,6x Stop 17mmx Comprimento 130mm

05.5.0291.116E	Broca Ø 1,6x Stop 17mmx Comprimento 135mm
05.5.0292.116E	Broca Ø 1,6x Stop 17mmx Comprimento 140mm
05.5.0293.116E	Broca Ø 1,6x Stop 17mmx Comprimento 145mm
05.5.0294.116E	Broca Ø 1,6x Stop 17mmx Comprimento 150mm
05.5.0295.116E	Broca Ø 1,6x Stop 18mmx Comprimento 50mm
05.5.0296.116E	Broca Ø 1,6x Stop 18mmx Comprimento 55mm
05.5.0297.116E	Broca Ø 1,6x Stop 18mmx Comprimento 60mm
05.5.0298.116E	Broca Ø 1,6x Stop 18mmx Comprimento 65mm
05.5.0299.116E	Broca Ø 1,6x Stop 18mmx Comprimento 70mm
05.5.0300.116E	Broca Ø 1,6x Stop 18mmx Comprimento 75mm
05.5.0301.116E	Broca Ø 1,6x Stop 18mmx Comprimento 80mm
05.5.0302.116E	Broca Ø 1,6x Stop 18mmx Comprimento 85mm
05.5.0303.116E	Broca Ø 1,6x Stop 18mmx Comprimento 90mm
05.5.0304.116E	Broca Ø 1,6x Stop 18mmx Comprimento 95mm
05.5.0305.116E	Broca Ø 1,6x Stop 18mmx Comprimento 100mm
05.5.0306.116E	Broca Ø 1,6x Stop 18mmx Comprimento 105mm
05.5.0307.116E	Broca Ø 1,6x Stop 18mmx Comprimento 110mm
05.5.0308.116E	Broca Ø 1,6x Stop 18mmx Comprimento 115mm
05.5.0309.116E	Broca Ø 1,6x Stop 18mmx Comprimento 120mm
05.5.0310.116E	Broca Ø 1,6x Stop 18mmx Comprimento 125mm
05.5.0311.116E	Broca Ø 1,6x Stop 18mmx Comprimento 130mm
05.5.0312.116E	Broca Ø 1,6x Stop 18mmx Comprimento 135mm
05.5.0313.116E	Broca Ø 1,6x Stop 18mmx Comprimento 140mm
05.5.0314.116E	Broca Ø 1,6x Stop 18mmx Comprimento 145mm
05.5.0315.116E	Broca Ø 1,6x Stop 18mmx Comprimento 150mm
05.5.0316.116E	Broca Ø 1,6x Stop 19mmx Comprimento 50mm
05.5.0317.116E	Broca Ø 1,6x Stop 19mmx Comprimento 55mm
05.5.0318.116E	Broca Ø 1,6x Stop 19mmx Comprimento 60mm
05.5.0319.116E	Broca Ø 1,6x Stop 19mmx Comprimento 65mm
05.5.0320.116E	Broca Ø 1,6x Stop 19mmx Comprimento 70mm
05.5.0321.116E	Broca Ø 1,6x Stop 19mmx Comprimento 75mm
05.5.0322.116E	Broca Ø 1,6x Stop 19mmx Comprimento 80mm
05.5.0323.116E	Broca Ø 1,6x Stop 19mmx Comprimento 85mm
05.5.0324.116E	Broca Ø 1,6x Stop 19mmx Comprimento 90mm
05.5.0325.116E	Broca Ø 1,6x Stop 19mmx Comprimento 95mm
05.5.0326.116E	Broca Ø 1,6x Stop 19mmx Comprimento 100mm
05.5.0327.116E	Broca Ø 1,6x Stop 19mmx Comprimento 105mm
05.5.0328.116E	Broca Ø 1,6x Stop 19mmx Comprimento 110mm
05.5.0329.116E	Broca Ø 1,6x Stop 19mmx Comprimento 115mm
05.5.0330.116E	Broca Ø 1,6x Stop 19mmx Comprimento 120mm
05.5.0331.116E	Broca Ø 1,6x Stop 19mmx Comprimento 125mm
05.5.0332.116E	Broca Ø 1,6x Stop 19mmx Comprimento 130mm
05.5.0333.116E	Broca Ø 1,6x Stop 19mmx Comprimento 135mm
05.5.0334.116E	Broca Ø 1,6x Stop 19mmx Comprimento 140mm
05.5.0335.116E	Broca Ø 1,6x Stop 19mmx Comprimento 145mm
05.5.0336.116E	Broca Ø 1,6x Stop 19mmx Comprimento 150mm
05.5.0337.116E	Broca Ø 1,6x Stop 20mmx Comprimento 50mm
05.5.0338.116E	Broca Ø 1,6x Stop 20mmx Comprimento 55mm
05.5.0339.116E	Broca Ø 1,6x Stop 20mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0340.116E	Broca Ø 1,6x Stop 20mmx Comprimento 65mm
05.5.0341.116E	Broca Ø 1,6x Stop 20mmx Comprimento 70mm
05.5.0342.116E	Broca Ø 1,6x Stop 20mmx Comprimento 75mm
05.5.0343.116E	Broca Ø 1,6x Stop 20mmx Comprimento 80mm
05.5.0344.116E	Broca Ø 1,6x Stop 20mmx Comprimento 85mm
05.5.0345.116E	Broca Ø 1,6x Stop 20mmx Comprimento 90mm
05.5.0346.116E	Broca Ø 1,6x Stop 20mmx Comprimento 95mm
05.5.0347.116E	Broca Ø 1,6x Stop 20mmx Comprimento 100mm
05.5.0348.116E	Broca Ø 1,6x Stop 20mmx Comprimento 105mm
05.5.0349.116E	Broca Ø 1,6x Stop 20mmx Comprimento 110mm
05.5.0350.116E	Broca Ø 1,6x Stop 20mmx Comprimento 115mm
05.5.0351.116E	Broca Ø 1,6x Stop 20mmx Comprimento 120mm
05.5.0352.116E	Broca Ø 1,6x Stop 20mmx Comprimento 125mm
05.5.0353.116E	Broca Ø 1,6x Stop 20mmx Comprimento 130mm
05.5.0354.116E	Broca Ø 1,6x Stop 20mmx Comprimento 135mm
05.5.0355.116E	Broca Ø 1,6x Stop 20mmx Comprimento 140mm
05.5.0356.116E	Broca Ø 1,6x Stop 20mmx Comprimento 145mm
05.5.0357.116E	Broca Ø 1,6x Stop 20mmx Comprimento 150mm
05.5.0358.116E	Broca Ø 1,6x Stop 21mmx Comprimento 50mm
05.5.0359.116E	Broca Ø 1,6x Stop 21mmx Comprimento 55mm
05.5.0360.116E	Broca Ø 1,6x Stop 21mmx Comprimento 60mm
05.5.0361.116E	Broca Ø 1,6x Stop 21mmx Comprimento 65mm
05.5.0362.116E	Broca Ø 1,6x Stop 21mmx Comprimento 70mm
05.5.0363.116E	Broca Ø 1,6x Stop 21mmx Comprimento 75mm
05.5.0364.116E	Broca Ø 1,6x Stop 21mmx Comprimento 80mm
05.5.0365.116E	Broca Ø 1,6x Stop 21mmx Comprimento 85mm
05.5.0366.116E	Broca Ø 1,6x Stop 21mmx Comprimento 90mm
05.5.0367.116E	Broca Ø 1,6x Stop 21mmx Comprimento 95mm
05.5.0368.116E	Broca Ø 1,6x Stop 21mmx Comprimento 100mm
05.5.0369.116E	Broca Ø 1,6x Stop 21mmx Comprimento 105mm
05.5.0370.116E	Broca Ø 1,6x Stop 21mmx Comprimento 110mm
05.5.0371.116E	Broca Ø 1,6x Stop 21mmx Comprimento 115mm
05.5.0372.116E	Broca Ø 1,6x Stop 21mmx Comprimento 120mm
05.5.0373.116E	Broca Ø 1,6x Stop 21mmx Comprimento 125mm
05.5.0374.116E	Broca Ø 1,6x Stop 21mmx Comprimento 130mm
05.5.0375.116E	Broca Ø 1,6x Stop 21mmx Comprimento 135mm
05.5.0376.116E	Broca Ø 1,6x Stop 21mmx Comprimento 140mm
05.5.0377.116E	Broca Ø 1,6x Stop 21mmx Comprimento 145mm
05.5.0378.116E	Broca Ø 1,6x Stop 21mmx Comprimento 150mm
05.5.0379.116E	Broca Ø 1,6x Stop 22mmx Comprimento 50mm
05.5.0380.116E	Broca Ø 1,6x Stop 22mmx Comprimento 55mm
05.5.0381.116E	Broca Ø 1,6x Stop 22mmx Comprimento 60mm
05.5.0382.116E	Broca Ø 1,6x Stop 22mmx Comprimento 65mm
05.5.0383.116E	Broca Ø 1,6x Stop 22mmx Comprimento 70mm
05.5.0384.116E	Broca Ø 1,6x Stop 22mmx Comprimento 75mm
05.5.0385.116E	Broca Ø 1,6x Stop 22mmx Comprimento 80mm
05.5.0386.116E	Broca Ø 1,6x Stop 22mmx Comprimento 85mm
05.5.0387.116E	Broca Ø 1,6x Stop 22mmx Comprimento 90mm
05.5.0388.116E	Broca Ø 1,6x Stop 22mmx Comprimento 95mm

05.5.0389.116E	Broca Ø 1,6x Stop 22mmx Comprimento 100mm
05.5.0390.116E	Broca Ø 1,6x Stop 22mmx Comprimento 105mm
05.5.0391.116E	Broca Ø 1,6x Stop 22mmx Comprimento 110mm
05.5.0392.116E	Broca Ø 1,6x Stop 22mmx Comprimento 115mm
05.5.0393.116E	Broca Ø 1,6x Stop 22mmx Comprimento 120mm
05.5.0394.116E	Broca Ø 1,6x Stop 22mmx Comprimento 125mm
05.5.0395.116E	Broca Ø 1,6x Stop 22mmx Comprimento 130mm
05.5.0396.116E	Broca Ø 1,6x Stop 22mmx Comprimento 135mm
05.5.0397.116E	Broca Ø 1,6x Stop 22mmx Comprimento 140mm
05.5.0398.116E	Broca Ø 1,6x Stop 22mmx Comprimento 145mm
05.5.0399.116E	Broca Ø 1,6x Stop 22mmx Comprimento 150mm
05.5.0400.116E	Broca Ø 1,6x Stop 23mmx Comprimento 50mm
05.5.0401.116E	Broca Ø 1,6x Stop 23mmx Comprimento 55mm
05.5.0402.116E	Broca Ø 1,6x Stop 23mmx Comprimento 60mm
05.5.0403.116E	Broca Ø 1,6x Stop 23mmx Comprimento 65mm
05.5.0404.116E	Broca Ø 1,6x Stop 23mmx Comprimento 70mm
05.5.0405.116E	Broca Ø 1,6x Stop 23mmx Comprimento 75mm
05.5.0406.116E	Broca Ø 1,6x Stop 23mmx Comprimento 80mm
05.5.0407.116E	Broca Ø 1,6x Stop 23mmx Comprimento 85mm
05.5.0408.116E	Broca Ø 1,6x Stop 23mmx Comprimento 90mm
05.5.0409.116E	Broca Ø 1,6x Stop 23mmx Comprimento 95mm
05.5.0410.116E	Broca Ø 1,6x Stop 23mmx Comprimento 100mm
05.5.0411.116E	Broca Ø 1,6x Stop 23mmx Comprimento 105mm
05.5.0412.116E	Broca Ø 1,6x Stop 23mmx Comprimento 110mm
05.5.0413.116E	Broca Ø 1,6x Stop 23mmx Comprimento 115mm
05.5.0414.116E	Broca Ø 1,6x Stop 23mmx Comprimento 120mm
05.5.0415.116E	Broca Ø 1,6x Stop 23mmx Comprimento 125mm
05.5.0416.116E	Broca Ø 1,6x Stop 23mmx Comprimento 130mm
05.5.0417.116E	Broca Ø 1,6x Stop 23mmx Comprimento 135mm
05.5.0418.116E	Broca Ø 1,6x Stop 23mmx Comprimento 140mm
05.5.0419.116E	Broca Ø 1,6x Stop 23mmx Comprimento 145mm
05.5.0420.116E	Broca Ø 1,6x Stop 23mmx Comprimento 150mm
05.5.0421.116E	Broca Ø 1,6x Stop 24mmx Comprimento 50mm
05.5.0422.116E	Broca Ø 1,6x Stop 24mmx Comprimento 55mm
05.5.0423.116E	Broca Ø 1,6x Stop 24mmx Comprimento 60mm
05.5.0424.116E	Broca Ø 1,6x Stop 24mmx Comprimento 65mm
05.5.0425.116E	Broca Ø 1,6x Stop 24mmx Comprimento 70mm
05.5.0426.116E	Broca Ø 1,6x Stop 24mmx Comprimento 75mm
05.5.0427.116E	Broca Ø 1,6x Stop 24mmx Comprimento 80mm
05.5.0428.116E	Broca Ø 1,6x Stop 24mmx Comprimento 85mm
05.5.0429.116E	Broca Ø 1,6x Stop 24mmx Comprimento 90mm
05.5.0430.116E	Broca Ø 1,6x Stop 24mmx Comprimento 95mm
05.5.0431.116E	Broca Ø 1,6x Stop 24mmx Comprimento 100mm
05.5.0432.116E	Broca Ø 1,6x Stop 24mmx Comprimento 105mm
05.5.0433.116E	Broca Ø 1,6x Stop 24mmx Comprimento 110mm
05.5.0434.116E	Broca Ø 1,6x Stop 24mmx Comprimento 115mm
05.5.0435.116E	Broca Ø 1,6x Stop 24mmx Comprimento 120mm
05.5.0436.116E	Broca Ø 1,6x Stop 24mmx Comprimento 125mm
05.5.0437.116E	Broca Ø 1,6x Stop 24mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0438.116E	Broca Ø 1,6x Stop 24mmx Comprimento 135mm
05.5.0439.116E	Broca Ø 1,6x Stop 24mmx Comprimento 140mm
05.5.0440.116E	Broca Ø 1,6x Stop 24mmx Comprimento 145mm
05.5.0441.116E	Broca Ø 1,6x Stop 24mmx Comprimento 150mm
05.5.0442.116E	Broca Ø 1,6x Stop 25mmx Comprimento 50mm
05.5.0443.116E	Broca Ø 1,6x Stop 25mmx Comprimento 55mm
05.5.0444.116E	Broca Ø 1,6x Stop 25mmx Comprimento 60mm
05.5.0445.116E	Broca Ø 1,6x Stop 25mmx Comprimento 65mm
05.5.0446.116E	Broca Ø 1,6x Stop 25mmx Comprimento 70mm
05.5.0447.116E	Broca Ø 1,6x Stop 25mmx Comprimento 75mm
05.5.0448.116E	Broca Ø 1,6x Stop 25mmx Comprimento 80mm
05.5.0449.116E	Broca Ø 1,6x Stop 25mmx Comprimento 85mm
05.5.0450.116E	Broca Ø 1,6x Stop 25mmx Comprimento 90mm
05.5.0451.116E	Broca Ø 1,6x Stop 25mmx Comprimento 95mm
05.5.0452.116E	Broca Ø 1,6x Stop 25mmx Comprimento 100mm
05.5.0453.116E	Broca Ø 1,6x Stop 25mmx Comprimento 105mm
05.5.0454.116E	Broca Ø 1,6x Stop 25mmx Comprimento 110mm
05.5.0455.116E	Broca Ø 1,6x Stop 25mmx Comprimento 115mm
05.5.0456.116E	Broca Ø 1,6x Stop 25mmx Comprimento 120mm
05.5.0457.116E	Broca Ø 1,6x Stop 25mmx Comprimento 125mm
05.5.0458.116E	Broca Ø 1,6x Stop 25mmx Comprimento 130mm
05.5.0459.116E	Broca Ø 1,6x Stop 25mmx Comprimento 135mm
05.5.0460.116E	Broca Ø 1,6x Stop 25mmx Comprimento 140mm
05.5.0461.116E	Broca Ø 1,6x Stop 25mmx Comprimento 145mm
05.5.0462.116E	Broca Ø 1,6x Stop 25mmx Comprimento 150mm
05.5.0463.116E	Broca Ø 1,6x Stop 30mmx Comprimento 50mm
05.5.0464.116E	Broca Ø 1,6x Stop 30mmx Comprimento 55mm
05.5.0465.116E	Broca Ø 1,6x Stop 30mmx Comprimento 60mm
05.5.0466.116E	Broca Ø 1,6x Stop 30mmx Comprimento 65mm
05.5.0467.116E	Broca Ø 1,6x Stop 30mmx Comprimento 70mm
05.5.0468.116E	Broca Ø 1,6x Stop 30mmx Comprimento 75mm
05.5.0469.116E	Broca Ø 1,6x Stop 30mmx Comprimento 80mm
05.5.0470.116E	Broca Ø 1,6x Stop 30mmx Comprimento 85mm
05.5.0471.116E	Broca Ø 1,6x Stop 30mmx Comprimento 90mm
05.5.0472.116E	Broca Ø 1,6x Stop 30mmx Comprimento 95mm
05.5.0473.116E	Broca Ø 1,6x Stop 30mmx Comprimento 100mm
05.5.0474.116E	Broca Ø 1,6x Stop 30mmx Comprimento 105mm
05.5.0475.116E	Broca Ø 1,6x Stop 30mmx Comprimento 110mm
05.5.0476.116E	Broca Ø 1,6x Stop 30mmx Comprimento 115mm
05.5.0477.116E	Broca Ø 1,6x Stop 30mmx Comprimento 120mm
05.5.0478.116E	Broca Ø 1,6x Stop 30mmx Comprimento 125mm
05.5.0479.116E	Broca Ø 1,6x Stop 30mmx Comprimento 130mm
05.5.0480.116E	Broca Ø 1,6x Stop 30mmx Comprimento 135mm
05.5.0481.116E	Broca Ø 1,6x Stop 30mmx Comprimento 140mm
05.5.0482.116E	Broca Ø 1,6x Stop 30mmx Comprimento 145mm
05.5.0483.116E	Broca Ø 1,6x Stop 30mmx Comprimento 150mm
05.5.0484.116E	Broca Ø 1,6x Stop 35mmx Comprimento 50mm
05.5.0485.116E	Broca Ø 1,6x Stop 35mmx Comprimento 55mm
05.5.0486.116E	Broca Ø 1,6x Stop 35mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0487.116E	Broca Ø 1,6x Stop 35mmx Comprimento 65mm
05.5.0488.116E	Broca Ø 1,6x Stop 35mmx Comprimento 70mm
05.5.0489.116E	Broca Ø 1,6x Stop 35mmx Comprimento 75mm
05.5.0490.116E	Broca Ø 1,6x Stop 35mmx Comprimento 80mm
05.5.0491.116E	Broca Ø 1,6x Stop 35mmx Comprimento 85mm
05.5.0492.116E	Broca Ø 1,6x Stop 35mmx Comprimento 90mm
05.5.0493.116E	Broca Ø 1,6x Stop 35mmx Comprimento 95mm
05.5.0494.116E	Broca Ø 1,6x Stop 35mmx Comprimento 100mm
05.5.0495.116E	Broca Ø 1,6x Stop 35mmx Comprimento 105mm
05.5.0496.116E	Broca Ø 1,6x Stop 35mmx Comprimento 110mm
05.5.0497.116E	Broca Ø 1,6x Stop 35mmx Comprimento 115mm
05.5.0498.116E	Broca Ø 1,6x Stop 35mmx Comprimento 120mm
05.5.0499.116E	Broca Ø 1,6x Stop 35mmx Comprimento 125mm
05.5.0500.116E	Broca Ø 1,6x Stop 35mmx Comprimento 130mm
05.5.0501.116E	Broca Ø 1,6x Stop 35mmx Comprimento 135mm
05.5.0502.116E	Broca Ø 1,6x Stop 35mmx Comprimento 140mm
05.5.0503.116E	Broca Ø 1,6x Stop 35mmx Comprimento 145mm
05.5.0504.116E	Broca Ø 1,6x Stop 35mmx Comprimento 150mm
05.5.0001.117E	Broca Ø 1,7x Stop 4mmx Comprimento 50mm
05.5.0002.117E	Broca Ø 1,7x Stop 4mmx Comprimento 55mm
05.5.0003.117E	Broca Ø 1,7x Stop 4mmx Comprimento 60mm
05.5.0004.117E	Broca Ø 1,7x Stop 4mmx Comprimento 65mm
05.5.0005.117E	Broca Ø 1,7x Stop 4mmx Comprimento 70mm
05.5.0006.117E	Broca Ø 1,7x Stop 4mmx Comprimento 75mm
05.5.0007.117E	Broca Ø 1,7x Stop 4mmx Comprimento 80mm
05.5.0008.117E	Broca Ø 1,7x Stop 4mmx Comprimento 85mm
05.5.0009.117E	Broca Ø 1,7x Stop 4mmx Comprimento 90mm
05.5.0010.117E	Broca Ø 1,7x Stop 4mmx Comprimento 95mm
05.5.0011.117E	Broca Ø 1,7x Stop 4mmx Comprimento 100mm
05.5.0012.117E	Broca Ø 1,7x Stop 4mmx Comprimento 105mm
05.5.0013.117E	Broca Ø 1,7x Stop 4mmx Comprimento 110mm
05.5.0014.117E	Broca Ø 1,7x Stop 4mmx Comprimento 115mm
05.5.0015.117E	Broca Ø 1,7x Stop 4mmx Comprimento 120mm
05.5.0016.117E	Broca Ø 1,7x Stop 4mmx Comprimento 125mm
05.5.0017.117E	Broca Ø 1,7x Stop 4mmx Comprimento 130mm
05.5.0018.117E	Broca Ø 1,7x Stop 4mmx Comprimento 135mm
05.5.0019.117E	Broca Ø 1,7x Stop 4mmx Comprimento 140mm
05.5.0020.117E	Broca Ø 1,7x Stop 4mmx Comprimento 145mm
05.5.0021.117E	Broca Ø 1,7x Stop 4mmx Comprimento 150mm
05.5.0022.117E	Broca Ø 1,7x Stop 5mmx Comprimento 50mm
05.5.0023.117E	Broca Ø 1,7x Stop 5mmx Comprimento 55mm
05.5.0024.117E	Broca Ø 1,7x Stop 5mmx Comprimento 60mm
05.5.0025.117E	Broca Ø 1,7x Stop 5mmx Comprimento 65mm
05.5.0026.117E	Broca Ø 1,7x Stop 5mmx Comprimento 70mm
05.5.0027.117E	Broca Ø 1,7x Stop 5mmx Comprimento 75mm
05.5.0028.117E	Broca Ø 1,7x Stop 5mmx Comprimento 80mm
05.5.0029.117E	Broca Ø 1,7x Stop 5mmx Comprimento 85mm
05.5.0030.117E	Broca Ø 1,7x Stop 5mmx Comprimento 90mm
05.5.0031.117E	Broca Ø 1,7x Stop 5mmx Comprimento 95mm

05.5.0032.117E	Broca Ø 1,7x Stop 5mmx Comprimento 100mm
05.5.0033.117E	Broca Ø 1,7x Stop 5mmx Comprimento 105mm
05.5.0034.117E	Broca Ø 1,7x Stop 5mmx Comprimento 110mm
05.5.0035.117E	Broca Ø 1,7x Stop 5mmx Comprimento 115mm
05.5.0036.117E	Broca Ø 1,7x Stop 5mmx Comprimento 120mm
05.5.0037.117E	Broca Ø 1,7x Stop 5mmx Comprimento 125mm
05.5.0038.117E	Broca Ø 1,7x Stop 5mmx Comprimento 130mm
05.5.0039.117E	Broca Ø 1,7x Stop 5mmx Comprimento 135mm
05.5.0040.117E	Broca Ø 1,7x Stop 5mmx Comprimento 140mm
05.5.0041.117E	Broca Ø 1,7x Stop 5mmx Comprimento 145mm
05.5.0042.117E	Broca Ø 1,7x Stop 5mmx Comprimento 150mm
05.5.0043.117E	Broca Ø 1,7x Stop 6mmx Comprimento 50mm
05.5.0044.117E	Broca Ø 1,7x Stop 6mmx Comprimento 55mm
05.5.0045.117E	Broca Ø 1,7x Stop 6mmx Comprimento 60mm
05.5.0046.117E	Broca Ø 1,7x Stop 6mmx Comprimento 65mm
05.5.0047.117E	Broca Ø 1,7x Stop 6mmx Comprimento 70mm
05.5.0048.117E	Broca Ø 1,7x Stop 6mmx Comprimento 75mm
05.5.0049.117E	Broca Ø 1,7x Stop 6mmx Comprimento 80mm
05.5.0050.117E	Broca Ø 1,7x Stop 6mmx Comprimento 85mm
05.5.0051.117E	Broca Ø 1,7x Stop 6mmx Comprimento 90mm
05.5.0052.117E	Broca Ø 1,7x Stop 6mmx Comprimento 95mm
05.5.0053.117E	Broca Ø 1,7x Stop 6mmx Comprimento 100mm
05.5.0054.117E	Broca Ø 1,7x Stop 6mmx Comprimento 105mm
05.5.0055.117E	Broca Ø 1,7x Stop 6mmx Comprimento 110mm
05.5.0056.117E	Broca Ø 1,7x Stop 6mmx Comprimento 115mm
05.5.0057.117E	Broca Ø 1,7x Stop 6mmx Comprimento 120mm
05.5.0058.117E	Broca Ø 1,7x Stop 6mmx Comprimento 125mm
05.5.0059.117E	Broca Ø 1,7x Stop 6mmx Comprimento 130mm
05.5.0060.117E	Broca Ø 1,7x Stop 6mmx Comprimento 135mm
05.5.0061.117E	Broca Ø 1,7x Stop 6mmx Comprimento 140mm
05.5.0062.117E	Broca Ø 1,7x Stop 6mmx Comprimento 145mm
05.5.0063.117E	Broca Ø 1,7x Stop 6mmx Comprimento 150mm
05.5.0064.117E	Broca Ø 1,7x Stop 7mmx Comprimento 50mm
05.5.0065.117E	Broca Ø 1,7x Stop 7mmx Comprimento 55mm
05.5.0066.117E	Broca Ø 1,7x Stop 7mmx Comprimento 60mm
05.5.0067.117E	Broca Ø 1,7x Stop 7mmx Comprimento 65mm
05.5.0068.117E	Broca Ø 1,7x Stop 7mmx Comprimento 70mm
05.5.0069.117E	Broca Ø 1,7x Stop 7mmx Comprimento 75mm
05.5.0070.117E	Broca Ø 1,7x Stop 7mmx Comprimento 80mm
05.5.0071.117E	Broca Ø 1,7x Stop 7mmx Comprimento 85mm
05.5.0072.117E	Broca Ø 1,7x Stop 7mmx Comprimento 90mm
05.5.0073.117E	Broca Ø 1,7x Stop 7mmx Comprimento 95mm
05.5.0074.117E	Broca Ø 1,7x Stop 7mmx Comprimento 100mm
05.5.0075.117E	Broca Ø 1,7x Stop 7mmx Comprimento 105mm
05.5.0076.117E	Broca Ø 1,7x Stop 7mmx Comprimento 110mm
05.5.0077.117E	Broca Ø 1,7x Stop 7mmx Comprimento 115mm
05.5.0078.117E	Broca Ø 1,7x Stop 7mmx Comprimento 120mm
05.5.0079.117E	Broca Ø 1,7x Stop 7mmx Comprimento 125mm
05.5.0080.117E	Broca Ø 1,7x Stop 7mmx Comprimento 130mm

05.5.0081.117E	Broca Ø 1,7x Stop 7mmx Comprimento 135mm
05.5.0082.117E	Broca Ø 1,7x Stop 7mmx Comprimento 140mm
05.5.0083.117E	Broca Ø 1,7x Stop 7mmx Comprimento 145mm
05.5.0084.117E	Broca Ø 1,7x Stop 7mmx Comprimento 150mm
05.5.0085.117E	Broca Ø 1,7x Stop 8mmx Comprimento 50mm
05.5.0086.117E	Broca Ø 1,7x Stop 8mmx Comprimento 55mm
05.5.0087.117E	Broca Ø 1,7x Stop 8mmx Comprimento 60mm
05.5.0088.117E	Broca Ø 1,7x Stop 8mmx Comprimento 65mm
05.5.0089.117E	Broca Ø 1,7x Stop 8mmx Comprimento 70mm
05.5.0090.117E	Broca Ø 1,7x Stop 8mmx Comprimento 75mm
05.5.0091.117E	Broca Ø 1,7x Stop 8mmx Comprimento 80mm
05.5.0092.117E	Broca Ø 1,7x Stop 8mmx Comprimento 85mm
05.5.0093.117E	Broca Ø 1,7x Stop 8mmx Comprimento 90mm
05.5.0094.117E	Broca Ø 1,7x Stop 8mmx Comprimento 95mm
05.5.0095.117E	Broca Ø 1,7x Stop 8mmx Comprimento 100mm
05.5.0096.117E	Broca Ø 1,7x Stop 8mmx Comprimento 105mm
05.5.0097.117E	Broca Ø 1,7x Stop 8mmx Comprimento 110mm
05.5.0098.117E	Broca Ø 1,7x Stop 8mmx Comprimento 115mm
05.5.0099.117E	Broca Ø 1,7x Stop 8mmx Comprimento 120mm
05.5.0100.117E	Broca Ø 1,7x Stop 8mmx Comprimento 125mm
05.5.0101.117E	Broca Ø 1,7x Stop 8mmx Comprimento 130mm
05.5.0102.117E	Broca Ø 1,7x Stop 8mmx Comprimento 135mm
05.5.0103.117E	Broca Ø 1,7x Stop 8mmx Comprimento 140mm
05.5.0104.117E	Broca Ø 1,7x Stop 8mmx Comprimento 145mm
05.5.0105.117E	Broca Ø 1,7x Stop 8mmx Comprimento 150mm
05.5.0106.117E	Broca Ø 1,7x Stop 9mmx Comprimento 50mm
05.5.0107.117E	Broca Ø 1,7x Stop 9mmx Comprimento 55mm
05.5.0108.117E	Broca Ø 1,7x Stop 9mmx Comprimento 60mm
05.5.0109.117E	Broca Ø 1,7x Stop 9mmx Comprimento 65mm
05.5.0110.117E	Broca Ø 1,7x Stop 9mmx Comprimento 70mm
05.5.0111.117E	Broca Ø 1,7x Stop 9mmx Comprimento 75mm
05.5.0112.117E	Broca Ø 1,7x Stop 9mmx Comprimento 80mm
05.5.0113.117E	Broca Ø 1,7x Stop 9mmx Comprimento 85mm
05.5.0114.117E	Broca Ø 1,7x Stop 9mmx Comprimento 90mm
05.5.0115.117E	Broca Ø 1,7x Stop 9mmx Comprimento 95mm
05.5.0116.117E	Broca Ø 1,7x Stop 9mmx Comprimento 100mm
05.5.0117.117E	Broca Ø 1,7x Stop 9mmx Comprimento 105mm
05.5.0118.117E	Broca Ø 1,7x Stop 9mmx Comprimento 110mm
05.5.0119.117E	Broca Ø 1,7x Stop 9mmx Comprimento 115mm
05.5.0120.117E	Broca Ø 1,7x Stop 9mmx Comprimento 120mm
05.5.0121.117E	Broca Ø 1,7x Stop 9mmx Comprimento 125mm
05.5.0122.117E	Broca Ø 1,7x Stop 9mmx Comprimento 130mm
05.5.0123.117E	Broca Ø 1,7x Stop 9mmx Comprimento 135mm
05.5.0124.117E	Broca Ø 1,7x Stop 9mmx Comprimento 140mm
05.5.0125.117E	Broca Ø 1,7x Stop 9mmx Comprimento 145mm
05.5.0126.117E	Broca Ø 1,7x Stop 9mmx Comprimento 150mm
05.5.0127.117E	Broca Ø 1,7x Stop 10mmx Comprimento 50mm
05.5.0128.117E	Broca Ø 1,7x Stop 10mmx Comprimento 55mm
05.5.0129.117E	Broca Ø 1,7x Stop 10mmx Comprimento 60mm

05.5.0130.117E	Broca Ø 1,7x Stop 10mmx Comprimento 65mm
05.5.0131.117E	Broca Ø 1,7x Stop 10mmx Comprimento 70mm
05.5.0132.117E	Broca Ø 1,7x Stop 10mmx Comprimento 75mm
05.5.0133.117E	Broca Ø 1,7x Stop 10mmx Comprimento 80mm
05.5.0134.117E	Broca Ø 1,7x Stop 10mmx Comprimento 85mm
05.5.0135.117E	Broca Ø 1,7x Stop 10mmx Comprimento 90mm
05.5.0136.117E	Broca Ø 1,7x Stop 10mmx Comprimento 95mm
05.5.0137.117E	Broca Ø 1,7x Stop 10mmx Comprimento 100mm
05.5.0138.117E	Broca Ø 1,7x Stop 10mmx Comprimento 105mm
05.5.0139.117E	Broca Ø 1,7x Stop 10mmx Comprimento 110mm
05.5.0140.117E	Broca Ø 1,7x Stop 10mmx Comprimento 115mm
05.5.0141.117E	Broca Ø 1,7x Stop 10mmx Comprimento 120mm
05.5.0142.117E	Broca Ø 1,7x Stop 10mmx Comprimento 125mm
05.5.0143.117E	Broca Ø 1,7x Stop 10mmx Comprimento 130mm
05.5.0144.117E	Broca Ø 1,7x Stop 10mmx Comprimento 135mm
05.5.0145.117E	Broca Ø 1,7x Stop 10mmx Comprimento 140mm
05.5.0146.117E	Broca Ø 1,7x Stop 10mmx Comprimento 145mm
05.5.0147.117E	Broca Ø 1,7x Stop 10mmx Comprimento 150mm
05.5.0148.117E	Broca Ø 1,7x Stop 11mmx Comprimento 50mm
05.5.0149.117E	Broca Ø 1,7x Stop 11mmx Comprimento 55mm
05.5.0150.117E	Broca Ø 1,7x Stop 11mmx Comprimento 60mm
05.5.0151.117E	Broca Ø 1,7x Stop 11mmx Comprimento 65mm
05.5.0152.117E	Broca Ø 1,7x Stop 11mmx Comprimento 70mm
05.5.0153.117E	Broca Ø 1,7x Stop 11mmx Comprimento 75mm
05.5.0154.117E	Broca Ø 1,7x Stop 11mmx Comprimento 80mm
05.5.0155.117E	Broca Ø 1,7x Stop 11mmx Comprimento 85mm
05.5.0156.117E	Broca Ø 1,7x Stop 11mmx Comprimento 90mm
05.5.0157.117E	Broca Ø 1,7x Stop 11mmx Comprimento 95mm
05.5.0158.117E	Broca Ø 1,7x Stop 11mmx Comprimento 100mm
05.5.0159.117E	Broca Ø 1,7x Stop 11mmx Comprimento 105mm
05.5.0160.117E	Broca Ø 1,7x Stop 11mmx Comprimento 110mm
05.5.0161.117E	Broca Ø 1,7x Stop 11mmx Comprimento 115mm
05.5.0162.117E	Broca Ø 1,7x Stop 11mmx Comprimento 120mm
05.5.0163.117E	Broca Ø 1,7x Stop 11mmx Comprimento 125mm
05.5.0164.117E	Broca Ø 1,7x Stop 11mmx Comprimento 130mm
05.5.0165.117E	Broca Ø 1,7x Stop 11mmx Comprimento 135mm
05.5.0166.117E	Broca Ø 1,7x Stop 11mmx Comprimento 140mm
05.5.0167.117E	Broca Ø 1,7x Stop 11mmx Comprimento 145mm
05.5.0168.117E	Broca Ø 1,7x Stop 11mmx Comprimento 150mm
05.5.0169.117E	Broca Ø 1,7x Stop 12mmx Comprimento 50mm
05.5.0170.117E	Broca Ø 1,7x Stop 12mmx Comprimento 55mm
05.5.0171.117E	Broca Ø 1,7x Stop 12mmx Comprimento 60mm
05.5.0172.117E	Broca Ø 1,7x Stop 12mmx Comprimento 65mm
05.5.0173.117E	Broca Ø 1,7x Stop 12mmx Comprimento 70mm
05.5.0174.117E	Broca Ø 1,7x Stop 12mmx Comprimento 75mm
05.5.0175.117E	Broca Ø 1,7x Stop 12mmx Comprimento 80mm
05.5.0176.117E	Broca Ø 1,7x Stop 12mmx Comprimento 85mm
05.5.0177.117E	Broca Ø 1,7x Stop 12mmx Comprimento 90mm
05.5.0178.117E	Broca Ø 1,7x Stop 12mmx Comprimento 95mm

05.5.0179.117E	Broca Ø 1,7x Stop 12mmx Comprimento 100mm
05.5.0180.117E	Broca Ø 1,7x Stop 12mmx Comprimento 105mm
05.5.0181.117E	Broca Ø 1,7x Stop 12mmx Comprimento 110mm
05.5.0182.117E	Broca Ø 1,7x Stop 12mmx Comprimento 115mm
05.5.0183.117E	Broca Ø 1,7x Stop 12mmx Comprimento 120mm
05.5.0184.117E	Broca Ø 1,7x Stop 12mmx Comprimento 125mm
05.5.0185.117E	Broca Ø 1,7x Stop 12mmx Comprimento 130mm
05.5.0186.117E	Broca Ø 1,7x Stop 12mmx Comprimento 135mm
05.5.0187.117E	Broca Ø 1,7x Stop 12mmx Comprimento 140mm
05.5.0188.117E	Broca Ø 1,7x Stop 12mmx Comprimento 145mm
05.5.0189.117E	Broca Ø 1,7x Stop 12mmx Comprimento 150mm
05.5.0190.117E	Broca Ø 1,7x Stop 13mmx Comprimento 50mm
05.5.0191.117E	Broca Ø 1,7x Stop 13mmx Comprimento 55mm
05.5.0192.117E	Broca Ø 1,7x Stop 13mmx Comprimento 60mm
05.5.0193.117E	Broca Ø 1,7x Stop 13mmx Comprimento 65mm
05.5.0194.117E	Broca Ø 1,7x Stop 13mmx Comprimento 70mm
05.5.0195.117E	Broca Ø 1,7x Stop 13mmx Comprimento 75mm
05.5.0196.117E	Broca Ø 1,7x Stop 13mmx Comprimento 80mm
05.5.0197.117E	Broca Ø 1,7x Stop 13mmx Comprimento 85mm
05.5.0198.117E	Broca Ø 1,7x Stop 13mmx Comprimento 90mm
05.5.0199.117E	Broca Ø 1,7x Stop 13mmx Comprimento 95mm
05.5.0200.117E	Broca Ø 1,7x Stop 13mmx Comprimento 100mm
05.5.0201.117E	Broca Ø 1,7x Stop 13mmx Comprimento 105mm
05.5.0202.117E	Broca Ø 1,7x Stop 13mmx Comprimento 110mm
05.5.0203.117E	Broca Ø 1,7x Stop 13mmx Comprimento 115mm
05.5.0204.117E	Broca Ø 1,7x Stop 13mmx Comprimento 120mm
05.5.0205.117E	Broca Ø 1,7x Stop 13mmx Comprimento 125mm
05.5.0206.117E	Broca Ø 1,7x Stop 13mmx Comprimento 130mm
05.5.0207.117E	Broca Ø 1,7x Stop 13mmx Comprimento 135mm
05.5.0208.117E	Broca Ø 1,7x Stop 13mmx Comprimento 140mm
05.5.0209.117E	Broca Ø 1,7x Stop 13mmx Comprimento 145mm
05.5.0210.117E	Broca Ø 1,7x Stop 13mmx Comprimento 150mm
05.5.0211.117E	Broca Ø 1,7x Stop 14mmx Comprimento 50mm
05.5.0212.117E	Broca Ø 1,7x Stop 14mmx Comprimento 55mm
05.5.0213.117E	Broca Ø 1,7x Stop 14mmx Comprimento 60mm
05.5.0214.117E	Broca Ø 1,7x Stop 14mmx Comprimento 65mm
05.5.0215.117E	Broca Ø 1,7x Stop 14mmx Comprimento 70mm
05.5.0216.117E	Broca Ø 1,7x Stop 14mmx Comprimento 75mm
05.5.0217.117E	Broca Ø 1,7x Stop 14mmx Comprimento 80mm
05.5.0218.117E	Broca Ø 1,7x Stop 14mmx Comprimento 85mm
05.5.0219.117E	Broca Ø 1,7x Stop 14mmx Comprimento 90mm
05.5.0220.117E	Broca Ø 1,7x Stop 14mmx Comprimento 95mm
05.5.0221.117E	Broca Ø 1,7x Stop 14mmx Comprimento 100mm
05.5.0222.117E	Broca Ø 1,7x Stop 14mmx Comprimento 105mm
05.5.0223.117E	Broca Ø 1,7x Stop 14mmx Comprimento 110mm
05.5.0224.117E	Broca Ø 1,7x Stop 14mmx Comprimento 115mm
05.5.0225.117E	Broca Ø 1,7x Stop 14mmx Comprimento 120mm
05.5.0226.117E	Broca Ø 1,7x Stop 14mmx Comprimento 125mm
05.5.0227.117E	Broca Ø 1,7x Stop 14mmx Comprimento 130mm

05.5.0228.117E	Broca Ø 1,7x Stop 14mmx Comprimento 135mm
05.5.0229.117E	Broca Ø 1,7x Stop 14mmx Comprimento 140mm
05.5.0230.117E	Broca Ø 1,7x Stop 14mmx Comprimento 145mm
05.5.0231.117E	Broca Ø 1,7x Stop 14mmx Comprimento 150mm
05.5.0232.117E	Broca Ø 1,7x Stop 15mmx Comprimento 50mm
05.5.0233.117E	Broca Ø 1,7x Stop 15mmx Comprimento 55mm
05.5.0234.117E	Broca Ø 1,7x Stop 15mmx Comprimento 60mm
05.5.0235.117E	Broca Ø 1,7x Stop 15mmx Comprimento 65mm
05.5.0236.117E	Broca Ø 1,7x Stop 15mmx Comprimento 70mm
05.5.0237.117E	Broca Ø 1,7x Stop 15mmx Comprimento 75mm
05.5.0238.117E	Broca Ø 1,7x Stop 15mmx Comprimento 80mm
05.5.0239.117E	Broca Ø 1,7x Stop 15mmx Comprimento 85mm
05.5.0240.117E	Broca Ø 1,7x Stop 15mmx Comprimento 90mm
05.5.0241.117E	Broca Ø 1,7x Stop 15mmx Comprimento 95mm
05.5.0242.117E	Broca Ø 1,7x Stop 15mmx Comprimento 100mm
05.5.0243.117E	Broca Ø 1,7x Stop 15mmx Comprimento 105mm
05.5.0244.117E	Broca Ø 1,7x Stop 15mmx Comprimento 110mm
05.5.0245.117E	Broca Ø 1,7x Stop 15mmx Comprimento 115mm
05.5.0246.117E	Broca Ø 1,7x Stop 15mmx Comprimento 120mm
05.5.0247.117E	Broca Ø 1,7x Stop 15mmx Comprimento 125mm
05.5.0248.117E	Broca Ø 1,7x Stop 15mmx Comprimento 130mm
05.5.0249.117E	Broca Ø 1,7x Stop 15mmx Comprimento 135mm
05.5.0250.117E	Broca Ø 1,7x Stop 15mmx Comprimento 140mm
05.5.0251.117E	Broca Ø 1,7x Stop 15mmx Comprimento 145mm
05.5.0252.117E	Broca Ø 1,7x Stop 15mmx Comprimento 150mm
05.5.0253.117E	Broca Ø 1,7x Stop 16mmx Comprimento 50mm
05.5.0254.117E	Broca Ø 1,7x Stop 16mmx Comprimento 55mm
05.5.0255.117E	Broca Ø 1,7x Stop 16mmx Comprimento 60mm
05.5.0256.117E	Broca Ø 1,7x Stop 16mmx Comprimento 65mm
05.5.0257.117E	Broca Ø 1,7x Stop 16mmx Comprimento 70mm
05.5.0258.117E	Broca Ø 1,7x Stop 16mmx Comprimento 75mm
05.5.0259.117E	Broca Ø 1,7x Stop 16mmx Comprimento 80mm
05.5.0260.117E	Broca Ø 1,7x Stop 16mmx Comprimento 85mm
05.5.0261.117E	Broca Ø 1,7x Stop 16mmx Comprimento 90mm
05.5.0262.117E	Broca Ø 1,7x Stop 16mmx Comprimento 95mm
05.5.0263.117E	Broca Ø 1,7x Stop 16mmx Comprimento 100mm
05.5.0264.117E	Broca Ø 1,7x Stop 16mmx Comprimento 105mm
05.5.0265.117E	Broca Ø 1,7x Stop 16mmx Comprimento 110mm
05.5.0266.117E	Broca Ø 1,7x Stop 16mmx Comprimento 115mm
05.5.0267.117E	Broca Ø 1,7x Stop 16mmx Comprimento 120mm
05.5.0268.117E	Broca Ø 1,7x Stop 16mmx Comprimento 125mm
05.5.0269.117E	Broca Ø 1,7x Stop 16mmx Comprimento 130mm
05.5.0270.117E	Broca Ø 1,7x Stop 16mmx Comprimento 135mm
05.5.0271.117E	Broca Ø 1,7x Stop 16mmx Comprimento 140mm
05.5.0272.117E	Broca Ø 1,7x Stop 16mmx Comprimento 145mm
05.5.0273.117E	Broca Ø 1,7x Stop 16mmx Comprimento 150mm
05.5.0274.117E	Broca Ø 1,7x Stop 17mmx Comprimento 50mm
05.5.0275.117E	Broca Ø 1,7x Stop 17mmx Comprimento 55mm
05.5.0276.117E	Broca Ø 1,7x Stop 17mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0277.117E	Broca Ø 1,7x Stop 17mmx Comprimento 65mm
05.5.0278.117E	Broca Ø 1,7x Stop 17mmx Comprimento 70mm
05.5.0279.117E	Broca Ø 1,7x Stop 17mmx Comprimento 75mm
05.5.0280.117E	Broca Ø 1,7x Stop 17mmx Comprimento 80mm
05.5.0281.117E	Broca Ø 1,7x Stop 17mmx Comprimento 85mm
05.5.0282.117E	Broca Ø 1,7x Stop 17mmx Comprimento 90mm
05.5.0283.117E	Broca Ø 1,7x Stop 17mmx Comprimento 95mm
05.5.0284.117E	Broca Ø 1,7x Stop 17mmx Comprimento 100mm
05.5.0285.117E	Broca Ø 1,7x Stop 17mmx Comprimento 105mm
05.5.0286.117E	Broca Ø 1,7x Stop 17mmx Comprimento 110mm
05.5.0287.117E	Broca Ø 1,7x Stop 17mmx Comprimento 115mm
05.5.0288.117E	Broca Ø 1,7x Stop 17mmx Comprimento 120mm
05.5.0289.117E	Broca Ø 1,7x Stop 17mmx Comprimento 125mm
05.5.0290.117E	Broca Ø 1,7x Stop 17mmx Comprimento 130mm
05.5.0291.117E	Broca Ø 1,7x Stop 17mmx Comprimento 135mm
05.5.0292.117E	Broca Ø 1,7x Stop 17mmx Comprimento 140mm
05.5.0293.117E	Broca Ø 1,7x Stop 17mmx Comprimento 145mm
05.5.0294.117E	Broca Ø 1,7x Stop 17mmx Comprimento 150mm
05.5.0295.117E	Broca Ø 1,7x Stop 18mmx Comprimento 50mm
05.5.0296.117E	Broca Ø 1,7x Stop 18mmx Comprimento 55mm
05.5.0297.117E	Broca Ø 1,7x Stop 18mmx Comprimento 60mm
05.5.0298.117E	Broca Ø 1,7x Stop 18mmx Comprimento 65mm
05.5.0299.117E	Broca Ø 1,7x Stop 18mmx Comprimento 70mm
05.5.0300.117E	Broca Ø 1,7x Stop 18mmx Comprimento 75mm
05.5.0301.117E	Broca Ø 1,7x Stop 18mmx Comprimento 80mm
05.5.0302.117E	Broca Ø 1,7x Stop 18mmx Comprimento 85mm
05.5.0303.117E	Broca Ø 1,7x Stop 18mmx Comprimento 90mm
05.5.0304.117E	Broca Ø 1,7x Stop 18mmx Comprimento 95mm
05.5.0305.117E	Broca Ø 1,7x Stop 18mmx Comprimento 100mm
05.5.0306.117E	Broca Ø 1,7x Stop 18mmx Comprimento 105mm
05.5.0307.117E	Broca Ø 1,7x Stop 18mmx Comprimento 110mm
05.5.0308.117E	Broca Ø 1,7x Stop 18mmx Comprimento 115mm
05.5.0309.117E	Broca Ø 1,7x Stop 18mmx Comprimento 120mm
05.5.0310.117E	Broca Ø 1,7x Stop 18mmx Comprimento 125mm
05.5.0311.117E	Broca Ø 1,7x Stop 18mmx Comprimento 130mm
05.5.0312.117E	Broca Ø 1,7x Stop 18mmx Comprimento 135mm
05.5.0313.117E	Broca Ø 1,7x Stop 18mmx Comprimento 140mm
05.5.0314.117E	Broca Ø 1,7x Stop 18mmx Comprimento 145mm
05.5.0315.117E	Broca Ø 1,7x Stop 18mmx Comprimento 150mm
05.5.0316.117E	Broca Ø 1,7x Stop 19mmx Comprimento 50mm
05.5.0317.117E	Broca Ø 1,7x Stop 19mmx Comprimento 55mm
05.5.0318.117E	Broca Ø 1,7x Stop 19mmx Comprimento 60mm
05.5.0319.117E	Broca Ø 1,7x Stop 19mmx Comprimento 65mm
05.5.0320.117E	Broca Ø 1,7x Stop 19mmx Comprimento 70mm
05.5.0321.117E	Broca Ø 1,7x Stop 19mmx Comprimento 75mm
05.5.0322.117E	Broca Ø 1,7x Stop 19mmx Comprimento 80mm
05.5.0323.117E	Broca Ø 1,7x Stop 19mmx Comprimento 85mm
05.5.0324.117E	Broca Ø 1,7x Stop 19mmx Comprimento 90mm
05.5.0325.117E	Broca Ø 1,7x Stop 19mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0326.117E	Broca Ø 1,7x Stop 19mmx Comprimento 100mm
05.5.0327.117E	Broca Ø 1,7x Stop 19mmx Comprimento 105mm
05.5.0328.117E	Broca Ø 1,7x Stop 19mmx Comprimento 110mm
05.5.0329.117E	Broca Ø 1,7x Stop 19mmx Comprimento 115mm
05.5.0330.117E	Broca Ø 1,7x Stop 19mmx Comprimento 120mm
05.5.0331.117E	Broca Ø 1,7x Stop 19mmx Comprimento 125mm
05.5.0332.117E	Broca Ø 1,7x Stop 19mmx Comprimento 130mm
05.5.0333.117E	Broca Ø 1,7x Stop 19mmx Comprimento 135mm
05.5.0334.117E	Broca Ø 1,7x Stop 19mmx Comprimento 140mm
05.5.0335.117E	Broca Ø 1,7x Stop 19mmx Comprimento 145mm
05.5.0336.117E	Broca Ø 1,7x Stop 19mmx Comprimento 150mm
05.5.0337.117E	Broca Ø 1,7x Stop 20mmx Comprimento 50mm
05.5.0338.117E	Broca Ø 1,7x Stop 20mmx Comprimento 55mm
05.5.0339.117E	Broca Ø 1,7x Stop 20mmx Comprimento 60mm
05.5.0340.117E	Broca Ø 1,7x Stop 20mmx Comprimento 65mm
05.5.0341.117E	Broca Ø 1,7x Stop 20mmx Comprimento 70mm
05.5.0342.117E	Broca Ø 1,7x Stop 20mmx Comprimento 75mm
05.5.0343.117E	Broca Ø 1,7x Stop 20mmx Comprimento 80mm
05.5.0344.117E	Broca Ø 1,7x Stop 20mmx Comprimento 85mm
05.5.0345.117E	Broca Ø 1,7x Stop 20mmx Comprimento 90mm
05.5.0346.117E	Broca Ø 1,7x Stop 20mmx Comprimento 95mm
05.5.0347.117E	Broca Ø 1,7x Stop 20mmx Comprimento 100mm
05.5.0348.117E	Broca Ø 1,7x Stop 20mmx Comprimento 105mm
05.5.0349.117E	Broca Ø 1,7x Stop 20mmx Comprimento 110mm
05.5.0350.117E	Broca Ø 1,7x Stop 20mmx Comprimento 115mm
05.5.0351.117E	Broca Ø 1,7x Stop 20mmx Comprimento 120mm
05.5.0352.117E	Broca Ø 1,7x Stop 20mmx Comprimento 125mm
05.5.0353.117E	Broca Ø 1,7x Stop 20mmx Comprimento 130mm
05.5.0354.117E	Broca Ø 1,7x Stop 20mmx Comprimento 135mm
05.5.0355.117E	Broca Ø 1,7x Stop 20mmx Comprimento 140mm
05.5.0356.117E	Broca Ø 1,7x Stop 20mmx Comprimento 145mm
05.5.0357.117E	Broca Ø 1,7x Stop 20mmx Comprimento 150mm
05.5.0358.117E	Broca Ø 1,7x Stop 21mmx Comprimento 50mm
05.5.0359.117E	Broca Ø 1,7x Stop 21mmx Comprimento 55mm
05.5.0360.117E	Broca Ø 1,7x Stop 21mmx Comprimento 60mm
05.5.0361.117E	Broca Ø 1,7x Stop 21mmx Comprimento 65mm
05.5.0362.117E	Broca Ø 1,7x Stop 21mmx Comprimento 70mm
05.5.0363.117E	Broca Ø 1,7x Stop 21mmx Comprimento 75mm
05.5.0364.117E	Broca Ø 1,7x Stop 21mmx Comprimento 80mm
05.5.0365.117E	Broca Ø 1,7x Stop 21mmx Comprimento 85mm
05.5.0366.117E	Broca Ø 1,7x Stop 21mmx Comprimento 90mm
05.5.0367.117E	Broca Ø 1,7x Stop 21mmx Comprimento 95mm
05.5.0368.117E	Broca Ø 1,7x Stop 21mmx Comprimento 100mm
05.5.0369.117E	Broca Ø 1,7x Stop 21mmx Comprimento 105mm
05.5.0370.117E	Broca Ø 1,7x Stop 21mmx Comprimento 110mm
05.5.0371.117E	Broca Ø 1,7x Stop 21mmx Comprimento 115mm
05.5.0372.117E	Broca Ø 1,7x Stop 21mmx Comprimento 120mm
05.5.0373.117E	Broca Ø 1,7x Stop 21mmx Comprimento 125mm
05.5.0374.117E	Broca Ø 1,7x Stop 21mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0375.117E	Broca Ø 1,7x Stop 21mmx Comprimento 135mm
05.5.0376.117E	Broca Ø 1,7x Stop 21mmx Comprimento 140mm
05.5.0377.117E	Broca Ø 1,7x Stop 21mmx Comprimento 145mm
05.5.0378.117E	Broca Ø 1,7x Stop 21mmx Comprimento 150mm
05.5.0379.117E	Broca Ø 1,7x Stop 22mmx Comprimento 50mm
05.5.0380.117E	Broca Ø 1,7x Stop 22mmx Comprimento 55mm
05.5.0381.117E	Broca Ø 1,7x Stop 22mmx Comprimento 60mm
05.5.0382.117E	Broca Ø 1,7x Stop 22mmx Comprimento 65mm
05.5.0383.117E	Broca Ø 1,7x Stop 22mmx Comprimento 70mm
05.5.0384.117E	Broca Ø 1,7x Stop 22mmx Comprimento 75mm
05.5.0385.117E	Broca Ø 1,7x Stop 22mmx Comprimento 80mm
05.5.0386.117E	Broca Ø 1,7x Stop 22mmx Comprimento 85mm
05.5.0387.117E	Broca Ø 1,7x Stop 22mmx Comprimento 90mm
05.5.0388.117E	Broca Ø 1,7x Stop 22mmx Comprimento 95mm
05.5.0389.117E	Broca Ø 1,7x Stop 22mmx Comprimento 100mm
05.5.0390.117E	Broca Ø 1,7x Stop 22mmx Comprimento 105mm
05.5.0391.117E	Broca Ø 1,7x Stop 22mmx Comprimento 110mm
05.5.0392.117E	Broca Ø 1,7x Stop 22mmx Comprimento 115mm
05.5.0393.117E	Broca Ø 1,7x Stop 22mmx Comprimento 120mm
05.5.0394.117E	Broca Ø 1,7x Stop 22mmx Comprimento 125mm
05.5.0395.117E	Broca Ø 1,7x Stop 22mmx Comprimento 130mm
05.5.0396.117E	Broca Ø 1,7x Stop 22mmx Comprimento 135mm
05.5.0397.117E	Broca Ø 1,7x Stop 22mmx Comprimento 140mm
05.5.0398.117E	Broca Ø 1,7x Stop 22mmx Comprimento 145mm
05.5.0399.117E	Broca Ø 1,7x Stop 22mmx Comprimento 150mm
05.5.0400.117E	Broca Ø 1,7x Stop 23mmx Comprimento 50mm
05.5.0401.117E	Broca Ø 1,7x Stop 23mmx Comprimento 55mm
05.5.0402.117E	Broca Ø 1,7x Stop 23mmx Comprimento 60mm
05.5.0403.117E	Broca Ø 1,7x Stop 23mmx Comprimento 65mm
05.5.0404.117E	Broca Ø 1,7x Stop 23mmx Comprimento 70mm
05.5.0405.117E	Broca Ø 1,7x Stop 23mmx Comprimento 75mm
05.5.0406.117E	Broca Ø 1,7x Stop 23mmx Comprimento 80mm
05.5.0407.117E	Broca Ø 1,7x Stop 23mmx Comprimento 85mm
05.5.0408.117E	Broca Ø 1,7x Stop 23mmx Comprimento 90mm
05.5.0409.117E	Broca Ø 1,7x Stop 23mmx Comprimento 95mm
05.5.0410.117E	Broca Ø 1,7x Stop 23mmx Comprimento 100mm
05.5.0411.117E	Broca Ø 1,7x Stop 23mmx Comprimento 105mm
05.5.0412.117E	Broca Ø 1,7x Stop 23mmx Comprimento 110mm
05.5.0413.117E	Broca Ø 1,7x Stop 23mmx Comprimento 115mm
05.5.0414.117E	Broca Ø 1,7x Stop 23mmx Comprimento 120mm
05.5.0415.117E	Broca Ø 1,7x Stop 23mmx Comprimento 125mm
05.5.0416.117E	Broca Ø 1,7x Stop 23mmx Comprimento 130mm
05.5.0417.117E	Broca Ø 1,7x Stop 23mmx Comprimento 135mm
05.5.0418.117E	Broca Ø 1,7x Stop 23mmx Comprimento 140mm
05.5.0419.117E	Broca Ø 1,7x Stop 23mmx Comprimento 145mm
05.5.0420.117E	Broca Ø 1,7x Stop 23mmx Comprimento 150mm
05.5.0421.117E	Broca Ø 1,7x Stop 24mmx Comprimento 50mm
05.5.0422.117E	Broca Ø 1,7x Stop 24mmx Comprimento 55mm
05.5.0423.117E	Broca Ø 1,7x Stop 24mmx Comprimento 60mm

05.5.0424.117E	Broca Ø 1,7x Stop 24mmx Comprimento 65mm
05.5.0425.117E	Broca Ø 1,7x Stop 24mmx Comprimento 70mm
05.5.0426.117E	Broca Ø 1,7x Stop 24mmx Comprimento 75mm
05.5.0427.117E	Broca Ø 1,7x Stop 24mmx Comprimento 80mm
05.5.0428.117E	Broca Ø 1,7x Stop 24mmx Comprimento 85mm
05.5.0429.117E	Broca Ø 1,7x Stop 24mmx Comprimento 90mm
05.5.0430.117E	Broca Ø 1,7x Stop 24mmx Comprimento 95mm
05.5.0431.117E	Broca Ø 1,7x Stop 24mmx Comprimento 100mm
05.5.0432.117E	Broca Ø 1,7x Stop 24mmx Comprimento 105mm
05.5.0433.117E	Broca Ø 1,7x Stop 24mmx Comprimento 110mm
05.5.0434.117E	Broca Ø 1,7x Stop 24mmx Comprimento 115mm
05.5.0435.117E	Broca Ø 1,7x Stop 24mmx Comprimento 120mm
05.5.0436.117E	Broca Ø 1,7x Stop 24mmx Comprimento 125mm
05.5.0437.117E	Broca Ø 1,7x Stop 24mmx Comprimento 130mm
05.5.0438.117E	Broca Ø 1,7x Stop 24mmx Comprimento 135mm
05.5.0439.117E	Broca Ø 1,7x Stop 24mmx Comprimento 140mm
05.5.0440.117E	Broca Ø 1,7x Stop 24mmx Comprimento 145mm
05.5.0441.117E	Broca Ø 1,7x Stop 24mmx Comprimento 150mm
05.5.0442.117E	Broca Ø 1,7x Stop 25mmx Comprimento 50mm
05.5.0443.117E	Broca Ø 1,7x Stop 25mmx Comprimento 55mm
05.5.0444.117E	Broca Ø 1,7x Stop 25mmx Comprimento 60mm
05.5.0445.117E	Broca Ø 1,7x Stop 25mmx Comprimento 65mm
05.5.0446.117E	Broca Ø 1,7x Stop 25mmx Comprimento 70mm
05.5.0447.117E	Broca Ø 1,7x Stop 25mmx Comprimento 75mm
05.5.0448.117E	Broca Ø 1,7x Stop 25mmx Comprimento 80mm
05.5.0449.117E	Broca Ø 1,7x Stop 25mmx Comprimento 85mm
05.5.0450.117E	Broca Ø 1,7x Stop 25mmx Comprimento 90mm
05.5.0451.117E	Broca Ø 1,7x Stop 25mmx Comprimento 95mm
05.5.0452.117E	Broca Ø 1,7x Stop 25mmx Comprimento 100mm
05.5.0453.117E	Broca Ø 1,7x Stop 25mmx Comprimento 105mm
05.5.0454.117E	Broca Ø 1,7x Stop 25mmx Comprimento 110mm
05.5.0455.117E	Broca Ø 1,7x Stop 25mmx Comprimento 115mm
05.5.0456.117E	Broca Ø 1,7x Stop 25mmx Comprimento 120mm
05.5.0457.117E	Broca Ø 1,7x Stop 25mmx Comprimento 125mm
05.5.0458.117E	Broca Ø 1,7x Stop 25mmx Comprimento 130mm
05.5.0459.117E	Broca Ø 1,7x Stop 25mmx Comprimento 135mm
05.5.0460.117E	Broca Ø 1,7x Stop 25mmx Comprimento 140mm
05.5.0461.117E	Broca Ø 1,7x Stop 25mmx Comprimento 145mm
05.5.0462.117E	Broca Ø 1,7x Stop 25mmx Comprimento 150mm
05.5.0463.117E	Broca Ø 1,7x Stop 30mmx Comprimento 50mm
05.5.0464.117E	Broca Ø 1,7x Stop 30mmx Comprimento 55mm
05.5.0465.117E	Broca Ø 1,7x Stop 30mmx Comprimento 60mm
05.5.0466.117E	Broca Ø 1,7x Stop 30mmx Comprimento 65mm
05.5.0467.117E	Broca Ø 1,7x Stop 30mmx Comprimento 70mm
05.5.0468.117E	Broca Ø 1,7x Stop 30mmx Comprimento 75mm
05.5.0469.117E	Broca Ø 1,7x Stop 30mmx Comprimento 80mm
05.5.0470.117E	Broca Ø 1,7x Stop 30mmx Comprimento 85mm
05.5.0471.117E	Broca Ø 1,7x Stop 30mmx Comprimento 90mm
05.5.0472.117E	Broca Ø 1,7x Stop 30mmx Comprimento 95mm

05.5.0473.117E	Broca Ø 1,7x Stop 30mmx Comprimento 100mm
05.5.0474.117E	Broca Ø 1,7x Stop 30mmx Comprimento 105mm
05.5.0475.117E	Broca Ø 1,7x Stop 30mmx Comprimento 110mm
05.5.0476.117E	Broca Ø 1,7x Stop 30mmx Comprimento 115mm
05.5.0477.117E	Broca Ø 1,7x Stop 30mmx Comprimento 120mm
05.5.0478.117E	Broca Ø 1,7x Stop 30mmx Comprimento 125mm
05.5.0479.117E	Broca Ø 1,7x Stop 30mmx Comprimento 130mm
05.5.0480.117E	Broca Ø 1,7x Stop 30mmx Comprimento 135mm
05.5.0481.117E	Broca Ø 1,7x Stop 30mmx Comprimento 140mm
05.5.0482.117E	Broca Ø 1,7x Stop 30mmx Comprimento 145mm
05.5.0483.117E	Broca Ø 1,7x Stop 30mmx Comprimento 150mm
05.5.0484.117E	Broca Ø 1,7x Stop 35mmx Comprimento 50mm
05.5.0485.117E	Broca Ø 1,7x Stop 35mmx Comprimento 55mm
05.5.0486.117E	Broca Ø 1,7x Stop 35mmx Comprimento 60mm
05.5.0487.117E	Broca Ø 1,7x Stop 35mmx Comprimento 65mm
05.5.0488.117E	Broca Ø 1,7x Stop 35mmx Comprimento 70mm
05.5.0489.117E	Broca Ø 1,7x Stop 35mmx Comprimento 75mm
05.5.0490.117E	Broca Ø 1,7x Stop 35mmx Comprimento 80mm
05.5.0491.117E	Broca Ø 1,7x Stop 35mmx Comprimento 85mm
05.5.0492.117E	Broca Ø 1,7x Stop 35mmx Comprimento 90mm
05.5.0493.117E	Broca Ø 1,7x Stop 35mmx Comprimento 95mm
05.5.0494.117E	Broca Ø 1,7x Stop 35mmx Comprimento 100mm
05.5.0495.117E	Broca Ø 1,7x Stop 35mmx Comprimento 105mm
05.5.0496.117E	Broca Ø 1,7x Stop 35mmx Comprimento 110mm
05.5.0497.117E	Broca Ø 1,7x Stop 35mmx Comprimento 115mm
05.5.0498.117E	Broca Ø 1,7x Stop 35mmx Comprimento 120mm
05.5.0499.117E	Broca Ø 1,7x Stop 35mmx Comprimento 125mm
05.5.0500.117E	Broca Ø 1,7x Stop 35mmx Comprimento 130mm
05.5.0501.117E	Broca Ø 1,7x Stop 35mmx Comprimento 135mm
05.5.0502.117E	Broca Ø 1,7x Stop 35mmx Comprimento 140mm
05.5.0503.117E	Broca Ø 1,7x Stop 35mmx Comprimento 145mm
05.5.0504.117E	Broca Ø 1,7x Stop 35mmx Comprimento 150mm
05.5.0001.118E	Broca Ø 1,8x Stop 4mmx Comprimento 50mm
05.5.0002.118E	Broca Ø 1,8x Stop 4mmx Comprimento 55mm
05.5.0003.118E	Broca Ø 1,8x Stop 4mmx Comprimento 60mm
05.5.0004.118E	Broca Ø 1,8x Stop 4mmx Comprimento 65mm
05.5.0005.118E	Broca Ø 1,8x Stop 4mmx Comprimento 70mm
05.5.0006.118E	Broca Ø 1,8x Stop 4mmx Comprimento 75mm
05.5.0007.118E	Broca Ø 1,8x Stop 4mmx Comprimento 80mm
05.5.0008.118E	Broca Ø 1,8x Stop 4mmx Comprimento 85mm
05.5.0009.118E	Broca Ø 1,8x Stop 4mmx Comprimento 90mm
05.5.0010.118E	Broca Ø 1,8x Stop 4mmx Comprimento 95mm
05.5.0011.118E	Broca Ø 1,8x Stop 4mmx Comprimento 100mm
05.5.0012.118E	Broca Ø 1,8x Stop 4mmx Comprimento 105mm
05.5.0013.118E	Broca Ø 1,8x Stop 4mmx Comprimento 110mm
05.5.0014.118E	Broca Ø 1,8x Stop 4mmx Comprimento 115mm
05.5.0015.118E	Broca Ø 1,8x Stop 4mmx Comprimento 120mm
05.5.0016.118E	Broca Ø 1,8x Stop 4mmx Comprimento 125mm
05.5.0017.118E	Broca Ø 1,8x Stop 4mmx Comprimento 130mm

05.5.0018.118E	Broca Ø 1,8x Stop 4mmx Comprimento 135mm
05.5.0019.118E	Broca Ø 1,8x Stop 4mmx Comprimento 140mm
05.5.0020.118E	Broca Ø 1,8x Stop 4mmx Comprimento 145mm
05.5.0021.118E	Broca Ø 1,8x Stop 4mmx Comprimento 150mm
05.5.0022.118E	Broca Ø 1,8x Stop 5mmx Comprimento 50mm
05.5.0023.118E	Broca Ø 1,8x Stop 5mmx Comprimento 55mm
05.5.0024.118E	Broca Ø 1,8x Stop 5mmx Comprimento 60mm
05.5.0025.118E	Broca Ø 1,8x Stop 5mmx Comprimento 65mm
05.5.0026.118E	Broca Ø 1,8x Stop 5mmx Comprimento 70mm
05.5.0027.118E	Broca Ø 1,8x Stop 5mmx Comprimento 75mm
05.5.0028.118E	Broca Ø 1,8x Stop 5mmx Comprimento 80mm
05.5.0029.118E	Broca Ø 1,8x Stop 5mmx Comprimento 85mm
05.5.0030.118E	Broca Ø 1,8x Stop 5mmx Comprimento 90mm
05.5.0031.118E	Broca Ø 1,8x Stop 5mmx Comprimento 95mm
05.5.0032.118E	Broca Ø 1,8x Stop 5mmx Comprimento 100mm
05.5.0033.118E	Broca Ø 1,8x Stop 5mmx Comprimento 105mm
05.5.0034.118E	Broca Ø 1,8x Stop 5mmx Comprimento 110mm
05.5.0035.118E	Broca Ø 1,8x Stop 5mmx Comprimento 115mm
05.5.0036.118E	Broca Ø 1,8x Stop 5mmx Comprimento 120mm
05.5.0037.118E	Broca Ø 1,8x Stop 5mmx Comprimento 125mm
05.5.0038.118E	Broca Ø 1,8x Stop 5mmx Comprimento 130mm
05.5.0039.118E	Broca Ø 1,8x Stop 5mmx Comprimento 135mm
05.5.0040.118E	Broca Ø 1,8x Stop 5mmx Comprimento 140mm
05.5.0041.118E	Broca Ø 1,8x Stop 5mmx Comprimento 145mm
05.5.0042.118E	Broca Ø 1,8x Stop 5mmx Comprimento 150mm
05.5.0043.118E	Broca Ø 1,8x Stop 6mmx Comprimento 50mm
05.5.0044.118E	Broca Ø 1,8x Stop 6mmx Comprimento 55mm
05.5.0045.118E	Broca Ø 1,8x Stop 6mmx Comprimento 60mm
05.5.0046.118E	Broca Ø 1,8x Stop 6mmx Comprimento 65mm
05.5.0047.118E	Broca Ø 1,8x Stop 6mmx Comprimento 70mm
05.5.0048.118E	Broca Ø 1,8x Stop 6mmx Comprimento 75mm
05.5.0049.118E	Broca Ø 1,8x Stop 6mmx Comprimento 80mm
05.5.0050.118E	Broca Ø 1,8x Stop 6mmx Comprimento 85mm
05.5.0051.118E	Broca Ø 1,8x Stop 6mmx Comprimento 90mm
05.5.0052.118E	Broca Ø 1,8x Stop 6mmx Comprimento 95mm
05.5.0053.118E	Broca Ø 1,8x Stop 6mmx Comprimento 100mm
05.5.0054.118E	Broca Ø 1,8x Stop 6mmx Comprimento 105mm
05.5.0055.118E	Broca Ø 1,8x Stop 6mmx Comprimento 110mm
05.5.0056.118E	Broca Ø 1,8x Stop 6mmx Comprimento 115mm
05.5.0057.118E	Broca Ø 1,8x Stop 6mmx Comprimento 120mm
05.5.0058.118E	Broca Ø 1,8x Stop 6mmx Comprimento 125mm
05.5.0059.118E	Broca Ø 1,8x Stop 6mmx Comprimento 130mm
05.5.0060.118E	Broca Ø 1,8x Stop 6mmx Comprimento 135mm
05.5.0061.118E	Broca Ø 1,8x Stop 6mmx Comprimento 140mm
05.5.0062.118E	Broca Ø 1,8x Stop 6mmx Comprimento 145mm
05.5.0063.118E	Broca Ø 1,8x Stop 6mmx Comprimento 150mm
05.5.0064.118E	Broca Ø 1,8x Stop 7mmx Comprimento 50mm
05.5.0065.118E	Broca Ø 1,8x Stop 7mmx Comprimento 55mm
05.5.0066.118E	Broca Ø 1,8x Stop 7mmx Comprimento 60mm

05.5.0067.118E	Broca Ø 1,8x Stop 7mmx Comprimento 65mm
05.5.0068.118E	Broca Ø 1,8x Stop 7mmx Comprimento 70mm
05.5.0069.118E	Broca Ø 1,8x Stop 7mmx Comprimento 75mm
05.5.0070.118E	Broca Ø 1,8x Stop 7mmx Comprimento 80mm
05.5.0071.118E	Broca Ø 1,8x Stop 7mmx Comprimento 85mm
05.5.0072.118E	Broca Ø 1,8x Stop 7mmx Comprimento 90mm
05.5.0073.118E	Broca Ø 1,8x Stop 7mmx Comprimento 95mm
05.5.0074.118E	Broca Ø 1,8x Stop 7mmx Comprimento 100mm
05.5.0075.118E	Broca Ø 1,8x Stop 7mmx Comprimento 105mm
05.5.0076.118E	Broca Ø 1,8x Stop 7mmx Comprimento 110mm
05.5.0077.118E	Broca Ø 1,8x Stop 7mmx Comprimento 115mm
05.5.0078.118E	Broca Ø 1,8x Stop 7mmx Comprimento 120mm
05.5.0079.118E	Broca Ø 1,8x Stop 7mmx Comprimento 125mm
05.5.0080.118E	Broca Ø 1,8x Stop 7mmx Comprimento 130mm
05.5.0081.118E	Broca Ø 1,8x Stop 7mmx Comprimento 135mm
05.5.0082.118E	Broca Ø 1,8x Stop 7mmx Comprimento 140mm
05.5.0083.118E	Broca Ø 1,8x Stop 7mmx Comprimento 145mm
05.5.0084.118E	Broca Ø 1,8x Stop 7mmx Comprimento 150mm
05.5.0085.118E	Broca Ø 1,8x Stop 8mmx Comprimento 50mm
05.5.0086.118E	Broca Ø 1,8x Stop 8mmx Comprimento 55mm
05.5.0087.118E	Broca Ø 1,8x Stop 8mmx Comprimento 60mm
05.5.0088.118E	Broca Ø 1,8x Stop 8mmx Comprimento 65mm
05.5.0089.118E	Broca Ø 1,8x Stop 8mmx Comprimento 70mm
05.5.0090.118E	Broca Ø 1,8x Stop 8mmx Comprimento 75mm
05.5.0091.118E	Broca Ø 1,8x Stop 8mmx Comprimento 80mm
05.5.0092.118E	Broca Ø 1,8x Stop 8mmx Comprimento 85mm
05.5.0093.118E	Broca Ø 1,8x Stop 8mmx Comprimento 90mm
05.5.0094.118E	Broca Ø 1,8x Stop 8mmx Comprimento 95mm
05.5.0095.118E	Broca Ø 1,8x Stop 8mmx Comprimento 100mm
05.5.0096.118E	Broca Ø 1,8x Stop 8mmx Comprimento 105mm
05.5.0097.118E	Broca Ø 1,8x Stop 8mmx Comprimento 110mm
05.5.0098.118E	Broca Ø 1,8x Stop 8mmx Comprimento 115mm
05.5.0099.118E	Broca Ø 1,8x Stop 8mmx Comprimento 120mm
05.5.0100.118E	Broca Ø 1,8x Stop 8mmx Comprimento 125mm
05.5.0101.118E	Broca Ø 1,8x Stop 8mmx Comprimento 130mm
05.5.0102.118E	Broca Ø 1,8x Stop 8mmx Comprimento 135mm
05.5.0103.118E	Broca Ø 1,8x Stop 8mmx Comprimento 140mm
05.5.0104.118E	Broca Ø 1,8x Stop 8mmx Comprimento 145mm
05.5.0105.118E	Broca Ø 1,8x Stop 8mmx Comprimento 150mm
05.5.0106.118E	Broca Ø 1,8x Stop 9mmx Comprimento 50mm
05.5.0107.118E	Broca Ø 1,8x Stop 9mmx Comprimento 55mm
05.5.0108.118E	Broca Ø 1,8x Stop 9mmx Comprimento 60mm
05.5.0109.118E	Broca Ø 1,8x Stop 9mmx Comprimento 65mm
05.5.0110.118E	Broca Ø 1,8x Stop 9mmx Comprimento 70mm
05.5.0111.118E	Broca Ø 1,8x Stop 9mmx Comprimento 75mm
05.5.0112.118E	Broca Ø 1,8x Stop 9mmx Comprimento 80mm
05.5.0113.118E	Broca Ø 1,8x Stop 9mmx Comprimento 85mm
05.5.0114.118E	Broca Ø 1,8x Stop 9mmx Comprimento 90mm
05.5.0115.118E	Broca Ø 1,8x Stop 9mmx Comprimento 95mm

05.5.0116.118E	Broca Ø 1,8x Stop 9mmx Comprimento 100mm
05.5.0117.118E	Broca Ø 1,8x Stop 9mmx Comprimento 105mm
05.5.0118.118E	Broca Ø 1,8x Stop 9mmx Comprimento 110mm
05.5.0119.118E	Broca Ø 1,8x Stop 9mmx Comprimento 115mm
05.5.0120.118E	Broca Ø 1,8x Stop 9mmx Comprimento 120mm
05.5.0121.118E	Broca Ø 1,8x Stop 9mmx Comprimento 125mm
05.5.0122.118E	Broca Ø 1,8x Stop 9mmx Comprimento 130mm
05.5.0123.118E	Broca Ø 1,8x Stop 9mmx Comprimento 135mm
05.5.0124.118E	Broca Ø 1,8x Stop 9mmx Comprimento 140mm
05.5.0125.118E	Broca Ø 1,8x Stop 9mmx Comprimento 145mm
05.5.0126.118E	Broca Ø 1,8x Stop 9mmx Comprimento 150mm
05.5.0127.118E	Broca Ø 1,8x Stop 10mmx Comprimento 50mm
05.5.0128.118E	Broca Ø 1,8x Stop 10mmx Comprimento 55mm
05.5.0129.118E	Broca Ø 1,8x Stop 10mmx Comprimento 60mm
05.5.0130.118E	Broca Ø 1,8x Stop 10mmx Comprimento 65mm
05.5.0131.118E	Broca Ø 1,8x Stop 10mmx Comprimento 70mm
05.5.0132.118E	Broca Ø 1,8x Stop 10mmx Comprimento 75mm
05.5.0133.118E	Broca Ø 1,8x Stop 10mmx Comprimento 80mm
05.5.0134.118E	Broca Ø 1,8x Stop 10mmx Comprimento 85mm
05.5.0135.118E	Broca Ø 1,8x Stop 10mmx Comprimento 90mm
05.5.0136.118E	Broca Ø 1,8x Stop 10mmx Comprimento 95mm
05.5.0137.118E	Broca Ø 1,8x Stop 10mmx Comprimento 100mm
05.5.0138.118E	Broca Ø 1,8x Stop 10mmx Comprimento 105mm
05.5.0139.118E	Broca Ø 1,8x Stop 10mmx Comprimento 110mm
05.5.0140.118E	Broca Ø 1,8x Stop 10mmx Comprimento 115mm
05.5.0141.118E	Broca Ø 1,8x Stop 10mmx Comprimento 120mm
05.5.0142.118E	Broca Ø 1,8x Stop 10mmx Comprimento 125mm
05.5.0143.118E	Broca Ø 1,8x Stop 10mmx Comprimento 130mm
05.5.0144.118E	Broca Ø 1,8x Stop 10mmx Comprimento 135mm
05.5.0145.118E	Broca Ø 1,8x Stop 10mmx Comprimento 140mm
05.5.0146.118E	Broca Ø 1,8x Stop 10mmx Comprimento 145mm
05.5.0147.118E	Broca Ø 1,8x Stop 10mmx Comprimento 150mm
05.5.0148.118E	Broca Ø 1,8x Stop 11mmx Comprimento 50mm
05.5.0149.118E	Broca Ø 1,8x Stop 11mmx Comprimento 55mm
05.5.0150.118E	Broca Ø 1,8x Stop 11mmx Comprimento 60mm
05.5.0151.118E	Broca Ø 1,8x Stop 11mmx Comprimento 65mm
05.5.0152.118E	Broca Ø 1,8x Stop 11mmx Comprimento 70mm
05.5.0153.118E	Broca Ø 1,8x Stop 11mmx Comprimento 75mm
05.5.0154.118E	Broca Ø 1,8x Stop 11mmx Comprimento 80mm
05.5.0155.118E	Broca Ø 1,8x Stop 11mmx Comprimento 85mm
05.5.0156.118E	Broca Ø 1,8x Stop 11mmx Comprimento 90mm
05.5.0157.118E	Broca Ø 1,8x Stop 11mmx Comprimento 95mm
05.5.0158.118E	Broca Ø 1,8x Stop 11mmx Comprimento 100mm
05.5.0159.118E	Broca Ø 1,8x Stop 11mmx Comprimento 105mm
05.5.0160.118E	Broca Ø 1,8x Stop 11mmx Comprimento 110mm
05.5.0161.118E	Broca Ø 1,8x Stop 11mmx Comprimento 115mm
05.5.0162.118E	Broca Ø 1,8x Stop 11mmx Comprimento 120mm
05.5.0163.118E	Broca Ø 1,8x Stop 11mmx Comprimento 125mm
05.5.0164.118E	Broca Ø 1,8x Stop 11mmx Comprimento 130mm

05.5.0165.118E	Broca Ø 1,8x Stop 11mmx Comprimento 135mm
05.5.0166.118E	Broca Ø 1,8x Stop 11mmx Comprimento 140mm
05.5.0167.118E	Broca Ø 1,8x Stop 11mmx Comprimento 145mm
05.5.0168.118E	Broca Ø 1,8x Stop 11mmx Comprimento 150mm
05.5.0169.118E	Broca Ø 1,8x Stop 12mmx Comprimento 50mm
05.5.0170.118E	Broca Ø 1,8x Stop 12mmx Comprimento 55mm
05.5.0171.118E	Broca Ø 1,8x Stop 12mmx Comprimento 60mm
05.5.0172.118E	Broca Ø 1,8x Stop 12mmx Comprimento 65mm
05.5.0173.118E	Broca Ø 1,8x Stop 12mmx Comprimento 70mm
05.5.0174.118E	Broca Ø 1,8x Stop 12mmx Comprimento 75mm
05.5.0175.118E	Broca Ø 1,8x Stop 12mmx Comprimento 80mm
05.5.0176.118E	Broca Ø 1,8x Stop 12mmx Comprimento 85mm
05.5.0177.118E	Broca Ø 1,8x Stop 12mmx Comprimento 90mm
05.5.0178.118E	Broca Ø 1,8x Stop 12mmx Comprimento 95mm
05.5.0179.118E	Broca Ø 1,8x Stop 12mmx Comprimento 100mm
05.5.0180.118E	Broca Ø 1,8x Stop 12mmx Comprimento 105mm
05.5.0181.118E	Broca Ø 1,8x Stop 12mmx Comprimento 110mm
05.5.0182.118E	Broca Ø 1,8x Stop 12mmx Comprimento 115mm
05.5.0183.118E	Broca Ø 1,8x Stop 12mmx Comprimento 120mm
05.5.0184.118E	Broca Ø 1,8x Stop 12mmx Comprimento 125mm
05.5.0185.118E	Broca Ø 1,8x Stop 12mmx Comprimento 130mm
05.5.0186.118E	Broca Ø 1,8x Stop 12mmx Comprimento 135mm
05.5.0187.118E	Broca Ø 1,8x Stop 12mmx Comprimento 140mm
05.5.0188.118E	Broca Ø 1,8x Stop 12mmx Comprimento 145mm
05.5.0189.118E	Broca Ø 1,8x Stop 12mmx Comprimento 150mm
05.5.0190.118E	Broca Ø 1,8x Stop 13mmx Comprimento 50mm
05.5.0191.118E	Broca Ø 1,8x Stop 13mmx Comprimento 55mm
05.5.0192.118E	Broca Ø 1,8x Stop 13mmx Comprimento 60mm
05.5.0193.118E	Broca Ø 1,8x Stop 13mmx Comprimento 65mm
05.5.0194.118E	Broca Ø 1,8x Stop 13mmx Comprimento 70mm
05.5.0195.118E	Broca Ø 1,8x Stop 13mmx Comprimento 75mm
05.5.0196.118E	Broca Ø 1,8x Stop 13mmx Comprimento 80mm
05.5.0197.118E	Broca Ø 1,8x Stop 13mmx Comprimento 85mm
05.5.0198.118E	Broca Ø 1,8x Stop 13mmx Comprimento 90mm
05.5.0199.118E	Broca Ø 1,8x Stop 13mmx Comprimento 95mm
05.5.0200.118E	Broca Ø 1,8x Stop 13mmx Comprimento 100mm
05.5.0201.118E	Broca Ø 1,8x Stop 13mmx Comprimento 105mm
05.5.0202.118E	Broca Ø 1,8x Stop 13mmx Comprimento 110mm
05.5.0203.118E	Broca Ø 1,8x Stop 13mmx Comprimento 115mm
05.5.0204.118E	Broca Ø 1,8x Stop 13mmx Comprimento 120mm
05.5.0205.118E	Broca Ø 1,8x Stop 13mmx Comprimento 125mm
05.5.0206.118E	Broca Ø 1,8x Stop 13mmx Comprimento 130mm
05.5.0207.118E	Broca Ø 1,8x Stop 13mmx Comprimento 135mm
05.5.0208.118E	Broca Ø 1,8x Stop 13mmx Comprimento 140mm
05.5.0209.118E	Broca Ø 1,8x Stop 13mmx Comprimento 145mm
05.5.0210.118E	Broca Ø 1,8x Stop 13mmx Comprimento 150mm
05.5.0211.118E	Broca Ø 1,8x Stop 14mmx Comprimento 50mm
05.5.0212.118E	Broca Ø 1,8x Stop 14mmx Comprimento 55mm
05.5.0213.118E	Broca Ø 1,8x Stop 14mmx Comprimento 60mm

05.5.0214.118E	Broca Ø 1,8x Stop 14mmx Comprimento 65mm
05.5.0215.118E	Broca Ø 1,8x Stop 14mmx Comprimento 70mm
05.5.0216.118E	Broca Ø 1,8x Stop 14mmx Comprimento 75mm
05.5.0217.118E	Broca Ø 1,8x Stop 14mmx Comprimento 80mm
05.5.0218.118E	Broca Ø 1,8x Stop 14mmx Comprimento 85mm
05.5.0219.118E	Broca Ø 1,8x Stop 14mmx Comprimento 90mm
05.5.0220.118E	Broca Ø 1,8x Stop 14mmx Comprimento 95mm
05.5.0221.118E	Broca Ø 1,8x Stop 14mmx Comprimento 100mm
05.5.0222.118E	Broca Ø 1,8x Stop 14mmx Comprimento 105mm
05.5.0223.118E	Broca Ø 1,8x Stop 14mmx Comprimento 110mm
05.5.0224.118E	Broca Ø 1,8x Stop 14mmx Comprimento 115mm
05.5.0225.118E	Broca Ø 1,8x Stop 14mmx Comprimento 120mm
05.5.0226.118E	Broca Ø 1,8x Stop 14mmx Comprimento 125mm
05.5.0227.118E	Broca Ø 1,8x Stop 14mmx Comprimento 130mm
05.5.0228.118E	Broca Ø 1,8x Stop 14mmx Comprimento 135mm
05.5.0229.118E	Broca Ø 1,8x Stop 14mmx Comprimento 140mm
05.5.0230.118E	Broca Ø 1,8x Stop 14mmx Comprimento 145mm
05.5.0231.118E	Broca Ø 1,8x Stop 14mmx Comprimento 150mm
05.5.0232.118E	Broca Ø 1,8x Stop 15mmx Comprimento 50mm
05.5.0233.118E	Broca Ø 1,8x Stop 15mmx Comprimento 55mm
05.5.0234.118E	Broca Ø 1,8x Stop 15mmx Comprimento 60mm
05.5.0235.118E	Broca Ø 1,8x Stop 15mmx Comprimento 65mm
05.5.0236.118E	Broca Ø 1,8x Stop 15mmx Comprimento 70mm
05.5.0237.118E	Broca Ø 1,8x Stop 15mmx Comprimento 75mm
05.5.0238.118E	Broca Ø 1,8x Stop 15mmx Comprimento 80mm
05.5.0239.118E	Broca Ø 1,8x Stop 15mmx Comprimento 85mm
05.5.0240.118E	Broca Ø 1,8x Stop 15mmx Comprimento 90mm
05.5.0241.118E	Broca Ø 1,8x Stop 15mmx Comprimento 95mm
05.5.0242.118E	Broca Ø 1,8x Stop 15mmx Comprimento 100mm
05.5.0243.118E	Broca Ø 1,8x Stop 15mmx Comprimento 105mm
05.5.0244.118E	Broca Ø 1,8x Stop 15mmx Comprimento 110mm
05.5.0245.118E	Broca Ø 1,8x Stop 15mmx Comprimento 115mm
05.5.0246.118E	Broca Ø 1,8x Stop 15mmx Comprimento 120mm
05.5.0247.118E	Broca Ø 1,8x Stop 15mmx Comprimento 125mm
05.5.0248.118E	Broca Ø 1,8x Stop 15mmx Comprimento 130mm
05.5.0249.118E	Broca Ø 1,8x Stop 15mmx Comprimento 135mm
05.5.0250.118E	Broca Ø 1,8x Stop 15mmx Comprimento 140mm
05.5.0251.118E	Broca Ø 1,8x Stop 15mmx Comprimento 145mm
05.5.0252.118E	Broca Ø 1,8x Stop 15mmx Comprimento 150mm
05.5.0253.118E	Broca Ø 1,8x Stop 16mmx Comprimento 50mm
05.5.0254.118E	Broca Ø 1,8x Stop 16mmx Comprimento 55mm
05.5.0255.118E	Broca Ø 1,8x Stop 16mmx Comprimento 60mm
05.5.0256.118E	Broca Ø 1,8x Stop 16mmx Comprimento 65mm
05.5.0257.118E	Broca Ø 1,8x Stop 16mmx Comprimento 70mm
05.5.0258.118E	Broca Ø 1,8x Stop 16mmx Comprimento 75mm
05.5.0259.118E	Broca Ø 1,8x Stop 16mmx Comprimento 80mm
05.5.0260.118E	Broca Ø 1,8x Stop 16mmx Comprimento 85mm
05.5.0261.118E	Broca Ø 1,8x Stop 16mmx Comprimento 90mm
05.5.0262.118E	Broca Ø 1,8x Stop 16mmx Comprimento 95mm

05.5.0263.118E	Broca Ø 1,8x Stop 16mmx Comprimento 100mm
05.5.0264.118E	Broca Ø 1,8x Stop 16mmx Comprimento 105mm
05.5.0265.118E	Broca Ø 1,8x Stop 16mmx Comprimento 110mm
05.5.0266.118E	Broca Ø 1,8x Stop 16mmx Comprimento 115mm
05.5.0267.118E	Broca Ø 1,8x Stop 16mmx Comprimento 120mm
05.5.0268.118E	Broca Ø 1,8x Stop 16mmx Comprimento 125mm
05.5.0269.118E	Broca Ø 1,8x Stop 16mmx Comprimento 130mm
05.5.0270.118E	Broca Ø 1,8x Stop 16mmx Comprimento 135mm
05.5.0271.118E	Broca Ø 1,8x Stop 16mmx Comprimento 140mm
05.5.0272.118E	Broca Ø 1,8x Stop 16mmx Comprimento 145mm
05.5.0273.118E	Broca Ø 1,8x Stop 16mmx Comprimento 150mm
05.5.0274.118E	Broca Ø 1,8x Stop 17mmx Comprimento 50mm
05.5.0275.118E	Broca Ø 1,8x Stop 17mmx Comprimento 55mm
05.5.0276.118E	Broca Ø 1,8x Stop 17mmx Comprimento 60mm
05.5.0277.118E	Broca Ø 1,8x Stop 17mmx Comprimento 65mm
05.5.0278.118E	Broca Ø 1,8x Stop 17mmx Comprimento 70mm
05.5.0279.118E	Broca Ø 1,8x Stop 17mmx Comprimento 75mm
05.5.0280.118E	Broca Ø 1,8x Stop 17mmx Comprimento 80mm
05.5.0281.118E	Broca Ø 1,8x Stop 17mmx Comprimento 85mm
05.5.0282.118E	Broca Ø 1,8x Stop 17mmx Comprimento 90mm
05.5.0283.118E	Broca Ø 1,8x Stop 17mmx Comprimento 95mm
05.5.0284.118E	Broca Ø 1,8x Stop 17mmx Comprimento 100mm
05.5.0285.118E	Broca Ø 1,8x Stop 17mmx Comprimento 105mm
05.5.0286.118E	Broca Ø 1,8x Stop 17mmx Comprimento 110mm
05.5.0287.118E	Broca Ø 1,8x Stop 17mmx Comprimento 115mm
05.5.0288.118E	Broca Ø 1,8x Stop 17mmx Comprimento 120mm
05.5.0289.118E	Broca Ø 1,8x Stop 17mmx Comprimento 125mm
05.5.0290.118E	Broca Ø 1,8x Stop 17mmx Comprimento 130mm
05.5.0291.118E	Broca Ø 1,8x Stop 17mmx Comprimento 135mm
05.5.0292.118E	Broca Ø 1,8x Stop 17mmx Comprimento 140mm
05.5.0293.118E	Broca Ø 1,8x Stop 17mmx Comprimento 145mm
05.5.0294.118E	Broca Ø 1,8x Stop 17mmx Comprimento 150mm
05.5.0295.118E	Broca Ø 1,8x Stop 18mmx Comprimento 50mm
05.5.0296.118E	Broca Ø 1,8x Stop 18mmx Comprimento 55mm
05.5.0297.118E	Broca Ø 1,8x Stop 18mmx Comprimento 60mm
05.5.0298.118E	Broca Ø 1,8x Stop 18mmx Comprimento 65mm
05.5.0299.118E	Broca Ø 1,8x Stop 18mmx Comprimento 70mm
05.5.0300.118E	Broca Ø 1,8x Stop 18mmx Comprimento 75mm
05.5.0301.118E	Broca Ø 1,8x Stop 18mmx Comprimento 80mm
05.5.0302.118E	Broca Ø 1,8x Stop 18mmx Comprimento 85mm
05.5.0303.118E	Broca Ø 1,8x Stop 18mmx Comprimento 90mm
05.5.0304.118E	Broca Ø 1,8x Stop 18mmx Comprimento 95mm
05.5.0305.118E	Broca Ø 1,8x Stop 18mmx Comprimento 100mm
05.5.0306.118E	Broca Ø 1,8x Stop 18mmx Comprimento 105mm
05.5.0307.118E	Broca Ø 1,8x Stop 18mmx Comprimento 110mm
05.5.0308.118E	Broca Ø 1,8x Stop 18mmx Comprimento 115mm
05.5.0309.118E	Broca Ø 1,8x Stop 18mmx Comprimento 120mm
05.5.0310.118E	Broca Ø 1,8x Stop 18mmx Comprimento 125mm
05.5.0311.118E	Broca Ø 1,8x Stop 18mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0312.118E	Broca Ø 1,8x Stop 18mmx Comprimento 135mm
05.5.0313.118E	Broca Ø 1,8x Stop 18mmx Comprimento 140mm
05.5.0314.118E	Broca Ø 1,8x Stop 18mmx Comprimento 145mm
05.5.0315.118E	Broca Ø 1,8x Stop 18mmx Comprimento 150mm
05.5.0316.118E	Broca Ø 1,8x Stop 19mmx Comprimento 50mm
05.5.0317.118E	Broca Ø 1,8x Stop 19mmx Comprimento 55mm
05.5.0318.118E	Broca Ø 1,8x Stop 19mmx Comprimento 60mm
05.5.0319.118E	Broca Ø 1,8x Stop 19mmx Comprimento 65mm
05.5.0320.118E	Broca Ø 1,8x Stop 19mmx Comprimento 70mm
05.5.0321.118E	Broca Ø 1,8x Stop 19mmx Comprimento 75mm
05.5.0322.118E	Broca Ø 1,8x Stop 19mmx Comprimento 80mm
05.5.0323.118E	Broca Ø 1,8x Stop 19mmx Comprimento 85mm
05.5.0324.118E	Broca Ø 1,8x Stop 19mmx Comprimento 90mm
05.5.0325.118E	Broca Ø 1,8x Stop 19mmx Comprimento 95mm
05.5.0326.118E	Broca Ø 1,8x Stop 19mmx Comprimento 100mm
05.5.0327.118E	Broca Ø 1,8x Stop 19mmx Comprimento 105mm
05.5.0328.118E	Broca Ø 1,8x Stop 19mmx Comprimento 110mm
05.5.0329.118E	Broca Ø 1,8x Stop 19mmx Comprimento 115mm
05.5.0330.118E	Broca Ø 1,8x Stop 19mmx Comprimento 120mm
05.5.0331.118E	Broca Ø 1,8x Stop 19mmx Comprimento 125mm
05.5.0332.118E	Broca Ø 1,8x Stop 19mmx Comprimento 130mm
05.5.0333.118E	Broca Ø 1,8x Stop 19mmx Comprimento 135mm
05.5.0334.118E	Broca Ø 1,8x Stop 19mmx Comprimento 140mm
05.5.0335.118E	Broca Ø 1,8x Stop 19mmx Comprimento 145mm
05.5.0336.118E	Broca Ø 1,8x Stop 19mmx Comprimento 150mm
05.5.0337.118E	Broca Ø 1,8x Stop 20mmx Comprimento 50mm
05.5.0338.118E	Broca Ø 1,8x Stop 20mmx Comprimento 55mm
05.5.0339.118E	Broca Ø 1,8x Stop 20mmx Comprimento 60mm
05.5.0340.118E	Broca Ø 1,8x Stop 20mmx Comprimento 65mm
05.5.0341.118E	Broca Ø 1,8x Stop 20mmx Comprimento 70mm
05.5.0342.118E	Broca Ø 1,8x Stop 20mmx Comprimento 75mm
05.5.0343.118E	Broca Ø 1,8x Stop 20mmx Comprimento 80mm
05.5.0344.118E	Broca Ø 1,8x Stop 20mmx Comprimento 85mm
05.5.0345.118E	Broca Ø 1,8x Stop 20mmx Comprimento 90mm
05.5.0346.118E	Broca Ø 1,8x Stop 20mmx Comprimento 95mm
05.5.0347.118E	Broca Ø 1,8x Stop 20mmx Comprimento 100mm
05.5.0348.118E	Broca Ø 1,8x Stop 20mmx Comprimento 105mm
05.5.0349.118E	Broca Ø 1,8x Stop 20mmx Comprimento 110mm
05.5.0350.118E	Broca Ø 1,8x Stop 20mmx Comprimento 115mm
05.5.0351.118E	Broca Ø 1,8x Stop 20mmx Comprimento 120mm
05.5.0352.118E	Broca Ø 1,8x Stop 20mmx Comprimento 125mm
05.5.0353.118E	Broca Ø 1,8x Stop 20mmx Comprimento 130mm
05.5.0354.118E	Broca Ø 1,8x Stop 20mmx Comprimento 135mm
05.5.0355.118E	Broca Ø 1,8x Stop 20mmx Comprimento 140mm
05.5.0356.118E	Broca Ø 1,8x Stop 20mmx Comprimento 145mm
05.5.0357.118E	Broca Ø 1,8x Stop 20mmx Comprimento 150mm
05.5.0358.118E	Broca Ø 1,8x Stop 21mmx Comprimento 50mm
05.5.0359.118E	Broca Ø 1,8x Stop 21mmx Comprimento 55mm
05.5.0360.118E	Broca Ø 1,8x Stop 21mmx Comprimento 60mm

05.5.0361.118E	Broca Ø 1,8x Stop 21mmx Comprimento 65mm
05.5.0362.118E	Broca Ø 1,8x Stop 21mmx Comprimento 70mm
05.5.0363.118E	Broca Ø 1,8x Stop 21mmx Comprimento 75mm
05.5.0364.118E	Broca Ø 1,8x Stop 21mmx Comprimento 80mm
05.5.0365.118E	Broca Ø 1,8x Stop 21mmx Comprimento 85mm
05.5.0366.118E	Broca Ø 1,8x Stop 21mmx Comprimento 90mm
05.5.0367.118E	Broca Ø 1,8x Stop 21mmx Comprimento 95mm
05.5.0368.118E	Broca Ø 1,8x Stop 21mmx Comprimento 100mm
05.5.0369.118E	Broca Ø 1,8x Stop 21mmx Comprimento 105mm
05.5.0370.118E	Broca Ø 1,8x Stop 21mmx Comprimento 110mm
05.5.0371.118E	Broca Ø 1,8x Stop 21mmx Comprimento 115mm
05.5.0372.118E	Broca Ø 1,8x Stop 21mmx Comprimento 120mm
05.5.0373.118E	Broca Ø 1,8x Stop 21mmx Comprimento 125mm
05.5.0374.118E	Broca Ø 1,8x Stop 21mmx Comprimento 130mm
05.5.0375.118E	Broca Ø 1,8x Stop 21mmx Comprimento 135mm
05.5.0376.118E	Broca Ø 1,8x Stop 21mmx Comprimento 140mm
05.5.0377.118E	Broca Ø 1,8x Stop 21mmx Comprimento 145mm
05.5.0378.118E	Broca Ø 1,8x Stop 21mmx Comprimento 150mm
05.5.0379.118E	Broca Ø 1,8x Stop 22mmx Comprimento 50mm
05.5.0380.118E	Broca Ø 1,8x Stop 22mmx Comprimento 55mm
05.5.0381.118E	Broca Ø 1,8x Stop 22mmx Comprimento 60mm
05.5.0382.118E	Broca Ø 1,8x Stop 22mmx Comprimento 65mm
05.5.0383.118E	Broca Ø 1,8x Stop 22mmx Comprimento 70mm
05.5.0384.118E	Broca Ø 1,8x Stop 22mmx Comprimento 75mm
05.5.0385.118E	Broca Ø 1,8x Stop 22mmx Comprimento 80mm
05.5.0386.118E	Broca Ø 1,8x Stop 22mmx Comprimento 85mm
05.5.0387.118E	Broca Ø 1,8x Stop 22mmx Comprimento 90mm
05.5.0388.118E	Broca Ø 1,8x Stop 22mmx Comprimento 95mm
05.5.0389.118E	Broca Ø 1,8x Stop 22mmx Comprimento 100mm
05.5.0390.118E	Broca Ø 1,8x Stop 22mmx Comprimento 105mm
05.5.0391.118E	Broca Ø 1,8x Stop 22mmx Comprimento 110mm
05.5.0392.118E	Broca Ø 1,8x Stop 22mmx Comprimento 115mm
05.5.0393.118E	Broca Ø 1,8x Stop 22mmx Comprimento 120mm
05.5.0394.118E	Broca Ø 1,8x Stop 22mmx Comprimento 125mm
05.5.0395.118E	Broca Ø 1,8x Stop 22mmx Comprimento 130mm
05.5.0396.118E	Broca Ø 1,8x Stop 22mmx Comprimento 135mm
05.5.0397.118E	Broca Ø 1,8x Stop 22mmx Comprimento 140mm
05.5.0398.118E	Broca Ø 1,8x Stop 22mmx Comprimento 145mm
05.5.0399.118E	Broca Ø 1,8x Stop 22mmx Comprimento 150mm
05.5.0400.118E	Broca Ø 1,8x Stop 23mmx Comprimento 50mm
05.5.0401.118E	Broca Ø 1,8x Stop 23mmx Comprimento 55mm
05.5.0402.118E	Broca Ø 1,8x Stop 23mmx Comprimento 60mm
05.5.0403.118E	Broca Ø 1,8x Stop 23mmx Comprimento 65mm
05.5.0404.118E	Broca Ø 1,8x Stop 23mmx Comprimento 70mm
05.5.0405.118E	Broca Ø 1,8x Stop 23mmx Comprimento 75mm
05.5.0406.118E	Broca Ø 1,8x Stop 23mmx Comprimento 80mm
05.5.0407.118E	Broca Ø 1,8x Stop 23mmx Comprimento 85mm
05.5.0408.118E	Broca Ø 1,8x Stop 23mmx Comprimento 90mm
05.5.0409.118E	Broca Ø 1,8x Stop 23mmx Comprimento 95mm

05.5.0410.118E	Broca Ø 1,8x Stop 23mmx Comprimento 100mm
05.5.0411.118E	Broca Ø 1,8x Stop 23mmx Comprimento 105mm
05.5.0412.118E	Broca Ø 1,8x Stop 23mmx Comprimento 110mm
05.5.0413.118E	Broca Ø 1,8x Stop 23mmx Comprimento 115mm
05.5.0414.118E	Broca Ø 1,8x Stop 23mmx Comprimento 120mm
05.5.0415.118E	Broca Ø 1,8x Stop 23mmx Comprimento 125mm
05.5.0416.118E	Broca Ø 1,8x Stop 23mmx Comprimento 130mm
05.5.0417.118E	Broca Ø 1,8x Stop 23mmx Comprimento 135mm
05.5.0418.118E	Broca Ø 1,8x Stop 23mmx Comprimento 140mm
05.5.0419.118E	Broca Ø 1,8x Stop 23mmx Comprimento 145mm
05.5.0420.118E	Broca Ø 1,8x Stop 23mmx Comprimento 150mm
05.5.0421.118E	Broca Ø 1,8x Stop 24mmx Comprimento 50mm
05.5.0422.118E	Broca Ø 1,8x Stop 24mmx Comprimento 55mm
05.5.0423.118E	Broca Ø 1,8x Stop 24mmx Comprimento 60mm
05.5.0424.118E	Broca Ø 1,8x Stop 24mmx Comprimento 65mm
05.5.0425.118E	Broca Ø 1,8x Stop 24mmx Comprimento 70mm
05.5.0426.118E	Broca Ø 1,8x Stop 24mmx Comprimento 75mm
05.5.0427.118E	Broca Ø 1,8x Stop 24mmx Comprimento 80mm
05.5.0428.118E	Broca Ø 1,8x Stop 24mmx Comprimento 85mm
05.5.0429.118E	Broca Ø 1,8x Stop 24mmx Comprimento 90mm
05.5.0430.118E	Broca Ø 1,8x Stop 24mmx Comprimento 95mm
05.5.0431.118E	Broca Ø 1,8x Stop 24mmx Comprimento 100mm
05.5.0432.118E	Broca Ø 1,8x Stop 24mmx Comprimento 105mm
05.5.0433.118E	Broca Ø 1,8x Stop 24mmx Comprimento 110mm
05.5.0434.118E	Broca Ø 1,8x Stop 24mmx Comprimento 115mm
05.5.0435.118E	Broca Ø 1,8x Stop 24mmx Comprimento 120mm
05.5.0436.118E	Broca Ø 1,8x Stop 24mmx Comprimento 125mm
05.5.0437.118E	Broca Ø 1,8x Stop 24mmx Comprimento 130mm
05.5.0438.118E	Broca Ø 1,8x Stop 24mmx Comprimento 135mm
05.5.0439.118E	Broca Ø 1,8x Stop 24mmx Comprimento 140mm
05.5.0440.118E	Broca Ø 1,8x Stop 24mmx Comprimento 145mm
05.5.0441.118E	Broca Ø 1,8x Stop 24mmx Comprimento 150mm
05.5.0442.118E	Broca Ø 1,8x Stop 25mmx Comprimento 50mm
05.5.0443.118E	Broca Ø 1,8x Stop 25mmx Comprimento 55mm
05.5.0444.118E	Broca Ø 1,8x Stop 25mmx Comprimento 60mm
05.5.0445.118E	Broca Ø 1,8x Stop 25mmx Comprimento 65mm
05.5.0446.118E	Broca Ø 1,8x Stop 25mmx Comprimento 70mm
05.5.0447.118E	Broca Ø 1,8x Stop 25mmx Comprimento 75mm
05.5.0448.118E	Broca Ø 1,8x Stop 25mmx Comprimento 80mm
05.5.0449.118E	Broca Ø 1,8x Stop 25mmx Comprimento 85mm
05.5.0450.118E	Broca Ø 1,8x Stop 25mmx Comprimento 90mm
05.5.0451.118E	Broca Ø 1,8x Stop 25mmx Comprimento 95mm
05.5.0452.118E	Broca Ø 1,8x Stop 25mmx Comprimento 100mm
05.5.0453.118E	Broca Ø 1,8x Stop 25mmx Comprimento 105mm
05.5.0454.118E	Broca Ø 1,8x Stop 25mmx Comprimento 110mm
05.5.0455.118E	Broca Ø 1,8x Stop 25mmx Comprimento 115mm
05.5.0456.118E	Broca Ø 1,8x Stop 25mmx Comprimento 120mm
05.5.0457.118E	Broca Ø 1,8x Stop 25mmx Comprimento 125mm
05.5.0458.118E	Broca Ø 1,8x Stop 25mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0459.118E	Broca Ø 1,8x Stop 25mmx Comprimento 135mm
05.5.0460.118E	Broca Ø 1,8x Stop 25mmx Comprimento 140mm
05.5.0461.118E	Broca Ø 1,8x Stop 25mmx Comprimento 145mm
05.5.0462.118E	Broca Ø 1,8x Stop 25mmx Comprimento 150mm
05.5.0463.118E	Broca Ø 1,8x Stop 30mmx Comprimento 50mm
05.5.0464.118E	Broca Ø 1,8x Stop 30mmx Comprimento 55mm
05.5.0465.118E	Broca Ø 1,8x Stop 30mmx Comprimento 60mm
05.5.0466.118E	Broca Ø 1,8x Stop 30mmx Comprimento 65mm
05.5.0467.118E	Broca Ø 1,8x Stop 30mmx Comprimento 70mm
05.5.0468.118E	Broca Ø 1,8x Stop 30mmx Comprimento 75mm
05.5.0469.118E	Broca Ø 1,8x Stop 30mmx Comprimento 80mm
05.5.0470.118E	Broca Ø 1,8x Stop 30mmx Comprimento 85mm
05.5.0471.118E	Broca Ø 1,8x Stop 30mmx Comprimento 90mm
05.5.0472.118E	Broca Ø 1,8x Stop 30mmx Comprimento 95mm
05.5.0473.118E	Broca Ø 1,8x Stop 30mmx Comprimento 100mm
05.5.0474.118E	Broca Ø 1,8x Stop 30mmx Comprimento 105mm
05.5.0475.118E	Broca Ø 1,8x Stop 30mmx Comprimento 110mm
05.5.0476.118E	Broca Ø 1,8x Stop 30mmx Comprimento 115mm
05.5.0477.118E	Broca Ø 1,8x Stop 30mmx Comprimento 120mm
05.5.0478.118E	Broca Ø 1,8x Stop 30mmx Comprimento 125mm
05.5.0479.118E	Broca Ø 1,8x Stop 30mmx Comprimento 130mm
05.5.0480.118E	Broca Ø 1,8x Stop 30mmx Comprimento 135mm
05.5.0481.118E	Broca Ø 1,8x Stop 30mmx Comprimento 140mm
05.5.0482.118E	Broca Ø 1,8x Stop 30mmx Comprimento 145mm
05.5.0483.118E	Broca Ø 1,8x Stop 30mmx Comprimento 150mm
05.5.0484.118E	Broca Ø 1,8x Stop 35mmx Comprimento 50mm
05.5.0485.118E	Broca Ø 1,8x Stop 35mmx Comprimento 55mm
05.5.0486.118E	Broca Ø 1,8x Stop 35mmx Comprimento 60mm
05.5.0487.118E	Broca Ø 1,8x Stop 35mmx Comprimento 65mm
05.5.0488.118E	Broca Ø 1,8x Stop 35mmx Comprimento 70mm
05.5.0489.118E	Broca Ø 1,8x Stop 35mmx Comprimento 75mm
05.5.0490.118E	Broca Ø 1,8x Stop 35mmx Comprimento 80mm
05.5.0491.118E	Broca Ø 1,8x Stop 35mmx Comprimento 85mm
05.5.0492.118E	Broca Ø 1,8x Stop 35mmx Comprimento 90mm
05.5.0493.118E	Broca Ø 1,8x Stop 35mmx Comprimento 95mm
05.5.0494.118E	Broca Ø 1,8x Stop 35mmx Comprimento 100mm
05.5.0495.118E	Broca Ø 1,8x Stop 35mmx Comprimento 105mm
05.5.0496.118E	Broca Ø 1,8x Stop 35mmx Comprimento 110mm
05.5.0497.118E	Broca Ø 1,8x Stop 35mmx Comprimento 115mm
05.5.0498.118E	Broca Ø 1,8x Stop 35mmx Comprimento 120mm
05.5.0499.118E	Broca Ø 1,8x Stop 35mmx Comprimento 125mm
05.5.0500.118E	Broca Ø 1,8x Stop 35mmx Comprimento 130mm
05.5.0501.118E	Broca Ø 1,8x Stop 35mmx Comprimento 135mm
05.5.0502.118E	Broca Ø 1,8x Stop 35mmx Comprimento 140mm
05.5.0503.118E	Broca Ø 1,8x Stop 35mmx Comprimento 145mm
05.5.0504.118E	Broca Ø 1,8x Stop 35mmx Comprimento 150mm
05.5.0001.119E	Broca Ø 1,9x Stop 4mmx Comprimento 50mm
05.5.0002.119E	Broca Ø 1,9x Stop 4mmx Comprimento 55mm
05.5.0003.119E	Broca Ø 1,9x Stop 4mmx Comprimento 60mm

05.5.0004.119E	Broca Ø 1,9x Stop 4mmx Comprimento 65mm
05.5.0005.119E	Broca Ø 1,9x Stop 4mmx Comprimento 70mm
05.5.0006.119E	Broca Ø 1,9x Stop 4mmx Comprimento 75mm
05.5.0007.119E	Broca Ø 1,9x Stop 4mmx Comprimento 80mm
05.5.0008.119E	Broca Ø 1,9x Stop 4mmx Comprimento 85mm
05.5.0009.119E	Broca Ø 1,9x Stop 4mmx Comprimento 90mm
05.5.0010.119E	Broca Ø 1,9x Stop 4mmx Comprimento 95mm
05.5.0011.119E	Broca Ø 1,9x Stop 4mmx Comprimento 100mm
05.5.0012.119E	Broca Ø 1,9x Stop 4mmx Comprimento 105mm
05.5.0013.119E	Broca Ø 1,9x Stop 4mmx Comprimento 110mm
05.5.0014.119E	Broca Ø 1,9x Stop 4mmx Comprimento 115mm
05.5.0015.119E	Broca Ø 1,9x Stop 4mmx Comprimento 120mm
05.5.0016.119E	Broca Ø 1,9x Stop 4mmx Comprimento 125mm
05.5.0017.119E	Broca Ø 1,9x Stop 4mmx Comprimento 130mm
05.5.0018.119E	Broca Ø 1,9x Stop 4mmx Comprimento 135mm
05.5.0019.119E	Broca Ø 1,9x Stop 4mmx Comprimento 140mm
05.5.0020.119E	Broca Ø 1,9x Stop 4mmx Comprimento 145mm
05.5.0021.119E	Broca Ø 1,9x Stop 4mmx Comprimento 150mm
05.5.0022.119E	Broca Ø 1,9x Stop 5mmx Comprimento 50mm
05.5.0023.119E	Broca Ø 1,9x Stop 5mmx Comprimento 55mm
05.5.0024.119E	Broca Ø 1,9x Stop 5mmx Comprimento 60mm
05.5.0025.119E	Broca Ø 1,9x Stop 5mmx Comprimento 65mm
05.5.0026.119E	Broca Ø 1,9x Stop 5mmx Comprimento 70mm
05.5.0027.119E	Broca Ø 1,9x Stop 5mmx Comprimento 75mm
05.5.0028.119E	Broca Ø 1,9x Stop 5mmx Comprimento 80mm
05.5.0029.119E	Broca Ø 1,9x Stop 5mmx Comprimento 85mm
05.5.0030.119E	Broca Ø 1,9x Stop 5mmx Comprimento 90mm
05.5.0031.119E	Broca Ø 1,9x Stop 5mmx Comprimento 95mm
05.5.0032.119E	Broca Ø 1,9x Stop 5mmx Comprimento 100mm
05.5.0033.119E	Broca Ø 1,9x Stop 5mmx Comprimento 105mm
05.5.0034.119E	Broca Ø 1,9x Stop 5mmx Comprimento 110mm
05.5.0035.119E	Broca Ø 1,9x Stop 5mmx Comprimento 115mm
05.5.0036.119E	Broca Ø 1,9x Stop 5mmx Comprimento 120mm
05.5.0037.119E	Broca Ø 1,9x Stop 5mmx Comprimento 125mm
05.5.0038.119E	Broca Ø 1,9x Stop 5mmx Comprimento 130mm
05.5.0039.119E	Broca Ø 1,9x Stop 5mmx Comprimento 135mm
05.5.0040.119E	Broca Ø 1,9x Stop 5mmx Comprimento 140mm
05.5.0041.119E	Broca Ø 1,9x Stop 5mmx Comprimento 145mm
05.5.0042.119E	Broca Ø 1,9x Stop 5mmx Comprimento 150mm
05.5.0043.119E	Broca Ø 1,9x Stop 6mmx Comprimento 50mm
05.5.0044.119E	Broca Ø 1,9x Stop 6mmx Comprimento 55mm
05.5.0045.119E	Broca Ø 1,9x Stop 6mmx Comprimento 60mm
05.5.0046.119E	Broca Ø 1,9x Stop 6mmx Comprimento 65mm
05.5.0047.119E	Broca Ø 1,9x Stop 6mmx Comprimento 70mm
05.5.0048.119E	Broca Ø 1,9x Stop 6mmx Comprimento 75mm
05.5.0049.119E	Broca Ø 1,9x Stop 6mmx Comprimento 80mm
05.5.0050.119E	Broca Ø 1,9x Stop 6mmx Comprimento 85mm
05.5.0051.119E	Broca Ø 1,9x Stop 6mmx Comprimento 90mm
05.5.0052.119E	Broca Ø 1,9x Stop 6mmx Comprimento 95mm

05.5.0053.119E	Broca Ø 1,9x Stop 6mmx Comprimento 100mm
05.5.0054.119E	Broca Ø 1,9x Stop 6mmx Comprimento 105mm
05.5.0055.119E	Broca Ø 1,9x Stop 6mmx Comprimento 110mm
05.5.0056.119E	Broca Ø 1,9x Stop 6mmx Comprimento 115mm
05.5.0057.119E	Broca Ø 1,9x Stop 6mmx Comprimento 120mm
05.5.0058.119E	Broca Ø 1,9x Stop 6mmx Comprimento 125mm
05.5.0059.119E	Broca Ø 1,9x Stop 6mmx Comprimento 130mm
05.5.0060.119E	Broca Ø 1,9x Stop 6mmx Comprimento 135mm
05.5.0061.119E	Broca Ø 1,9x Stop 6mmx Comprimento 140mm
05.5.0062.119E	Broca Ø 1,9x Stop 6mmx Comprimento 145mm
05.5.0063.119E	Broca Ø 1,9x Stop 6mmx Comprimento 150mm
05.5.0064.119E	Broca Ø 1,9x Stop 7mmx Comprimento 50mm
05.5.0065.119E	Broca Ø 1,9x Stop 7mmx Comprimento 55mm
05.5.0066.119E	Broca Ø 1,9x Stop 7mmx Comprimento 60mm
05.5.0067.119E	Broca Ø 1,9x Stop 7mmx Comprimento 65mm
05.5.0068.119E	Broca Ø 1,9x Stop 7mmx Comprimento 70mm
05.5.0069.119E	Broca Ø 1,9x Stop 7mmx Comprimento 75mm
05.5.0070.119E	Broca Ø 1,9x Stop 7mmx Comprimento 80mm
05.5.0071.119E	Broca Ø 1,9x Stop 7mmx Comprimento 85mm
05.5.0072.119E	Broca Ø 1,9x Stop 7mmx Comprimento 90mm
05.5.0073.119E	Broca Ø 1,9x Stop 7mmx Comprimento 95mm
05.5.0074.119E	Broca Ø 1,9x Stop 7mmx Comprimento 100mm
05.5.0075.119E	Broca Ø 1,9x Stop 7mmx Comprimento 105mm
05.5.0076.119E	Broca Ø 1,9x Stop 7mmx Comprimento 110mm
05.5.0077.119E	Broca Ø 1,9x Stop 7mmx Comprimento 115mm
05.5.0078.119E	Broca Ø 1,9x Stop 7mmx Comprimento 120mm
05.5.0079.119E	Broca Ø 1,9x Stop 7mmx Comprimento 125mm
05.5.0080.119E	Broca Ø 1,9x Stop 7mmx Comprimento 130mm
05.5.0081.119E	Broca Ø 1,9x Stop 7mmx Comprimento 135mm
05.5.0082.119E	Broca Ø 1,9x Stop 7mmx Comprimento 140mm
05.5.0083.119E	Broca Ø 1,9x Stop 7mmx Comprimento 145mm
05.5.0084.119E	Broca Ø 1,9x Stop 7mmx Comprimento 150mm
05.5.0085.119E	Broca Ø 1,9x Stop 8mmx Comprimento 50mm
05.5.0086.119E	Broca Ø 1,9x Stop 8mmx Comprimento 55mm
05.5.0087.119E	Broca Ø 1,9x Stop 8mmx Comprimento 60mm
05.5.0088.119E	Broca Ø 1,9x Stop 8mmx Comprimento 65mm
05.5.0089.119E	Broca Ø 1,9x Stop 8mmx Comprimento 70mm
05.5.0090.119E	Broca Ø 1,9x Stop 8mmx Comprimento 75mm
05.5.0091.119E	Broca Ø 1,9x Stop 8mmx Comprimento 80mm
05.5.0092.119E	Broca Ø 1,9x Stop 8mmx Comprimento 85mm
05.5.0093.119E	Broca Ø 1,9x Stop 8mmx Comprimento 90mm
05.5.0094.119E	Broca Ø 1,9x Stop 8mmx Comprimento 95mm
05.5.0095.119E	Broca Ø 1,9x Stop 8mmx Comprimento 100mm
05.5.0096.119E	Broca Ø 1,9x Stop 8mmx Comprimento 105mm
05.5.0097.119E	Broca Ø 1,9x Stop 8mmx Comprimento 110mm
05.5.0098.119E	Broca Ø 1,9x Stop 8mmx Comprimento 115mm
05.5.0099.119E	Broca Ø 1,9x Stop 8mmx Comprimento 120mm
05.5.0100.119E	Broca Ø 1,9x Stop 8mmx Comprimento 125mm
05.5.0101.119E	Broca Ø 1,9x Stop 8mmx Comprimento 130mm

05.5.0102.119E	Broca Ø 1,9x Stop 8mmx Comprimento 135mm
05.5.0103.119E	Broca Ø 1,9x Stop 8mmx Comprimento 140mm
05.5.0104.119E	Broca Ø 1,9x Stop 8mmx Comprimento 145mm
05.5.0105.119E	Broca Ø 1,9x Stop 8mmx Comprimento 150mm
05.5.0106.119E	Broca Ø 1,9x Stop 9mmx Comprimento 50mm
05.5.0107.119E	Broca Ø 1,9x Stop 9mmx Comprimento 55mm
05.5.0108.119E	Broca Ø 1,9x Stop 9mmx Comprimento 60mm
05.5.0109.119E	Broca Ø 1,9x Stop 9mmx Comprimento 65mm
05.5.0110.119E	Broca Ø 1,9x Stop 9mmx Comprimento 70mm
05.5.0111.119E	Broca Ø 1,9x Stop 9mmx Comprimento 75mm
05.5.0112.119E	Broca Ø 1,9x Stop 9mmx Comprimento 80mm
05.5.0113.119E	Broca Ø 1,9x Stop 9mmx Comprimento 85mm
05.5.0114.119E	Broca Ø 1,9x Stop 9mmx Comprimento 90mm
05.5.0115.119E	Broca Ø 1,9x Stop 9mmx Comprimento 95mm
05.5.0116.119E	Broca Ø 1,9x Stop 9mmx Comprimento 100mm
05.5.0117.119E	Broca Ø 1,9x Stop 9mmx Comprimento 105mm
05.5.0118.119E	Broca Ø 1,9x Stop 9mmx Comprimento 110mm
05.5.0119.119E	Broca Ø 1,9x Stop 9mmx Comprimento 115mm
05.5.0120.119E	Broca Ø 1,9x Stop 9mmx Comprimento 120mm
05.5.0121.119E	Broca Ø 1,9x Stop 9mmx Comprimento 125mm
05.5.0122.119E	Broca Ø 1,9x Stop 9mmx Comprimento 130mm
05.5.0123.119E	Broca Ø 1,9x Stop 9mmx Comprimento 135mm
05.5.0124.119E	Broca Ø 1,9x Stop 9mmx Comprimento 140mm
05.5.0125.119E	Broca Ø 1,9x Stop 9mmx Comprimento 145mm
05.5.0126.119E	Broca Ø 1,9x Stop 9mmx Comprimento 150mm
05.5.0127.119E	Broca Ø 1,9x Stop 10mmx Comprimento 50mm
05.5.0128.119E	Broca Ø 1,9x Stop 10mmx Comprimento 55mm
05.5.0129.119E	Broca Ø 1,9x Stop 10mmx Comprimento 60mm
05.5.0130.119E	Broca Ø 1,9x Stop 10mmx Comprimento 65mm
05.5.0131.119E	Broca Ø 1,9x Stop 10mmx Comprimento 70mm
05.5.0132.119E	Broca Ø 1,9x Stop 10mmx Comprimento 75mm
05.5.0133.119E	Broca Ø 1,9x Stop 10mmx Comprimento 80mm
05.5.0134.119E	Broca Ø 1,9x Stop 10mmx Comprimento 85mm
05.5.0135.119E	Broca Ø 1,9x Stop 10mmx Comprimento 90mm
05.5.0136.119E	Broca Ø 1,9x Stop 10mmx Comprimento 95mm
05.5.0137.119E	Broca Ø 1,9x Stop 10mmx Comprimento 100mm
05.5.0138.119E	Broca Ø 1,9x Stop 10mmx Comprimento 105mm
05.5.0139.119E	Broca Ø 1,9x Stop 10mmx Comprimento 110mm
05.5.0140.119E	Broca Ø 1,9x Stop 10mmx Comprimento 115mm
05.5.0141.119E	Broca Ø 1,9x Stop 10mmx Comprimento 120mm
05.5.0142.119E	Broca Ø 1,9x Stop 10mmx Comprimento 125mm
05.5.0143.119E	Broca Ø 1,9x Stop 10mmx Comprimento 130mm
05.5.0144.119E	Broca Ø 1,9x Stop 10mmx Comprimento 135mm
05.5.0145.119E	Broca Ø 1,9x Stop 10mmx Comprimento 140mm
05.5.0146.119E	Broca Ø 1,9x Stop 10mmx Comprimento 145mm
05.5.0147.119E	Broca Ø 1,9x Stop 10mmx Comprimento 150mm
05.5.0148.119E	Broca Ø 1,9x Stop 11mmx Comprimento 50mm
05.5.0149.119E	Broca Ø 1,9x Stop 11mmx Comprimento 55mm
05.5.0150.119E	Broca Ø 1,9x Stop 11mmx Comprimento 60mm

05.5.0151.119E	Broca Ø 1,9x Stop 11mmx Comprimento 65mm
05.5.0152.119E	Broca Ø 1,9x Stop 11mmx Comprimento 70mm
05.5.0153.119E	Broca Ø 1,9x Stop 11mmx Comprimento 75mm
05.5.0154.119E	Broca Ø 1,9x Stop 11mmx Comprimento 80mm
05.5.0155.119E	Broca Ø 1,9x Stop 11mmx Comprimento 85mm
05.5.0156.119E	Broca Ø 1,9x Stop 11mmx Comprimento 90mm
05.5.0157.119E	Broca Ø 1,9x Stop 11mmx Comprimento 95mm
05.5.0158.119E	Broca Ø 1,9x Stop 11mmx Comprimento 100mm
05.5.0159.119E	Broca Ø 1,9x Stop 11mmx Comprimento 105mm
05.5.0160.119E	Broca Ø 1,9x Stop 11mmx Comprimento 110mm
05.5.0161.119E	Broca Ø 1,9x Stop 11mmx Comprimento 115mm
05.5.0162.119E	Broca Ø 1,9x Stop 11mmx Comprimento 120mm
05.5.0163.119E	Broca Ø 1,9x Stop 11mmx Comprimento 125mm
05.5.0164.119E	Broca Ø 1,9x Stop 11mmx Comprimento 130mm
05.5.0165.119E	Broca Ø 1,9x Stop 11mmx Comprimento 135mm
05.5.0166.119E	Broca Ø 1,9x Stop 11mmx Comprimento 140mm
05.5.0167.119E	Broca Ø 1,9x Stop 11mmx Comprimento 145mm
05.5.0168.119E	Broca Ø 1,9x Stop 11mmx Comprimento 150mm
05.5.0169.119E	Broca Ø 1,9x Stop 12mmx Comprimento 50mm
05.5.0170.119E	Broca Ø 1,9x Stop 12mmx Comprimento 55mm
05.5.0171.119E	Broca Ø 1,9x Stop 12mmx Comprimento 60mm
05.5.0172.119E	Broca Ø 1,9x Stop 12mmx Comprimento 65mm
05.5.0173.119E	Broca Ø 1,9x Stop 12mmx Comprimento 70mm
05.5.0174.119E	Broca Ø 1,9x Stop 12mmx Comprimento 75mm
05.5.0175.119E	Broca Ø 1,9x Stop 12mmx Comprimento 80mm
05.5.0176.119E	Broca Ø 1,9x Stop 12mmx Comprimento 85mm
05.5.0177.119E	Broca Ø 1,9x Stop 12mmx Comprimento 90mm
05.5.0178.119E	Broca Ø 1,9x Stop 12mmx Comprimento 95mm
05.5.0179.119E	Broca Ø 1,9x Stop 12mmx Comprimento 100mm
05.5.0180.119E	Broca Ø 1,9x Stop 12mmx Comprimento 105mm
05.5.0181.119E	Broca Ø 1,9x Stop 12mmx Comprimento 110mm
05.5.0182.119E	Broca Ø 1,9x Stop 12mmx Comprimento 115mm
05.5.0183.119E	Broca Ø 1,9x Stop 12mmx Comprimento 120mm
05.5.0184.119E	Broca Ø 1,9x Stop 12mmx Comprimento 125mm
05.5.0185.119E	Broca Ø 1,9x Stop 12mmx Comprimento 130mm
05.5.0186.119E	Broca Ø 1,9x Stop 12mmx Comprimento 135mm
05.5.0187.119E	Broca Ø 1,9x Stop 12mmx Comprimento 140mm
05.5.0188.119E	Broca Ø 1,9x Stop 12mmx Comprimento 145mm
05.5.0189.119E	Broca Ø 1,9x Stop 12mmx Comprimento 150mm
05.5.0190.119E	Broca Ø 1,9x Stop 13mmx Comprimento 50mm
05.5.0191.119E	Broca Ø 1,9x Stop 13mmx Comprimento 55mm
05.5.0192.119E	Broca Ø 1,9x Stop 13mmx Comprimento 60mm
05.5.0193.119E	Broca Ø 1,9x Stop 13mmx Comprimento 65mm
05.5.0194.119E	Broca Ø 1,9x Stop 13mmx Comprimento 70mm
05.5.0195.119E	Broca Ø 1,9x Stop 13mmx Comprimento 75mm
05.5.0196.119E	Broca Ø 1,9x Stop 13mmx Comprimento 80mm
05.5.0197.119E	Broca Ø 1,9x Stop 13mmx Comprimento 85mm
05.5.0198.119E	Broca Ø 1,9x Stop 13mmx Comprimento 90mm
05.5.0199.119E	Broca Ø 1,9x Stop 13mmx Comprimento 95mm

05.5.0200.119E	Broca Ø 1,9x Stop 13mmx Comprimento 100mm
05.5.0201.119E	Broca Ø 1,9x Stop 13mmx Comprimento 105mm
05.5.0202.119E	Broca Ø 1,9x Stop 13mmx Comprimento 110mm
05.5.0203.119E	Broca Ø 1,9x Stop 13mmx Comprimento 115mm
05.5.0204.119E	Broca Ø 1,9x Stop 13mmx Comprimento 120mm
05.5.0205.119E	Broca Ø 1,9x Stop 13mmx Comprimento 125mm
05.5.0206.119E	Broca Ø 1,9x Stop 13mmx Comprimento 130mm
05.5.0207.119E	Broca Ø 1,9x Stop 13mmx Comprimento 135mm
05.5.0208.119E	Broca Ø 1,9x Stop 13mmx Comprimento 140mm
05.5.0209.119E	Broca Ø 1,9x Stop 13mmx Comprimento 145mm
05.5.0210.119E	Broca Ø 1,9x Stop 13mmx Comprimento 150mm
05.5.0211.119E	Broca Ø 1,9x Stop 14mmx Comprimento 50mm
05.5.0212.119E	Broca Ø 1,9x Stop 14mmx Comprimento 55mm
05.5.0213.119E	Broca Ø 1,9x Stop 14mmx Comprimento 60mm
05.5.0214.119E	Broca Ø 1,9x Stop 14mmx Comprimento 65mm
05.5.0215.119E	Broca Ø 1,9x Stop 14mmx Comprimento 70mm
05.5.0216.119E	Broca Ø 1,9x Stop 14mmx Comprimento 75mm
05.5.0217.119E	Broca Ø 1,9x Stop 14mmx Comprimento 80mm
05.5.0218.119E	Broca Ø 1,9x Stop 14mmx Comprimento 85mm
05.5.0219.119E	Broca Ø 1,9x Stop 14mmx Comprimento 90mm
05.5.0220.119E	Broca Ø 1,9x Stop 14mmx Comprimento 95mm
05.5.0221.119E	Broca Ø 1,9x Stop 14mmx Comprimento 100mm
05.5.0222.119E	Broca Ø 1,9x Stop 14mmx Comprimento 105mm
05.5.0223.119E	Broca Ø 1,9x Stop 14mmx Comprimento 110mm
05.5.0224.119E	Broca Ø 1,9x Stop 14mmx Comprimento 115mm
05.5.0225.119E	Broca Ø 1,9x Stop 14mmx Comprimento 120mm
05.5.0226.119E	Broca Ø 1,9x Stop 14mmx Comprimento 125mm
05.5.0227.119E	Broca Ø 1,9x Stop 14mmx Comprimento 130mm
05.5.0228.119E	Broca Ø 1,9x Stop 14mmx Comprimento 135mm
05.5.0229.119E	Broca Ø 1,9x Stop 14mmx Comprimento 140mm
05.5.0230.119E	Broca Ø 1,9x Stop 14mmx Comprimento 145mm
05.5.0231.119E	Broca Ø 1,9x Stop 14mmx Comprimento 150mm
05.5.0232.119E	Broca Ø 1,9x Stop 15mmx Comprimento 50mm
05.5.0233.119E	Broca Ø 1,9x Stop 15mmx Comprimento 55mm
05.5.0234.119E	Broca Ø 1,9x Stop 15mmx Comprimento 60mm
05.5.0235.119E	Broca Ø 1,9x Stop 15mmx Comprimento 65mm
05.5.0236.119E	Broca Ø 1,9x Stop 15mmx Comprimento 70mm
05.5.0237.119E	Broca Ø 1,9x Stop 15mmx Comprimento 75mm
05.5.0238.119E	Broca Ø 1,9x Stop 15mmx Comprimento 80mm
05.5.0239.119E	Broca Ø 1,9x Stop 15mmx Comprimento 85mm
05.5.0240.119E	Broca Ø 1,9x Stop 15mmx Comprimento 90mm
05.5.0241.119E	Broca Ø 1,9x Stop 15mmx Comprimento 95mm
05.5.0242.119E	Broca Ø 1,9x Stop 15mmx Comprimento 100mm
05.5.0243.119E	Broca Ø 1,9x Stop 15mmx Comprimento 105mm
05.5.0244.119E	Broca Ø 1,9x Stop 15mmx Comprimento 110mm
05.5.0245.119E	Broca Ø 1,9x Stop 15mmx Comprimento 115mm
05.5.0246.119E	Broca Ø 1,9x Stop 15mmx Comprimento 120mm
05.5.0247.119E	Broca Ø 1,9x Stop 15mmx Comprimento 125mm
05.5.0248.119E	Broca Ø 1,9x Stop 15mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0249.119E	Broca Ø 1,9x Stop 15mmx Comprimento 135mm
05.5.0250.119E	Broca Ø 1,9x Stop 15mmx Comprimento 140mm
05.5.0251.119E	Broca Ø 1,9x Stop 15mmx Comprimento 145mm
05.5.0252.119E	Broca Ø 1,9x Stop 15mmx Comprimento 150mm
05.5.0253.119E	Broca Ø 1,9x Stop 16mmx Comprimento 50mm
05.5.0254.119E	Broca Ø 1,9x Stop 16mmx Comprimento 55mm
05.5.0255.119E	Broca Ø 1,9x Stop 16mmx Comprimento 60mm
05.5.0256.119E	Broca Ø 1,9x Stop 16mmx Comprimento 65mm
05.5.0257.119E	Broca Ø 1,9x Stop 16mmx Comprimento 70mm
05.5.0258.119E	Broca Ø 1,9x Stop 16mmx Comprimento 75mm
05.5.0259.119E	Broca Ø 1,9x Stop 16mmx Comprimento 80mm
05.5.0260.119E	Broca Ø 1,9x Stop 16mmx Comprimento 85mm
05.5.0261.119E	Broca Ø 1,9x Stop 16mmx Comprimento 90mm
05.5.0262.119E	Broca Ø 1,9x Stop 16mmx Comprimento 95mm
05.5.0263.119E	Broca Ø 1,9x Stop 16mmx Comprimento 100mm
05.5.0264.119E	Broca Ø 1,9x Stop 16mmx Comprimento 105mm
05.5.0265.119E	Broca Ø 1,9x Stop 16mmx Comprimento 110mm
05.5.0266.119E	Broca Ø 1,9x Stop 16mmx Comprimento 115mm
05.5.0267.119E	Broca Ø 1,9x Stop 16mmx Comprimento 120mm
05.5.0268.119E	Broca Ø 1,9x Stop 16mmx Comprimento 125mm
05.5.0269.119E	Broca Ø 1,9x Stop 16mmx Comprimento 130mm
05.5.0270.119E	Broca Ø 1,9x Stop 16mmx Comprimento 135mm
05.5.0271.119E	Broca Ø 1,9x Stop 16mmx Comprimento 140mm
05.5.0272.119E	Broca Ø 1,9x Stop 16mmx Comprimento 145mm
05.5.0273.119E	Broca Ø 1,9x Stop 16mmx Comprimento 150mm
05.5.0274.119E	Broca Ø 1,9x Stop 17mmx Comprimento 50mm
05.5.0275.119E	Broca Ø 1,9x Stop 17mmx Comprimento 55mm
05.5.0276.119E	Broca Ø 1,9x Stop 17mmx Comprimento 60mm
05.5.0277.119E	Broca Ø 1,9x Stop 17mmx Comprimento 65mm
05.5.0278.119E	Broca Ø 1,9x Stop 17mmx Comprimento 70mm
05.5.0279.119E	Broca Ø 1,9x Stop 17mmx Comprimento 75mm
05.5.0280.119E	Broca Ø 1,9x Stop 17mmx Comprimento 80mm
05.5.0281.119E	Broca Ø 1,9x Stop 17mmx Comprimento 85mm
05.5.0282.119E	Broca Ø 1,9x Stop 17mmx Comprimento 90mm
05.5.0283.119E	Broca Ø 1,9x Stop 17mmx Comprimento 95mm
05.5.0284.119E	Broca Ø 1,9x Stop 17mmx Comprimento 100mm
05.5.0285.119E	Broca Ø 1,9x Stop 17mmx Comprimento 105mm
05.5.0286.119E	Broca Ø 1,9x Stop 17mmx Comprimento 110mm
05.5.0287.119E	Broca Ø 1,9x Stop 17mmx Comprimento 115mm
05.5.0288.119E	Broca Ø 1,9x Stop 17mmx Comprimento 120mm
05.5.0289.119E	Broca Ø 1,9x Stop 17mmx Comprimento 125mm
05.5.0290.119E	Broca Ø 1,9x Stop 17mmx Comprimento 130mm
05.5.0291.119E	Broca Ø 1,9x Stop 17mmx Comprimento 135mm
05.5.0292.119E	Broca Ø 1,9x Stop 17mmx Comprimento 140mm
05.5.0293.119E	Broca Ø 1,9x Stop 17mmx Comprimento 145mm
05.5.0294.119E	Broca Ø 1,9x Stop 17mmx Comprimento 150mm
05.5.0295.119E	Broca Ø 1,9x Stop 18mmx Comprimento 50mm
05.5.0296.119E	Broca Ø 1,9x Stop 18mmx Comprimento 55mm
05.5.0297.119E	Broca Ø 1,9x Stop 18mmx Comprimento 60mm

05.5.0298.119E	Broca Ø 1,9x Stop 18mmx Comprimento 65mm
05.5.0299.119E	Broca Ø 1,9x Stop 18mmx Comprimento 70mm
05.5.0300.119E	Broca Ø 1,9x Stop 18mmx Comprimento 75mm
05.5.0301.119E	Broca Ø 1,9x Stop 18mmx Comprimento 80mm
05.5.0302.119E	Broca Ø 1,9x Stop 18mmx Comprimento 85mm
05.5.0303.119E	Broca Ø 1,9x Stop 18mmx Comprimento 90mm
05.5.0304.119E	Broca Ø 1,9x Stop 18mmx Comprimento 95mm
05.5.0305.119E	Broca Ø 1,9x Stop 18mmx Comprimento 100mm
05.5.0306.119E	Broca Ø 1,9x Stop 18mmx Comprimento 105mm
05.5.0307.119E	Broca Ø 1,9x Stop 18mmx Comprimento 110mm
05.5.0308.119E	Broca Ø 1,9x Stop 18mmx Comprimento 115mm
05.5.0309.119E	Broca Ø 1,9x Stop 18mmx Comprimento 120mm
05.5.0310.119E	Broca Ø 1,9x Stop 18mmx Comprimento 125mm
05.5.0311.119E	Broca Ø 1,9x Stop 18mmx Comprimento 130mm
05.5.0312.119E	Broca Ø 1,9x Stop 18mmx Comprimento 135mm
05.5.0313.119E	Broca Ø 1,9x Stop 18mmx Comprimento 140mm
05.5.0314.119E	Broca Ø 1,9x Stop 18mmx Comprimento 145mm
05.5.0315.119E	Broca Ø 1,9x Stop 18mmx Comprimento 150mm
05.5.0316.119E	Broca Ø 1,9x Stop 19mmx Comprimento 50mm
05.5.0317.119E	Broca Ø 1,9x Stop 19mmx Comprimento 55mm
05.5.0318.119E	Broca Ø 1,9x Stop 19mmx Comprimento 60mm
05.5.0319.119E	Broca Ø 1,9x Stop 19mmx Comprimento 65mm
05.5.0320.119E	Broca Ø 1,9x Stop 19mmx Comprimento 70mm
05.5.0321.119E	Broca Ø 1,9x Stop 19mmx Comprimento 75mm
05.5.0322.119E	Broca Ø 1,9x Stop 19mmx Comprimento 80mm
05.5.0323.119E	Broca Ø 1,9x Stop 19mmx Comprimento 85mm
05.5.0324.119E	Broca Ø 1,9x Stop 19mmx Comprimento 90mm
05.5.0325.119E	Broca Ø 1,9x Stop 19mmx Comprimento 95mm
05.5.0326.119E	Broca Ø 1,9x Stop 19mmx Comprimento 100mm
05.5.0327.119E	Broca Ø 1,9x Stop 19mmx Comprimento 105mm
05.5.0328.119E	Broca Ø 1,9x Stop 19mmx Comprimento 110mm
05.5.0329.119E	Broca Ø 1,9x Stop 19mmx Comprimento 115mm
05.5.0330.119E	Broca Ø 1,9x Stop 19mmx Comprimento 120mm
05.5.0331.119E	Broca Ø 1,9x Stop 19mmx Comprimento 125mm
05.5.0332.119E	Broca Ø 1,9x Stop 19mmx Comprimento 130mm
05.5.0333.119E	Broca Ø 1,9x Stop 19mmx Comprimento 135mm
05.5.0334.119E	Broca Ø 1,9x Stop 19mmx Comprimento 140mm
05.5.0335.119E	Broca Ø 1,9x Stop 19mmx Comprimento 145mm
05.5.0336.119E	Broca Ø 1,9x Stop 19mmx Comprimento 150mm
05.5.0337.119E	Broca Ø 1,9x Stop 20mmx Comprimento 50mm
05.5.0338.119E	Broca Ø 1,9x Stop 20mmx Comprimento 55mm
05.5.0339.119E	Broca Ø 1,9x Stop 20mmx Comprimento 60mm
05.5.0340.119E	Broca Ø 1,9x Stop 20mmx Comprimento 65mm
05.5.0341.119E	Broca Ø 1,9x Stop 20mmx Comprimento 70mm
05.5.0342.119E	Broca Ø 1,9x Stop 20mmx Comprimento 75mm
05.5.0343.119E	Broca Ø 1,9x Stop 20mmx Comprimento 80mm
05.5.0344.119E	Broca Ø 1,9x Stop 20mmx Comprimento 85mm
05.5.0345.119E	Broca Ø 1,9x Stop 20mmx Comprimento 90mm
05.5.0346.119E	Broca Ø 1,9x Stop 20mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0347.119E	Broca Ø 1,9x Stop 20mmx Comprimento 100mm
05.5.0348.119E	Broca Ø 1,9x Stop 20mmx Comprimento 105mm
05.5.0349.119E	Broca Ø 1,9x Stop 20mmx Comprimento 110mm
05.5.0350.119E	Broca Ø 1,9x Stop 20mmx Comprimento 115mm
05.5.0351.119E	Broca Ø 1,9x Stop 20mmx Comprimento 120mm
05.5.0352.119E	Broca Ø 1,9x Stop 20mmx Comprimento 125mm
05.5.0353.119E	Broca Ø 1,9x Stop 20mmx Comprimento 130mm
05.5.0354.119E	Broca Ø 1,9x Stop 20mmx Comprimento 135mm
05.5.0355.119E	Broca Ø 1,9x Stop 20mmx Comprimento 140mm
05.5.0356.119E	Broca Ø 1,9x Stop 20mmx Comprimento 145mm
05.5.0357.119E	Broca Ø 1,9x Stop 20mmx Comprimento 150mm
05.5.0358.119E	Broca Ø 1,9x Stop 21mmx Comprimento 50mm
05.5.0359.119E	Broca Ø 1,9x Stop 21mmx Comprimento 55mm
05.5.0360.119E	Broca Ø 1,9x Stop 21mmx Comprimento 60mm
05.5.0361.119E	Broca Ø 1,9x Stop 21mmx Comprimento 65mm
05.5.0362.119E	Broca Ø 1,9x Stop 21mmx Comprimento 70mm
05.5.0363.119E	Broca Ø 1,9x Stop 21mmx Comprimento 75mm
05.5.0364.119E	Broca Ø 1,9x Stop 21mmx Comprimento 80mm
05.5.0365.119E	Broca Ø 1,9x Stop 21mmx Comprimento 85mm
05.5.0366.119E	Broca Ø 1,9x Stop 21mmx Comprimento 90mm
05.5.0367.119E	Broca Ø 1,9x Stop 21mmx Comprimento 95mm
05.5.0368.119E	Broca Ø 1,9x Stop 21mmx Comprimento 100mm
05.5.0369.119E	Broca Ø 1,9x Stop 21mmx Comprimento 105mm
05.5.0370.119E	Broca Ø 1,9x Stop 21mmx Comprimento 110mm
05.5.0371.119E	Broca Ø 1,9x Stop 21mmx Comprimento 115mm
05.5.0372.119E	Broca Ø 1,9x Stop 21mmx Comprimento 120mm
05.5.0373.119E	Broca Ø 1,9x Stop 21mmx Comprimento 125mm
05.5.0374.119E	Broca Ø 1,9x Stop 21mmx Comprimento 130mm
05.5.0375.119E	Broca Ø 1,9x Stop 21mmx Comprimento 135mm
05.5.0376.119E	Broca Ø 1,9x Stop 21mmx Comprimento 140mm
05.5.0377.119E	Broca Ø 1,9x Stop 21mmx Comprimento 145mm
05.5.0378.119E	Broca Ø 1,9x Stop 21mmx Comprimento 150mm
05.5.0379.119E	Broca Ø 1,9x Stop 22mmx Comprimento 50mm
05.5.0380.119E	Broca Ø 1,9x Stop 22mmx Comprimento 55mm
05.5.0381.119E	Broca Ø 1,9x Stop 22mmx Comprimento 60mm
05.5.0382.119E	Broca Ø 1,9x Stop 22mmx Comprimento 65mm
05.5.0383.119E	Broca Ø 1,9x Stop 22mmx Comprimento 70mm
05.5.0384.119E	Broca Ø 1,9x Stop 22mmx Comprimento 75mm
05.5.0385.119E	Broca Ø 1,9x Stop 22mmx Comprimento 80mm
05.5.0386.119E	Broca Ø 1,9x Stop 22mmx Comprimento 85mm
05.5.0387.119E	Broca Ø 1,9x Stop 22mmx Comprimento 90mm
05.5.0388.119E	Broca Ø 1,9x Stop 22mmx Comprimento 95mm
05.5.0389.119E	Broca Ø 1,9x Stop 22mmx Comprimento 100mm
05.5.0390.119E	Broca Ø 1,9x Stop 22mmx Comprimento 105mm
05.5.0391.119E	Broca Ø 1,9x Stop 22mmx Comprimento 110mm
05.5.0392.119E	Broca Ø 1,9x Stop 22mmx Comprimento 115mm
05.5.0393.119E	Broca Ø 1,9x Stop 22mmx Comprimento 120mm
05.5.0394.119E	Broca Ø 1,9x Stop 22mmx Comprimento 125mm
05.5.0395.119E	Broca Ø 1,9x Stop 22mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0396.119E	Broca Ø 1,9x Stop 22mmx Comprimento 135mm
05.5.0397.119E	Broca Ø 1,9x Stop 22mmx Comprimento 140mm
05.5.0398.119E	Broca Ø 1,9x Stop 22mmx Comprimento 145mm
05.5.0399.119E	Broca Ø 1,9x Stop 22mmx Comprimento 150mm
05.5.0400.119E	Broca Ø 1,9x Stop 23mmx Comprimento 50mm
05.5.0401.119E	Broca Ø 1,9x Stop 23mmx Comprimento 55mm
05.5.0402.119E	Broca Ø 1,9x Stop 23mmx Comprimento 60mm
05.5.0403.119E	Broca Ø 1,9x Stop 23mmx Comprimento 65mm
05.5.0404.119E	Broca Ø 1,9x Stop 23mmx Comprimento 70mm
05.5.0405.119E	Broca Ø 1,9x Stop 23mmx Comprimento 75mm
05.5.0406.119E	Broca Ø 1,9x Stop 23mmx Comprimento 80mm
05.5.0407.119E	Broca Ø 1,9x Stop 23mmx Comprimento 85mm
05.5.0408.119E	Broca Ø 1,9x Stop 23mmx Comprimento 90mm
05.5.0409.119E	Broca Ø 1,9x Stop 23mmx Comprimento 95mm
05.5.0410.119E	Broca Ø 1,9x Stop 23mmx Comprimento 100mm
05.5.0411.119E	Broca Ø 1,9x Stop 23mmx Comprimento 105mm
05.5.0412.119E	Broca Ø 1,9x Stop 23mmx Comprimento 110mm
05.5.0413.119E	Broca Ø 1,9x Stop 23mmx Comprimento 115mm
05.5.0414.119E	Broca Ø 1,9x Stop 23mmx Comprimento 120mm
05.5.0415.119E	Broca Ø 1,9x Stop 23mmx Comprimento 125mm
05.5.0416.119E	Broca Ø 1,9x Stop 23mmx Comprimento 130mm
05.5.0417.119E	Broca Ø 1,9x Stop 23mmx Comprimento 135mm
05.5.0418.119E	Broca Ø 1,9x Stop 23mmx Comprimento 140mm
05.5.0419.119E	Broca Ø 1,9x Stop 23mmx Comprimento 145mm
05.5.0420.119E	Broca Ø 1,9x Stop 23mmx Comprimento 150mm
05.5.0421.119E	Broca Ø 1,9x Stop 24mmx Comprimento 50mm
05.5.0422.119E	Broca Ø 1,9x Stop 24mmx Comprimento 55mm
05.5.0423.119E	Broca Ø 1,9x Stop 24mmx Comprimento 60mm
05.5.0424.119E	Broca Ø 1,9x Stop 24mmx Comprimento 65mm
05.5.0425.119E	Broca Ø 1,9x Stop 24mmx Comprimento 70mm
05.5.0426.119E	Broca Ø 1,9x Stop 24mmx Comprimento 75mm
05.5.0427.119E	Broca Ø 1,9x Stop 24mmx Comprimento 80mm
05.5.0428.119E	Broca Ø 1,9x Stop 24mmx Comprimento 85mm
05.5.0429.119E	Broca Ø 1,9x Stop 24mmx Comprimento 90mm
05.5.0430.119E	Broca Ø 1,9x Stop 24mmx Comprimento 95mm
05.5.0431.119E	Broca Ø 1,9x Stop 24mmx Comprimento 100mm
05.5.0432.119E	Broca Ø 1,9x Stop 24mmx Comprimento 105mm
05.5.0433.119E	Broca Ø 1,9x Stop 24mmx Comprimento 110mm
05.5.0434.119E	Broca Ø 1,9x Stop 24mmx Comprimento 115mm
05.5.0435.119E	Broca Ø 1,9x Stop 24mmx Comprimento 120mm
05.5.0436.119E	Broca Ø 1,9x Stop 24mmx Comprimento 125mm
05.5.0437.119E	Broca Ø 1,9x Stop 24mmx Comprimento 130mm
05.5.0438.119E	Broca Ø 1,9x Stop 24mmx Comprimento 135mm
05.5.0439.119E	Broca Ø 1,9x Stop 24mmx Comprimento 140mm
05.5.0440.119E	Broca Ø 1,9x Stop 24mmx Comprimento 145mm
05.5.0441.119E	Broca Ø 1,9x Stop 24mmx Comprimento 150mm
05.5.0442.119E	Broca Ø 1,9x Stop 25mmx Comprimento 50mm
05.5.0443.119E	Broca Ø 1,9x Stop 25mmx Comprimento 55mm
05.5.0444.119E	Broca Ø 1,9x Stop 25mmx Comprimento 60mm

05.5.0445.119E	Broca Ø 1,9x Stop 25mmx Comprimento 65mm
05.5.0446.119E	Broca Ø 1,9x Stop 25mmx Comprimento 70mm
05.5.0447.119E	Broca Ø 1,9x Stop 25mmx Comprimento 75mm
05.5.0448.119E	Broca Ø 1,9x Stop 25mmx Comprimento 80mm
05.5.0449.119E	Broca Ø 1,9x Stop 25mmx Comprimento 85mm
05.5.0450.119E	Broca Ø 1,9x Stop 25mmx Comprimento 90mm
05.5.0451.119E	Broca Ø 1,9x Stop 25mmx Comprimento 95mm
05.5.0452.119E	Broca Ø 1,9x Stop 25mmx Comprimento 100mm
05.5.0453.119E	Broca Ø 1,9x Stop 25mmx Comprimento 105mm
05.5.0454.119E	Broca Ø 1,9x Stop 25mmx Comprimento 110mm
05.5.0455.119E	Broca Ø 1,9x Stop 25mmx Comprimento 115mm
05.5.0456.119E	Broca Ø 1,9x Stop 25mmx Comprimento 120mm
05.5.0457.119E	Broca Ø 1,9x Stop 25mmx Comprimento 125mm
05.5.0458.119E	Broca Ø 1,9x Stop 25mmx Comprimento 130mm
05.5.0459.119E	Broca Ø 1,9x Stop 25mmx Comprimento 135mm
05.5.0460.119E	Broca Ø 1,9x Stop 25mmx Comprimento 140mm
05.5.0461.119E	Broca Ø 1,9x Stop 25mmx Comprimento 145mm
05.5.0462.119E	Broca Ø 1,9x Stop 25mmx Comprimento 150mm
05.5.0463.119E	Broca Ø 1,9x Stop 30mmx Comprimento 50mm
05.5.0464.119E	Broca Ø 1,9x Stop 30mmx Comprimento 55mm
05.5.0465.119E	Broca Ø 1,9x Stop 30mmx Comprimento 60mm
05.5.0466.119E	Broca Ø 1,9x Stop 30mmx Comprimento 65mm
05.5.0467.119E	Broca Ø 1,9x Stop 30mmx Comprimento 70mm
05.5.0468.119E	Broca Ø 1,9x Stop 30mmx Comprimento 75mm
05.5.0469.119E	Broca Ø 1,9x Stop 30mmx Comprimento 80mm
05.5.0470.119E	Broca Ø 1,9x Stop 30mmx Comprimento 85mm
05.5.0471.119E	Broca Ø 1,9x Stop 30mmx Comprimento 90mm
05.5.0472.119E	Broca Ø 1,9x Stop 30mmx Comprimento 95mm
05.5.0473.119E	Broca Ø 1,9x Stop 30mmx Comprimento 100mm
05.5.0474.119E	Broca Ø 1,9x Stop 30mmx Comprimento 105mm
05.5.0475.119E	Broca Ø 1,9x Stop 30mmx Comprimento 110mm
05.5.0476.119E	Broca Ø 1,9x Stop 30mmx Comprimento 115mm
05.5.0477.119E	Broca Ø 1,9x Stop 30mmx Comprimento 120mm
05.5.0478.119E	Broca Ø 1,9x Stop 30mmx Comprimento 125mm
05.5.0479.119E	Broca Ø 1,9x Stop 30mmx Comprimento 130mm
05.5.0480.119E	Broca Ø 1,9x Stop 30mmx Comprimento 135mm
05.5.0481.119E	Broca Ø 1,9x Stop 30mmx Comprimento 140mm
05.5.0482.119E	Broca Ø 1,9x Stop 30mmx Comprimento 145mm
05.5.0483.119E	Broca Ø 1,9x Stop 30mmx Comprimento 150mm
05.5.0484.119E	Broca Ø 1,9x Stop 35mmx Comprimento 50mm
05.5.0485.119E	Broca Ø 1,9x Stop 35mmx Comprimento 55mm
05.5.0486.119E	Broca Ø 1,9x Stop 35mmx Comprimento 60mm
05.5.0487.119E	Broca Ø 1,9x Stop 35mmx Comprimento 65mm
05.5.0488.119E	Broca Ø 1,9x Stop 35mmx Comprimento 70mm
05.5.0489.119E	Broca Ø 1,9x Stop 35mmx Comprimento 75mm
05.5.0490.119E	Broca Ø 1,9x Stop 35mmx Comprimento 80mm
05.5.0491.119E	Broca Ø 1,9x Stop 35mmx Comprimento 85mm
05.5.0492.119E	Broca Ø 1,9x Stop 35mmx Comprimento 90mm
05.5.0493.119E	Broca Ø 1,9x Stop 35mmx Comprimento 95mm

05.5.0494.119E	Broca Ø 1,9x Stop 35mmx Comprimento 100mm
05.5.0495.119E	Broca Ø 1,9x Stop 35mmx Comprimento 105mm
05.5.0496.119E	Broca Ø 1,9x Stop 35mmx Comprimento 110mm
05.5.0497.119E	Broca Ø 1,9x Stop 35mmx Comprimento 115mm
05.5.0498.119E	Broca Ø 1,9x Stop 35mmx Comprimento 120mm
05.5.0499.119E	Broca Ø 1,9x Stop 35mmx Comprimento 125mm
05.5.0500.119E	Broca Ø 1,9x Stop 35mmx Comprimento 130mm
05.5.0501.119E	Broca Ø 1,9x Stop 35mmx Comprimento 135mm
05.5.0502.119E	Broca Ø 1,9x Stop 35mmx Comprimento 140mm
05.5.0503.119E	Broca Ø 1,9x Stop 35mmx Comprimento 145mm
05.5.0504.119E	Broca Ø 1,9x Stop 35mmx Comprimento 150mm
05.5.0001.120E	Broca Ø 2,0x Stop 4mmx Comprimento 50mm
05.5.0002.120E	Broca Ø 2,0x Stop 4mmx Comprimento 55mm
05.5.0003.120E	Broca Ø 2,0x Stop 4mmx Comprimento 60mm
05.5.0004.120E	Broca Ø 2,0x Stop 4mmx Comprimento 65mm
05.5.0005.120E	Broca Ø 2,0x Stop 4mmx Comprimento 70mm
05.5.0006.120E	Broca Ø 2,0x Stop 4mmx Comprimento 75mm
05.5.0007.120E	Broca Ø 2,0x Stop 4mmx Comprimento 80mm
05.5.0008.120E	Broca Ø 2,0x Stop 4mmx Comprimento 85mm
05.5.0009.120E	Broca Ø 2,0x Stop 4mmx Comprimento 90mm
05.5.0010.120E	Broca Ø 2,0x Stop 4mmx Comprimento 95mm
05.5.0011.120E	Broca Ø 2,0x Stop 4mmx Comprimento 100mm
05.5.0012.120E	Broca Ø 2,0x Stop 4mmx Comprimento 105mm
05.5.0013.120E	Broca Ø 2,0x Stop 4mmx Comprimento 110mm
05.5.0014.120E	Broca Ø 2,0x Stop 4mmx Comprimento 115mm
05.5.0015.120E	Broca Ø 2,0x Stop 4mmx Comprimento 120mm
05.5.0016.120E	Broca Ø 2,0x Stop 4mmx Comprimento 125mm
05.5.0017.120E	Broca Ø 2,0x Stop 4mmx Comprimento 130mm
05.5.0018.120E	Broca Ø 2,0x Stop 4mmx Comprimento 135mm
05.5.0019.120E	Broca Ø 2,0x Stop 4mmx Comprimento 140mm
05.5.0020.120E	Broca Ø 2,0x Stop 4mmx Comprimento 145mm
05.5.0021.120E	Broca Ø 2,0x Stop 4mmx Comprimento 150mm
05.5.0022.120E	Broca Ø 2,0x Stop 5mmx Comprimento 50mm
05.5.0023.120E	Broca Ø 2,0x Stop 5mmx Comprimento 55mm
05.5.0024.120E	Broca Ø 2,0x Stop 5mmx Comprimento 60mm
05.5.0025.120E	Broca Ø 2,0x Stop 5mmx Comprimento 65mm
05.5.0026.120E	Broca Ø 2,0x Stop 5mmx Comprimento 70mm
05.5.0027.120E	Broca Ø 2,0x Stop 5mmx Comprimento 75mm
05.5.0028.120E	Broca Ø 2,0x Stop 5mmx Comprimento 80mm
05.5.0029.120E	Broca Ø 2,0x Stop 5mmx Comprimento 85mm
05.5.0030.120E	Broca Ø 2,0x Stop 5mmx Comprimento 90mm
05.5.0031.120E	Broca Ø 2,0x Stop 5mmx Comprimento 95mm
05.5.0032.120E	Broca Ø 2,0x Stop 5mmx Comprimento 100mm
05.5.0033.120E	Broca Ø 2,0x Stop 5mmx Comprimento 105mm
05.5.0034.120E	Broca Ø 2,0x Stop 5mmx Comprimento 110mm
05.5.0035.120E	Broca Ø 2,0x Stop 5mmx Comprimento 115mm
05.5.0036.120E	Broca Ø 2,0x Stop 5mmx Comprimento 120mm
05.5.0037.120E	Broca Ø 2,0x Stop 5mmx Comprimento 125mm
05.5.0038.120E	Broca Ø 2,0x Stop 5mmx Comprimento 130mm

05.5.0039.120E	Broca Ø 2,0x Stop 5mmx Comprimento 135mm
05.5.0040.120E	Broca Ø 2,0x Stop 5mmx Comprimento 140mm
05.5.0041.120E	Broca Ø 2,0x Stop 5mmx Comprimento 145mm
05.5.0042.120E	Broca Ø 2,0x Stop 5mmx Comprimento 150mm
05.5.0043.120E	Broca Ø 2,0x Stop 6mmx Comprimento 50mm
05.5.0044.120E	Broca Ø 2,0x Stop 6mmx Comprimento 55mm
05.5.0045.120E	Broca Ø 2,0x Stop 6mmx Comprimento 60mm
05.5.0046.120E	Broca Ø 2,0x Stop 6mmx Comprimento 65mm
05.5.0047.120E	Broca Ø 2,0x Stop 6mmx Comprimento 70mm
05.5.0048.120E	Broca Ø 2,0x Stop 6mmx Comprimento 75mm
05.5.0049.120E	Broca Ø 2,0x Stop 6mmx Comprimento 80mm
05.5.0050.120E	Broca Ø 2,0x Stop 6mmx Comprimento 85mm
05.5.0051.120E	Broca Ø 2,0x Stop 6mmx Comprimento 90mm
05.5.0052.120E	Broca Ø 2,0x Stop 6mmx Comprimento 95mm
05.5.0053.120E	Broca Ø 2,0x Stop 6mmx Comprimento 100mm
05.5.0054.120E	Broca Ø 2,0x Stop 6mmx Comprimento 105mm
05.5.0055.120E	Broca Ø 2,0x Stop 6mmx Comprimento 110mm
05.5.0056.120E	Broca Ø 2,0x Stop 6mmx Comprimento 115mm
05.5.0057.120E	Broca Ø 2,0x Stop 6mmx Comprimento 120mm
05.5.0058.120E	Broca Ø 2,0x Stop 6mmx Comprimento 125mm
05.5.0059.120E	Broca Ø 2,0x Stop 6mmx Comprimento 130mm
05.5.0060.120E	Broca Ø 2,0x Stop 6mmx Comprimento 135mm
05.5.0061.120E	Broca Ø 2,0x Stop 6mmx Comprimento 140mm
05.5.0062.120E	Broca Ø 2,0x Stop 6mmx Comprimento 145mm
05.5.0063.120E	Broca Ø 2,0x Stop 6mmx Comprimento 150mm
05.5.0064.120E	Broca Ø 2,0x Stop 7mmx Comprimento 50mm
05.5.0065.120E	Broca Ø 2,0x Stop 7mmx Comprimento 55mm
05.5.0066.120E	Broca Ø 2,0x Stop 7mmx Comprimento 60mm
05.5.0067.120E	Broca Ø 2,0x Stop 7mmx Comprimento 65mm
05.5.0068.120E	Broca Ø 2,0x Stop 7mmx Comprimento 70mm
05.5.0069.120E	Broca Ø 2,0x Stop 7mmx Comprimento 75mm
05.5.0070.120E	Broca Ø 2,0x Stop 7mmx Comprimento 80mm
05.5.0071.120E	Broca Ø 2,0x Stop 7mmx Comprimento 85mm
05.5.0072.120E	Broca Ø 2,0x Stop 7mmx Comprimento 90mm
05.5.0073.120E	Broca Ø 2,0x Stop 7mmx Comprimento 95mm
05.5.0074.120E	Broca Ø 2,0x Stop 7mmx Comprimento 100mm
05.5.0075.120E	Broca Ø 2,0x Stop 7mmx Comprimento 105mm
05.5.0076.120E	Broca Ø 2,0x Stop 7mmx Comprimento 110mm
05.5.0077.120E	Broca Ø 2,0x Stop 7mmx Comprimento 115mm
05.5.0078.120E	Broca Ø 2,0x Stop 7mmx Comprimento 120mm
05.5.0079.120E	Broca Ø 2,0x Stop 7mmx Comprimento 125mm
05.5.0080.120E	Broca Ø 2,0x Stop 7mmx Comprimento 130mm
05.5.0081.120E	Broca Ø 2,0x Stop 7mmx Comprimento 135mm
05.5.0082.120E	Broca Ø 2,0x Stop 7mmx Comprimento 140mm
05.5.0083.120E	Broca Ø 2,0x Stop 7mmx Comprimento 145mm
05.5.0084.120E	Broca Ø 2,0x Stop 7mmx Comprimento 150mm
05.5.0085.120E	Broca Ø 2,0x Stop 8mmx Comprimento 50mm
05.5.0086.120E	Broca Ø 2,0x Stop 8mmx Comprimento 55mm
05.5.0087.120E	Broca Ø 2,0x Stop 8mmx Comprimento 60mm

05.5.0088.120E	Broca Ø 2,0x Stop 8mmx Comprimento 65mm
05.5.0089.120E	Broca Ø 2,0x Stop 8mmx Comprimento 70mm
05.5.0090.120E	Broca Ø 2,0x Stop 8mmx Comprimento 75mm
05.5.0091.120E	Broca Ø 2,0x Stop 8mmx Comprimento 80mm
05.5.0092.120E	Broca Ø 2,0x Stop 8mmx Comprimento 85mm
05.5.0093.120E	Broca Ø 2,0x Stop 8mmx Comprimento 90mm
05.5.0094.120E	Broca Ø 2,0x Stop 8mmx Comprimento 95mm
05.5.0095.120E	Broca Ø 2,0x Stop 8mmx Comprimento 100mm
05.5.0096.120E	Broca Ø 2,0x Stop 8mmx Comprimento 105mm
05.5.0097.120E	Broca Ø 2,0x Stop 8mmx Comprimento 110mm
05.5.0098.120E	Broca Ø 2,0x Stop 8mmx Comprimento 115mm
05.5.0099.120E	Broca Ø 2,0x Stop 8mmx Comprimento 120mm
05.5.0100.120E	Broca Ø 2,0x Stop 8mmx Comprimento 125mm
05.5.0101.120E	Broca Ø 2,0x Stop 8mmx Comprimento 130mm
05.5.0102.120E	Broca Ø 2,0x Stop 8mmx Comprimento 135mm
05.5.0103.120E	Broca Ø 2,0x Stop 8mmx Comprimento 140mm
05.5.0104.120E	Broca Ø 2,0x Stop 8mmx Comprimento 145mm
05.5.0105.120E	Broca Ø 2,0x Stop 8mmx Comprimento 150mm
05.5.0106.120E	Broca Ø 2,0x Stop 9mmx Comprimento 50mm
05.5.0107.120E	Broca Ø 2,0x Stop 9mmx Comprimento 55mm
05.5.0108.120E	Broca Ø 2,0x Stop 9mmx Comprimento 60mm
05.5.0109.120E	Broca Ø 2,0x Stop 9mmx Comprimento 65mm
05.5.0110.120E	Broca Ø 2,0x Stop 9mmx Comprimento 70mm
05.5.0111.120E	Broca Ø 2,0x Stop 9mmx Comprimento 75mm
05.5.0112.120E	Broca Ø 2,0x Stop 9mmx Comprimento 80mm
05.5.0113.120E	Broca Ø 2,0x Stop 9mmx Comprimento 85mm
05.5.0114.120E	Broca Ø 2,0x Stop 9mmx Comprimento 90mm
05.5.0115.120E	Broca Ø 2,0x Stop 9mmx Comprimento 95mm
05.5.0116.120E	Broca Ø 2,0x Stop 9mmx Comprimento 100mm
05.5.0117.120E	Broca Ø 2,0x Stop 9mmx Comprimento 105mm
05.5.0118.120E	Broca Ø 2,0x Stop 9mmx Comprimento 110mm
05.5.0119.120E	Broca Ø 2,0x Stop 9mmx Comprimento 115mm
05.5.0120.120E	Broca Ø 2,0x Stop 9mmx Comprimento 120mm
05.5.0121.120E	Broca Ø 2,0x Stop 9mmx Comprimento 125mm
05.5.0122.120E	Broca Ø 2,0x Stop 9mmx Comprimento 130mm
05.5.0123.120E	Broca Ø 2,0x Stop 9mmx Comprimento 135mm
05.5.0124.120E	Broca Ø 2,0x Stop 9mmx Comprimento 140mm
05.5.0125.120E	Broca Ø 2,0x Stop 9mmx Comprimento 145mm
05.5.0126.120E	Broca Ø 2,0x Stop 9mmx Comprimento 150mm
05.5.0127.120E	Broca Ø 2,0x Stop 10mmx Comprimento 50mm
05.5.0128.120E	Broca Ø 2,0x Stop 10mmx Comprimento 55mm
05.5.0129.120E	Broca Ø 2,0x Stop 10mmx Comprimento 60mm
05.5.0130.120E	Broca Ø 2,0x Stop 10mmx Comprimento 65mm
05.5.0131.120E	Broca Ø 2,0x Stop 10mmx Comprimento 70mm
05.5.0132.120E	Broca Ø 2,0x Stop 10mmx Comprimento 75mm
05.5.0133.120E	Broca Ø 2,0x Stop 10mmx Comprimento 80mm
05.5.0134.120E	Broca Ø 2,0x Stop 10mmx Comprimento 85mm
05.5.0135.120E	Broca Ø 2,0x Stop 10mmx Comprimento 90mm
05.5.0136.120E	Broca Ø 2,0x Stop 10mmx Comprimento 95mm

05.5.0137.120E	Broca Ø 2,0x Stop 10mmx Comprimento 100mm
05.5.0138.120E	Broca Ø 2,0x Stop 10mmx Comprimento 105mm
05.5.0139.120E	Broca Ø 2,0x Stop 10mmx Comprimento 110mm
05.5.0140.120E	Broca Ø 2,0x Stop 10mmx Comprimento 115mm
05.5.0141.120E	Broca Ø 2,0x Stop 10mmx Comprimento 120mm
05.5.0142.120E	Broca Ø 2,0x Stop 10mmx Comprimento 125mm
05.5.0143.120E	Broca Ø 2,0x Stop 10mmx Comprimento 130mm
05.5.0144.120E	Broca Ø 2,0x Stop 10mmx Comprimento 135mm
05.5.0145.120E	Broca Ø 2,0x Stop 10mmx Comprimento 140mm
05.5.0146.120E	Broca Ø 2,0x Stop 10mmx Comprimento 145mm
05.5.0147.120E	Broca Ø 2,0x Stop 10mmx Comprimento 150mm
05.5.0148.120E	Broca Ø 2,0x Stop 11mmx Comprimento 50mm
05.5.0149.120E	Broca Ø 2,0x Stop 11mmx Comprimento 55mm
05.5.0150.120E	Broca Ø 2,0x Stop 11mmx Comprimento 60mm
05.5.0151.120E	Broca Ø 2,0x Stop 11mmx Comprimento 65mm
05.5.0152.120E	Broca Ø 2,0x Stop 11mmx Comprimento 70mm
05.5.0153.120E	Broca Ø 2,0x Stop 11mmx Comprimento 75mm
05.5.0154.120E	Broca Ø 2,0x Stop 11mmx Comprimento 80mm
05.5.0155.120E	Broca Ø 2,0x Stop 11mmx Comprimento 85mm
05.5.0156.120E	Broca Ø 2,0x Stop 11mmx Comprimento 90mm
05.5.0157.120E	Broca Ø 2,0x Stop 11mmx Comprimento 95mm
05.5.0158.120E	Broca Ø 2,0x Stop 11mmx Comprimento 100mm
05.5.0159.120E	Broca Ø 2,0x Stop 11mmx Comprimento 105mm
05.5.0160.120E	Broca Ø 2,0x Stop 11mmx Comprimento 110mm
05.5.0161.120E	Broca Ø 2,0x Stop 11mmx Comprimento 115mm
05.5.0162.120E	Broca Ø 2,0x Stop 11mmx Comprimento 120mm
05.5.0163.120E	Broca Ø 2,0x Stop 11mmx Comprimento 125mm
05.5.0164.120E	Broca Ø 2,0x Stop 11mmx Comprimento 130mm
05.5.0165.120E	Broca Ø 2,0x Stop 11mmx Comprimento 135mm
05.5.0166.120E	Broca Ø 2,0x Stop 11mmx Comprimento 140mm
05.5.0167.120E	Broca Ø 2,0x Stop 11mmx Comprimento 145mm
05.5.0168.120E	Broca Ø 2,0x Stop 11mmx Comprimento 150mm
05.5.0169.120E	Broca Ø 2,0x Stop 12mmx Comprimento 50mm
05.5.0170.120E	Broca Ø 2,0x Stop 12mmx Comprimento 55mm
05.5.0171.120E	Broca Ø 2,0x Stop 12mmx Comprimento 60mm
05.5.0172.120E	Broca Ø 2,0x Stop 12mmx Comprimento 65mm
05.5.0173.120E	Broca Ø 2,0x Stop 12mmx Comprimento 70mm
05.5.0174.120E	Broca Ø 2,0x Stop 12mmx Comprimento 75mm
05.5.0175.120E	Broca Ø 2,0x Stop 12mmx Comprimento 80mm
05.5.0176.120E	Broca Ø 2,0x Stop 12mmx Comprimento 85mm
05.5.0177.120E	Broca Ø 2,0x Stop 12mmx Comprimento 90mm
05.5.0178.120E	Broca Ø 2,0x Stop 12mmx Comprimento 95mm
05.5.0179.120E	Broca Ø 2,0x Stop 12mmx Comprimento 100mm
05.5.0180.120E	Broca Ø 2,0x Stop 12mmx Comprimento 105mm
05.5.0181.120E	Broca Ø 2,0x Stop 12mmx Comprimento 110mm
05.5.0182.120E	Broca Ø 2,0x Stop 12mmx Comprimento 115mm
05.5.0183.120E	Broca Ø 2,0x Stop 12mmx Comprimento 120mm
05.5.0184.120E	Broca Ø 2,0x Stop 12mmx Comprimento 125mm
05.5.0185.120E	Broca Ø 2,0x Stop 12mmx Comprimento 130mm

05.5.0186.120E	Broca Ø 2,0x Stop 12mmx Comprimento 135mm
05.5.0187.120E	Broca Ø 2,0x Stop 12mmx Comprimento 140mm
05.5.0188.120E	Broca Ø 2,0x Stop 12mmx Comprimento 145mm
05.5.0189.120E	Broca Ø 2,0x Stop 12mmx Comprimento 150mm
05.5.0190.120E	Broca Ø 2,0x Stop 13mmx Comprimento 50mm
05.5.0191.120E	Broca Ø 2,0x Stop 13mmx Comprimento 55mm
05.5.0192.120E	Broca Ø 2,0x Stop 13mmx Comprimento 60mm
05.5.0193.120E	Broca Ø 2,0x Stop 13mmx Comprimento 65mm
05.5.0194.120E	Broca Ø 2,0x Stop 13mmx Comprimento 70mm
05.5.0195.120E	Broca Ø 2,0x Stop 13mmx Comprimento 75mm
05.5.0196.120E	Broca Ø 2,0x Stop 13mmx Comprimento 80mm
05.5.0197.120E	Broca Ø 2,0x Stop 13mmx Comprimento 85mm
05.5.0198.120E	Broca Ø 2,0x Stop 13mmx Comprimento 90mm
05.5.0199.120E	Broca Ø 2,0x Stop 13mmx Comprimento 95mm
05.5.0200.120E	Broca Ø 2,0x Stop 13mmx Comprimento 100mm
05.5.0201.120E	Broca Ø 2,0x Stop 13mmx Comprimento 105mm
05.5.0202.120E	Broca Ø 2,0x Stop 13mmx Comprimento 110mm
05.5.0203.120E	Broca Ø 2,0x Stop 13mmx Comprimento 115mm
05.5.0204.120E	Broca Ø 2,0x Stop 13mmx Comprimento 120mm
05.5.0205.120E	Broca Ø 2,0x Stop 13mmx Comprimento 125mm
05.5.0206.120E	Broca Ø 2,0x Stop 13mmx Comprimento 130mm
05.5.0207.120E	Broca Ø 2,0x Stop 13mmx Comprimento 135mm
05.5.0208.120E	Broca Ø 2,0x Stop 13mmx Comprimento 140mm
05.5.0209.120E	Broca Ø 2,0x Stop 13mmx Comprimento 145mm
05.5.0210.120E	Broca Ø 2,0x Stop 13mmx Comprimento 150mm
05.5.0211.120E	Broca Ø 2,0x Stop 14mmx Comprimento 50mm
05.5.0212.120E	Broca Ø 2,0x Stop 14mmx Comprimento 55mm
05.5.0213.120E	Broca Ø 2,0x Stop 14mmx Comprimento 60mm
05.5.0214.120E	Broca Ø 2,0x Stop 14mmx Comprimento 65mm
05.5.0215.120E	Broca Ø 2,0x Stop 14mmx Comprimento 70mm
05.5.0216.120E	Broca Ø 2,0x Stop 14mmx Comprimento 75mm
05.5.0217.120E	Broca Ø 2,0x Stop 14mmx Comprimento 80mm
05.5.0218.120E	Broca Ø 2,0x Stop 14mmx Comprimento 85mm
05.5.0219.120E	Broca Ø 2,0x Stop 14mmx Comprimento 90mm
05.5.0220.120E	Broca Ø 2,0x Stop 14mmx Comprimento 95mm
05.5.0221.120E	Broca Ø 2,0x Stop 14mmx Comprimento 100mm
05.5.0222.120E	Broca Ø 2,0x Stop 14mmx Comprimento 105mm
05.5.0223.120E	Broca Ø 2,0x Stop 14mmx Comprimento 110mm
05.5.0224.120E	Broca Ø 2,0x Stop 14mmx Comprimento 115mm
05.5.0225.120E	Broca Ø 2,0x Stop 14mmx Comprimento 120mm
05.5.0226.120E	Broca Ø 2,0x Stop 14mmx Comprimento 125mm
05.5.0227.120E	Broca Ø 2,0x Stop 14mmx Comprimento 130mm
05.5.0228.120E	Broca Ø 2,0x Stop 14mmx Comprimento 135mm
05.5.0229.120E	Broca Ø 2,0x Stop 14mmx Comprimento 140mm
05.5.0230.120E	Broca Ø 2,0x Stop 14mmx Comprimento 145mm
05.5.0231.120E	Broca Ø 2,0x Stop 14mmx Comprimento 150mm
05.5.0232.120E	Broca Ø 2,0x Stop 15mmx Comprimento 50mm
05.5.0233.120E	Broca Ø 2,0x Stop 15mmx Comprimento 55mm
05.5.0234.120E	Broca Ø 2,0x Stop 15mmx Comprimento 60mm

05.5.0235.120E	Broca Ø 2,0x Stop 15mmx Comprimento 65mm
05.5.0236.120E	Broca Ø 2,0x Stop 15mmx Comprimento 70mm
05.5.0237.120E	Broca Ø 2,0x Stop 15mmx Comprimento 75mm
05.5.0238.120E	Broca Ø 2,0x Stop 15mmx Comprimento 80mm
05.5.0239.120E	Broca Ø 2,0x Stop 15mmx Comprimento 85mm
05.5.0240.120E	Broca Ø 2,0x Stop 15mmx Comprimento 90mm
05.5.0241.120E	Broca Ø 2,0x Stop 15mmx Comprimento 95mm
05.5.0242.120E	Broca Ø 2,0x Stop 15mmx Comprimento 100mm
05.5.0243.120E	Broca Ø 2,0x Stop 15mmx Comprimento 105mm
05.5.0244.120E	Broca Ø 2,0x Stop 15mmx Comprimento 110mm
05.5.0245.120E	Broca Ø 2,0x Stop 15mmx Comprimento 115mm
05.5.0246.120E	Broca Ø 2,0x Stop 15mmx Comprimento 120mm
05.5.0247.120E	Broca Ø 2,0x Stop 15mmx Comprimento 125mm
05.5.0248.120E	Broca Ø 2,0x Stop 15mmx Comprimento 130mm
05.5.0249.120E	Broca Ø 2,0x Stop 15mmx Comprimento 135mm
05.5.0250.120E	Broca Ø 2,0x Stop 15mmx Comprimento 140mm
05.5.0251.120E	Broca Ø 2,0x Stop 15mmx Comprimento 145mm
05.5.0252.120E	Broca Ø 2,0x Stop 15mmx Comprimento 150mm
05.5.0253.120E	Broca Ø 2,0x Stop 16mmx Comprimento 50mm
05.5.0254.120E	Broca Ø 2,0x Stop 16mmx Comprimento 55mm
05.5.0255.120E	Broca Ø 2,0x Stop 16mmx Comprimento 60mm
05.5.0256.120E	Broca Ø 2,0x Stop 16mmx Comprimento 65mm
05.5.0257.120E	Broca Ø 2,0x Stop 16mmx Comprimento 70mm
05.5.0258.120E	Broca Ø 2,0x Stop 16mmx Comprimento 75mm
05.5.0259.120E	Broca Ø 2,0x Stop 16mmx Comprimento 80mm
05.5.0260.120E	Broca Ø 2,0x Stop 16mmx Comprimento 85mm
05.5.0261.120E	Broca Ø 2,0x Stop 16mmx Comprimento 90mm
05.5.0262.120E	Broca Ø 2,0x Stop 16mmx Comprimento 95mm
05.5.0263.120E	Broca Ø 2,0x Stop 16mmx Comprimento 100mm
05.5.0264.120E	Broca Ø 2,0x Stop 16mmx Comprimento 105mm
05.5.0265.120E	Broca Ø 2,0x Stop 16mmx Comprimento 110mm
05.5.0266.120E	Broca Ø 2,0x Stop 16mmx Comprimento 115mm
05.5.0267.120E	Broca Ø 2,0x Stop 16mmx Comprimento 120mm
05.5.0268.120E	Broca Ø 2,0x Stop 16mmx Comprimento 125mm
05.5.0269.120E	Broca Ø 2,0x Stop 16mmx Comprimento 130mm
05.5.0270.120E	Broca Ø 2,0x Stop 16mmx Comprimento 135mm
05.5.0271.120E	Broca Ø 2,0x Stop 16mmx Comprimento 140mm
05.5.0272.120E	Broca Ø 2,0x Stop 16mmx Comprimento 145mm
05.5.0273.120E	Broca Ø 2,0x Stop 16mmx Comprimento 150mm
05.5.0274.120E	Broca Ø 2,0x Stop 17mmx Comprimento 50mm
05.5.0275.120E	Broca Ø 2,0x Stop 17mmx Comprimento 55mm
05.5.0276.120E	Broca Ø 2,0x Stop 17mmx Comprimento 60mm
05.5.0277.120E	Broca Ø 2,0x Stop 17mmx Comprimento 65mm
05.5.0278.120E	Broca Ø 2,0x Stop 17mmx Comprimento 70mm
05.5.0279.120E	Broca Ø 2,0x Stop 17mmx Comprimento 75mm
05.5.0280.120E	Broca Ø 2,0x Stop 17mmx Comprimento 80mm
05.5.0281.120E	Broca Ø 2,0x Stop 17mmx Comprimento 85mm
05.5.0282.120E	Broca Ø 2,0x Stop 17mmx Comprimento 90mm
05.5.0283.120E	Broca Ø 2,0x Stop 17mmx Comprimento 95mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0284.120E	Broca Ø 2,0x Stop 17mmx Comprimento 100mm
05.5.0285.120E	Broca Ø 2,0x Stop 17mmx Comprimento 105mm
05.5.0286.120E	Broca Ø 2,0x Stop 17mmx Comprimento 110mm
05.5.0287.120E	Broca Ø 2,0x Stop 17mmx Comprimento 115mm
05.5.0288.120E	Broca Ø 2,0x Stop 17mmx Comprimento 120mm
05.5.0289.120E	Broca Ø 2,0x Stop 17mmx Comprimento 125mm
05.5.0290.120E	Broca Ø 2,0x Stop 17mmx Comprimento 130mm
05.5.0291.120E	Broca Ø 2,0x Stop 17mmx Comprimento 135mm
05.5.0292.120E	Broca Ø 2,0x Stop 17mmx Comprimento 140mm
05.5.0293.120E	Broca Ø 2,0x Stop 17mmx Comprimento 145mm
05.5.0294.120E	Broca Ø 2,0x Stop 17mmx Comprimento 150mm
05.5.0295.120E	Broca Ø 2,0x Stop 18mmx Comprimento 50mm
05.5.0296.120E	Broca Ø 2,0x Stop 18mmx Comprimento 55mm
05.5.0297.120E	Broca Ø 2,0x Stop 18mmx Comprimento 60mm
05.5.0298.120E	Broca Ø 2,0x Stop 18mmx Comprimento 65mm
05.5.0299.120E	Broca Ø 2,0x Stop 18mmx Comprimento 70mm
05.5.0300.120E	Broca Ø 2,0x Stop 18mmx Comprimento 75mm
05.5.0301.120E	Broca Ø 2,0x Stop 18mmx Comprimento 80mm
05.5.0302.120E	Broca Ø 2,0x Stop 18mmx Comprimento 85mm
05.5.0303.120E	Broca Ø 2,0x Stop 18mmx Comprimento 90mm
05.5.0304.120E	Broca Ø 2,0x Stop 18mmx Comprimento 95mm
05.5.0305.120E	Broca Ø 2,0x Stop 18mmx Comprimento 100mm
05.5.0306.120E	Broca Ø 2,0x Stop 18mmx Comprimento 105mm
05.5.0307.120E	Broca Ø 2,0x Stop 18mmx Comprimento 110mm
05.5.0308.120E	Broca Ø 2,0x Stop 18mmx Comprimento 115mm
05.5.0309.120E	Broca Ø 2,0x Stop 18mmx Comprimento 120mm
05.5.0310.120E	Broca Ø 2,0x Stop 18mmx Comprimento 125mm
05.5.0311.120E	Broca Ø 2,0x Stop 18mmx Comprimento 130mm
05.5.0312.120E	Broca Ø 2,0x Stop 18mmx Comprimento 135mm
05.5.0313.120E	Broca Ø 2,0x Stop 18mmx Comprimento 140mm
05.5.0314.120E	Broca Ø 2,0x Stop 18mmx Comprimento 145mm
05.5.0315.120E	Broca Ø 2,0x Stop 18mmx Comprimento 150mm
05.5.0316.120E	Broca Ø 2,0x Stop 19mmx Comprimento 50mm
05.5.0317.120E	Broca Ø 2,0x Stop 19mmx Comprimento 55mm
05.5.0318.120E	Broca Ø 2,0x Stop 19mmx Comprimento 60mm
05.5.0319.120E	Broca Ø 2,0x Stop 19mmx Comprimento 65mm
05.5.0320.120E	Broca Ø 2,0x Stop 19mmx Comprimento 70mm
05.5.0321.120E	Broca Ø 2,0x Stop 19mmx Comprimento 75mm
05.5.0322.120E	Broca Ø 2,0x Stop 19mmx Comprimento 80mm
05.5.0323.120E	Broca Ø 2,0x Stop 19mmx Comprimento 85mm
05.5.0324.120E	Broca Ø 2,0x Stop 19mmx Comprimento 90mm
05.5.0325.120E	Broca Ø 2,0x Stop 19mmx Comprimento 95mm
05.5.0326.120E	Broca Ø 2,0x Stop 19mmx Comprimento 100mm
05.5.0327.120E	Broca Ø 2,0x Stop 19mmx Comprimento 105mm
05.5.0328.120E	Broca Ø 2,0x Stop 19mmx Comprimento 110mm
05.5.0329.120E	Broca Ø 2,0x Stop 19mmx Comprimento 115mm
05.5.0330.120E	Broca Ø 2,0x Stop 19mmx Comprimento 120mm
05.5.0331.120E	Broca Ø 2,0x Stop 19mmx Comprimento 125mm
05.5.0332.120E	Broca Ø 2,0x Stop 19mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0333.120E	Broca Ø 2,0x Stop 19mmx Comprimento 135mm
05.5.0334.120E	Broca Ø 2,0x Stop 19mmx Comprimento 140mm
05.5.0335.120E	Broca Ø 2,0x Stop 19mmx Comprimento 145mm
05.5.0336.120E	Broca Ø 2,0x Stop 19mmx Comprimento 150mm
05.5.0337.120E	Broca Ø 2,0x Stop 20mmx Comprimento 50mm
05.5.0338.120E	Broca Ø 2,0x Stop 20mmx Comprimento 55mm
05.5.0339.120E	Broca Ø 2,0x Stop 20mmx Comprimento 60mm
05.5.0340.120E	Broca Ø 2,0x Stop 20mmx Comprimento 65mm
05.5.0341.120E	Broca Ø 2,0x Stop 20mmx Comprimento 70mm
05.5.0342.120E	Broca Ø 2,0x Stop 20mmx Comprimento 75mm
05.5.0343.120E	Broca Ø 2,0x Stop 20mmx Comprimento 80mm
05.5.0344.120E	Broca Ø 2,0x Stop 20mmx Comprimento 85mm
05.5.0345.120E	Broca Ø 2,0x Stop 20mmx Comprimento 90mm
05.5.0346.120E	Broca Ø 2,0x Stop 20mmx Comprimento 95mm
05.5.0347.120E	Broca Ø 2,0x Stop 20mmx Comprimento 100mm
05.5.0348.120E	Broca Ø 2,0x Stop 20mmx Comprimento 105mm
05.5.0349.120E	Broca Ø 2,0x Stop 20mmx Comprimento 110mm
05.5.0350.120E	Broca Ø 2,0x Stop 20mmx Comprimento 115mm
05.5.0351.120E	Broca Ø 2,0x Stop 20mmx Comprimento 120mm
05.5.0352.120E	Broca Ø 2,0x Stop 20mmx Comprimento 125mm
05.5.0353.120E	Broca Ø 2,0x Stop 20mmx Comprimento 130mm
05.5.0354.120E	Broca Ø 2,0x Stop 20mmx Comprimento 135mm
05.5.0355.120E	Broca Ø 2,0x Stop 20mmx Comprimento 140mm
05.5.0356.120E	Broca Ø 2,0x Stop 20mmx Comprimento 145mm
05.5.0357.120E	Broca Ø 2,0x Stop 20mmx Comprimento 150mm
05.5.0358.120E	Broca Ø 2,0x Stop 21mmx Comprimento 50mm
05.5.0359.120E	Broca Ø 2,0x Stop 21mmx Comprimento 55mm
05.5.0360.120E	Broca Ø 2,0x Stop 21mmx Comprimento 60mm
05.5.0361.120E	Broca Ø 2,0x Stop 21mmx Comprimento 65mm
05.5.0362.120E	Broca Ø 2,0x Stop 21mmx Comprimento 70mm
05.5.0363.120E	Broca Ø 2,0x Stop 21mmx Comprimento 75mm
05.5.0364.120E	Broca Ø 2,0x Stop 21mmx Comprimento 80mm
05.5.0365.120E	Broca Ø 2,0x Stop 21mmx Comprimento 85mm
05.5.0366.120E	Broca Ø 2,0x Stop 21mmx Comprimento 90mm
05.5.0367.120E	Broca Ø 2,0x Stop 21mmx Comprimento 95mm
05.5.0368.120E	Broca Ø 2,0x Stop 21mmx Comprimento 100mm
05.5.0369.120E	Broca Ø 2,0x Stop 21mmx Comprimento 105mm
05.5.0370.120E	Broca Ø 2,0x Stop 21mmx Comprimento 110mm
05.5.0371.120E	Broca Ø 2,0x Stop 21mmx Comprimento 115mm
05.5.0372.120E	Broca Ø 2,0x Stop 21mmx Comprimento 120mm
05.5.0373.120E	Broca Ø 2,0x Stop 21mmx Comprimento 125mm
05.5.0374.120E	Broca Ø 2,0x Stop 21mmx Comprimento 130mm
05.5.0375.120E	Broca Ø 2,0x Stop 21mmx Comprimento 135mm
05.5.0376.120E	Broca Ø 2,0x Stop 21mmx Comprimento 140mm
05.5.0377.120E	Broca Ø 2,0x Stop 21mmx Comprimento 145mm
05.5.0378.120E	Broca Ø 2,0x Stop 21mmx Comprimento 150mm
05.5.0379.120E	Broca Ø 2,0x Stop 22mmx Comprimento 50mm
05.5.0380.120E	Broca Ø 2,0x Stop 22mmx Comprimento 55mm
05.5.0381.120E	Broca Ø 2,0x Stop 22mmx Comprimento 60mm

05.5.0382.120E	Broca Ø 2,0x Stop 22mmx Comprimento 65mm
05.5.0383.120E	Broca Ø 2,0x Stop 22mmx Comprimento 70mm
05.5.0384.120E	Broca Ø 2,0x Stop 22mmx Comprimento 75mm
05.5.0385.120E	Broca Ø 2,0x Stop 22mmx Comprimento 80mm
05.5.0386.120E	Broca Ø 2,0x Stop 22mmx Comprimento 85mm
05.5.0387.120E	Broca Ø 2,0x Stop 22mmx Comprimento 90mm
05.5.0388.120E	Broca Ø 2,0x Stop 22mmx Comprimento 95mm
05.5.0389.120E	Broca Ø 2,0x Stop 22mmx Comprimento 100mm
05.5.0390.120E	Broca Ø 2,0x Stop 22mmx Comprimento 105mm
05.5.0391.120E	Broca Ø 2,0x Stop 22mmx Comprimento 110mm
05.5.0392.120E	Broca Ø 2,0x Stop 22mmx Comprimento 115mm
05.5.0393.120E	Broca Ø 2,0x Stop 22mmx Comprimento 120mm
05.5.0394.120E	Broca Ø 2,0x Stop 22mmx Comprimento 125mm
05.5.0395.120E	Broca Ø 2,0x Stop 22mmx Comprimento 130mm
05.5.0396.120E	Broca Ø 2,0x Stop 22mmx Comprimento 135mm
05.5.0397.120E	Broca Ø 2,0x Stop 22mmx Comprimento 140mm
05.5.0398.120E	Broca Ø 2,0x Stop 22mmx Comprimento 145mm
05.5.0399.120E	Broca Ø 2,0x Stop 22mmx Comprimento 150mm
05.5.0400.120E	Broca Ø 2,0x Stop 23mmx Comprimento 50mm
05.5.0401.120E	Broca Ø 2,0x Stop 23mmx Comprimento 55mm
05.5.0402.120E	Broca Ø 2,0x Stop 23mmx Comprimento 60mm
05.5.0403.120E	Broca Ø 2,0x Stop 23mmx Comprimento 65mm
05.5.0404.120E	Broca Ø 2,0x Stop 23mmx Comprimento 70mm
05.5.0405.120E	Broca Ø 2,0x Stop 23mmx Comprimento 75mm
05.5.0406.120E	Broca Ø 2,0x Stop 23mmx Comprimento 80mm
05.5.0407.120E	Broca Ø 2,0x Stop 23mmx Comprimento 85mm
05.5.0408.120E	Broca Ø 2,0x Stop 23mmx Comprimento 90mm
05.5.0409.120E	Broca Ø 2,0x Stop 23mmx Comprimento 95mm
05.5.0410.120E	Broca Ø 2,0x Stop 23mmx Comprimento 100mm
05.5.0411.120E	Broca Ø 2,0x Stop 23mmx Comprimento 105mm
05.5.0412.120E	Broca Ø 2,0x Stop 23mmx Comprimento 110mm
05.5.0413.120E	Broca Ø 2,0x Stop 23mmx Comprimento 115mm
05.5.0414.120E	Broca Ø 2,0x Stop 23mmx Comprimento 120mm
05.5.0415.120E	Broca Ø 2,0x Stop 23mmx Comprimento 125mm
05.5.0416.120E	Broca Ø 2,0x Stop 23mmx Comprimento 130mm
05.5.0417.120E	Broca Ø 2,0x Stop 23mmx Comprimento 135mm
05.5.0418.120E	Broca Ø 2,0x Stop 23mmx Comprimento 140mm
05.5.0419.120E	Broca Ø 2,0x Stop 23mmx Comprimento 145mm
05.5.0420.120E	Broca Ø 2,0x Stop 23mmx Comprimento 150mm
05.5.0421.120E	Broca Ø 2,0x Stop 24mmx Comprimento 50mm
05.5.0422.120E	Broca Ø 2,0x Stop 24mmx Comprimento 55mm
05.5.0423.120E	Broca Ø 2,0x Stop 24mmx Comprimento 60mm
05.5.0424.120E	Broca Ø 2,0x Stop 24mmx Comprimento 65mm
05.5.0425.120E	Broca Ø 2,0x Stop 24mmx Comprimento 70mm
05.5.0426.120E	Broca Ø 2,0x Stop 24mmx Comprimento 75mm
05.5.0427.120E	Broca Ø 2,0x Stop 24mmx Comprimento 80mm
05.5.0428.120E	Broca Ø 2,0x Stop 24mmx Comprimento 85mm
05.5.0429.120E	Broca Ø 2,0x Stop 24mmx Comprimento 90mm
05.5.0430.120E	Broca Ø 2,0x Stop 24mmx Comprimento 95mm

05.5.0431.120E	Broca Ø 2,0x Stop 24mmx Comprimento 100mm
05.5.0432.120E	Broca Ø 2,0x Stop 24mmx Comprimento 105mm
05.5.0433.120E	Broca Ø 2,0x Stop 24mmx Comprimento 110mm
05.5.0434.120E	Broca Ø 2,0x Stop 24mmx Comprimento 115mm
05.5.0435.120E	Broca Ø 2,0x Stop 24mmx Comprimento 120mm
05.5.0436.120E	Broca Ø 2,0x Stop 24mmx Comprimento 125mm
05.5.0437.120E	Broca Ø 2,0x Stop 24mmx Comprimento 130mm
05.5.0438.120E	Broca Ø 2,0x Stop 24mmx Comprimento 135mm
05.5.0439.120E	Broca Ø 2,0x Stop 24mmx Comprimento 140mm
05.5.0440.120E	Broca Ø 2,0x Stop 24mmx Comprimento 145mm
05.5.0441.120E	Broca Ø 2,0x Stop 24mmx Comprimento 150mm
05.5.0442.120E	Broca Ø 2,0x Stop 25mmx Comprimento 50mm
05.5.0443.120E	Broca Ø 2,0x Stop 25mmx Comprimento 55mm
05.5.0444.120E	Broca Ø 2,0x Stop 25mmx Comprimento 60mm
05.5.0445.120E	Broca Ø 2,0x Stop 25mmx Comprimento 65mm
05.5.0446.120E	Broca Ø 2,0x Stop 25mmx Comprimento 70mm
05.5.0447.120E	Broca Ø 2,0x Stop 25mmx Comprimento 75mm
05.5.0448.120E	Broca Ø 2,0x Stop 25mmx Comprimento 80mm
05.5.0449.120E	Broca Ø 2,0x Stop 25mmx Comprimento 85mm
05.5.0450.120E	Broca Ø 2,0x Stop 25mmx Comprimento 90mm
05.5.0451.120E	Broca Ø 2,0x Stop 25mmx Comprimento 95mm
05.5.0452.120E	Broca Ø 2,0x Stop 25mmx Comprimento 100mm
05.5.0453.120E	Broca Ø 2,0x Stop 25mmx Comprimento 105mm
05.5.0454.120E	Broca Ø 2,0x Stop 25mmx Comprimento 110mm
05.5.0455.120E	Broca Ø 2,0x Stop 25mmx Comprimento 115mm
05.5.0456.120E	Broca Ø 2,0x Stop 25mmx Comprimento 120mm
05.5.0457.120E	Broca Ø 2,0x Stop 25mmx Comprimento 125mm
05.5.0458.120E	Broca Ø 2,0x Stop 25mmx Comprimento 130mm
05.5.0459.120E	Broca Ø 2,0x Stop 25mmx Comprimento 135mm
05.5.0460.120E	Broca Ø 2,0x Stop 25mmx Comprimento 140mm
05.5.0461.120E	Broca Ø 2,0x Stop 25mmx Comprimento 145mm
05.5.0462.120E	Broca Ø 2,0x Stop 25mmx Comprimento 150mm
05.5.0463.120E	Broca Ø 2,0x Stop 30mmx Comprimento 50mm
05.5.0464.120E	Broca Ø 2,0x Stop 30mmx Comprimento 55mm
05.5.0465.120E	Broca Ø 2,0x Stop 30mmx Comprimento 60mm
05.5.0466.120E	Broca Ø 2,0x Stop 30mmx Comprimento 65mm
05.5.0467.120E	Broca Ø 2,0x Stop 30mmx Comprimento 70mm
05.5.0468.120E	Broca Ø 2,0x Stop 30mmx Comprimento 75mm
05.5.0469.120E	Broca Ø 2,0x Stop 30mmx Comprimento 80mm
05.5.0470.120E	Broca Ø 2,0x Stop 30mmx Comprimento 85mm
05.5.0471.120E	Broca Ø 2,0x Stop 30mmx Comprimento 90mm
05.5.0472.120E	Broca Ø 2,0x Stop 30mmx Comprimento 95mm
05.5.0473.120E	Broca Ø 2,0x Stop 30mmx Comprimento 100mm
05.5.0474.120E	Broca Ø 2,0x Stop 30mmx Comprimento 105mm
05.5.0475.120E	Broca Ø 2,0x Stop 30mmx Comprimento 110mm
05.5.0476.120E	Broca Ø 2,0x Stop 30mmx Comprimento 115mm
05.5.0477.120E	Broca Ø 2,0x Stop 30mmx Comprimento 120mm
05.5.0478.120E	Broca Ø 2,0x Stop 30mmx Comprimento 125mm
05.5.0479.120E	Broca Ø 2,0x Stop 30mmx Comprimento 130mm

05.5.0480.120E	Broca Ø 2,0x Stop 30mmx Comprimento 135mm
05.5.0481.120E	Broca Ø 2,0x Stop 30mmx Comprimento 140mm
05.5.0482.120E	Broca Ø 2,0x Stop 30mmx Comprimento 145mm
05.5.0483.120E	Broca Ø 2,0x Stop 30mmx Comprimento 150mm
05.5.0484.120E	Broca Ø 2,0x Stop 35mmx Comprimento 50mm
05.5.0485.120E	Broca Ø 2,0x Stop 35mmx Comprimento 55mm
05.5.0486.120E	Broca Ø 2,0x Stop 35mmx Comprimento 60mm
05.5.0487.120E	Broca Ø 2,0x Stop 35mmx Comprimento 65mm
05.5.0488.120E	Broca Ø 2,0x Stop 35mmx Comprimento 70mm
05.5.0489.120E	Broca Ø 2,0x Stop 35mmx Comprimento 75mm
05.5.0490.120E	Broca Ø 2,0x Stop 35mmx Comprimento 80mm
05.5.0491.120E	Broca Ø 2,0x Stop 35mmx Comprimento 85mm
05.5.0492.120E	Broca Ø 2,0x Stop 35mmx Comprimento 90mm
05.5.0493.120E	Broca Ø 2,0x Stop 35mmx Comprimento 95mm
05.5.0494.120E	Broca Ø 2,0x Stop 35mmx Comprimento 100mm
05.5.0495.120E	Broca Ø 2,0x Stop 35mmx Comprimento 105mm
05.5.0496.120E	Broca Ø 2,0x Stop 35mmx Comprimento 110mm
05.5.0497.120E	Broca Ø 2,0x Stop 35mmx Comprimento 115mm
05.5.0498.120E	Broca Ø 2,0x Stop 35mmx Comprimento 120mm
05.5.0499.120E	Broca Ø 2,0x Stop 35mmx Comprimento 125mm
05.5.0500.120E	Broca Ø 2,0x Stop 35mmx Comprimento 130mm
05.5.0501.120E	Broca Ø 2,0x Stop 35mmx Comprimento 135mm
05.5.0502.120E	Broca Ø 2,0x Stop 35mmx Comprimento 140mm
05.5.0503.120E	Broca Ø 2,0x Stop 35mmx Comprimento 145mm
05.5.0504.120E	Broca Ø 2,0x Stop 35mmx Comprimento 150mm
05.5.0001.121E	Broca Ø 2,1x Stop 4mmx Comprimento 50mm
05.5.0002.121E	Broca Ø 2,1x Stop 4mmx Comprimento 55mm
05.5.0003.121E	Broca Ø 2,1x Stop 4mmx Comprimento 60mm
05.5.0004.121E	Broca Ø 2,1x Stop 4mmx Comprimento 65mm
05.5.0005.121E	Broca Ø 2,1x Stop 4mmx Comprimento 70mm
05.5.0006.121E	Broca Ø 2,1x Stop 4mmx Comprimento 75mm
05.5.0007.121E	Broca Ø 2,1x Stop 4mmx Comprimento 80mm
05.5.0008.121E	Broca Ø 2,1x Stop 4mmx Comprimento 85mm
05.5.0009.121E	Broca Ø 2,1x Stop 4mmx Comprimento 90mm
05.5.0010.121E	Broca Ø 2,1x Stop 4mmx Comprimento 95mm
05.5.0011.121E	Broca Ø 2,1x Stop 4mmx Comprimento 100mm
05.5.0012.121E	Broca Ø 2,1x Stop 4mmx Comprimento 105mm
05.5.0013.121E	Broca Ø 2,1x Stop 4mmx Comprimento 110mm
05.5.0014.121E	Broca Ø 2,1x Stop 4mmx Comprimento 115mm
05.5.0015.121E	Broca Ø 2,1x Stop 4mmx Comprimento 120mm
05.5.0016.121E	Broca Ø 2,1x Stop 4mmx Comprimento 125mm
05.5.0017.121E	Broca Ø 2,1x Stop 4mmx Comprimento 130mm
05.5.0018.121E	Broca Ø 2,1x Stop 4mmx Comprimento 135mm
05.5.0019.121E	Broca Ø 2,1x Stop 4mmx Comprimento 140mm
05.5.0020.121E	Broca Ø 2,1x Stop 4mmx Comprimento 145mm
05.5.0021.121E	Broca Ø 2,1x Stop 4mmx Comprimento 150mm
05.5.0022.121E	Broca Ø 2,1x Stop 5mmx Comprimento 50mm
05.5.0023.121E	Broca Ø 2,1x Stop 5mmx Comprimento 55mm
05.5.0024.121E	Broca Ø 2,1x Stop 5mmx Comprimento 60mm

05.5.0025.121E	Broca Ø 2,1x Stop 5mmx Comprimento 65mm
05.5.0026.121E	Broca Ø 2,1x Stop 5mmx Comprimento 70mm
05.5.0027.121E	Broca Ø 2,1x Stop 5mmx Comprimento 75mm
05.5.0028.121E	Broca Ø 2,1x Stop 5mmx Comprimento 80mm
05.5.0029.121E	Broca Ø 2,1x Stop 5mmx Comprimento 85mm
05.5.0030.121E	Broca Ø 2,1x Stop 5mmx Comprimento 90mm
05.5.0031.121E	Broca Ø 2,1x Stop 5mmx Comprimento 95mm
05.5.0032.121E	Broca Ø 2,1x Stop 5mmx Comprimento 100mm
05.5.0033.121E	Broca Ø 2,1x Stop 5mmx Comprimento 105mm
05.5.0034.121E	Broca Ø 2,1x Stop 5mmx Comprimento 110mm
05.5.0035.121E	Broca Ø 2,1x Stop 5mmx Comprimento 115mm
05.5.0036.121E	Broca Ø 2,1x Stop 5mmx Comprimento 120mm
05.5.0037.121E	Broca Ø 2,1x Stop 5mmx Comprimento 125mm
05.5.0038.121E	Broca Ø 2,1x Stop 5mmx Comprimento 130mm
05.5.0039.121E	Broca Ø 2,1x Stop 5mmx Comprimento 135mm
05.5.0040.121E	Broca Ø 2,1x Stop 5mmx Comprimento 140mm
05.5.0041.121E	Broca Ø 2,1x Stop 5mmx Comprimento 145mm
05.5.0042.121E	Broca Ø 2,1x Stop 5mmx Comprimento 150mm
05.5.0043.121E	Broca Ø 2,1x Stop 6mmx Comprimento 50mm
05.5.0044.121E	Broca Ø 2,1x Stop 6mmx Comprimento 55mm
05.5.0045.121E	Broca Ø 2,1x Stop 6mmx Comprimento 60mm
05.5.0046.121E	Broca Ø 2,1x Stop 6mmx Comprimento 65mm
05.5.0047.121E	Broca Ø 2,1x Stop 6mmx Comprimento 70mm
05.5.0048.121E	Broca Ø 2,1x Stop 6mmx Comprimento 75mm
05.5.0049.121E	Broca Ø 2,1x Stop 6mmx Comprimento 80mm
05.5.0050.121E	Broca Ø 2,1x Stop 6mmx Comprimento 85mm
05.5.0051.121E	Broca Ø 2,1x Stop 6mmx Comprimento 90mm
05.5.0052.121E	Broca Ø 2,1x Stop 6mmx Comprimento 95mm
05.5.0053.121E	Broca Ø 2,1x Stop 6mmx Comprimento 100mm
05.5.0054.121E	Broca Ø 2,1x Stop 6mmx Comprimento 105mm
05.5.0055.121E	Broca Ø 2,1x Stop 6mmx Comprimento 110mm
05.5.0056.121E	Broca Ø 2,1x Stop 6mmx Comprimento 115mm
05.5.0057.121E	Broca Ø 2,1x Stop 6mmx Comprimento 120mm
05.5.0058.121E	Broca Ø 2,1x Stop 6mmx Comprimento 125mm
05.5.0059.121E	Broca Ø 2,1x Stop 6mmx Comprimento 130mm
05.5.0060.121E	Broca Ø 2,1x Stop 6mmx Comprimento 135mm
05.5.0061.121E	Broca Ø 2,1x Stop 6mmx Comprimento 140mm
05.5.0062.121E	Broca Ø 2,1x Stop 6mmx Comprimento 145mm
05.5.0063.121E	Broca Ø 2,1x Stop 6mmx Comprimento 150mm
05.5.0064.121E	Broca Ø 2,1x Stop 7mmx Comprimento 50mm
05.5.0065.121E	Broca Ø 2,1x Stop 7mmx Comprimento 55mm
05.5.0066.121E	Broca Ø 2,1x Stop 7mmx Comprimento 60mm
05.5.0067.121E	Broca Ø 2,1x Stop 7mmx Comprimento 65mm
05.5.0068.121E	Broca Ø 2,1x Stop 7mmx Comprimento 70mm
05.5.0069.121E	Broca Ø 2,1x Stop 7mmx Comprimento 75mm
05.5.0070.121E	Broca Ø 2,1x Stop 7mmx Comprimento 80mm
05.5.0071.121E	Broca Ø 2,1x Stop 7mmx Comprimento 85mm
05.5.0072.121E	Broca Ø 2,1x Stop 7mmx Comprimento 90mm
05.5.0073.121E	Broca Ø 2,1x Stop 7mmx Comprimento 95mm

05.5.0074.121E	Broca Ø 2,1x Stop 7mmx Comprimento 100mm
05.5.0075.121E	Broca Ø 2,1x Stop 7mmx Comprimento 105mm
05.5.0076.121E	Broca Ø 2,1x Stop 7mmx Comprimento 110mm
05.5.0077.121E	Broca Ø 2,1x Stop 7mmx Comprimento 115mm
05.5.0078.121E	Broca Ø 2,1x Stop 7mmx Comprimento 120mm
05.5.0079.121E	Broca Ø 2,1x Stop 7mmx Comprimento 125mm
05.5.0080.121E	Broca Ø 2,1x Stop 7mmx Comprimento 130mm
05.5.0081.121E	Broca Ø 2,1x Stop 7mmx Comprimento 135mm
05.5.0082.121E	Broca Ø 2,1x Stop 7mmx Comprimento 140mm
05.5.0083.121E	Broca Ø 2,1x Stop 7mmx Comprimento 145mm
05.5.0084.121E	Broca Ø 2,1x Stop 7mmx Comprimento 150mm
05.5.0085.121E	Broca Ø 2,1x Stop 8mmx Comprimento 50mm
05.5.0086.121E	Broca Ø 2,1x Stop 8mmx Comprimento 55mm
05.5.0087.121E	Broca Ø 2,1x Stop 8mmx Comprimento 60mm
05.5.0088.121E	Broca Ø 2,1x Stop 8mmx Comprimento 65mm
05.5.0089.121E	Broca Ø 2,1x Stop 8mmx Comprimento 70mm
05.5.0090.121E	Broca Ø 2,1x Stop 8mmx Comprimento 75mm
05.5.0091.121E	Broca Ø 2,1x Stop 8mmx Comprimento 80mm
05.5.0092.121E	Broca Ø 2,1x Stop 8mmx Comprimento 85mm
05.5.0093.121E	Broca Ø 2,1x Stop 8mmx Comprimento 90mm
05.5.0094.121E	Broca Ø 2,1x Stop 8mmx Comprimento 95mm
05.5.0095.121E	Broca Ø 2,1x Stop 8mmx Comprimento 100mm
05.5.0096.121E	Broca Ø 2,1x Stop 8mmx Comprimento 105mm
05.5.0097.121E	Broca Ø 2,1x Stop 8mmx Comprimento 110mm
05.5.0098.121E	Broca Ø 2,1x Stop 8mmx Comprimento 115mm
05.5.0099.121E	Broca Ø 2,1x Stop 8mmx Comprimento 120mm
05.5.0100.121E	Broca Ø 2,1x Stop 8mmx Comprimento 125mm
05.5.0101.121E	Broca Ø 2,1x Stop 8mmx Comprimento 130mm
05.5.0102.121E	Broca Ø 2,1x Stop 8mmx Comprimento 135mm
05.5.0103.121E	Broca Ø 2,1x Stop 8mmx Comprimento 140mm
05.5.0104.121E	Broca Ø 2,1x Stop 8mmx Comprimento 145mm
05.5.0105.121E	Broca Ø 2,1x Stop 8mmx Comprimento 150mm
05.5.0106.121E	Broca Ø 2,1x Stop 9mmx Comprimento 50mm
05.5.0107.121E	Broca Ø 2,1x Stop 9mmx Comprimento 55mm
05.5.0108.121E	Broca Ø 2,1x Stop 9mmx Comprimento 60mm
05.5.0109.121E	Broca Ø 2,1x Stop 9mmx Comprimento 65mm
05.5.0110.121E	Broca Ø 2,1x Stop 9mmx Comprimento 70mm
05.5.0111.121E	Broca Ø 2,1x Stop 9mmx Comprimento 75mm
05.5.0112.121E	Broca Ø 2,1x Stop 9mmx Comprimento 80mm
05.5.0113.121E	Broca Ø 2,1x Stop 9mmx Comprimento 85mm
05.5.0114.121E	Broca Ø 2,1x Stop 9mmx Comprimento 90mm
05.5.0115.121E	Broca Ø 2,1x Stop 9mmx Comprimento 95mm
05.5.0116.121E	Broca Ø 2,1x Stop 9mmx Comprimento 100mm
05.5.0117.121E	Broca Ø 2,1x Stop 9mmx Comprimento 105mm
05.5.0118.121E	Broca Ø 2,1x Stop 9mmx Comprimento 110mm
05.5.0119.121E	Broca Ø 2,1x Stop 9mmx Comprimento 115mm
05.5.0120.121E	Broca Ø 2,1x Stop 9mmx Comprimento 120mm
05.5.0121.121E	Broca Ø 2,1x Stop 9mmx Comprimento 125mm
05.5.0122.121E	Broca Ø 2,1x Stop 9mmx Comprimento 130mm

05.5.0123.121E	Broca Ø 2,1x Stop 9mmx Comprimento 135mm
05.5.0124.121E	Broca Ø 2,1x Stop 9mmx Comprimento 140mm
05.5.0125.121E	Broca Ø 2,1x Stop 9mmx Comprimento 145mm
05.5.0126.121E	Broca Ø 2,1x Stop 9mmx Comprimento 150mm
05.5.0127.121E	Broca Ø 2,1x Stop 10mmx Comprimento 50mm
05.5.0128.121E	Broca Ø 2,1x Stop 10mmx Comprimento 55mm
05.5.0129.121E	Broca Ø 2,1x Stop 10mmx Comprimento 60mm
05.5.0130.121E	Broca Ø 2,1x Stop 10mmx Comprimento 65mm
05.5.0131.121E	Broca Ø 2,1x Stop 10mmx Comprimento 70mm
05.5.0132.121E	Broca Ø 2,1x Stop 10mmx Comprimento 75mm
05.5.0133.121E	Broca Ø 2,1x Stop 10mmx Comprimento 80mm
05.5.0134.121E	Broca Ø 2,1x Stop 10mmx Comprimento 85mm
05.5.0135.121E	Broca Ø 2,1x Stop 10mmx Comprimento 90mm
05.5.0136.121E	Broca Ø 2,1x Stop 10mmx Comprimento 95mm
05.5.0137.121E	Broca Ø 2,1x Stop 10mmx Comprimento 100mm
05.5.0138.121E	Broca Ø 2,1x Stop 10mmx Comprimento 105mm
05.5.0139.121E	Broca Ø 2,1x Stop 10mmx Comprimento 110mm
05.5.0140.121E	Broca Ø 2,1x Stop 10mmx Comprimento 115mm
05.5.0141.121E	Broca Ø 2,1x Stop 10mmx Comprimento 120mm
05.5.0142.121E	Broca Ø 2,1x Stop 10mmx Comprimento 125mm
05.5.0143.121E	Broca Ø 2,1x Stop 10mmx Comprimento 130mm
05.5.0144.121E	Broca Ø 2,1x Stop 10mmx Comprimento 135mm
05.5.0145.121E	Broca Ø 2,1x Stop 10mmx Comprimento 140mm
05.5.0146.121E	Broca Ø 2,1x Stop 10mmx Comprimento 145mm
05.5.0147.121E	Broca Ø 2,1x Stop 10mmx Comprimento 150mm
05.5.0148.121E	Broca Ø 2,1x Stop 11mmx Comprimento 50mm
05.5.0149.121E	Broca Ø 2,1x Stop 11mmx Comprimento 55mm
05.5.0150.121E	Broca Ø 2,1x Stop 11mmx Comprimento 60mm
05.5.0151.121E	Broca Ø 2,1x Stop 11mmx Comprimento 65mm
05.5.0152.121E	Broca Ø 2,1x Stop 11mmx Comprimento 70mm
05.5.0153.121E	Broca Ø 2,1x Stop 11mmx Comprimento 75mm
05.5.0154.121E	Broca Ø 2,1x Stop 11mmx Comprimento 80mm
05.5.0155.121E	Broca Ø 2,1x Stop 11mmx Comprimento 85mm
05.5.0156.121E	Broca Ø 2,1x Stop 11mmx Comprimento 90mm
05.5.0157.121E	Broca Ø 2,1x Stop 11mmx Comprimento 95mm
05.5.0158.121E	Broca Ø 2,1x Stop 11mmx Comprimento 100mm
05.5.0159.121E	Broca Ø 2,1x Stop 11mmx Comprimento 105mm
05.5.0160.121E	Broca Ø 2,1x Stop 11mmx Comprimento 110mm
05.5.0161.121E	Broca Ø 2,1x Stop 11mmx Comprimento 115mm
05.5.0162.121E	Broca Ø 2,1x Stop 11mmx Comprimento 120mm
05.5.0163.121E	Broca Ø 2,1x Stop 11mmx Comprimento 125mm
05.5.0164.121E	Broca Ø 2,1x Stop 11mmx Comprimento 130mm
05.5.0165.121E	Broca Ø 2,1x Stop 11mmx Comprimento 135mm
05.5.0166.121E	Broca Ø 2,1x Stop 11mmx Comprimento 140mm
05.5.0167.121E	Broca Ø 2,1x Stop 11mmx Comprimento 145mm
05.5.0168.121E	Broca Ø 2,1x Stop 11mmx Comprimento 150mm
05.5.0169.121E	Broca Ø 2,1x Stop 12mmx Comprimento 50mm
05.5.0170.121E	Broca Ø 2,1x Stop 12mmx Comprimento 55mm
05.5.0171.121E	Broca Ø 2,1x Stop 12mmx Comprimento 60mm

05.5.0172.121E	Broca Ø 2,1x Stop 12mmx Comprimento 65mm
05.5.0173.121E	Broca Ø 2,1x Stop 12mmx Comprimento 70mm
05.5.0174.121E	Broca Ø 2,1x Stop 12mmx Comprimento 75mm
05.5.0175.121E	Broca Ø 2,1x Stop 12mmx Comprimento 80mm
05.5.0176.121E	Broca Ø 2,1x Stop 12mmx Comprimento 85mm
05.5.0177.121E	Broca Ø 2,1x Stop 12mmx Comprimento 90mm
05.5.0178.121E	Broca Ø 2,1x Stop 12mmx Comprimento 95mm
05.5.0179.121E	Broca Ø 2,1x Stop 12mmx Comprimento 100mm
05.5.0180.121E	Broca Ø 2,1x Stop 12mmx Comprimento 105mm
05.5.0181.121E	Broca Ø 2,1x Stop 12mmx Comprimento 110mm
05.5.0182.121E	Broca Ø 2,1x Stop 12mmx Comprimento 115mm
05.5.0183.121E	Broca Ø 2,1x Stop 12mmx Comprimento 120mm
05.5.0184.121E	Broca Ø 2,1x Stop 12mmx Comprimento 125mm
05.5.0185.121E	Broca Ø 2,1x Stop 12mmx Comprimento 130mm
05.5.0186.121E	Broca Ø 2,1x Stop 12mmx Comprimento 135mm
05.5.0187.121E	Broca Ø 2,1x Stop 12mmx Comprimento 140mm
05.5.0188.121E	Broca Ø 2,1x Stop 12mmx Comprimento 145mm
05.5.0189.121E	Broca Ø 2,1x Stop 12mmx Comprimento 150mm
05.5.0190.121E	Broca Ø 2,1x Stop 13mmx Comprimento 50mm
05.5.0191.121E	Broca Ø 2,1x Stop 13mmx Comprimento 55mm
05.5.0192.121E	Broca Ø 2,1x Stop 13mmx Comprimento 60mm
05.5.0193.121E	Broca Ø 2,1x Stop 13mmx Comprimento 65mm
05.5.0194.121E	Broca Ø 2,1x Stop 13mmx Comprimento 70mm
05.5.0195.121E	Broca Ø 2,1x Stop 13mmx Comprimento 75mm
05.5.0196.121E	Broca Ø 2,1x Stop 13mmx Comprimento 80mm
05.5.0197.121E	Broca Ø 2,1x Stop 13mmx Comprimento 85mm
05.5.0198.121E	Broca Ø 2,1x Stop 13mmx Comprimento 90mm
05.5.0199.121E	Broca Ø 2,1x Stop 13mmx Comprimento 95mm
05.5.0200.121E	Broca Ø 2,1x Stop 13mmx Comprimento 100mm
05.5.0201.121E	Broca Ø 2,1x Stop 13mmx Comprimento 105mm
05.5.0202.121E	Broca Ø 2,1x Stop 13mmx Comprimento 110mm
05.5.0203.121E	Broca Ø 2,1x Stop 13mmx Comprimento 115mm
05.5.0204.121E	Broca Ø 2,1x Stop 13mmx Comprimento 120mm
05.5.0205.121E	Broca Ø 2,1x Stop 13mmx Comprimento 125mm
05.5.0206.121E	Broca Ø 2,1x Stop 13mmx Comprimento 130mm
05.5.0207.121E	Broca Ø 2,1x Stop 13mmx Comprimento 135mm
05.5.0208.121E	Broca Ø 2,1x Stop 13mmx Comprimento 140mm
05.5.0209.121E	Broca Ø 2,1x Stop 13mmx Comprimento 145mm
05.5.0210.121E	Broca Ø 2,1x Stop 13mmx Comprimento 150mm
05.5.0211.121E	Broca Ø 2,1x Stop 14mmx Comprimento 50mm
05.5.0212.121E	Broca Ø 2,1x Stop 14mmx Comprimento 55mm
05.5.0213.121E	Broca Ø 2,1x Stop 14mmx Comprimento 60mm
05.5.0214.121E	Broca Ø 2,1x Stop 14mmx Comprimento 65mm
05.5.0215.121E	Broca Ø 2,1x Stop 14mmx Comprimento 70mm
05.5.0216.121E	Broca Ø 2,1x Stop 14mmx Comprimento 75mm
05.5.0217.121E	Broca Ø 2,1x Stop 14mmx Comprimento 80mm
05.5.0218.121E	Broca Ø 2,1x Stop 14mmx Comprimento 85mm
05.5.0219.121E	Broca Ø 2,1x Stop 14mmx Comprimento 90mm
05.5.0220.121E	Broca Ø 2,1x Stop 14mmx Comprimento 95mm

05.5.0221.121E	Broca Ø 2,1x Stop 14mmx Comprimento 100mm
05.5.0222.121E	Broca Ø 2,1x Stop 14mmx Comprimento 105mm
05.5.0223.121E	Broca Ø 2,1x Stop 14mmx Comprimento 110mm
05.5.0224.121E	Broca Ø 2,1x Stop 14mmx Comprimento 115mm
05.5.0225.121E	Broca Ø 2,1x Stop 14mmx Comprimento 120mm
05.5.0226.121E	Broca Ø 2,1x Stop 14mmx Comprimento 125mm
05.5.0227.121E	Broca Ø 2,1x Stop 14mmx Comprimento 130mm
05.5.0228.121E	Broca Ø 2,1x Stop 14mmx Comprimento 135mm
05.5.0229.121E	Broca Ø 2,1x Stop 14mmx Comprimento 140mm
05.5.0230.121E	Broca Ø 2,1x Stop 14mmx Comprimento 145mm
05.5.0231.121E	Broca Ø 2,1x Stop 14mmx Comprimento 150mm
05.5.0232.121E	Broca Ø 2,1x Stop 15mmx Comprimento 50mm
05.5.0233.121E	Broca Ø 2,1x Stop 15mmx Comprimento 55mm
05.5.0234.121E	Broca Ø 2,1x Stop 15mmx Comprimento 60mm
05.5.0235.121E	Broca Ø 2,1x Stop 15mmx Comprimento 65mm
05.5.0236.121E	Broca Ø 2,1x Stop 15mmx Comprimento 70mm
05.5.0237.121E	Broca Ø 2,1x Stop 15mmx Comprimento 75mm
05.5.0238.121E	Broca Ø 2,1x Stop 15mmx Comprimento 80mm
05.5.0239.121E	Broca Ø 2,1x Stop 15mmx Comprimento 85mm
05.5.0240.121E	Broca Ø 2,1x Stop 15mmx Comprimento 90mm
05.5.0241.121E	Broca Ø 2,1x Stop 15mmx Comprimento 95mm
05.5.0242.121E	Broca Ø 2,1x Stop 15mmx Comprimento 100mm
05.5.0243.121E	Broca Ø 2,1x Stop 15mmx Comprimento 105mm
05.5.0244.121E	Broca Ø 2,1x Stop 15mmx Comprimento 110mm
05.5.0245.121E	Broca Ø 2,1x Stop 15mmx Comprimento 115mm
05.5.0246.121E	Broca Ø 2,1x Stop 15mmx Comprimento 120mm
05.5.0247.121E	Broca Ø 2,1x Stop 15mmx Comprimento 125mm
05.5.0248.121E	Broca Ø 2,1x Stop 15mmx Comprimento 130mm
05.5.0249.121E	Broca Ø 2,1x Stop 15mmx Comprimento 135mm
05.5.0250.121E	Broca Ø 2,1x Stop 15mmx Comprimento 140mm
05.5.0251.121E	Broca Ø 2,1x Stop 15mmx Comprimento 145mm
05.5.0252.121E	Broca Ø 2,1x Stop 15mmx Comprimento 150mm
05.5.0253.121E	Broca Ø 2,1x Stop 16mmx Comprimento 50mm
05.5.0254.121E	Broca Ø 2,1x Stop 16mmx Comprimento 55mm
05.5.0255.121E	Broca Ø 2,1x Stop 16mmx Comprimento 60mm
05.5.0256.121E	Broca Ø 2,1x Stop 16mmx Comprimento 65mm
05.5.0257.121E	Broca Ø 2,1x Stop 16mmx Comprimento 70mm
05.5.0258.121E	Broca Ø 2,1x Stop 16mmx Comprimento 75mm
05.5.0259.121E	Broca Ø 2,1x Stop 16mmx Comprimento 80mm
05.5.0260.121E	Broca Ø 2,1x Stop 16mmx Comprimento 85mm
05.5.0261.121E	Broca Ø 2,1x Stop 16mmx Comprimento 90mm
05.5.0262.121E	Broca Ø 2,1x Stop 16mmx Comprimento 95mm
05.5.0263.121E	Broca Ø 2,1x Stop 16mmx Comprimento 100mm
05.5.0264.121E	Broca Ø 2,1x Stop 16mmx Comprimento 105mm
05.5.0265.121E	Broca Ø 2,1x Stop 16mmx Comprimento 110mm
05.5.0266.121E	Broca Ø 2,1x Stop 16mmx Comprimento 115mm
05.5.0267.121E	Broca Ø 2,1x Stop 16mmx Comprimento 120mm
05.5.0268.121E	Broca Ø 2,1x Stop 16mmx Comprimento 125mm
05.5.0269.121E	Broca Ø 2,1x Stop 16mmx Comprimento 130mm

05.5.0270.121E	Broca Ø 2,1x Stop 16mmx Comprimento 135mm
05.5.0271.121E	Broca Ø 2,1x Stop 16mmx Comprimento 140mm
05.5.0272.121E	Broca Ø 2,1x Stop 16mmx Comprimento 145mm
05.5.0273.121E	Broca Ø 2,1x Stop 16mmx Comprimento 150mm
05.5.0274.121E	Broca Ø 2,1x Stop 17mmx Comprimento 50mm
05.5.0275.121E	Broca Ø 2,1x Stop 17mmx Comprimento 55mm
05.5.0276.121E	Broca Ø 2,1x Stop 17mmx Comprimento 60mm
05.5.0277.121E	Broca Ø 2,1x Stop 17mmx Comprimento 65mm
05.5.0278.121E	Broca Ø 2,1x Stop 17mmx Comprimento 70mm
05.5.0279.121E	Broca Ø 2,1x Stop 17mmx Comprimento 75mm
05.5.0280.121E	Broca Ø 2,1x Stop 17mmx Comprimento 80mm
05.5.0281.121E	Broca Ø 2,1x Stop 17mmx Comprimento 85mm
05.5.0282.121E	Broca Ø 2,1x Stop 17mmx Comprimento 90mm
05.5.0283.121E	Broca Ø 2,1x Stop 17mmx Comprimento 95mm
05.5.0284.121E	Broca Ø 2,1x Stop 17mmx Comprimento 100mm
05.5.0285.121E	Broca Ø 2,1x Stop 17mmx Comprimento 105mm
05.5.0286.121E	Broca Ø 2,1x Stop 17mmx Comprimento 110mm
05.5.0287.121E	Broca Ø 2,1x Stop 17mmx Comprimento 115mm
05.5.0288.121E	Broca Ø 2,1x Stop 17mmx Comprimento 120mm
05.5.0289.121E	Broca Ø 2,1x Stop 17mmx Comprimento 125mm
05.5.0290.121E	Broca Ø 2,1x Stop 17mmx Comprimento 130mm
05.5.0291.121E	Broca Ø 2,1x Stop 17mmx Comprimento 135mm
05.5.0292.121E	Broca Ø 2,1x Stop 17mmx Comprimento 140mm
05.5.0293.121E	Broca Ø 2,1x Stop 17mmx Comprimento 145mm
05.5.0294.121E	Broca Ø 2,1x Stop 17mmx Comprimento 150mm
05.5.0295.121E	Broca Ø 2,1x Stop 18mmx Comprimento 50mm
05.5.0296.121E	Broca Ø 2,1x Stop 18mmx Comprimento 55mm
05.5.0297.121E	Broca Ø 2,1x Stop 18mmx Comprimento 60mm
05.5.0298.121E	Broca Ø 2,1x Stop 18mmx Comprimento 65mm
05.5.0299.121E	Broca Ø 2,1x Stop 18mmx Comprimento 70mm
05.5.0300.121E	Broca Ø 2,1x Stop 18mmx Comprimento 75mm
05.5.0301.121E	Broca Ø 2,1x Stop 18mmx Comprimento 80mm
05.5.0302.121E	Broca Ø 2,1x Stop 18mmx Comprimento 85mm
05.5.0303.121E	Broca Ø 2,1x Stop 18mmx Comprimento 90mm
05.5.0304.121E	Broca Ø 2,1x Stop 18mmx Comprimento 95mm
05.5.0305.121E	Broca Ø 2,1x Stop 18mmx Comprimento 100mm
05.5.0306.121E	Broca Ø 2,1x Stop 18mmx Comprimento 105mm
05.5.0307.121E	Broca Ø 2,1x Stop 18mmx Comprimento 110mm
05.5.0308.121E	Broca Ø 2,1x Stop 18mmx Comprimento 115mm
05.5.0309.121E	Broca Ø 2,1x Stop 18mmx Comprimento 120mm
05.5.0310.121E	Broca Ø 2,1x Stop 18mmx Comprimento 125mm
05.5.0311.121E	Broca Ø 2,1x Stop 18mmx Comprimento 130mm
05.5.0312.121E	Broca Ø 2,1x Stop 18mmx Comprimento 135mm
05.5.0313.121E	Broca Ø 2,1x Stop 18mmx Comprimento 140mm
05.5.0314.121E	Broca Ø 2,1x Stop 18mmx Comprimento 145mm
05.5.0315.121E	Broca Ø 2,1x Stop 18mmx Comprimento 150mm
05.5.0316.121E	Broca Ø 2,1x Stop 19mmx Comprimento 50mm
05.5.0317.121E	Broca Ø 2,1x Stop 19mmx Comprimento 55mm
05.5.0318.121E	Broca Ø 2,1x Stop 19mmx Comprimento 60mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0319.121E	Broca Ø 2,1x Stop 19mmx Comprimento 65mm
05.5.0320.121E	Broca Ø 2,1x Stop 19mmx Comprimento 70mm
05.5.0321.121E	Broca Ø 2,1x Stop 19mmx Comprimento 75mm
05.5.0322.121E	Broca Ø 2,1x Stop 19mmx Comprimento 80mm
05.5.0323.121E	Broca Ø 2,1x Stop 19mmx Comprimento 85mm
05.5.0324.121E	Broca Ø 2,1x Stop 19mmx Comprimento 90mm
05.5.0325.121E	Broca Ø 2,1x Stop 19mmx Comprimento 95mm
05.5.0326.121E	Broca Ø 2,1x Stop 19mmx Comprimento 100mm
05.5.0327.121E	Broca Ø 2,1x Stop 19mmx Comprimento 105mm
05.5.0328.121E	Broca Ø 2,1x Stop 19mmx Comprimento 110mm
05.5.0329.121E	Broca Ø 2,1x Stop 19mmx Comprimento 115mm
05.5.0330.121E	Broca Ø 2,1x Stop 19mmx Comprimento 120mm
05.5.0331.121E	Broca Ø 2,1x Stop 19mmx Comprimento 125mm
05.5.0332.121E	Broca Ø 2,1x Stop 19mmx Comprimento 130mm
05.5.0333.121E	Broca Ø 2,1x Stop 19mmx Comprimento 135mm
05.5.0334.121E	Broca Ø 2,1x Stop 19mmx Comprimento 140mm
05.5.0335.121E	Broca Ø 2,1x Stop 19mmx Comprimento 145mm
05.5.0336.121E	Broca Ø 2,1x Stop 19mmx Comprimento 150mm
05.5.0337.121E	Broca Ø 2,1x Stop 20mmx Comprimento 50mm
05.5.0338.121E	Broca Ø 2,1x Stop 20mmx Comprimento 55mm
05.5.0339.121E	Broca Ø 2,1x Stop 20mmx Comprimento 60mm
05.5.0340.121E	Broca Ø 2,1x Stop 20mmx Comprimento 65mm
05.5.0341.121E	Broca Ø 2,1x Stop 20mmx Comprimento 70mm
05.5.0342.121E	Broca Ø 2,1x Stop 20mmx Comprimento 75mm
05.5.0343.121E	Broca Ø 2,1x Stop 20mmx Comprimento 80mm
05.5.0344.121E	Broca Ø 2,1x Stop 20mmx Comprimento 85mm
05.5.0345.121E	Broca Ø 2,1x Stop 20mmx Comprimento 90mm
05.5.0346.121E	Broca Ø 2,1x Stop 20mmx Comprimento 95mm
05.5.0347.121E	Broca Ø 2,1x Stop 20mmx Comprimento 100mm
05.5.0348.121E	Broca Ø 2,1x Stop 20mmx Comprimento 105mm
05.5.0349.121E	Broca Ø 2,1x Stop 20mmx Comprimento 110mm
05.5.0350.121E	Broca Ø 2,1x Stop 20mmx Comprimento 115mm
05.5.0351.121E	Broca Ø 2,1x Stop 20mmx Comprimento 120mm
05.5.0352.121E	Broca Ø 2,1x Stop 20mmx Comprimento 125mm
05.5.0353.121E	Broca Ø 2,1x Stop 20mmx Comprimento 130mm
05.5.0354.121E	Broca Ø 2,1x Stop 20mmx Comprimento 135mm
05.5.0355.121E	Broca Ø 2,1x Stop 20mmx Comprimento 140mm
05.5.0356.121E	Broca Ø 2,1x Stop 20mmx Comprimento 145mm
05.5.0357.121E	Broca Ø 2,1x Stop 20mmx Comprimento 150mm
05.5.0358.121E	Broca Ø 2,1x Stop 21mmx Comprimento 50mm
05.5.0359.121E	Broca Ø 2,1x Stop 21mmx Comprimento 55mm
05.5.0360.121E	Broca Ø 2,1x Stop 21mmx Comprimento 60mm
05.5.0361.121E	Broca Ø 2,1x Stop 21mmx Comprimento 65mm
05.5.0362.121E	Broca Ø 2,1x Stop 21mmx Comprimento 70mm
05.5.0363.121E	Broca Ø 2,1x Stop 21mmx Comprimento 75mm
05.5.0364.121E	Broca Ø 2,1x Stop 21mmx Comprimento 80mm
05.5.0365.121E	Broca Ø 2,1x Stop 21mmx Comprimento 85mm
05.5.0366.121E	Broca Ø 2,1x Stop 21mmx Comprimento 90mm
05.5.0367.121E	Broca Ø 2,1x Stop 21mmx Comprimento 95mm

05.5.0368.121E	Broca Ø 2,1x Stop 21mmx Comprimento 100mm
05.5.0369.121E	Broca Ø 2,1x Stop 21mmx Comprimento 105mm
05.5.0370.121E	Broca Ø 2,1x Stop 21mmx Comprimento 110mm
05.5.0371.121E	Broca Ø 2,1x Stop 21mmx Comprimento 115mm
05.5.0372.121E	Broca Ø 2,1x Stop 21mmx Comprimento 120mm
05.5.0373.121E	Broca Ø 2,1x Stop 21mmx Comprimento 125mm
05.5.0374.121E	Broca Ø 2,1x Stop 21mmx Comprimento 130mm
05.5.0375.121E	Broca Ø 2,1x Stop 21mmx Comprimento 135mm
05.5.0376.121E	Broca Ø 2,1x Stop 21mmx Comprimento 140mm
05.5.0377.121E	Broca Ø 2,1x Stop 21mmx Comprimento 145mm
05.5.0378.121E	Broca Ø 2,1x Stop 21mmx Comprimento 150mm
05.5.0379.121E	Broca Ø 2,1x Stop 22mmx Comprimento 50mm
05.5.0380.121E	Broca Ø 2,1x Stop 22mmx Comprimento 55mm
05.5.0381.121E	Broca Ø 2,1x Stop 22mmx Comprimento 60mm
05.5.0382.121E	Broca Ø 2,1x Stop 22mmx Comprimento 65mm
05.5.0383.121E	Broca Ø 2,1x Stop 22mmx Comprimento 70mm
05.5.0384.121E	Broca Ø 2,1x Stop 22mmx Comprimento 75mm
05.5.0385.121E	Broca Ø 2,1x Stop 22mmx Comprimento 80mm
05.5.0386.121E	Broca Ø 2,1x Stop 22mmx Comprimento 85mm
05.5.0387.121E	Broca Ø 2,1x Stop 22mmx Comprimento 90mm
05.5.0388.121E	Broca Ø 2,1x Stop 22mmx Comprimento 95mm
05.5.0389.121E	Broca Ø 2,1x Stop 22mmx Comprimento 100mm
05.5.0390.121E	Broca Ø 2,1x Stop 22mmx Comprimento 105mm
05.5.0391.121E	Broca Ø 2,1x Stop 22mmx Comprimento 110mm
05.5.0392.121E	Broca Ø 2,1x Stop 22mmx Comprimento 115mm
05.5.0393.121E	Broca Ø 2,1x Stop 22mmx Comprimento 120mm
05.5.0394.121E	Broca Ø 2,1x Stop 22mmx Comprimento 125mm
05.5.0395.121E	Broca Ø 2,1x Stop 22mmx Comprimento 130mm
05.5.0396.121E	Broca Ø 2,1x Stop 22mmx Comprimento 135mm
05.5.0397.121E	Broca Ø 2,1x Stop 22mmx Comprimento 140mm
05.5.0398.121E	Broca Ø 2,1x Stop 22mmx Comprimento 145mm
05.5.0399.121E	Broca Ø 2,1x Stop 22mmx Comprimento 150mm
05.5.0400.121E	Broca Ø 2,1x Stop 23mmx Comprimento 50mm
05.5.0401.121E	Broca Ø 2,1x Stop 23mmx Comprimento 55mm
05.5.0402.121E	Broca Ø 2,1x Stop 23mmx Comprimento 60mm
05.5.0403.121E	Broca Ø 2,1x Stop 23mmx Comprimento 65mm
05.5.0404.121E	Broca Ø 2,1x Stop 23mmx Comprimento 70mm
05.5.0405.121E	Broca Ø 2,1x Stop 23mmx Comprimento 75mm
05.5.0406.121E	Broca Ø 2,1x Stop 23mmx Comprimento 80mm
05.5.0407.121E	Broca Ø 2,1x Stop 23mmx Comprimento 85mm
05.5.0408.121E	Broca Ø 2,1x Stop 23mmx Comprimento 90mm
05.5.0409.121E	Broca Ø 2,1x Stop 23mmx Comprimento 95mm
05.5.0410.121E	Broca Ø 2,1x Stop 23mmx Comprimento 100mm
05.5.0411.121E	Broca Ø 2,1x Stop 23mmx Comprimento 105mm
05.5.0412.121E	Broca Ø 2,1x Stop 23mmx Comprimento 110mm
05.5.0413.121E	Broca Ø 2,1x Stop 23mmx Comprimento 115mm
05.5.0414.121E	Broca Ø 2,1x Stop 23mmx Comprimento 120mm
05.5.0415.121E	Broca Ø 2,1x Stop 23mmx Comprimento 125mm
05.5.0416.121E	Broca Ø 2,1x Stop 23mmx Comprimento 130mm

INSTRUÇÃO DE USO



Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003

05.5.0417.121E	Broca Ø 2,1x Stop 23mmx Comprimento 135mm
05.5.0418.121E	Broca Ø 2,1x Stop 23mmx Comprimento 140mm
05.5.0419.121E	Broca Ø 2,1x Stop 23mmx Comprimento 145mm
05.5.0420.121E	Broca Ø 2,1x Stop 23mmx Comprimento 150mm
05.5.0421.121E	Broca Ø 2,1x Stop 24mmx Comprimento 50mm
05.5.0422.121E	Broca Ø 2,1x Stop 24mmx Comprimento 55mm
05.5.0423.121E	Broca Ø 2,1x Stop 24mmx Comprimento 60mm
05.5.0424.121E	Broca Ø 2,1x Stop 24mmx Comprimento 65mm
05.5.0425.121E	Broca Ø 2,1x Stop 24mmx Comprimento 70mm
05.5.0426.121E	Broca Ø 2,1x Stop 24mmx Comprimento 75mm
05.5.0427.121E	Broca Ø 2,1x Stop 24mmx Comprimento 80mm
05.5.0428.121E	Broca Ø 2,1x Stop 24mmx Comprimento 85mm
05.5.0429.121E	Broca Ø 2,1x Stop 24mmx Comprimento 90mm
05.5.0430.121E	Broca Ø 2,1x Stop 24mmx Comprimento 95mm
05.5.0431.121E	Broca Ø 2,1x Stop 24mmx Comprimento 100mm
05.5.0432.121E	Broca Ø 2,1x Stop 24mmx Comprimento 105mm
05.5.0433.121E	Broca Ø 2,1x Stop 24mmx Comprimento 110mm
05.5.0434.121E	Broca Ø 2,1x Stop 24mmx Comprimento 115mm
05.5.0435.121E	Broca Ø 2,1x Stop 24mmx Comprimento 120mm
05.5.0436.121E	Broca Ø 2,1x Stop 24mmx Comprimento 125mm
05.5.0437.121E	Broca Ø 2,1x Stop 24mmx Comprimento 130mm
05.5.0438.121E	Broca Ø 2,1x Stop 24mmx Comprimento 135mm
05.5.0439.121E	Broca Ø 2,1x Stop 24mmx Comprimento 140mm
05.5.0440.121E	Broca Ø 2,1x Stop 24mmx Comprimento 145mm
05.5.0441.121E	Broca Ø 2,1x Stop 24mmx Comprimento 150mm
05.5.0442.121E	Broca Ø 2,1x Stop 25mmx Comprimento 50mm
05.5.0443.121E	Broca Ø 2,1x Stop 25mmx Comprimento 55mm
05.5.0444.121E	Broca Ø 2,1x Stop 25mmx Comprimento 60mm
05.5.0445.121E	Broca Ø 2,1x Stop 25mmx Comprimento 65mm
05.5.0446.121E	Broca Ø 2,1x Stop 25mmx Comprimento 70mm
05.5.0447.121E	Broca Ø 2,1x Stop 25mmx Comprimento 75mm
05.5.0448.121E	Broca Ø 2,1x Stop 25mmx Comprimento 80mm
05.5.0449.121E	Broca Ø 2,1x Stop 25mmx Comprimento 85mm
05.5.0450.121E	Broca Ø 2,1x Stop 25mmx Comprimento 90mm
05.5.0451.121E	Broca Ø 2,1x Stop 25mmx Comprimento 95mm
05.5.0452.121E	Broca Ø 2,1x Stop 25mmx Comprimento 100mm
05.5.0453.121E	Broca Ø 2,1x Stop 25mmx Comprimento 105mm
05.5.0454.121E	Broca Ø 2,1x Stop 25mmx Comprimento 110mm
05.5.0455.121E	Broca Ø 2,1x Stop 25mmx Comprimento 115mm
05.5.0456.121E	Broca Ø 2,1x Stop 25mmx Comprimento 120mm
05.5.0457.121E	Broca Ø 2,1x Stop 25mmx Comprimento 125mm
05.5.0458.121E	Broca Ø 2,1x Stop 25mmx Comprimento 130mm
05.5.0459.121E	Broca Ø 2,1x Stop 25mmx Comprimento 135mm
05.5.0460.121E	Broca Ø 2,1x Stop 25mmx Comprimento 140mm
05.5.0461.121E	Broca Ø 2,1x Stop 25mmx Comprimento 145mm
05.5.0462.121E	Broca Ø 2,1x Stop 25mmx Comprimento 150mm
05.5.0463.121E	Broca Ø 2,1x Stop 30mmx Comprimento 50mm
05.5.0464.121E	Broca Ø 2,1x Stop 30mmx Comprimento 55mm
05.5.0465.121E	Broca Ø 2,1x Stop 30mmx Comprimento 60mm

05.5.0466.121E	Broca Ø 2,1x Stop 30mmx Comprimento 65mm
05.5.0467.121E	Broca Ø 2,1x Stop 30mmx Comprimento 70mm
05.5.0468.121E	Broca Ø 2,1x Stop 30mmx Comprimento 75mm
05.5.0469.121E	Broca Ø 2,1x Stop 30mmx Comprimento 80mm
05.5.0470.121E	Broca Ø 2,1x Stop 30mmx Comprimento 85mm
05.5.0471.121E	Broca Ø 2,1x Stop 30mmx Comprimento 90mm
05.5.0472.121E	Broca Ø 2,1x Stop 30mmx Comprimento 95mm
05.5.0473.121E	Broca Ø 2,1x Stop 30mmx Comprimento 100mm
05.5.0474.121E	Broca Ø 2,1x Stop 30mmx Comprimento 105mm
05.5.0475.121E	Broca Ø 2,1x Stop 30mmx Comprimento 110mm
05.5.0476.121E	Broca Ø 2,1x Stop 30mmx Comprimento 115mm
05.5.0477.121E	Broca Ø 2,1x Stop 30mmx Comprimento 120mm
05.5.0478.121E	Broca Ø 2,1x Stop 30mmx Comprimento 125mm
05.5.0479.121E	Broca Ø 2,1x Stop 30mmx Comprimento 130mm
05.5.0480.121E	Broca Ø 2,1x Stop 30mmx Comprimento 135mm
05.5.0481.121E	Broca Ø 2,1x Stop 30mmx Comprimento 140mm
05.5.0482.121E	Broca Ø 2,1x Stop 30mmx Comprimento 145mm
05.5.0483.121E	Broca Ø 2,1x Stop 30mmx Comprimento 150mm
05.5.0484.121E	Broca Ø 2,1x Stop 35mmx Comprimento 50mm
05.5.0485.121E	Broca Ø 2,1x Stop 35mmx Comprimento 55mm
05.5.0486.121E	Broca Ø 2,1x Stop 35mmx Comprimento 60mm
05.5.0487.121E	Broca Ø 2,1x Stop 35mmx Comprimento 65mm
05.5.0488.121E	Broca Ø 2,1x Stop 35mmx Comprimento 70mm
05.5.0489.121E	Broca Ø 2,1x Stop 35mmx Comprimento 75mm
05.5.0490.121E	Broca Ø 2,1x Stop 35mmx Comprimento 80mm
05.5.0491.121E	Broca Ø 2,1x Stop 35mmx Comprimento 85mm
05.5.0492.121E	Broca Ø 2,1x Stop 35mmx Comprimento 90mm
05.5.0493.121E	Broca Ø 2,1x Stop 35mmx Comprimento 95mm
05.5.0494.121E	Broca Ø 2,1x Stop 35mmx Comprimento 100mm
05.5.0495.121E	Broca Ø 2,1x Stop 35mmx Comprimento 105mm
05.5.0496.121E	Broca Ø 2,1x Stop 35mmx Comprimento 110mm
05.5.0497.121E	Broca Ø 2,1x Stop 35mmx Comprimento 115mm
05.5.0498.121E	Broca Ø 2,1x Stop 35mmx Comprimento 120mm
05.5.0499.121E	Broca Ø 2,1x Stop 35mmx Comprimento 125mm
05.5.0500.121E	Broca Ø 2,1x Stop 35mmx Comprimento 130mm
05.5.0501.121E	Broca Ø 2,1x Stop 35mmx Comprimento 135mm
05.5.0502.121E	Broca Ø 2,1x Stop 35mmx Comprimento 140mm
05.5.0503.121E	Broca Ø 2,1x Stop 35mmx Comprimento 145mm
05.5.0504.121E	Broca Ø 2,1x Stop 35mmx Comprimento 150mm

IMPORTANTE

Os instrumentais BIOMETAL devem ser conectadas motor cirúrgico, fornecendo ao cirurgião as respectivas modalidades de uso.

Durante seu uso em procedimento único, o desgaste habitual, o exercício de forças excessivas e o uso fora do especificado, podem prejudicar a evolução do procedimento e danos aos instrumentos e/ou paciente.

As instrumentais cortantes são componentes metálicos que estão sujeitas a importantes solicitações mecânicas durante o uso continuado por prazo variável e indefinido, sendo necessário a inspeção e revisão das condições de uso durante o procedimento cirúrgico.

Esteja alerta a caso ocorra variações no desempenho pretendido, como: perda de precisão, instabilidade ou falta de corte.

Caso ocorra estas negativas variações no desempenho, o material deve imediatamente ser substituído por outra unidade.

Assegure-se de que seu instrumental cirúrgico esteja devidamente encaixado e travado na peça de mão antes de acionar o motor cirúrgico.

ADVERTÊNCIAS E PRECAUÇÕES

Os instrumentais BIOMETAL somente devem ser utilizadas por profissionais capacitados e após uma análise detalhada do procedimento cirúrgico a ser adotado e da leitura completa destas INSTRUÇÕES DE USO do produto.

Sempre verificar se a esterilização se encontra dentro do prazo de validade, conforme data de validade fixada no rótulo do produto. Caso esta data não esteja dentro do período de validade, não utilizar.

A esterilização será garantida somente se a embalagem não estiver danificada, em caso de danos não utilizar o produto.

A vida útil dos instrumentais está diretamente relacionada ao número de vezes que são utilizados e da maneira como são conservados, para isso são necessários alguns cuidados conforme descritos abaixo:

- ✓ *Os instrumentos cirúrgicos devem ser conectados motor cirúrgico, fornecendo ao cirurgião as respectivas modalidades de uso – por isso realize o planejamento prévio antes do procedimento cirúrgico.*
- ✓ *Os instrumentos cirúrgicos devem ser utilizados especificamente para os fins propostos;*
- ✓ *Os instrumentos cirúrgicos estão sujeitos a desgaste ao longo de sua utilização normal, portanto o uso inapropriado, abuso ou força excessiva dos instrumentos pode causar a quebra deles em procedimentos;*
- ✓ *O instrumental deve ser manipulado cuidadosamente, em pequenas porções, evitando-se batidas ou quedas;*

- ✓ O uso de irrigação durante a operação reduzirá a possibilidade à necrose térmica;
- ✓ O fabricante se isenta de responsabilidade por danos causados pelo uso incorreto ou inadequado do material;
- ✓ Deve-se evitar inclinar ou efetuar movimento de avalanche durante a cirurgia;
- ✓ Manipule com cuidado;
- ✓ A má seleção, colocação, posicionamento e acoplamento dos instrumentos cirúrgicos em motores cirúrgicos podem acarretar resultados indesejados. O cirurgião responsável deve se familiarizar com o produto, sua técnica de manuseio em cirurgia antes de sua utilização.
- ✓ Este é um produto não biocompatível, não absorvente, que tem como finalidade auxiliar procedimentos e/ou a implantação cirúrgica de implantes;
- ✓ Para manuseio e descarte dos instrumentos após uso em cirurgia, seguir os procedimentos legais locais do país para manuseio de produtos potencialmente contaminantes.

MANIPULAÇÃO

Manipular os instrumentos cirúrgicos esterilizados de acordo com procedimentos assépticos adequados.

Proteger os itens de contato acidental, quebras, sobreposição com outros materiais – visando evitar danos a superfície que podem gerar desvios no momento do uso.

Inspeccionar visualmente todo instrumento cirúrgico antes da sua utilização clínica.

TRANSPORTE

Os Instrumentais devem ser transportados com cuidado, evitando-se danos ao sistema de embalagem estéril, quedas e atritos que possam causar a perda da condição estéril, e até defeitos no acabamento da superfície da peça.

Manter o instrumental sempre nas embalagens originais, até o momento de sua utilização, sob a responsabilidade da equipe médico/hospitalar designada para este fim, lembrando-se que essa deve sempre se atentar para a integridade da embalagem.

Data de fabricação, prazo de validade e nº do lote: vide rótulo.

ARMAZENAMENTO

Para o armazenamento, recomenda-se local seco e arejado, sem exposição à incidência de luz, à umidade ou a substâncias contaminantes.

Os instrumentos cirúrgicos não devem ser armazenados diretamente ao chão (altura mínima = 20 cm), assim, recomenda-se a utilização de prateleiras.

O armazenamento em prateleiras altas e próximas de lâmpadas não deve ser realizado (para não ressecar a embalagem ou apagar o rótulo), além de não poderem ser armazenados em áreas onde sejam utilizadas substâncias contaminantes como inseticidas, pesticidas ou materiais de limpeza.

Cuidados no recebimento, estocagem e transporte devem ser adotados em conjunto com as Boas Práticas de Armazenamento e Distribuição de Produtos Médicos estabelecidos pela ANVISA.

INSTRUÇÕES DE USO

Os instrumentais BIOMETAL são material para saúde destinados ao uso temporário.

Remover o instrumental de sua embalagem somente no momento da sua utilização clínica e manusear em campos estéreis para assegurar a esterilidade.

Após abrir a embalagem, inserir o INSTRUMENTAL CIRÚRGICO somente em peça de mão que atendam aos requisitos das normas aplicáveis para produtos para saúde.

Durante o uso, movimentar o INSTRUMENTAL CIRÚRGICO continuamente de modo a evitar o aquecimento localizado e consequente quebra e riscos ao paciente.

Não há contra indicações e nem efeitos colaterais.

Pode ocorrer quebra do instrumental devido ao excesso de força ou rotação empregada, assim como partículas metálicas podem se soltar. Em caso de reações alérgicas, o paciente deve passar por acompanhamento clínico periódico para verificar possíveis efeitos adversos.

RESTRIÇÕES DE USO

Não utilizar o INSTRUMENTAL CIRÚRGICO se a embalagem estiver aberta ou danificada.

Não utilizar o INSTRUMENTAL CIRÚRGICO se houver sinais de danos.

O produto somente deve ser utilizado por profissional habilitado, recomenda-se médico cirurgião.

Quaisquer complicações devidas a indicações erradas, a uma seleção incorreta da técnica operatória, assim como complicações devidas a limites do método de tratamento aplicado ou a falta de assepsia clínica não são responsabilidade desta fabricante.

DESCARTE DO PRODUTO

A BIOMETAL recomenda o descarte após o uso.

Para manuseio dos instrumentos após uso em cirurgia, seguir os procedimentos legais vigentes do país para manuseio de produtos potencialmente contaminantes.

RASTREABILIDADE

Para garantir a rastreabilidade, recomenda-se ao cirurgião ou sua Equipe o registro no prontuário do paciente das informações acerca do INSTRUMENTAL CIRÚRGICO estéril utilizado.

As informações necessárias para rastreabilidade são as relativas ao produto utilizado, cirurgia e paciente, conforme abaixo:

- ✓ Nome do paciente;
- ✓ Nome do cirurgião;
- ✓ Nome do Hospital;
- ✓ Nome do Fabricante;
- ✓ Nome do Fornecedor;
- ✓ Data da cirurgia;
- ✓ Código do produto;
- ✓ Número de lote do produto;
- ✓ Quantidades utilizadas;
- ✓ Registro do produto na ANVISA

As informações necessárias para a rastreabilidade do produto, a seguir, encontram-se gravadas na peça e/ou podem ser obtidas através da ETIQUETA DE RASTREABILIDADE e RÓTULO FIXADO NA EMBALAGEM DO PRODUTO.

A BIOMETAL disponibiliza de forma adicional 3 VIAS da ETIQUETA DE RASTREABILIDADES autoadesiva e avulsas, com as mesmas informações do rótulo aprovado pela ANVISA, destinados à rastreabilidade do produto.

INSTRUÇÃO DE USO

Instrumental não articulado cortante – CADASTRO ANVISA nº 80482490003



Os rótulos autoadesivos devem ser devidamente distribuídos de forma a garantir a rastreabilidade do produto.

Caso seja observado algum EVENTO ADVERSO ou haja necessidade de realização de QUEIXA TÉCNICA deve-se proceder à notificação no Sistema Nacional de Notificações para a Vigilância Sanitária – NOTIVISA, que pode ser encontrado no site da Agência Nacional de Vigilância Sanitária – ANVISA no endereço www.anvisa.gov.br, link NOTIVISA.

Nota: Se algum instrumento objeto de desvios de qualidade necessitar e for devolvido ao fabricante BIOMETAL e/ou as Empresas Distribuidoras, este item deve ser descontaminado e esterilizado, e ser acompanhado com a evidência relevante destes processos.

FABRICADO E DISTRIBUÍDO POR:

BIOMETAL INDUSTRIA E COMERCIO EIRELI

CNPJ: 07.165.928/0001-75

Endereço: V LUIZ VARGA 2010 - GLEBA BEATRIZ

Cidade/UF: LIMEIRA/SP

CEP: 13.486-606

Fones: (19) 3495-2485 / (19) 3495 3296

RESPONSÁVEL TÉCNICO: RENATA VICELLI – CREA/SP 5063105420

CADASTRO ANVISA Nº: 80482490003

REVISÃO: 01